

The factory setting is $^{\circ}\mathbf{F}$. If you wish to change it to $^{\circ}\mathbf{C}$, follow the steps below:

- | | | | |
|---|---|---|---|
| 0 | 7 | 0 | 7 |
| 7 | 7 | 7 | 7 |
| 7 | 7 | 7 | 7 |

Why take temperature measurements at the skin surface over the temporal artery?

Why do I get a different reading with each scan?

Multiple scans in rapid succession will change the temperature of the skin and result in variability of the readings. It is important to wait a full minute before repeating the scan to allow the skin over the temporal artery to recover from the probe temperature of the probe head. The probe head is at room temperature and cooler temperature of the forehead will be correct 98% of the time when there is no hand (or a kiss on the forehead) will be correct 98% of the time when there is a fever, but wrong half of the time when the individual actually has a fever. The reason is that body heat is released or retained to keep our body temperature in the normal range. This contributes to the erroneous assessment as felt by the hand or lips. The core body temperature stays normal by the release of heat (like opening the windows when the house is too hot) or by retention of the heat (closing the windows when the house is too cold).

the windows when the house is too cold)

Display Messages

Sen	hi
-----	----

5	2
---	---

LG	H1.A
----	------

11. M	12. M
13. M	14. M

Err	CU. M
-----	-------

Blinking battery icon with 2 dashes: not enough energy in the battery to measure correct temperature. Replace battery.

Sweat will cause low temperature readings: If the individual is sweaty

Clinical Accuracy: Meets ASTM E1965-98 and EN60601-1 standards for electronic and radiation thermometers to the extent applicable to thermometers


Five and Three Year Warranty

misuse or abuse, lack of reasonable care, the affixing of any attachment not provided with the product or loss of parts or subjecting the product to any but the specified battery.* Use of unauthorized replacement parts will void this warranty.

Used and Recommended by Millions of Doctors and Nurses



“The Carapax Thermoplastic is for more conventional than measuring a body temperature orally – who can resist?”
— *Annex News*

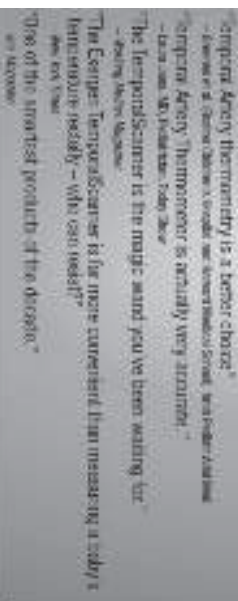

EXERGEN Corporation

medical@exergen.com, www.exergen.com

P/N 818688 Rev 3

CEC
This symbol on the product's
compliance means it is listed

Underwriters' Laboratories, Inc.

EXERGEN Corporation

Clinical studies,
multiple language
videos, & user
manuals:
exergien.com/s

