Product Name: Ellipse By Legxercise
Model No.:LEGX-013
Distributed By :IntelliBrands, LLC.
Customer service address:9990 NW 14th Street, Suite
109,Miami, FL 33172, USA

If you have any questions, contact Customer Service at
Support@intellibrands.tv Add:

Purchases made directly from Intellibrands LLC for Legxercise or Ellipse branded products, from our online sites or our call centers, include a 60 day money back guarantee. Products purchased pursuant to such offers may be returned within 60 days of receipt of the product for a full refund of the purchase price (not including any original processing and handling fees nor return processing fees). If you purchased through a Retailer or Amazon, the money back guarantee timeframe may vary. Please contact the Retailer for details on the timeframe for returns. Please contact customer service at Support@Intellibrands.tv for any questions pertaining to the money back guarantee, 1 year manufacturers warranty or any product specific questions. You can call us at 786-460-8394

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Instructions
Instrucciones


Ellipse by Legxercise is the new revolutionary elliptical machine that keeps your legs in movement while you are sitting, working at your office desk, watching TV, anywhere you can sit and pedal. Ellipse effortlessly keeps your legs moving, anytime, anywhere.

- Moves forward and reverse.
- 3 speed levels in each direction.

Speed one: 105 steps per minute

Speed two: 140 steps per minute


Speed three: 165 steps per minute

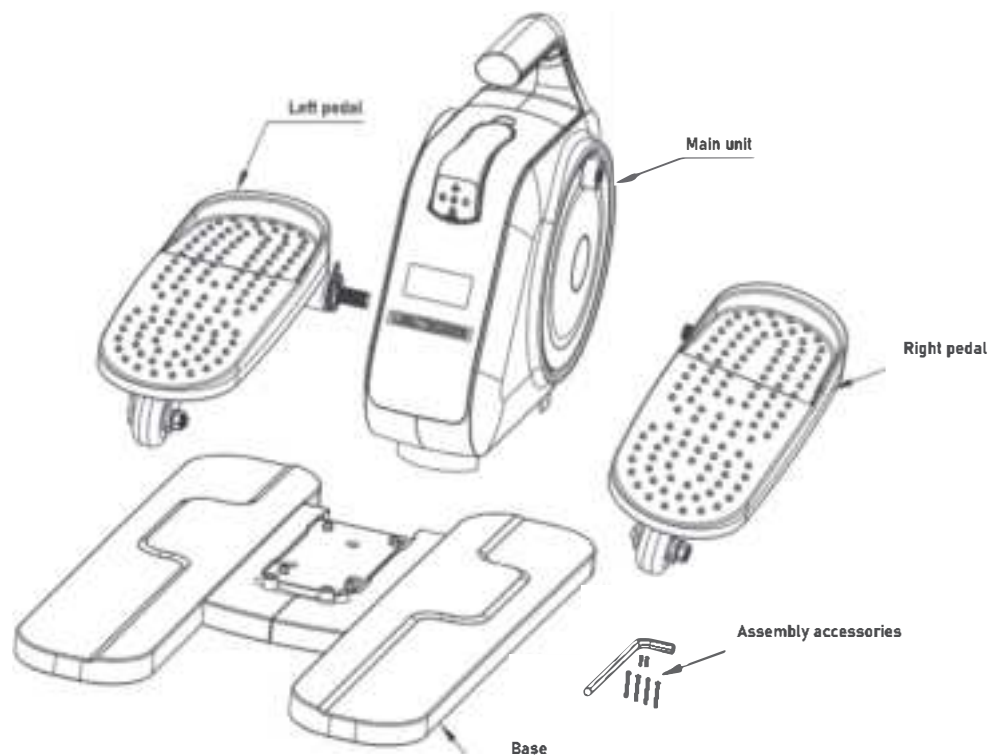Assembly Instruction


Step one:


Please take out everything within the Ellipse By Legxercise box. Make sure nothing is left inside. The assembly kit with 6 screws and allen wrench is taped to the top of the pedals.

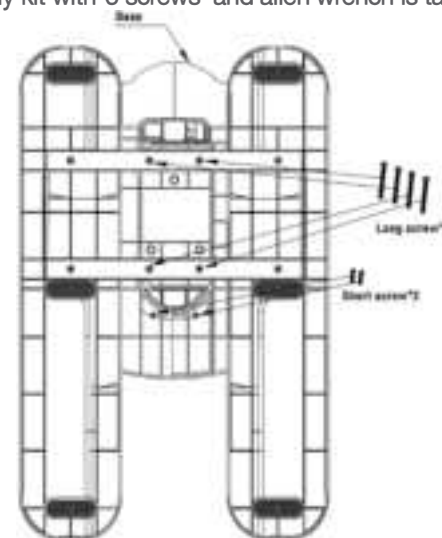Step two:

Secure the main unit and the base together by using the six screws provided in the assembly kit.

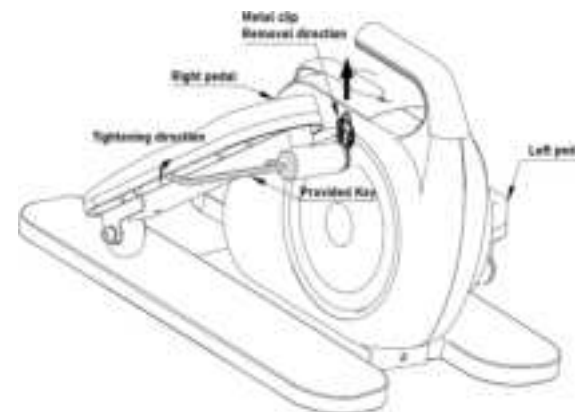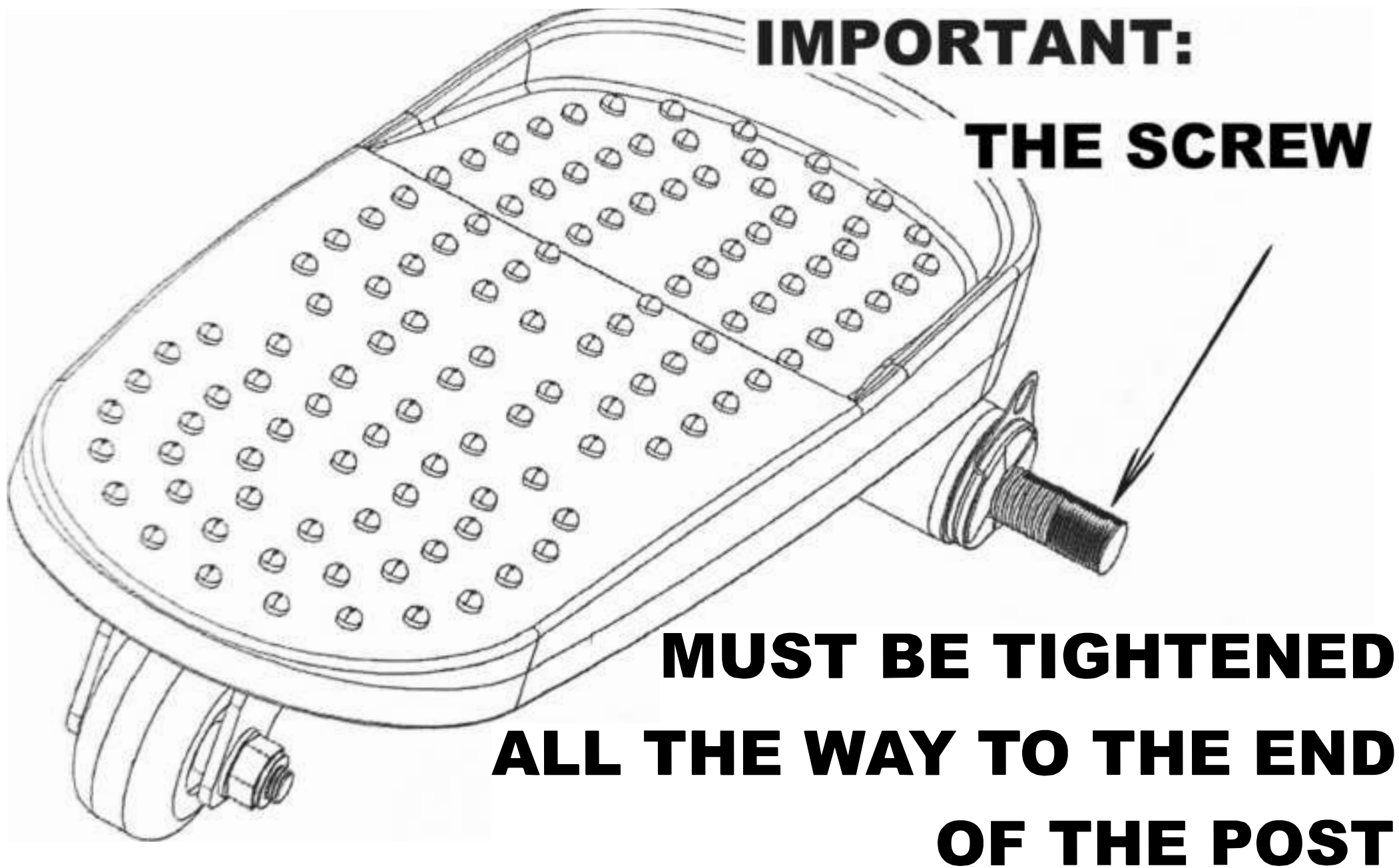Step three:

The pedals must be tightened so that the bolt is not visible. Use the Allen wrench to tighten the pedals onto the main unit. It is easiest to attach the pedals if they are lying upside down, flat on the floor, toward the front of the unit. If attached in this manner, you do not need to hold onto the pedal while tightening the bolt. Remove the metal clip from the pedal before the final tightening of the pedal and make sure the washers are aligned as you finish tightening.

Right Pedal


Left Pedal

Place the pedals flat (UPSIDE DOWN)

ATTENTION: REMOVE THE METAL CLIP BEFORE FINAL TIGHTENING AND MAKE SURE THE WASHERS STAY ALIGNED AS YOU FINISH TIGHTENING.

Next, find the remote control located in the upper part of the machine, 2 x AAA batteries are required, look for the lid compartment in the handle of the remote control, batteries are not included.

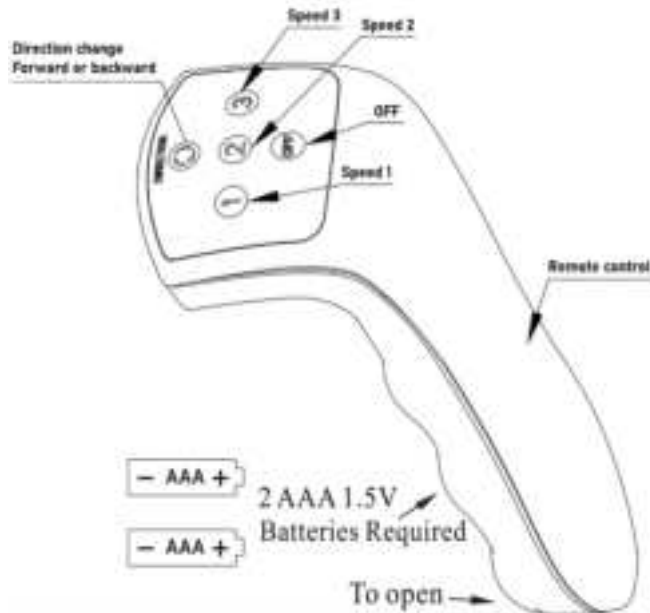How to Operate:

Connect the power adapter to the unit and plug into an outlet


The power adapter works with a 100V to 240V power source.

The remote has 3 speeds as well as the ability to change the direction (forward/backward) and Off.

Place your feet onto the pedals of the Ellipse, select your desired speed level using the remote control, sit back and relax.

You can choose to pedal backwards by clicking the Direction Change button on the remote.


Make sure you are seated in a manner that allows you to pedal comfortably and without interference from surrounding objects. You must be able to easily dismount the Ellipse without tripping over any other furniture

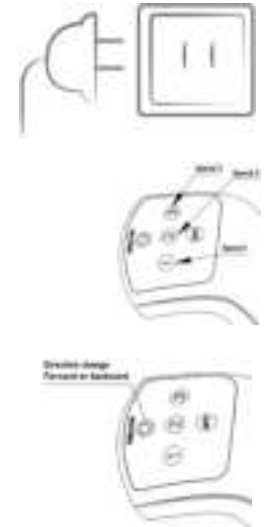During operation, the Ellipse will count each rotation of the foot plates, once you reach 9,999 steps, the unit will stop automatically. You can reset the Ellipse counter by choosing a speed, and pedaling. The step counter will reset to zero and begin counting again.

Manual Push Mode:



Ellipse will also count your steps during manual pedal mode operation (forward only). The step counter display will turn off once you stop pedaling. Make sure the unit is unplugged from a power source.

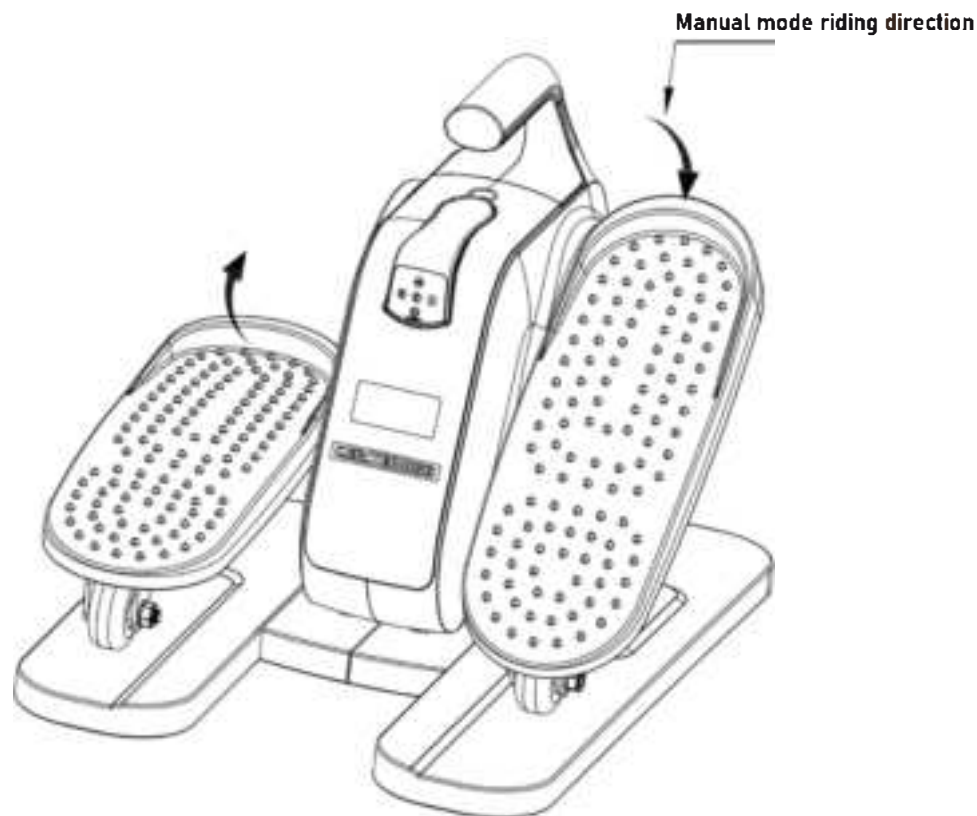For best results you can use both modes, Manual and assisted

WARNING:

PLEASE KEEP FINGERS AWAY FROM ROTATION MECHANISM AT ALL TIMES.
PLEASE UNPLUG IF NOT IN USE. WARNING:

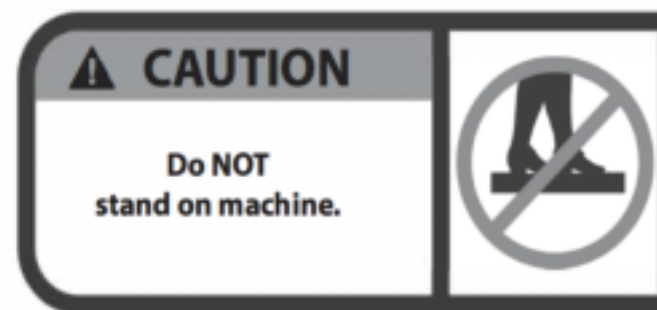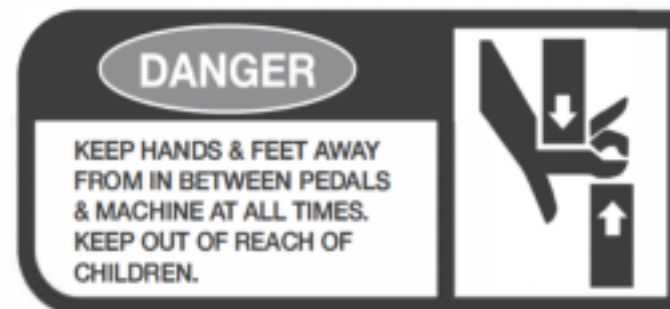DO NOT STAND on unit.

Ellipse by Legxercise is designed for seated use only. Do not stand on the machine. You may injure yourself or damage the rotation mechanism.

The Ellipse should be used indoors only.

The unit should not be exposed to water. Do not place it in water or near water. If liquids are spilled on it, turn it off and unplug it. Wipe it dry and allow to completely air dry before use.

The maximum weight of legs and feet should not exceed 100 Lbs.

Keep children and Pets away from the unit.

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If you have any questions, contact Customer Service at
Support@intellibrands.tv

Purchases made directly from Intellibrands LLC for Legxercise or Ellipse branded products, from our online sites or our call centers, include a 60 day money back guarantee. Products purchased pursuant to such offers may be returned within 60 days of receipt of the product for a full refund of the purchase price (not including any original processing and handling fees nor return processing fees). If you purchased through a Retailer or Amazon, the money back guarantee timeframe may vary. Please contact the Retailer for details on the timeframe for returns. Please contact customer service at Support@Intellibrands.tv for any questions pertaining to the money back guarantee, 1 year manufacturers warranty or any product specific questions. You can call us at 786-460-8394

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