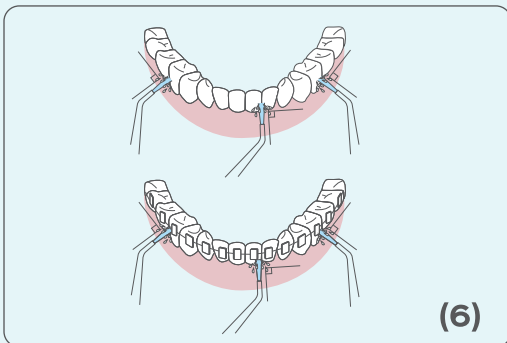
(6)

EN 6. Press "On/Off" button to start irrigating. (6)
1. Direct water flow at 90° to the teeth and gums
2. Direct water flow slowly along the teeth
3. Direct water flow along the gum line with nozzle tip at 90°
4. Clean between the teeth and around any braces, crowns or bridgework thoroughly

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
• After Using (8)

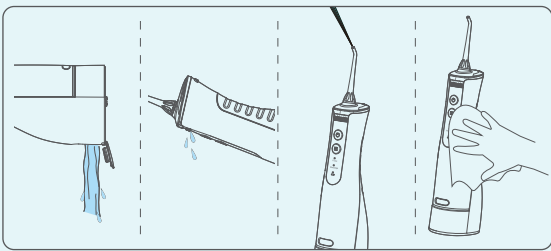EN Note: To prevent bacteria breeding, make sure there is no remaining water inside. If you're not going to use the appliance for a long time, be sure to wipe and dry it before storing.

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
CLEANING TIP INSTRUCTION

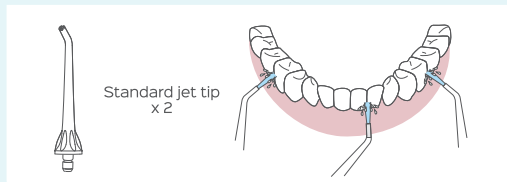
EN **Standard Jet Tip:** daily clean teeth, gum line, fissure, remove food residue, stain, plaque, etc.

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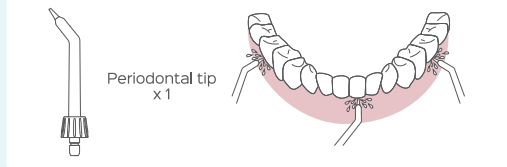EN **Periodontal Tip:** clean gingival sulcus, especially suitable for those who are suffering from periodontal diseases. When cleaning periodontal pocket, please use the lowest water pressure. Direct water flow along the gum line with nozzle tip at 45° and turn it on.

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
EN Note: It is recommended to replace cleaning tip every 6 months for better cleaning experience.

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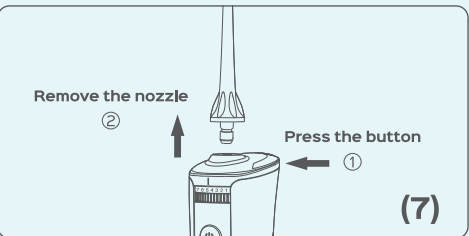
EN 7. Press "On/Off" button to turn it off, then press nozzle release button to replace other nozzle you need. (7)

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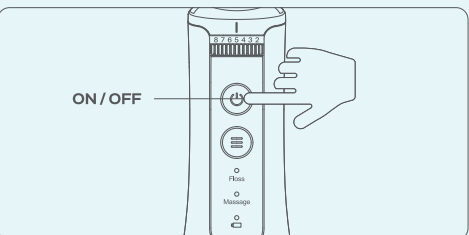EN 8. After irrigating is finished, press "On/Off" button again to turn off.

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
EN 1. Open water tank cap and dispose of water inside the tank.
2. Tilt and shake the main unit to drain the water from the speed control knob.
3. Press "On/Off" button to turn on and drain water inside the appliance.
4. Wipe the appliance with dry cloth.

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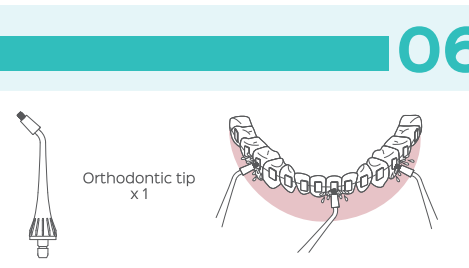
EN **Orthodontic Tip:** clean bridge and teeth with braces, suitable for those who wear braces or have dental implant.

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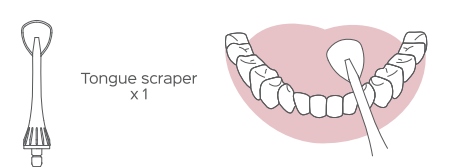EN **Tongue Scraper:** clean tongue surface, remove bad breath, food residue and thick tongue fur. When cleaning tongue, please use the soft mode. Lean forward, put scraper on your tongue and turn it on, then move gently from back to front. You can adjust water pressure according to your need.

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MAINTENANCE

EN Use only water or neutral cleaning agent to clean the appliance. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may cause damage to the appliance. Do not use water over 40°C. Keep the appliance away from high temperature and avoid direct sunlight.

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EN

- **Main unit**
 1. Wipe the dirt by cloth with soapy water.
 2. Do not submerge the main unit in water to clean.
 3. Dry after washing with water.
- **Water tank cap**
 1. Wipe off any grime and dirt with a cloth.
- **Nozzle**
 1. Wash with water and wipe with soft cloth.
 2. Replace the nozzle every 6 months.
- **Water suction tube**
 1. Wash with water and wipe with soft cloth.
 2. Do not bend, pull or twist the water suction tube.
- **Tank**
 1. Wash with water and wipe with soft cloth.
 2. Dry the tank when not using for a long period.

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08 TROUBLESHOOTING EN


Problem	Cause	Solution
Water does not come out	Negative pressure caused by not using for a long time.	Solution 1: Full fill the tank with water, invert the device and shake it to eliminate negative pressure.
	You've just purchased the appliance or it has been left unused for a period.	Solution 2: Install the nozzle tip on the main unit, suck out the air from nozzle tip with mouth to eliminate negative pressure.
It does not work	You've just purchased the appliance or it has been left unused for more than 3 months.	Pull out the nozzle, then flush nozzle hole on the body.
	The battery has reached the end of its service life.	Charge the appliance for at least 4 hours.
It only works for a few minutes even after charging	Charging time was insufficient.	Contact service center for support.
	There is dirt stuck in the nozzle.	Charge the unit for at least 4 hours.
Water pressure is very low	The nozzle is deformed.	Replace the nozzle.
	There is dirt stuck in the filter.	Clean the filter.

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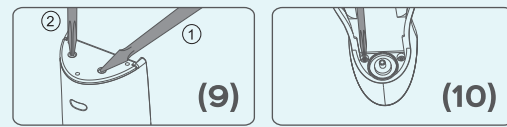EN 1. Remove 2pcs rubber feet on the bottom of the appliance. (9)
2. Use screwdriver counterclockwise to loosen bottom cap. (9)
3. Slide water tank down off the main unit and unscrew two screws beside nozzle suction tube. (10)
4. Open top cover and put inner part out of the main unit, then unscrew all screws of the inner part. (11)
5. Cut off battery wire and pull out batteries from the inner. (12)

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Warm Tip: If you have any question about recycling, please contact your local waste management office.

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