

SCAN QR CODE

To access the full
user manual



Or visit: www.hairmax.com/usermanual

You can also call customer service free of charge to request a hard copy

CUSTOMER SERVICE

Direct line +1 561-417-0200

United States and Canada (toll-free) 1-866-527-3726

Authorized Representative in Europe (Regulatory Affairs Only)



EMERGO EUROPE
Westervoortsedijk 60
6827 AT Arnhem
The Netherlands

Tel: +31 (0) 70 345-8570 Fax: +31 70 346-7299

Australia Sponsor



EMERGO AUSTRALIA
201 Sussex Street
Darling Park, Tower II, Level 20
Sydney NSW 2000
Tel: +61 2 9006 1662

Tel: +61 2 9006 1662

Brazil Sponsor



EMERGO BRAZIL
Avenida Francisco Matarazzo 1.752, Salas 502/503
Água Branca
São Paulo - SP CEP - 05001-200
Tel: +55 11 4800 9226

Tel: +55 11 4800 9226

United Kingdom Sponsor



EMERGO UNITED KINGDOM
Emergo Consulting (UK) Limited c/o Cr360 - UL International
Compass House, Vision Park Histon
Cambridge CB24 9BZ
Tel: +44 (0) 1223 214 044

Tel: +44 (0) 1223 214 044



LEXINGTON INTERNATIONAL, LLC

1040 Holland Drive | Boca Raton, Florida | 33487 | USA

www.hairmax.com

©2024 Lexington Intl., LLC. All Rights Reserved

LQH012CLEDXXX 00 0422



hairmax®
TRANSFORMATIVE HAIR REGROWTH

QUICK START GUIDE

ULTIMA 12 CLASSIC LASERCOMB

Hair Growth **Laser Light Device**

QUICK START GUIDE

Ultima 12 Classic LaserComb




3 times
per week

FOR MAXIMUM BENEFIT


- Scalp should be freshly washed (hair can be wet, dry or damp)
- Use the laser device prior to applying styling products (these can build up and impede the laser light)
- Regularly clean the device with a lint-free cloth. Lightly wipe the surfaces with 70% Isopropyl alcohol for 2 minutes, allowing it to dry before starting the treatment. Do not use harsh detergents or cleaning products on any part of the device.

CHARGING OPTIONS

You must charge the device for 2-3 hours before initial use.

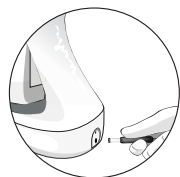


OPTION 1 Charge the laser device using the charging stand and power cord. Prior to charging the laser device, make sure the device is turned off and the teeth are facing away from you. Plug the power cord into the charging cradle/stand and the other end into an electrical wall outlet.

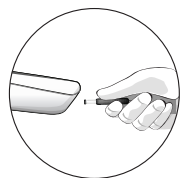

OPTION 2 Charge the laser device by placing the power cord directly into the laser device and the other end into an electrical wall outlet.

The indicator light on the power button will flash during the charging process. Once charging is complete, the light will stop flashing.

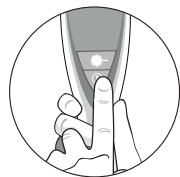

NOTE – Do not use the LaserComb while charging.

Turn on the laser device by pressing the power button (lower button) and hold until the power button light turns on. Next, press the laser button (upper button) until it lights up. That means the device is ready for use.

If the laser device needs to be recharged during your treatment, allow it to fully charge before continuing.

Prior to treatment, test your skin for sensitivity by shining the light on the inside of your forearm for 2 minutes. If you see any reaction, do not use your LaserComb and contact customer service.

TREATMENT METHOD



To start your treatment, place the laser device flat on your scalp so that both rows of teeth are touching the scalp.



Start at your hairline. Move the laser device from front to back of the scalp in a slow combing motion, pausing every 1/2 inch (1 cm) for 4 seconds. Do not lift the device off your scalp while combing from front to back.



Next, lift the laser device off your scalp and start at a new position at the hair line. Continue in this manner from front to back, and then side to side until you cover your entire scalp with multiple passes.

When the device is on it will beep and vibrate every four (4) seconds. This means that it is time to move it to the next area of the scalp. Repeat these passes for the total treatment time of 8 minutes.

After treatment, turn off the device by pressing and holding the power button (lower button) until the indicator light turns off. Store the device in the charging cradle for future use. You may leave the power cord plugged into the charging cradle at all times.

NOTE: The Ultima 12 Classic LaserComb will turn off automatically after 12 minutes.

Based on the durability of the lasers, electronic components and charging cradle, the laser device is designed to last for 5 years of regular use.

For detailed information about troubleshooting, proper use, warnings and cautions, access the full user manual by scanning the QR code located on the back panel. Please read the entire user manual before operating.

The Hairmax laser device has a limited warranty against defects in materials or workmanship from Lexington International, LLC. Items that are damaged or abused are not covered under the warranty and will be charged at their regular repair rates. Any attempt to self-service the device, unless instructed by LEXINGTON INTERNATIONAL CUSTOMER SERVICE, will void your warranty. Please refer to our website for complete warranty terms.

To validate your warranty, register online at
www.hairmax.com/warranty

Laser Light - Avoid Direct Eye Exposure

Hairmax® laser devices are indicated to treat Androgenetic Alopecia, and promote hair growth in males who have Norwood- Hamilton Classifications of IIa to V and in females who have Ludwig (Savin) I-4, II-1, II-2, or frontal patterns of hair loss and who both have Fitzpatrick Skin Types I to IV.