

## User Manual



*Enjoy it.*

**EYE MASSAGER**

HYH-03

## Catalogue


Note of caution.....2	Maintenance .....10
Taboo people .....2	<b>5</b> Cleaning method .....10
<b>1</b> Use precautions .....3	Abandonedmode .....10
Precautions for use of the site ...4	
Maintenance precautions .....4	<b>6</b> Fault diagnosis .....11
<b>2</b> Product function.....5	<b>7</b> Main technical parameters.....11
<b>3</b> Product component.....5	<b>8</b> Environmental statement .....12
Usage method.....6	<b>9</b> Packing list.....12
<b>4</b> Power switch.....6	Note: The product pictures in this manual are for reference only. All contents have been reviewed and proofread, if there are printing errors and technical improvements, will be included in the new manual!
Music keys .....8	
Indicator status .....8	
Charging method .....9	


**This product is a household product, do not use for commercial purposes**

## Note of caution

### Description

To avoid personal injury or property damage , observe the following.  
 • Use the following illustrations to indicate what you want to comply with.

 Express prohibition

 Indicates that you need to follow the instructions

### Taboo group

 The following groups/situations should be used with caution

If you are being treated in a medical facility, or if you have any of the following conditions, be sure to consult a doctor or physical therapy before using the doctor consulted.

- Pregnant women or women giving birth
- People with body temperature above 38 , People with loss of skin temperature sensation
- People with hypothermia burns
- People with abnormal body sensation
- For contact lens wearers, please remove contact lenses before use
- Keep this product away from water sources, fire sources and corrosive environments
- For use in specific environments, please comply with applicable regulations
- Sensitive skin, facial allergic dermatitis, facial parts with cosmetic plastic materials  
Please use this product under the guidance of your physician
- This product is not a medical device and is not intended for diagnostic, therapeutic or other commercial use
- Please charge only the charging equipment configured by our company

- Please pay attention to safe use of electricity
- The appliance has a hot surface, and people who are not sensitive to heat must pay attention when using it

### **The following people/situations are prohibited**







Doctors advise people not to exercise or stretch






For example: thrombosis, severe aneurysms, acute venous tumors, various dermatitis, various skin infections,

Because it is easy to cause deterioration, it is prohibited to use .






- Do not use this product while driving
- Do not use this product when drunk or unwell
- It should not be used in patients with eye surgery and heart disease
- Do not use when there are symptoms of glaucoma, cataract, retinal detachment, pink eye and other eye diseases
- Do not use it for pregnant women and menstruating women
- Do not massage other parts of the body than the eyes
- Do not wear this product during sleep
- It is prohibited for children under 6 years old, and for young people between 6 and 14 years old, please use it under the guidance of an adult guardian

## **Precautions for use**







-  Wear glasses, hair accessories and other things, should be taken off before use before use.
-  If you do not use it again for a long time, be sure to confirm whether the switch is normal in advance.
-  Please close your eyes when using the product. Blurred vision will occur when removing the product after use. Close It can be relieved in 1-3 minutes.
-  Power off immediately if any abnormality occurs .
-  This product is not intended for diagnostic or therapeutic use .
-  When the product does not have waterproof performance, be careful not to let water enter the product.

-  Do not deliberately put too much effort on the product.
-  Do not use non-specified voltage power supply, easy to damage the product circuit.
-  Do not drop the product, nor strong vibration, easy to lead to product damage .
-  Never break down the product, transform or repair the product yourself.
-  Do not pull the product hard, and do not use a knife or sharp objects to stab the product.

## **Precautions for use of the site**

-  Do not go near the stove, or the electric heating is used inside and on electric heating appliances such as electric blankets, easy cause fire or malfunction.
-  Do not cover the product with a towel or cloth.
-  Do not use this product near bathtubs, showers, washbasins or other high wet utensils.
-  Do not use in bathrooms or places with high humidity .
-  Do not use in places above room temperature 35° C or in direct sunlight.



## **Maintenance precautions**

-  The product is turned off during maintenance.
-  Be sure to unplug the charger when not in use for a long time.
-  Do not drop the body Do not use two people at the same time.
-  Do not drizzle or wash the product.
-  Please do not use volatile oil, diluent, strong alkaline detergent , bleach and other cleaning agents.
-  When not in use for a long time, the power should be turned off and kept charged once every 1 month.

## Product function

- Two buttons (on key, music key) + three indicators (on key, music key, status indicator)
- Constant temperature compress under eye contact
- Bluetooth Music + Bone conduction speakers (2 speakers)
- 4 massage modes
- Voice broadcast
- Built-in HD ambient white noise
- 15 minutes intelligent timed shutdown

## Product component

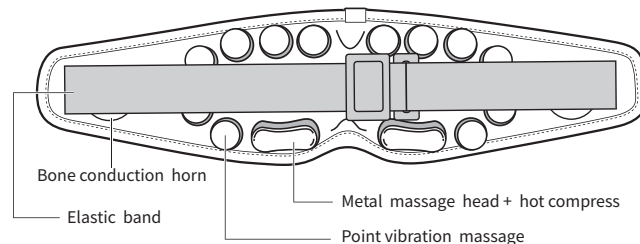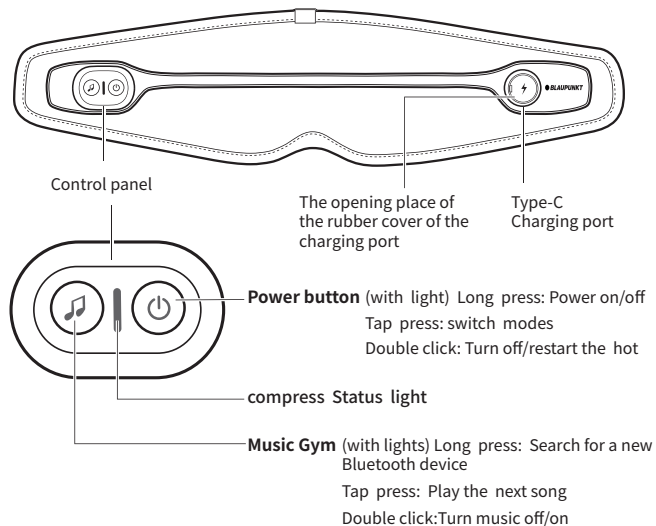

## Usage method

### 1. Power switch

#### 1. On/standby:

- (1) Long press the power button to start/standby, and the voice will announce "Welcome to use Lan bao Eye massage" "/" message "End", the power button display light is displayed according to the current operating mode.
- (2) After the boot, the default into the "vitality mode, hot compress open, automatically play built-in music, voice broadcast "Vitality Mode". Press the power button to switch the massage combination mode.

### 2. Switch music:

Press the music button and play a song.

### 3. Turn music off/on:

Double click the music button to turn off/turn on music, and the voice will announce "Turn off music" "/" Play music ".

### 4. Default Bluetooth connection

After powering on, the background automatically searches for the Bluetooth device (the search time is 1 minute, and the indicator is over (stro bosco pic blue), automatically connect the paired Bluetooth device, after the pairing is successful, the indicator will change from blue stro bosco pic to cyan breathing light, the voice will announce "Bluetooth connected", and play the audio of the paired device.

If the pairing is not successful within 1 minute, the blue strobe turns to cyan steady on, and the built-in music continues to play. No voice broadcast.

## Four massage combination modes

### 【 Vitality mode 】

Vibratory massage technique 1+ heating low temperature (38-40 °C), voice broadcast "vitality mode", Background music "Waves" single loop will automatically run for 15 minutes when switched to this mode, After 15 minutes, it will enter standby mode and the voice will announce "Massage is over".

### 【 Eye protection mode 】

Vibration massage technique 2+ heating high temperature (41-43 °C), voice broadcast "eye protection mode", Background music "Bird Call" single loop will automatically run for 15 minutes when switched to this mode, After 15 minutes, it will enter standby mode and the voice will announce "Massage is over".

### 【 Sleep mode 】

Vibration massage 3+ heating low temperature (38-40°C), voice broadcast "sleep mode", Background music "River 1" single loop, switching to this mode will automatically run for 15 minutes. After 15 minutes, the system enters the standby state. No voice prompt is displayed.

### 【 Music mode 】

Enjoy music, voice broadcast "music mode", background music "River 2" single cycle After switching to this mode, it will automatically run for 15 minutes, and enter the standby state after 15 minutes. Wait There is no voice prompt.

## 2. Music keys

1、Built- in music mode (turn on music cyan light , turn off music without light) :

Tap to —— play the next song

Double click —— turn music off/on (voice announcement "Turn music off" / " Play Music ")

Long press —— search for a new Bluetooth device (music stops playing , music light flashes blue , voice announcement "Search for new Bluetooth devices ") .

2、Bluetooth music mode (music light is cyan breathing light) :

Bluetooth Name: HB

Short press —— play song

Double click —— turn music off/on (voice announcement "Turn music off"/ " Play Music ")

Long press —— disconnect the currently connected Bluetooth device , search for a new Bluetooth device (music stops playing , sound Music lights flashing blue, voice announcement "Search for new Bluetooth device")

3、The new Bluetooth device search time is 1 minute, long press again to stop searching, back to the previous music mode such as pairing failure, will play the built-in music, the indicator light from the blue strobe to the cyan light such as pairing success, will play the paired device audio, the indicator light from the blue strobe to the cyan breathing light, voice announcements "Bluetooth is connected".

## 3. Indicator status

### 1. Power lamp

Blue breathing light :Vitality model

Green breathing light :Eye protection mode

Cyan breathing light :Sleep pattern

Blue light :Musical pattern

## 2. Music lights


Blue long light: Bluetooth music mode is not connected

Cyan breath light: Bluetooth music mode

Blue cheek flashing light: Bluetooth pairing

## 4. Charging method

1. Connect the USB charging cable to the Type-C charging port.
2. The indicator light shows an orange breathing state during charging, and the indicator light is blue after charging is completed.
3. During the charging process, the product cannot be started to work.
4. The host adopts built-in rechargeable battery, the charging time is about 3.5 hours.
5. In the process of use, the voice prompts the voice alarm of low battery. Please suspend the use until it is fully charged.
6. Support a charging mode computer USB/mobile power /5V -1 A Power adapter.

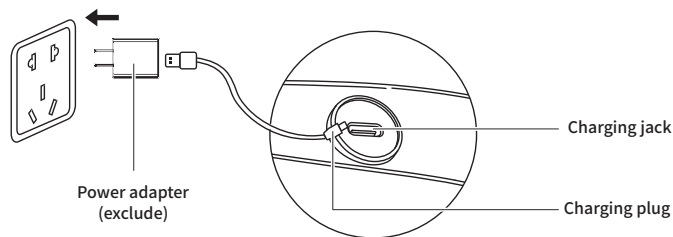

Note: This machine is not equipped with a power adapter (charging), please choose a qualified USB output port Power adapter with DC 5V voltage and 1000mA output current.

## Maintenance

### Cleaning method

1. The surface dirt can be removed with a clean rag
2. Wet a cloth with a neutral detergent for deep dirt and gently wipe it off
3. Do not let water enter the product;
4. Let dry naturally

- ❗ Do not store in direct sunlight or humid areas after use. If long term If not in use, after removing the dirt, put it in a storage bag and do not catch dust.
- ❗ When maintaining, please remove the charging cable first, but do not use wet hands to remove the wire.

### Abandonment mode

The lithium-ion battery used in the product is a valuable resource that can be recycled and recycled. Please take out the lithium-ion battery when it is used and discarded, and hand it to a professional institution that can recycle and recycle the rechargeable battery for processing.

### Battery removal steps:

1. Loosen the screws on the circuit board and take out the circuit board
2. Cut the battery cable and wrap the cable with insulation tape

- ❗ The battery cannot be replaced by the user;
- ❗ The battery must be removed from the appliance before discarding it;
- ❗ The appliance must be powered off when the battery is removed;
- ❗ Batteries should be disposed of safely.

## Fault diagnosis

This chapter lists the most common questions about this product.

If you can't find what you see in it, answers to your questions, please contact the service hotline or an authorized service center in your area.

Fault description	Possible causes	Basic exclusion method
Unable to boot	The lithium battery is low	Fully charge and then turn it on
Bluetooth playback music sound is low	Mobile phone playback volume is too low	Adjust the phone playback volume
No built-in music	The product is connected to Bluetooth	Disconnect Bluetooth device
No music	The music is off or connected Bluetooth's not playing music	Turn on the music

Attention:

If the fault cannot be removed or the product still cannot operate normally, contact the after-sales service personnel Line disassembly or attempt repair.

## Main technical parameters

Product name: EYE MASSAGER	Battery capacity: 1200mAh
Product model: HYH-03	Interface Type: Type-c
Rated voltage: 3.7V ---	Product net weight: 242g
Rated power: 4W	Product size: 310x82x90mm
Rated working time: 15 minutes	Executive standard :GB4706.1-2005
Charging voltage: 5V --- 1A	GB4706.10-2008

## The name and content of harmful substances in the product

1. This product meets the requirements of the "Management Method for the Restricted Use of Harmful Substances in Electrical and Electronic Products";

2. Within the environmental protection use period, consumers will not leak harmful substances during normal use; Precipitation and other issues affecting consumer health can be used with confidence;

3. The environmental protection life of the company's products is 10 years, only under the normal circumstances described in this manual; The "Environmental Life" is valid only when using this product.

Part name	The name and content of harmful substances in the product					
	(Pb)	(Hg)	(Cd)	(Cr(VI))	(PBB)	(PBDE)
Plastic component	○	○	○	○	○	○
Circuit board (electronic) components	×	○	○	○	○	○
Power cord	○	○	○	○	○	○
Hardware components	×	○	○	○	○	○
Battery	○	○	○	○	○	○

This form is prepared in accordance with SJ/T 11364.

○: It means that the content of the hazardous substance in all homogeneous materials of the component is below the limit stipulated in GB/T26572.  
The content of the hazardous substance in all homogeneous materials of the part is below the limit specified in GB/T26572.

×: Indicates that the content of the hazardous substance in at least one homogeneous material of the component exceeds the limit requirements stipulated in GB/T26572. Limit requirement.

(Note: Some of the components marked with "X" above contain hazardous substances exceeding the standard due to the limitations of the current technical level of the industry and the temporary inability to achieve substitution or reduction.)



## Packing list

Host machine	1×	Type-C Charging cable	1×
		User Manual	1×
		Certificate of Conformity	1×

## **FCC Requirement**

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.