



7.c

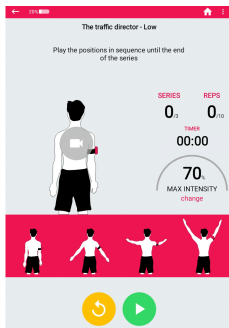



**Start** the exercise and move as the avatar

7.d

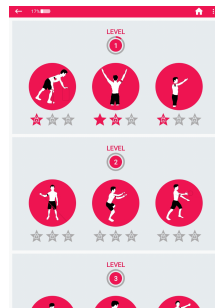


**Continue** your rehabilitation program

8

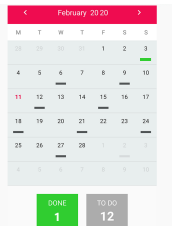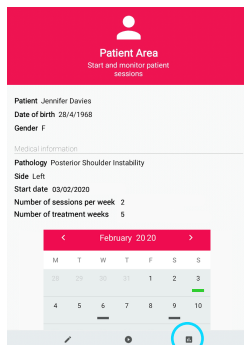


View the updated session calendar and your **history**

9

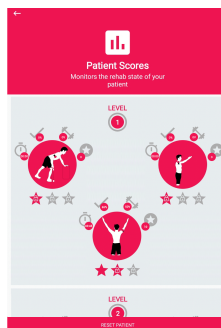

View your **performance**

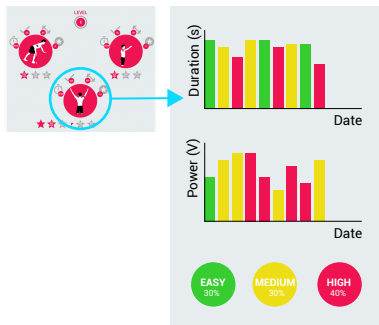

View **details** for each exercise

## **WARNING**

Make sure the device battery is fully charged before each session.  
MySPM@Home requires internet connection and active bluetooth.  
MySPM@Home requires confirmation of location and storage permissions.

**MySPM@Home**  
Shoulder Pacemaker™ App

**REF** **SPM01A2**



NCS Lab srl  
Via Pola Esterna 4/12  
41012 Carpi (MO) - Italy  
Ph: +39 059 669813  
info@ncs-company.com



**www.shoulderpacemaker.com**  
**info@shoulderpacemaker.com**