human touch®


humantouch.com

800.355.2762

©2022 Human Touch®, LLC. Patents pending and intellectual property rights enforced.
No medical claims are warranted or implied by the use of this product.
Human Touch is a trademark of Human Touch®, LLC.
HT-ReflexPopUp-NA-A0

Register your product **ONLINE** today

www.humantouch.com/register

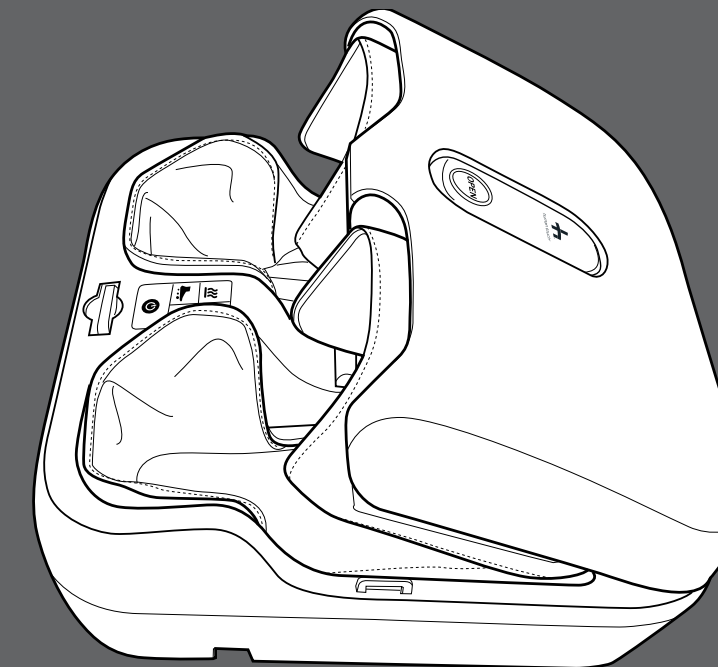HUMAN TOUCH®
**REFLEX
POPUP**
FOOT MASSAGER

human touch®
Your best. Every day.™

use & care manual

Register your product **ONLINE** today

www.humantouch.com/register

ENG

USE & CARE MANUAL

FRC

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For **more than four decades**, Human Touch has established itself as an industry leader that is changing the lives of its customers with products and services that help them **feel, perform, and live their best**. Developed in cooperation with the medical community and industry experts, we deliver products that blend patented technology, craftsmanship and intelligent ergonomics with superior design to enable therapeutic relief for everyone.

Our industry-leading massage chairs, Perfect Chair® recliners and evolving set of innovative solutions share attributes in their ease-of-use, attention to performance and quality, and design to adapt to most any lifestyle. These attributes are highly valued by our customers and garner broad industry award and recognition. As a result, Human Touch empowers people to enjoy greater health and vitality in their professions, pursuits and everyday lives. And with a population increasingly aware of the benefits of wellness, Human Touch extends its brand leadership and available **solutions for the benefit of everyone**.

Proudly endorsed by:

Register your product at www.humantouch.com to ensure your product is in our system to simplify the warranty process should you ever need it.

humantouch.com

HUMAN TOUCH®
**REFLEX
POPUP**
FOOT MASSAGER

welcome to **human touch**®

Your Human Touch® Reflex PopUp Foot Massager is the perfect bridge between state-of-the-art technology and the ancient healing art of massage, emulating the same techniques used by massage professionals. It's hard to believe it's not human! It's a key ingredient for a balanced life—and it's right in your own home!

Our Targeted Relief products deliver therapy to some of the body's most stressed muscles and joints. Beyond delivering unparalleled relaxation and rejuvenation, these amazing devices provide measurable health and wellness benefits that will keep you performing at your best.

The Reflex PopUp Foot Massager provides a fully-encompassing air compression massage that rejuvenates tired feet by forcing blood toward the core of the body to improve circulation. The Reflex PopUp's soothing benefits are complemented by holistic reflexology massage to help improve general health and relieve stress and pain, so you can hit the ground running—every day.

To ensure that you don't overlook any of your foot and calf massager's features, and that you use it safely, read this manual carefully. Then sit back, relax, and enjoy!

keep it safe

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this appliance.

DANGER: To reduce the risk of electrical shock

Always unplug this product from the electrical outlet immediately after using and before cleaning.

WARNING: To reduce the risk of burns, fire, electric shock or injury to persons:

- Do not put hands or fingers near the massage mechanism or rollers while this product is running. Touching the mechanism can result in injury from the rollers squeezing together.
- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Close supervision is necessary when this product is used by, on, or near children, invalids, or disabled persons.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Do not operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return this product to a service center for examination and repair.
- Always keep the power cord away from heated surfaces.
- If your feet feel uncomfortable, discontinue use immediately.
- Never operate this product with the air openings blocked. Keep the air openings free of lint, hair, and the like.

- Never drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the OFF position, then remove plug from outlet.
- Do not massage any area of the body that is swollen, inflamed or covered with skin eruptions.
- Do not use this product when standing on a damp floor, or when any part of the body is in contact with plumbing or any similar ground.
- Do not stand on or in appliance. Use only while seated.
- Anyone with a medical condition, including those who are or may be pregnant, should consult a physician before using any massage product.
- Use heated surfaces carefully. May cause serious burns. Do not use over sensitive skin area or in the presence or poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

Safety Precautions

- Comfortable clothing is recommended while using this product.
- Do not sit on the controller.
- Do not pull on the AC cord.
- When using this product after it has been in storage, check to see that the product is operating properly prior to use.

SAVE THESE INSTRUCTIONS

Do not attempt to repair this product yourself.

- This product is designed for personal in-home use. Commercial use will void warranty.
- To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wide than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not change the plug in any way.

Care and Maintenance

The product is specially designed and constructed to be maintenance free. No periodic lubrication or servicing is required. Your product should only be serviced by an authorized service center. Foreign objects should not be inserted between the massage nodes or the unit housing. Care should be taken not to overload the unit.

STORAGE:

Your product has been carefully designed to give you many years of trouble-free service. Always unplug the unit when not in use. Additionally, if not in use for a long period of time, we recommend that the power cord be coiled and that the unit be placed in an environment free of dust and moisture.

Do not store near heat or open flame. Do not leave this product exposed to direct sunlight for long periods of time. This may cause fading or damage.


WARNING:

Be sure to test any cleaning solution on a small, inconspicuous part of the product to ensure that the color does not fade prior to using the solution on the entire product surface.


WARNING: This appliance is equipped with a thermal resetting protection device. This is a safety feature to protect the appliance from overheating. If the appliance should suddenly stop and will not start, turn off the main switch and pull the plug off the socket outlet. Do not operate the appliance for at least 30 minutes. Failure to turn the appliance off may result in the appliance starting unexpectedly when the device cools.


get ready

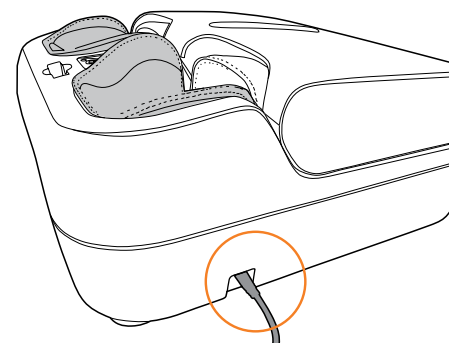
Unpack the box

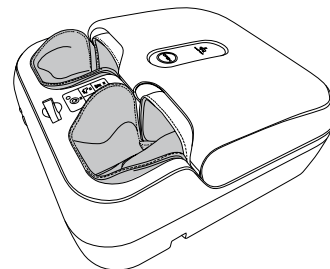Unit

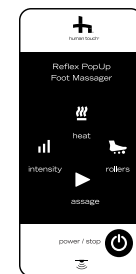Remote Control

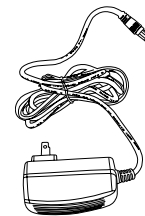AC Power Adapter

Plug It In

Plug the power cord into the right-side of the foot massager and the AC Adapter into an appropriate power source.



Position Your Foot Massager in Front of Your Favorite Chair

take a seat

Open The Cover

Open the cover by pressing down on the **OPEN** button. The cover will pop open, ready for your feet.

OPEN button

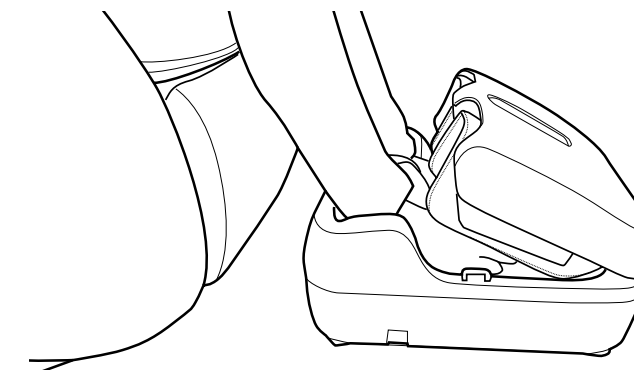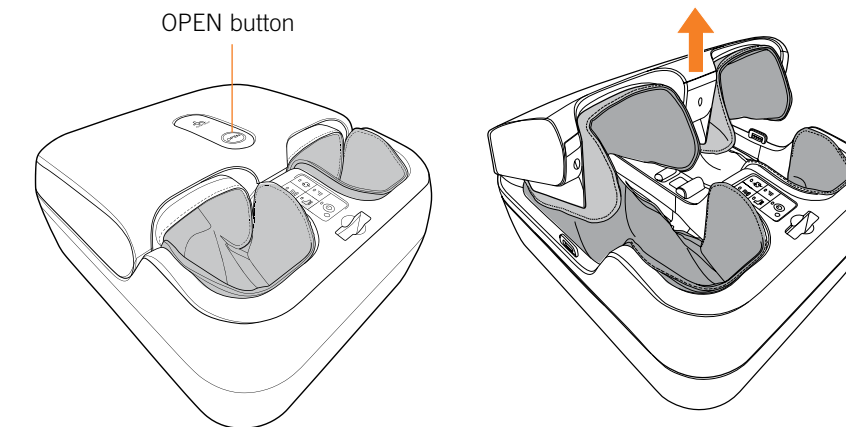
Place Your Feet Into the Foot Massage Wells

Press down on the cover to close. You will hear the cover click into place when fully closed.

Your foot massager comes fully assembled. Simply plug it in, take a seat, pop your feet in, close the cover, and you're set to go!

take control

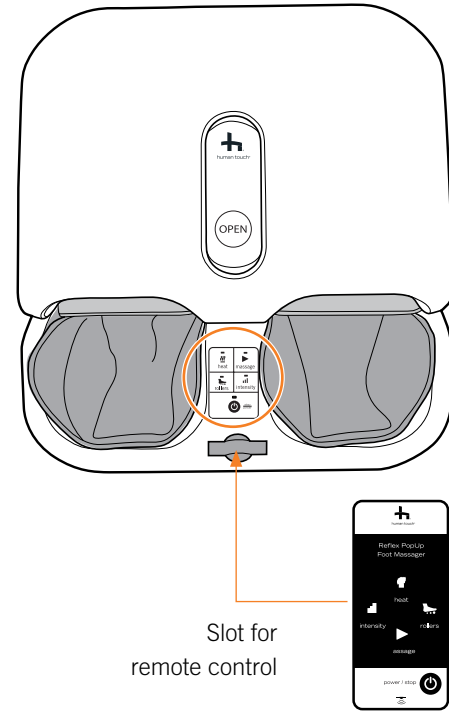Slot for remote control

Heat – Press to add soothing heat to relax your feet and enhance your experience

Rollers – Press to toggle between three underfoot massage speeds: Slow, Medium, or Fast

Power/Stop – Stops all motion


Infrared Sensor – Keep clear to ensure your remote control can communicate with the massager

Massage – Press to toggle between three unique auto-wellness massage experiences

Intensity – Press to toggle between three massage intensities: Low, Medium, High

Power Indicator – Light flashes when unit is on and is steady when massage is in use

⚠ If your feet feel uncomfortable, discontinue use immediately.

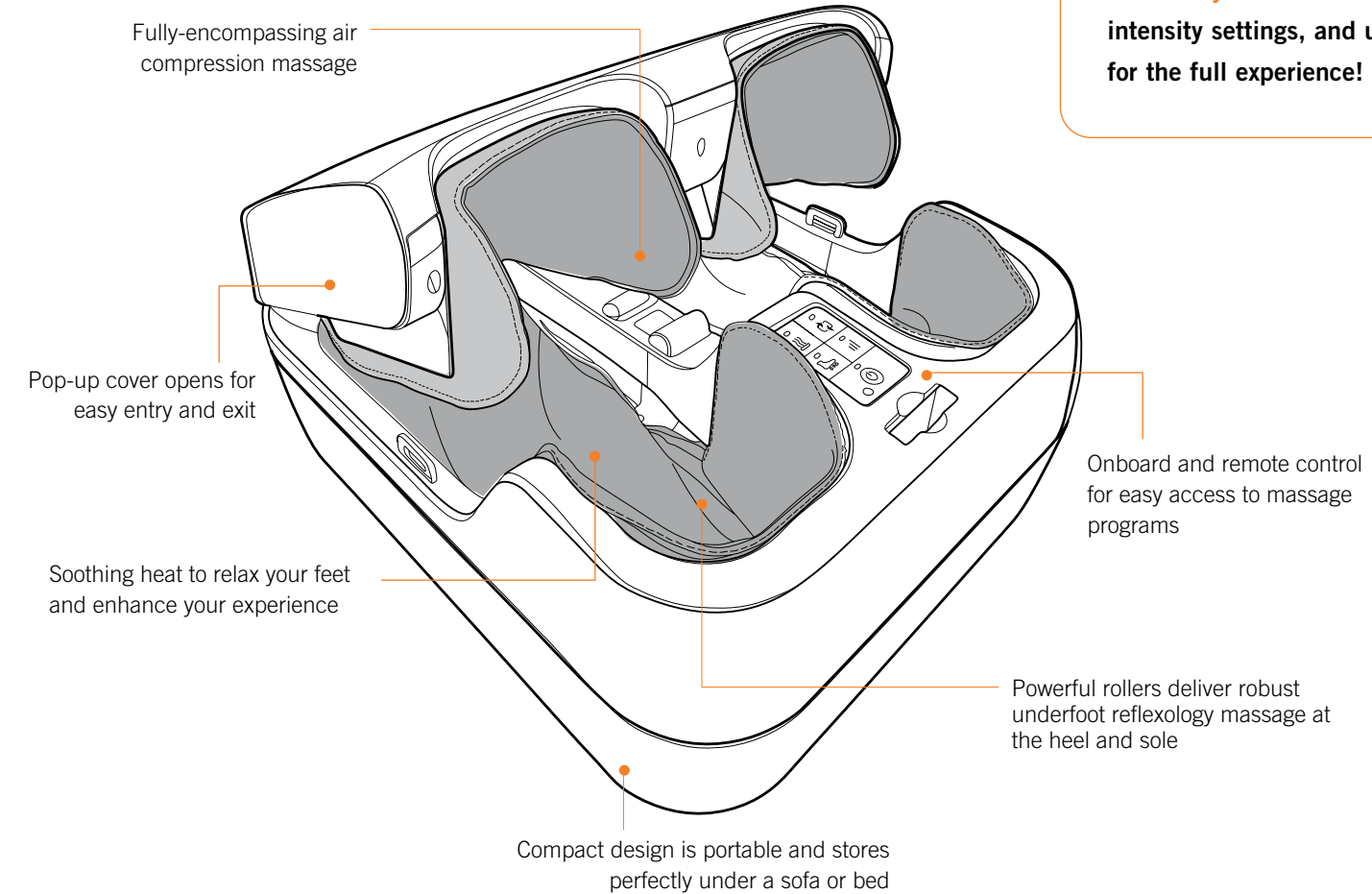Fully-encompassing air compression massage

Pop-up cover opens for easy entry and exit

Soothing heat to relax your feet and enhance your experience

Your foot massager offers a customizable experience to best suit your needs. Choose from three massage programs, intensity settings, and underfoot roller speeds, and add heat for the full experience!

Onboard and remote control for easy access to massage programs

Powerful rollers deliver robust underfoot reflexology massage at the heel and sole

Compact design is portable and stores perfectly under a sofa or bed

Specifications*

Rating: AC 100-240 V, 50-60 Hz

Current: 0.8 A

Product Weight: 11 lb.

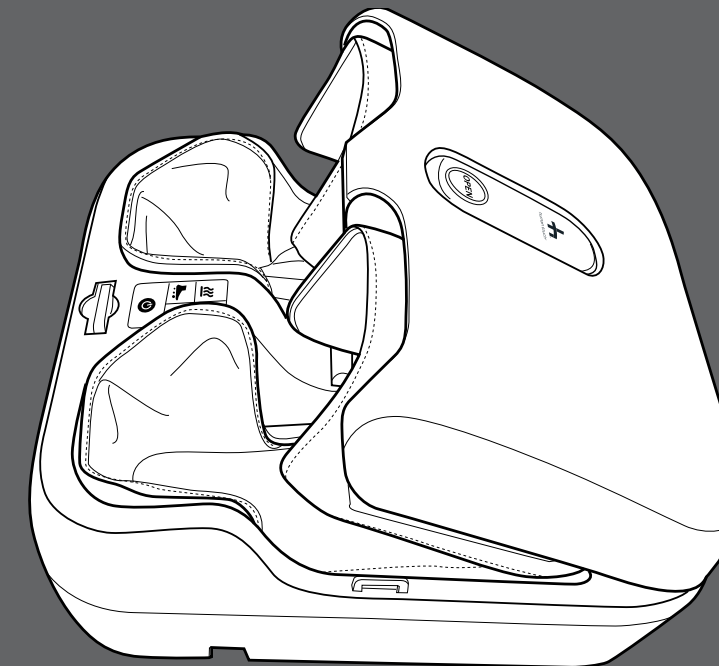
*Specifications are subject to change without prior notice.

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