

DE Fußmassagegerät

GB Foot Massager

FR Masseur de pieds

IT Massaggiatore per piedi

ES Masajeador de pies

SLY-527



Gebrauchsanweisung
Instruction Manual
Mode d'emploi
Istruzioni per l'uso
Instrucciones de manejo

Bitte sorgfältig lesen!
Please read carefully!
A lire attentivement s.v.p.!
Da leggere con attenzione!
¡Por favor lea con cuidado!



IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

Explanation of symbols



This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.



WARNING

These warning notes must be observed to prevent any injury to the user.



CAUTION

These notes must be observed to prevent any damage to the device.



NOTE

These notes give you useful additional information on the installation or operation.



Protection category II

LOT

Lot number



Manufacturer

1 Safety Instructions



Safety Instructions

about the power supply

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only connect the mains plug to the wall socket when the device is switched off.
- Only use the mains adapter that is provided with the foot massager.
- Keep the power supply, the cable, the control unit and the foot massager away from heat, hot surfaces, moisture and liquids. Never get hold of the power supply with wet hands.
- Do not take hold of a device if it falls into water. Pull out the power plug and the power supply unit from the socket immediately.
- Always unplug the equipment immediately after use.
- To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
- Do not carry, pull or twist the foot massager on the cable, the power supply unit.
- You must not continue to use the unit if the lead or the mains adaptor are damaged. For safety reasons, these parts may only be replaced by authorized service personnel. In order to prevent hazards, always send the unit in to the service centre.
- Ensure that the cable will not cause tripping. The cable must not be kinked, trapped or twisted.
- Do not use this product for medical treatment.

specific groups of people

- The surface of the appliance is hot. People who are insensitive to heat must be careful when using the appliance.
- Children must be supervised to ensure that they not play with the device.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- You should refrain from using the FOOT MASSAGE or should con-sult your doctor before using it, if
 - you are pregnant,
 - you have a cardiac pacemaker, artificial joints or electronic implants,
 - you suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Never leave the unit unattended if it is connected to the power supply.
- Take care when using the heater. Do not use if you have circulatory problems or on areas not sensitive to heat. Risk of burning!
- Consult your doctor if you have any questions concerning the therapeutic use of foot massager.
- If you have unexplained pain, are receiving treatment from a doctor and/or use medical equipment, consult your doctor before using the foot massager.

1 Safety Instructions

before using the unit

- Always carefully inspect the power supply, the cable, the control unit and the foot massager for damage before use. Do not attempt to use a defective unit.
- Do not use the unit if there are visible signs of damage on the unit or the cable components, the unit does not work properly or the massager or the control unit have been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.

operating the unit

- The foot massager must only be used for its intended purpose as described in the instruction manual.
- Using the device for any other purpose invalidates the warranty.
- Never leave the unit unattended if it is connected to the power supply.
- The maximum operating time for one session is 15 minutes.
- Do not allow the foot massager to come into contact with pointed or sharp objects.
- Never put the unit down next to an electric stove or other sources of heat.
- Do not stand on the unit.
- Do not touch the power supply unit when you are standing in water and always take hold of the plug with dry hands.



The device is not to be used for commercial use or medical purposes.

If you have health concerns, consult your doctor before using the foot massager.

**The driver may not use the foot massager when driving.
Only use the foot massager in a parked vehicle.**



Only use the foot massager indoors.

Do not use the foot massager in damp rooms (e.g. when bathing or showering).

maintenance and cleaning

- The unit is maintenance free. If a fault does occur, just check whether the mains power adapter is properly connected.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all guarantee rights. Ask your specialist dealer and have repairs carried out by the service.

1 Safety Instructions

- Do not use the unit if there are visible signs of damage on the unit or the cable components, the unit does not work properly or foot massager or the control unit have been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.
- Do not immerse the device in water or other fluids.
- If, nonetheless, liquid should penetrate the device, unplug the mains cable immediately.
- Cleaning and user maintenance shall not be made by children without supervision.



The whole device, Do not wash!

The whole device, Do not dry clean!



The detachable foot sleeve is machine washable


2 Useful Information

Thank you very much

Thank you for purchasing SNAILAX foot massager.
With normal care and proper treatment, it will provide years of reliable service.

Please read all instructions carefully before using this product.

2.1 Contents

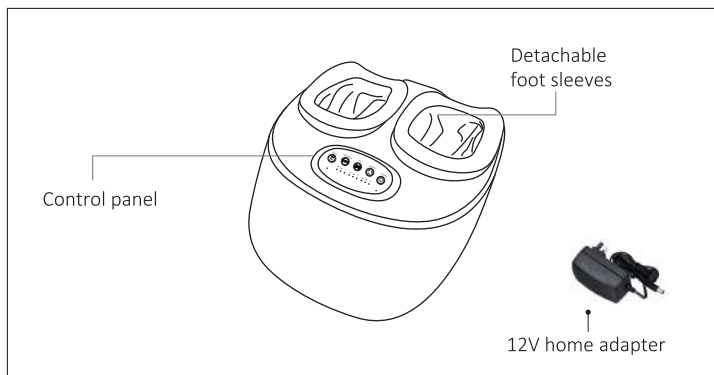

2.2 Technical Data

Name and model	:	Foot massager SLY-527
Voltage supply	:	100-240 V~50/60 Hz
Power supply unit	:	12V \approx 2000 mA
Power consumption	:	24W
Timer	:	15 minutes
Operating conditions	:	only use in dry rooms
Storage conditions	:	in a cool, dry place
Dimensions	:	approx. 34.8×32.5×17.3 cm
Weight	:	approx. 3.1 kg
Adapter cable length	:	approx. 1.8 m




WARNING

Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

3 Operating

1. Place the massage machine on the flat floor.
2. Connect the adapter cable to the corresponding port on the massager

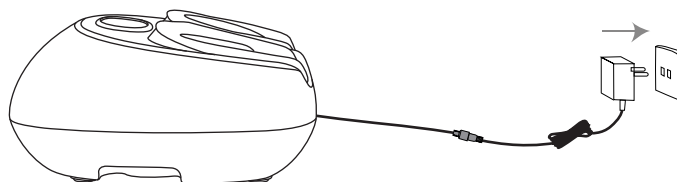


3. Plug the home adapter into electrical outlet.
4. Turn on the device through operating control panel (instructions on page 17).



- Before plug ,please check and make sure the power cord is undamaged.
- Do not use the wet hand to plug in the adapter.
- Please don't use the product by standing posture or put any heavy object on the product to avoid any broken caused.
- Please don't stand or seat on the product .
- When moving the product, please don't pull the power cord or drag the product directly to avoid any damaged caused.
- Please don't let children and pets play the products for preventing any accident happened.
- Do not use the wet hand to pull out the power cord in order to avoid any electric shock.
- Do not use this product in damp environment to avoid electric shock.

Setup and Operation

- Please choose a chair with desired height and be seated, then put your feet in the massager.
- Operate below control panel to get foot massages.

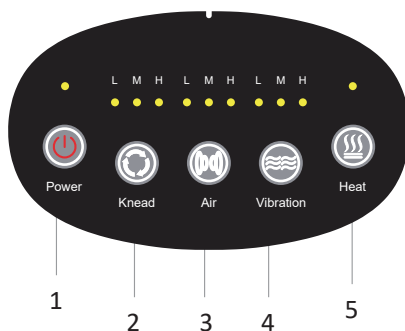

Setup and Operation

1. **Power button:** Press this button to turn on or off the massager;

Turn on the massager and it automatically will work in the default mode:

i.e. The rolling +compression +vibration+heating functions will work at Low speed

2. **Rolling/Kneading button:**

After power on, press this button to turn on kneading massage function, continuously press to cycle through Low/Medium/High Speed or Off setting, the corresponding indicator L/M/H will light up.

3 **Compression button:**

After power on, press this button to turn on compression function, continuously press to cycle through Low/Medium/High level of intensity of Off setting, the corresponding indicator L/M/H will light up.

4. **Vibration Button:**

After power on, press this button to turn on Vibration function, continuously press to cycle through Low/Medium/High level of intensity of Off setting, the corresponding indicator L/M/H will light up.

5. **Heat button:**

After power on, press this button to turn on Heat function, continuously press to choose Low/High heat level or Off setting.

Low level: the temperature is 45-50 °C (113-122 °F), the indicator light is blue;

High Level: the temperature is 50-55 °C (122-131 °F), the indicator light is red.



- Please don't put water or other liquid on the control panel.
- Please don't press the button too long time at a time to avoid the program confused.
- When power failure happens during massage, please turn off the power and pull out of the plug from socket.
- Please do not let your hand or fingers close to and inserted into the gap of products for preventing any injuries caused.
- Please do not let children and pets play the products for preventing any accident happened.

4 Miscellaneous

4.1

Cleaning and maintenance

- Before cleaning, make sure the power supply is turned off and the plug has been pulled out from the socket.
- It is forbidden to dip the whole product in water or other liquid when cleaning.
- Use a soft cloth to wipe the surface of the product and control panel.
- Use a piece of cloth that has been dipped with cleaning agents to wipe the main unit of the foot Massager, and then use a piece of dry and soft cloth to wipe it.
- Do not pull the power plug by damp hands.
- Do not use water washing product, do not spray the water or other liquids to the interior of the product.
- Do not use corrosive cleaning agents like benzene or thinner to clean the Foot Massager, otherwise, the color of upholstery will fade or the Foot Massager will crack or be eroded.
- Do not iron the upholstery.
- When product is wet, let the product dry naturally and thoroughly.

Removable & Washable Cover

- Detach the foot sleeves and wash them in machine or by hand. Allow to air dry.

Maintenance

- Make sure the power supply is turned off and the plug has been pulled out from the socket.
- Sweep off the dust. Keep the foot massager clean and away from damp places.
- Store the massager in a cool, dry place out of reach of children.
- Avoid contact with sharp edges or pointed objects which might cut or puncture the cover of massager.

4.2

Disposal



This product must not be disposed of together with domestic waste.

All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

Consult your municipal authority or your dealer for information about disposal.



In accordance with our policy of continual product improvement, we reserve the right to make technical and optical changes without notice.

5 Warranty

If the product(s) is damaged or defective within a reasonable time, Snailax offers the warranty as below.

Period	Return for refund	Return for replacement	Who paid shipping?
Within 30 Days	Yes	Yes if undamaged	Snailax
31- 90 Days	Yes after being justified	Yes after being justified	Snailax
91 Days - 1 Year	No	Yes after being justified	Buyer
1 Year - 3 Years	No	Only to buyers who have received extended warranty	Buyer

Snailax warrants the product to be free from defects in workmanship and materials, under normal use, for a period of 12 months from the original invoice date, except as noted below. Snailax products can be repaired or replaced within warranty period or returned for refund within 90 days if the products break down in the period of proper use.

The warranty extends only to consumers, NOT to any retailers or resellers. This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. Warranty claims for product defects, unless extended warranty is registered, expire 12 months from the purchase day. Warranty claims CANNOT be processed for items that have exceeded the warranty period. The warranty does NOT cover the replaced, refurbished product. The warranty does have other limitations as below "Warranty Limitations". To obtain warranty service on Snailax products, follow "Warranty Claim Procedures" below.

Warranty Claim Procedures:

1. Any product the customer believes is defective and seeks for warranty, the customer must provide to Customer Support:
 - a) Valid proof of purchase (see below the definition of valid proof of purchase);
 - b) Pictures and/or description are required to claim defective merchandise;
 - c) Model number and Series Number label attached to each product
2. If the warranty claim is justified by Snailax Customer Support, customer(s) will obtain from Snailax a Warranty/Return Authorization Number (RMA Number)
3. Snailax reserves the right to specify that items be returned to the designated warehouse for inspection or be inspected by our representative in the field or be retained by customer.
4. Any claim for defective merchandise returns must be packed in original packaging
5. The item(s) will be repaired or replaced or a refund credit will be issued once the item(s) are delivered. When returning items with RMA# and a prepaid shipping label provided by Snailax, Snailax takes responsibility for any damage or loss incurred in transit.

Valid proof of purchase under this Warranty:

Order numbers from online purchases made through Snailax official store (www.snailax.com) or Snailax's direct-authorized resellers, currently only at Amazon/eBay/Walmart. For product(s) obtained from other Channels, please contact the customer service of the Channel directly for appropriate solution.

Warranty

Warranty Limitations(this warranty does not include):

Normal wear and tear; Any condition resulting from other than ordinary residential wear or any use for which the product was not intended, such as use in rental or contract trade or commercial use; Any condition resulting from incorrect or inadequate maintenance or care; Damage resulting from misuse, abuse, negligence, accidents or shipping damage, exposed to environmental conditions temperatures, water, improper operation, improper installation, improper use of electrical supply; transportation damage; theft, vandalism or; Damage resulting from the attachment or modification of any unauthorized accessories; Fail to show the valid proof of purchase; Lost or stolen products; Items that have expired their warranty period; No quality-related issues (after 30 days of purchase); Free product/ Gift offer by Snailax; unauthorized repairs or modification; Replaced or refurbished product. Dissatisfaction due to buyer's remorse (after 30 days of purchase).

Refund Procedures:

(Refer to <https://www.snailax.com/pages/refund-policy> for detail instructions)

Warranty Disclaimer

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. SNAILAX SHALL HAVE NO LIABILITY FOR, AND EXPRESSLY DISCLAIMS ANY WARRANTY OR AFFIRMATION OF FACT, EXPRESS OR IMPLIED, OTHER THAN AS SET FORTH IN THIS WARRANTY STATEMENT, INCLUDING, WITHOUT LIMITATION

- (1) THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE;
- (2) ANY WARRANTY OR AFFIRMATION OF FACT RELATED TO MISUSE, IMPROPER SELECTION, RECOMMENDATION, OR MISAPPLICATION OF ANY PRODUCT; AND
- (3) ANY WARRANTY OR AFFIRMATION OF FACT THAT THE CATALOGS, LITERATURE AND WEBSITES IT PROVIDES ACCURATELY ILLUSTRATE AND DESCRIBE PRODUCTS.

Snailax reserves all the right for final interpretation of these Terms of Warranty and Return Policy above and the right to change, modify, add or remove portions of these terms at any time without prior notification.

Extend Warranty for Free

1) Extend your warranty to 3 years by visiting our Snailax online service center at <http://support.snailax.com> or by scanning the QR code below:



2) Alternatively, like our Snailax Facebook page at <https://facebook.com/snailax> or scan the QR code below, to message us with your purchase proof



Your application for a free warranty extension must be completed within 60 days of purchase.



Email: supporteu@snailax.com

Internet: www.snailax.com

UK	REP
-----------	------------

TANMET INT'L BUSINESS LTD
9 Pantygraigwen Road, Pontypridd,
Mid Glamorgan, CF37 2RR, UK
tanmetbiz@outlook.com