



Dimensions

Size:18*16*13cm

3.7V 3000MAH 18650 liion battery (include)

Before using this product, please read it carefully and keep it properly.



Vibration



Heating



Infrared



laser

Structure



Features

1:Indicator light

3:Mesh cloth

5:Laser infrared

7:Type-C power jack

2:LED light touch screen


4:Adjusting button

6:Vibration Area


8:Velcro

Button panel introduction

- LED light touch screen, clearly visible and easy to operate.
- Ergonomic design, fully wrapped to fit the knee joint.
- Vibration massage automatically shutdown within 2 minutes if no any operates.
- Built-in 3 infrared laser lights, carbon fiber hot compress.
- The buzzer prompts button is touched.
- 10 minutes timer.
- Low-grade heating 40 degrees, high-grade heating 45 degrees.

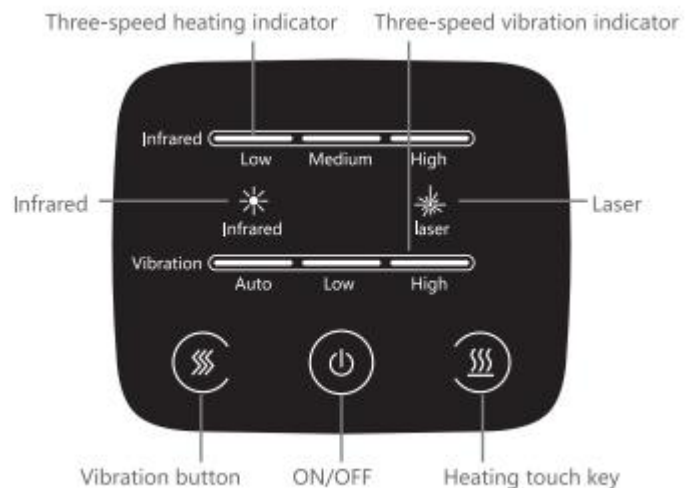

Instructions

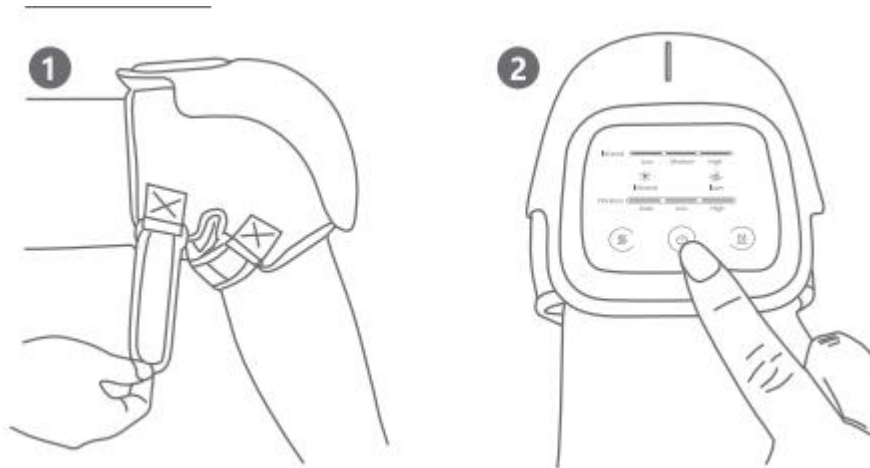

1.Put the product on your knees,fasten the product on your knees with a fixing strap,and pull the product against your kness near the fixed belt.

2.Long press "ON/OFF" to turn on the Auto made (vibration mode,low-grade heating,infrared,and laser)

3.When powered on:Short press "vibration button"to

adjust the vibration mode (auto-low-high-turn off the vibration cycle in turn...),Light turn white.
Short press "Heating touch key"to adjust the herting temperature (low-medium-high-close heating cycle in turn...),Light turn white.

Short press "ON/OFF"to turn on and off the laser and infrared (turn on infrared-turn on the laser-ture off the infrared -ture off the laser cycle in ture ...),Light turn white.

4.Long press"ON/OFF"to turn off the device.

Disposal

- Coopers of Stortford use recyclable or recycled packaging where possible.
- Please dispose of all packaging, paper, cartons, packaging in accordance with your local recycling regulations.
- Plastics, polybags-the handle is made from the following recvclable plastic.

Code & Symbol	 OTHER
Type of Plastic	ABS
Commonly used for	CD's, crisp packets, various flexible packaging, baby bottles, sunglass lenses
Notes	Not often or widely recycled

- At the end of the product's lifespan please check with your local council authorised household waste recycling centre for disposal.



*Calls cost 7 pence per minute plus your telephone company's network access charge

Head Office, 11 Bridge Street,
Bishop's Stortford CM23 2JU
Tel: 0844 482 4400*
Fax: 01279 756 595
www.coopersofstortford.co.uk

Coopers
OF STORTFORD