

# Beautify® Micro-Infusion Treatment

## Quick Guide



< 10 Minutes Total



Every 2 Weeks

**Video Treatment Instructions**

**Scan QR Code:**



**KEEP MICRONEEDLE TIP & INFUSION VIALS  
OUT OF THE REACH OF CHILDREN & PETS;  
THEY CAN BE A CHOKING HAZARD.**

**IMPORTANT:**

**PLEASE READ ALL INSTRUCTIONS  
FOR SAFETY & BEST RESULTS**

**Take Before & After Pics!**

**Be Sure to Take Before, During & After Pics on Day 1 - Day 15 & Day 30.  
This allows you to better see your actual progress & results!**

***Set Your Calendar Reminder.***

## PRE-CARE INSTRUCTIONS:

- **No Botox or Injectables 2 weeks** before and after a Micro-Infusion Treatment.
- **Do not lay out in the Sun** three (3) days before, or seven (7) days (*minimum*) after your Micro-Infusion treatment. If you are sunburnt, *do not do a treatment*.
- **If you use any Glycolic/Lactic Acid or Prescription Retin-A Products;** discontinue use three (3) days before, *and for seven (7) days after treatment*.
- **If you use Non-Prescription Retin-A Products;** discontinue use at least one (1) day, preferably three (3) days, before treatment and discontinue use for three (3) days after treatment.
- **Cold Sores or Herpes Simplex:** If you are highly prone to cold sores, a micro-infusion treatment can sometimes cause a break-out, especially if we are doing a treatment on your lips. We recommend that you take an anti-viral drug such as; Zovirax, Famvir, Valtrex or L-lysine as prescribed by your Physician. *This is at your discretion*. Take the medication as prescribed for 1 week prior and one week after the lip or full face Micro-Infusion Treatment.

## PRE-TREATMENT - IMPORTANT

1. Tie your hair back and keep it out of your face.
2. **Three Cleanses:** First; double cleanse your Face & Neck (*Two Cleanses*) to ensure all makeup and products are fully removed on skin. *Use makeup remover if necessary.*  
**Cleanse #3:** Use 70% Isopropyl alcohol on Face & Neck with a makeup cotton round. **It's very important to remove every bit of oil, dirt, and bad bacteria to prevent irritation or adverse reactions during treatment as you will have (invisible) open wounds in your skin.**
3. **Make sure your hands and the area you're using are clean.** If you touch anything besides your skin or the device, like a bottle or brush, it could cross-contaminate. Stick to just touching your face and the tool during the treatment.

## STEP 1: PREPARING THE SERUM VIAL FOR TREATMENT

1. **Open Sterile Pack Roller Tip:** Remove the Roller tip from its sterile package - (*Leave the Roller cap on that covers the micro-needle Roller*).
2. **Open Sterile Treatment Serum Vial:** Your sealed glass vial is already filled with our Million Dollar Growth Factor Serum. Hold the glass vial firmly: Keep a steady grip on the body of the vial with one hand. Look for the two arrows on the Gold cap...**Break the Gold Cap seal:** Simply lift up on the clear plastic piece that covers the Gold cap from where the two arrows point. Once the seal is broken, the plastic piece and Gold cap should easily pop off. Gently remove the remaining Gold seal from the top of the vial and dispose of it.
3. **Triple Infusion Treatment:** Pour a generous amount of serum into a clean hand and quickly apply it to your face & neck and even your lips. *Then quickly do Step 4 so that the serum does not dry.* (This is first step of your triple infusion).
4. **Insert Roller Tip into Serum Vial:** Remove the rubber stopper and *Keep Rubber Stopper – Do Not Discard*. Leave Roller Tip Cap on and Insert the Roller tip by firmly pressing it down into the bottle and snap it into the top of the glass Serum vial...**it should be snug and seated straight (not at an angle).**

## STEP 2: PRIMING THE MICRO-NEEDLE ROLLER TIP

- Remove the clear cap on the Microneedle Tip Roller. Turn the Roller Bottle upside down, and *gently* press the Roller needle head onto a clean paper towel into a paper towel on a table or Counter. Then, roll it gently on the paper towel – This will help prime the Serum. (*There is a spring mechanism that moves up and down that allows the serum to flow out*). **Do Not press the Roller down too hard as it may bend the needles. Do Not Discard the Clear Cap.**
- 2. **Look for bubbles:** You should see bubbles forming in the vial, which shows that the serum is beginning to flow. Keep pressing until you see the serum flowing normally on to the paper towel. While rolling; gently push the vial up & down to pump the serum; this allows it to flow out more easily. **Note: Serum does not flow out excessively.**

## STEP 3: TREATMENT – MICRO-NEEDLING THE SERUM INTO YOUR SKIN

### 1. Treatment Tips for Best Results:

- **Use a Wall mirror.** You can lay on your bed with a hand-held mirror as well for better serum flow. Total actual treatment time is about 3 to 5 minutes.
- **Keep the treatment microneedle vial upside down (tip down) when treating your skin.**
- **TO IMPROVE SERUM FLOW:** *While rolling your skin; push the vial 'up & down' into the skin to help pump out the serum, this allows the serum to flow out more easily.*
- **Note: Serum does not flow out excessively.**
- **Push Microneedle Tip down firmly into the skin - not too hard or roll too fast - do a slow speed or a medium-slow speed.**
- **One pass is all you need...if you have a problem area; you can do a second pass over it.**
- **DO NOT REUSE TIP OR INFUSION VIAL.** *Bacteria growth would cause treatment issues if re-used and could cause serious skin issues!*
- **Keep glass Serum vial & Microneedle Roller Tip Out of the reach of Children & Pets. Discard safely.**

### 2. Treatment Directions (Keep Serum bottle tipped up)

- **It's best to start on the Forehead;** Do vertical passes; up & down – You can do horizontal passes as well. You can do double passes if so desired. You can also do cross-hatching – X's.
- **Temples;** Next do your Temple Areas.
- **Face;** Do the entire right side of your face; from your Upper Cheek to the jawline.
- Then do the left side of the face; from your Upper Cheek to the jawline.
- **Under Eyes; Be very careful around your eyes; Roll Horizontally ONLY.** A slower, more controlled pace is recommended around the eye area. You can do the upper brow area. **DO NOT TREAT EYELIDS!** *You can use a magnifying mirror for your under eye area, if you prefer.*
- **Chin & Upper Lip;** Do your Chin area & upper lip area. You can treat your Lips as well; *It's not painful.*
- **Neck; Important -** Be sure to treat your Neck; From your Jawline to the bottom of your Neck.
- *If so desired; you can treat your upper Décolleté as well.*
- Put the clear cap back on the Roller tip and carefully pull (remove) the Roller tip out of the bottle by wiggling it back and forth and dispose of it properly - **Do Not Flush. Keep away from Children & Pets.**
- **Left-Over Serum:** Once your treatment is completed; you will have some serum left in the vial. Immediately pour some of the remaining serum into your (Clean) hand and gently apply the serum evenly to your face & neck. **Do not press too hard. No need to use all the Serum.**
- You still will have some serum is still left, put the gray rubber cap back into the serum bottle. Apply any remaining serum **8 to 12 hours later** on a clean face & neck. *And again 8 to 12 hours later, if there is any Serum still left.*
- Properly dispose of the glass Serum vial & microneedle Tip in the trash when done – **Do Not Flush.** The vial and Roller tip are designed for one-use only. **DO NOT REUSE TIP OR INFUSION VIAL.** *Bacteria growth would cause treatment issues if re-used and could cause serious skin issues.*

3. **Optimal Collagen Growth Results:** Only perform one treatment every two weeks. **Do not treat every week.** Your Body will begin Collagen production immediately and your skin will improve over time. **Perform your second treatment two weeks later:** Collagen growth peaks at Day 14, so performing another treatment on Day 14 will help maintain that Collagen peak, producing maximum results! For even better results; use the Beautify Million Dollar Growth Factor Daily Serum, twice a day on your face & neck – **Save \$111! MillionDollarSerum.com Use QR Code on next page.**


## POST-CARE INSTRUCTIONS:

- **Sleep on a Clean Pillow Case; Change your pillow case.** Do not touch the treated area as it may be more prone to infection; *this includes children & kissing pets.* It's important to keep it as clean as possible. If your skin itches; lightly tap the itch with a new Q-tip.
- **Do Not Wash Your Face after Your Micro-Infusion Treatment Before Going to Sleep.** The following day after your treatment, you may gently cleanse the treated area. Be sure to clean (sanitize) your hands prior. **DO NOT SCRUB YOUR FACE.**
- **Do not apply ANY make-up cosmetics, or skin care products for a minimum of 24 hours** following your Micro-Infusion treatment. We recommend that you do not use any make-up cosmetics on the treated area for 48 hours after your Micro-Infusion treatment as make-up may contain bacteria. **Clean your make-up brushes with Alcohol. It is best if you use Mineral make-up.**
- **DO NOT lay in direct sunlight for at least one (1) week.** If you are going to be exposed to the Sun, wear a hat & apply an SPF 30 or higher Natural Sun Screen starting on the 2nd day, after your Micro-Infusion treatment. Additionally, it is important to avoid Tanning Bed exposure for at least one week post treatment. **Only use a Mineral-based Sunscreen – Not Chemical Based.**
- **For lasting results, reduce the amount of direct sun exposure and tanning beds as part of your healthy skin care regimen.**
- **Do not use any Glycolic/Lactic Acid or Retin-A Products for a minimum of three (3) days after your Micro-Infusion Treatment.** Beautify has the Worl's Best Mineral-based Sunscreen. See QR Code Below:
- **Do not sweat heavily, workout, swim, use a steam room or sauna, or do hot yoga for one day** after your Micro-Infusion treatment.
- **Do not get a Massage, a Chiropractic adjustment, do Physical Therapy or go to the Dentist for a minimum of one day after your Treatment.** It's best if you do not travel for 24 hours after treatment. Wait 3 days before doing an LED Skin Rejuvenation Treatment.
- **Do not or have any other facial therapy performed for a minimum of (7) seven days. Wait two (2) weeks to get Botox or Injectables after your Micro-Infusion Treatment.**

**Thank You for Your Purchase – We Hope You Enjoy Your MicroInfusion Results!**

If you have any questions or issues, please contact us at: [Support@JustBeautify.com](mailto:Support@JustBeautify.com)

**If You Like Your Treatment Results;  
Please Write Us A Review!**

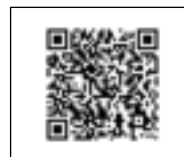

**You Can Purchase Beautify® Million Dollar Growth Factor  
1 oz. Daily Serum - The Gold Standard of Serums.  
Save \$111!**



**You Can Purchase Beautify® World's Best  
Natural Sunscreen. 2 oz. & 4 oz.:  
Our Sunscreen Anti-Aging Formula Contains 7 Anti-Oxidants  
to Fight Free-Radicals and it even contains a Copper Peptide,  
which helps grow Collagen! SPF 30+ Save up to 50%!**



**MillionDollarSerum.com**

Beautify® Skin Care 2025

Beautify I AM LLC

# Receive 3 Beautify® Skin Care Products **100% FREE!**

## Pick Any Three! Up to \$580 Value!

### It's Easy! Just Use our Beautify® Million Dollar Growth Factor Micro-Infusion System for 90 Days:

**Do 1 Treatment Every 2 Weeks: Only 6 Treatments Total.**

**Take a Selfie Pic On:**


- Day 1
- Day 45
- Day 90 (3 pics total)

Upload your pics, write a review and if approved, we will send you a one-time Coupon Code to receive 100% Off on Three Full-Size Beautify® Skin Care Products, *including Free Shipping!* Retail Value up to \$580.00!

**Choose Any of the Following Products – 100% FREE!**

- ◆ Moisturizers
- ◆ Serums
- ◆ Cleansers
- ◆ Eye Treatment Crème
- ◆ Natural Sunscreen
- ◆ Scrub
- ◆ Mask

**USE QR CODE TO REGISTER  
FOR YOUR 3 FREE PRODUCTS:**

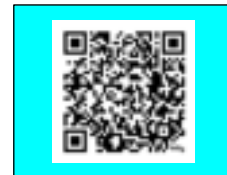

**NOTE: You can choose any 3 Beautify® Skin Care Products!\***  
**\*(Each Product must be different - you cannot choose Million Dollar Serum).**

## BEAUTIFUL, MORE YOUTHFUL LOOKING SKIN!

### MILLION DOLLAR GROWTH FACTOR MICRO-INFUSION SYSTEM

**MillionDollarSerum.com**  
Beautify® Skin Care 2025  
Beautify I AM LLC