CG10F

User Manual

EN

Multifunctional Scale for Body Weight

1- 5

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FR

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6- 10

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ES

Báscula multifuncional para peso corporal

11- 15

Precautions


Read carefully before use:

- ⚠ Do not use the scale on a wet surface. This could lead to slipping.
- ⚠ Dispose of the scale properly following your local regulations.
- ⚠ Do not jump on the scale. This scale has a weight limit of 400 lbs.

Scale Overview


This smart scale can automatically calculate and display the difference between two consecutive weight measurements.

BMI (Body mass index) is a person's weight in kilograms divided by the square of height in meters. This indicator will be calculated and shown in Fit Profile App.

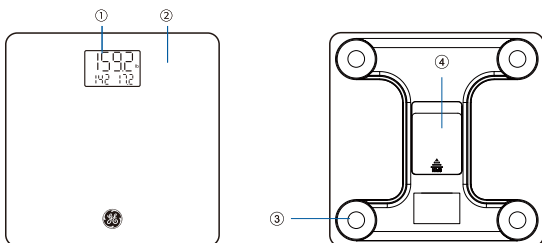Scale Layout

- ① LED Display
- ② Tempered Glass Surface
- ③ Scale Feet
- ④ Battery Box/Compartment

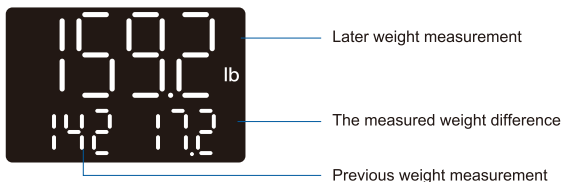Scale Overview

Model:	CG10F
Dimensions:	11.0" x 11.0" x 1.0" / 280 x 280 x 25 mm
Weight:	2.6 lb / 1.2 kg
Weighing Range:	0.45 lb - 400 lb / 0.2 kg - 180 kg
Power:	3 x 1.5V AAA battery

Get Started

1. Download the app and create an account

- (1) Download the app. Go to the App Store for iOS or Google Play Store for Android, and search for "Fit Profile", or scan the QR code below.
- (2) Sign up and complete your profile, including height, etc.

Scan QR code to download the app

Sign up

- (3) For Android, turn on "Location". For iPhone, allow "Fit Profile" to use Bluetooth.
 - a) Open the "Fit Profile" app.
 - b) Tap "+" Pair the Device".
 - c) Turn on the scale by briefly stepping on and off it. On the app, tap the digital scale icon, then select "CG10F" and click "Connect"

Your scale has now been successfully paired with your smartphone and is ready to show more data.

2. Setup your scale


- (1) Before putting your digital scale into operation, please open the battery compartment on the underside of the device and insert the three 1.5 V batteries (type AAA) included.
- (2) Place the scale on a hard flat floor surface.
- (3) Turn on the scale by briefly stepping on and off it. Then, wait for the scale to be set to 0.00.
- (4) Step onto the scale barefoot and stand still until the numbers flash twice and stabilize.
(A single measurement is completed, and you can track all you need using the app for routine use.)
If you need track your BMI and more detailed values, please measure your weight after completing app registration.

(5) Difference Function Usage

- a) Second Weighing: Step off, then quickly step back on before the screen turns off to complete the second measurement.
- b) View Results: The weight difference between the two measurements will be displayed in the lower right corner of the screen.

3. Take measurements

- (1) Open the "Fit Profile" app, and then keep the app open on the "Home" screen.
- (2) Place the scale on a hard flat floor surface, do not place it on carpet or any other soft surfaces.

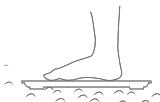- (3) Turn on the scale by briefly stepping on and off it. Then, wait for the scale to set to 0.
- (4) Step on the scale while barefoot till the numbers display on the scale. Your measurement results will also display on the app.

Recommendations

To have accurate, consistent, and comparable measurement results, please follow the recommendations below.

1. DO' s

- Ensure the personal information in your profile is correct.
- Place the scale on a hard flat floor surface when you weigh.
- Make sure the scale displays a 0 before stepping on to measure. While barefoot, stand still on the scale until the measurement is completed.
- Weigh yourself under similar conditions and around the same time every day.
- Empty your bladder before weighing.

2. DON' Ts

- Do not place the scale on a carpet or a soft surface when you weigh.
- Do not measure with socks or shoes on.
- Do not measure when your body is wet. For example, immediately after a shower, bath, sauna, or strenuous exercise.
- Do not measure after drinking a lot of water or alcohol.
- Do not measure if you have a cold or fever.

Troubleshooting

 Normal working	No display Auto-off, or no power	 No wireless connection	---- Try Again
 Change battery	 Overload	 Unstable weighing	 The scale needs to be calibrated

Messages Displayed

Symptom	Cause	Solution
Lo	Low power	Use new 1.5V AAA batteries
EEr	Overload	Use the scale within 400 lb.
Err1	Unstable measurement	Place the scale on the flat, hard surface and measure again.
----	Unable to calibrate	Place the scale on the flat, hard surface and let the scale calibrate.
No display	Auto-off, or no power	Auto-off: Step on then off to turn on the scale ; No power: replace new batteries.
Display Weight Only	Not completed profile data in the app	No User Profile: Create an app profile and weigh once with Bluetooth connected. No Electrical Current: Need to stand on the scale barefoot. Users: Not suitable for children under 10. Other Issues: Contact Customer Support for assistance.
Unable to Sync with App	Lost connection, or scale turns off	Reconnect the scale: For iPhones, authorize bluetooth; for Android, enable location. Click the icon on the top-right corner and tap "My Device" ->"+ Add device" on the app.

Warranty

This product is warranted by the manufacturer against defects in materials and workmanship for one (1) year from the original date of purchase. Proof of purchase is required.

This warranty does not cover accidental damage, mistreatment, or unauthorized repair. Detailed warranty policy: www.fit-profile.com/warranty

If you have any issues or questions, please contact Global eFit Inc.

Contact: Global eFit Inc.

cs@fit-profile.com

(866) 786-5988 (Mon – Fri. 9:00 am – 5:00 pm, US Central Time)

FCC Compliance Statements

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, under Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user authority to operate the equipment.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction

ISED Compliance Statements

This device contains license-exempt transmitter(s)/receiver(s) that comply with Innovation, Science, and Economic Development Canada's license-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

This equipment complies with IC RSS-102 radiation exposure limits set forth for an uncontrolled environment.

Global eFit Inc.
17304 Preston Road, Suite 800, Dallas, TX 75252, USA
www.fit-profile.com
cs@fit-profile.com
1-866-786-5988
(toll-free US & Canada; M-F 9:00am - 5:00pm Central Time)

FCC ID: 2ANDX-CS20M2
IC ID: 26148-CS20N

GE is a registered trademark of General Electric Company.
Manufactured under trademark license.

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Si tiene algún problema o pregunta, comuníquese con Global eFit Inc.