DE

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EN

Operating instructions

IT

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SAFETY INSTRUCTIONS

- A temperature range of 36.0°C to 38.0°C is recommended for bathing water for a baby.
- **CAUTION:** Water temperatures above 40 °C can injure the baby's skin.
- Please supervise your baby or toddler when the bath thermometer is in the bath water. This device is not a toy. It contains small parts (battery, etc.) that could be swallowed by children. In the unlikely event that a battery or other part is swallowed, seek medical attention immediately.
- Batteries must not be recharged or reactivated by any other means, disassembled, thrown into a fire, or short-circuited.
- **DANGER OF ELECTRIC SHOCK!** This device contains electronic components powered by a power source (battery). Never leave children unattended when handling the device! Only use the device as described in the instruction manual, otherwise there is a **RISK** of an **ELECTRIC SHOCK!**
- **FIRE/EXPLOSION HAZARD!** Do not expose the device to high temperatures. Use only the recommended batteries. Do not short-circuit or throw the device and batteries into a fire! Excessive heat and improper handling can cause short circuits, fires, and even explosions!
- **RISK OF BURNS!** Batteries must be kept out of children's reach! Ensure correct polarity when inserting the battery. Leaking or damaged batteries will cause chemical burns if they come into contact with the skin. Use suitable protective gloves if necessary.
- Leaking battery acid can cause chemical burns! Avoid contact of battery acid with skin, eyes, and mucous membranes. In case of contact with the acid, immediately rinse the affected areas with plenty of clear water and seek medical attention.
- Do not disassemble the device! In the event of a defect, please contact your dealer. The dealer will contact the service centre and send the device for repair if necessary.
- This device often requires the use of sharp-edged and pointed tools. Therefore, keep this device, all accessories, and tools out of the reach of children. There is a **DANGER OF INJURY!**
- Keep the instructions and packaging as they contain important information.

DOWNLOAD INSTRUCTIONS:

www.bresser.de/download/9810402

INTENDED USE

This bath thermometer is used to measure the water temperature

On/off switch

LCD display

Rear cover

Sensor

Temperature

Low temperature display

High temperature display

Battery low indicator

when the baby takes a bath. It has a digital display and can be used in the household. It is not a medical device; children should use it under adult supervision.

FEATURES

- Easy-to-read digital LCD display.
- Compact, accurate, and durable device.
- The device automatically turns off after about 30 minutes if not manually turned off earlier.
- The temperature display is updated every 10 seconds.
- Low battery indicator
- The device is waterproof.

OPERATION

Turn the device on with the on/off switch.

The device displays:

The display now shows L°C with °C blinking, and then shows the water temperature, which is continuously updated every 10 seconds.

To save battery life, turn the device off after the bath. Otherwise, the device will automatically turn off.

Bath temperature reminder: If the measured value is above or below the specified temperature, the high or low temperature indication starts to flash.

NOTE:

36.0°C ~ 38.0°C is recommended for baby bath water.

TECHNICAL DATA

Measurement range:	0.0°C ~ 60.0°C
Resolution:	0.1°C
Measurement accuracy:	30°C ≤ T ≤ 50°C: +/-0.5°C, Other: +/- 1.0°C
Power supply:	1 x LR41
Battery life:	More than 24 months (30 minutes/day)
Low temperature indication	If the Temperature is lower than 35.0°C,  flashes.
High temperature indication	If the Temperature is higher than 38.0°C,  flashes.
Automatic Shutdown	after about 30 minutes
Dimensions	90 mm x 47 mm x 30 mm
Weight	approx. 40 g

CLEANING

Clean the bath thermometer if necessary with lukewarm water and a standard detergent. Do not clean in a dishwasher. Clean the sensor before storage.

BATTERY REPLACEMENT

If the symbol [] appears on the LCD display, a battery replacement is required. Proceed as follows:

1. Turn the device over.


2. Unscrew the screw.

3. Remove the cover.

4. Turn the device over and press on the display.

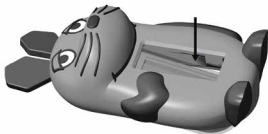5. Turn the device over and remove the module.

6. Pull the cap from the module.

7. Remove the old battery with a toothpick or other non-metallic tool and then insert a new battery (LR 41) according to the polarity. ("+" is on top, "-" is at the bottom).

8. Attach the cap to the module.

9. Reinsert the module. Once the main module is securely in place, a clicking sound will be heard.

10. Attach the cover.

11. Screw in the screw.

WARNING!

Dispose of the old battery carefully and out of the reach of small children. Swallowing the battery can be fatal. If the battery is swallowed, contact a hospital immediately to have it removed. Do not throw the battery into a fire. It may explode.

TROUBLESHOOTING

Err

The automatic detection is faulty. Please contact the dealer.

H

Please repeat the measurements.

L

Please repeat the measurements.

The device detects a low battery, please replace the battery.

EXPLANATION OF SYMBOLS

Follow the instruction manual

Manufacturer

MMYY: Batch number consisting of month (M) and year (Y) of manufacture

Type BF: Protection against electric shock, type of applied part, Body Floating.

REACH symbol: Registration, Evaluation, Authorisation, and Restriction of Chemicals (REACH) is a European Union regulation.

DISPOSAL

Dispose of the packaging materials properly according to their type, such as paper or cardboard. Contact your local waste disposal service or environmental authority for information on proper disposal.

Do not dispose of electronic devices in household waste! According to the European Directive 2012/19/EU on Waste Electrical and Electronic Equipment and its transposition into national law, used electrical equipment must be collected separately and recycled in an environmentally sound manner. Discharged old batteries and accumulators must be disposed of by consumers in battery collection containers. For information on the disposal of old devices or batteries produced after 01.06.2006, please contact your municipal waste disposal service provider or environmental agency. Batteries and rechargeable batteries may not be disposed of in

household waste, but you are legally obliged to return used batteries and rechargeable batteries. You can return the batteries after use either at our sales point or in the immediate vicinity (e.g., in retail or at municipal collection points) free of charge.

Batteries are marked with a crossed-out waste bin symbol and the chemical symbol of the contaminant: "Cd" stands for cadmium, "Hg" stands for mercury and "Pb" stands for lead.

Cd¹

Hg²

Pb³

¹ Battery contains cadmium

² Battery contains mercury

³ Battery contains lead

EC DECLARATION OF CONFORMITY

A "Declaration of conformity" in accordance with the applicable directives and corresponding standards has been prepared by Bresser GmbH. The full text of the EC declaration of conformity is available at the following Internet address:

www.bresser.de/download/9810402/CE/9810402_CE.pdf

WARRANTY & SERVICE

The regular warranty period is 2 years and begins on the day of purchase. To benefit from an extended voluntary warranty period as stated on the gift box, registration on our website is required. The full warranty conditions and information on warranty extension and service can be found at www.bresser.de/warranty_terms

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© I. Schmitt-Menzel
WDR mediagroup GmbH