Contents	
Englis	01-05
Deutsch	06-10
Francia	11-15
İspanol	16-20
İtalyan	21-25

Shenzhen Beiship Technology Co., Ltd.
Web: www.naipo.com
Email: beiship@naipo.com
Add: 303 Mantao Building 729 Minshi Avenue Minshi
Street Longhua District, Shenzhen, Guangdong 518001

CE FCC RoHS Made in China

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10. Keep the cord away from heated surface, otherwise it may cause serious burns. Do not use over sensitive skin areas or in the presence of poor circulation.

11. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

12. The appliance must only be supplied at safety extra low voltage (as per the marking on the appliance).

13. The appliance is only to be used with the power supply unit provided with the appliance.

14. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or persons lacking experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

15. Special attention should be given to components most susceptible to wear and it is advised to replace defective components immediately and/or keep the equipment out of use until fixed.

16. Use this appliance only for its intended use as described in this manual.

17. Never operate the appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if it has dropped into water. Return the appliance to a service center for examination and repair.


SAVE THESE INSTRUCTIONS

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using the product.

Only use indoors!

After using this product for 15 minutes, rest for 30 minutes, then use again.

Do not wash!
Do not dry clean!
After using this product for 15 minutes, rest for 30 minutes, then use again.

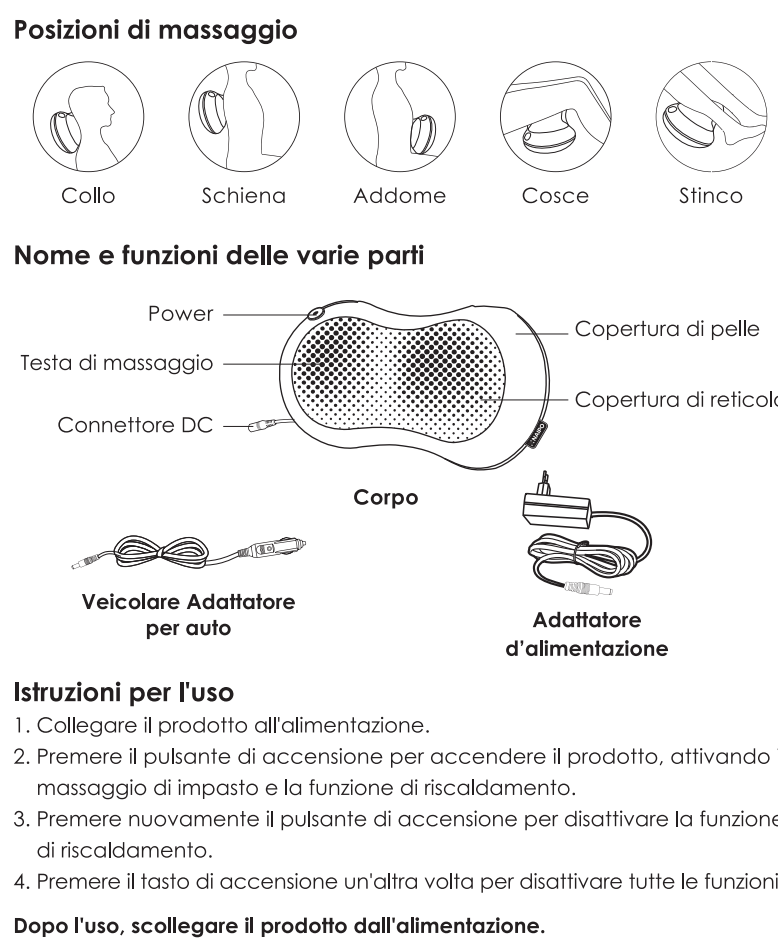Important Safety Instruction

When using the product, instructions and basic precautions should always be followed, including the following:

DANGER

To reduce the risk of electric shock:
Always unplug the product from the electrical outlet immediately after using and before cleaning.

WARNING

To reduce the risk of burns, fire, electric shock, or injury to persons:
01. The product should never be left unattended when plugged in. Unplug from outlet when it is not used and before putting on or taking off parts.
02. Do not operate under a blanket or a pillow. Excessive heat can occur and cause fire, electric shock, or injury to persons.
03. Do not carry this product by its supply cord or use the cord as a handle.
04. Never operate the product with the air openings blocked. Keep the air openings free of hair, dirt, and the like.
05. Do not drop or insert any object into any opening.
06. Do not use outdoors.
07. Do not operate where aerosol (spray) products are being used or where oxygen is administered.
08. To disconnect, switch off all controls and remove the plug from outlet.
09. Close supervision is necessary when this appliance is used by, or near children, invalids, or disabled persons to ensure that they do not misuse the product.

-11-

Product Maintenance

1. Before cleaning, please make sure that the product is off and that it's disconnected from the power supply. Please clean using dry cloth.
2. Keep the product clean and dry away from humid environment.

Attention

Please make sure that the product is not cleaned using benzene, thinner, or other aggressive detergents, otherwise the product color might fade.

Contact Us

If you need any assistance or support, please contact us. We will do our best to provide you with a viable solution.
Web: www.naipo.com
Email: support@naipo.com

Technology parameter

Nome: Massage Pillow
Modelo: MGP-129G
Rated Voltage: DC 12V ± 2A
Rated Power: 24W

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