

Warranty card

Warranty Card (Receipt)

Customer name: _____

Customer telephone: _____

Customer Address: _____

Purchase Machine Type and Number: _____

Maintenance record _____

Thank you for your suggestions and comments
(once adopted)

Company receipt

Product Qualification Certificate

Product name: _____

Product type: _____

Date of production: _____

Warranty card

Thank you very much for buying our products!
In order to protect your interests, our company promises as follows:This product can enjoy one-year quality guarantee from the date of purchase. The warranty period for power supply is three months.All of them can enjoy lifelong maintenance if they meet the following terms.

1. Damage caused by non-human factors, natural disasters, etc.

2. Machine damage caused by non-error operation.

3. The whole machine has not been disassembled, modified and assembled before the product is sent to the factory for maintenance.

4. Maintenance by customer with completed warranty card.

5. Once the warranty card is altered, the warranty will be invalid.

6. Normal wear and tear of non-core parts and loss of function.

Customer name: . _____

Date of purchase: _____

Product Type and Number: _____


Distributors: _____

Contact number: _____

NEATCELL Picosecond pen



Name of the relevant operating parts

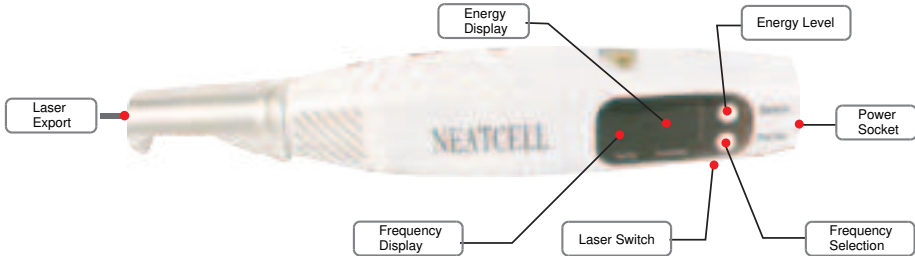

Light picosecond

Operation method

General usage and precautions

- 1.Plug in the power, and the machine is energized, and the digital tube displays
- 2.Press the frequency mode select button to select the frequency, then set the output energy size
- 3.After setting, press the laser emission switch to output energy
- 4.During operation, wear protective glasses
- 5.Keep the machine clean and avoid contact with other contaminants
- 6.No unauthorized operation of the hand
- 7.After use, please keep it properly after cleaning

Red light picosecond

After the normal display of the digital tube, press the frequency button, the current selection has three frequency modes,1.2.3 light mode,Through the energy increase button to adjust the size, press the switch after the setting, the machine will run out of light and operate outAfter that, the machine can be turned off

Blue light picosecond

After the power of the digital tube, the frequency button will be pressed, and the blue ray will have nine frequency modes1 to 9 out of the light mode ,Press the energy increment button to adjust the energy size, finish the setting and output the laser, and operate the closing machine, The blue light emits a laser every time.