Manufacturer: ZHONGSHAN CAMRY ELECTRONIC CO., LTD.

Address: Baishawan Industrial Park, Qiwan Road, East District,
ZHONGSHAN, GUANGDONG, CHINA, 528400
Email: am-contact@camry.com.cn

E-CrossStu GmbH
Mainzer Landstr.69, Frankfurt am Main, Germany 60329
Tel: +49 6933296767
Email: e-crossstu@outlook.com


YH Consulting Limited
C/O YH Consulting Limited Office 147, Centurion House,
London Road, Staines-upon-Thames, Surrey, London,
TW18 4AX
Tel: +44 07514-677868 Email: H2YHUK@gmail.com

This device complies with part 15 of the FCC Rules. Operation is subject
to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including
interference that may cause undesired operation

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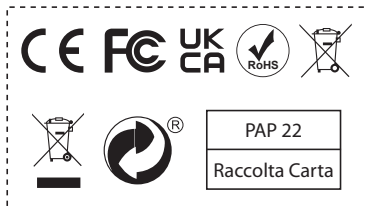Important Safety Instructions

SAVE THESE INSTRUCTIONS

READ ALL INSTRUCTIONS CAREFULLY BEFORE USE.

- 1.If the scale malfunctions, first check the battery power. Replace the batteries if necessary.
- 2.Check the scale before each use. Do not use the scale if damaged. The continuous use of a damaged scale may cause inaccurate results or injury.
- 3.Please keep the scale dry. Take caution when using the scale on wet and slippery surfaces.
- 4.The scale is fragile. Please handle it with care.
- 5.This scale is a high-precision measuring device. Never jump on, stomp on, or disassemble the scale. Handle it carefully to avoid breakage.
- 6.The scale is not waterproof; never immerse the scale in water. To clean the surface, use a damp cloth or glass cleaner. Do not use soap or other chemical cleaning agents.
- 7.Do not expose the scale or batteries to excessive temperatures.

BATTERIES

- 1.Clean the battery contacts and also those of the scale prior to battery installation.
- 2.Always select the correct size and type of battery most suitable for the intended use. Do not mix new and used batteries or batteries of different types or brands.
- 3.Always insert batteries correctly with regard to polarity (+ and -) marked on the battery and scale.
- 4.If "Lo" appears on the display, please replace the batteries and make sure to replace all batteries of a set at the same time.
- 5.Do not short-circuit, charge or force discharge the batteries.
- 6.Remove batteries from the scale if it is not to be used for an extended period of time unless it is for emergency purposes.
- 7.Keep batteries out of the reach of children. In case of ingestion of a battery, seek medical assistance promptly.
- 8.Exhausted batteries should be immediately removed from the scale and properly disposed of. When discharged batteries are kept in the scale for a long time, electrolyte leakage can occur causing damage to the scale and/or personal injury.
- 9.Do not heat, deform, dismantle, encapsulate and/or modify batteries.
- 10.Do not weld or solder directly to batteries. Do not dispose of batteries in fire.
- 11.Store unused batteries in their original packaging away from metal objects. If already unpacked, do not mix or jumble batteries.

This product is RoHS compliant. This product is in compliance with Directive 2011/65/EU, and its amendments, on the restriction of the use of certain hazardous substances in electrical and electronic equipment.

This symbol means the product must not be discarded as household waste, and should be delivered to an appropriate collection facility for recycling. Proper disposal and recycling helps protect natural resources, human health and the environment. For more information on disposal and recycling of this product, contact your local municipality, disposal service, or the shop where you bought this product.

Instructions for Users on Removal, Recycling and Disposal of Used Batteries

To remove the batteries from your equipment or remote control, reverse the procedure described in the owner's manual for inserting batteries. For products with a built-in battery that lasts for the lifetime of the product removal may not be possible for the user. In this case, recycling or recovery centers handle the dismantling of the product and the removal of the battery. If, for any reason, it becomes necessary to replace such a battery, this procedure must be performed by authorized service centers. In the European Union and other locations, it is illegal to dispose of any battery with household trash. All batteries must be disposed of in an environmentally sound manner. Contact your local waste-management officials for information regarding the environmentally sound collection, recycling and disposal of used batteries.

WARNING: Danger of explosion if battery is incorrectly replaced. To reduce risk of fire, explosion or leakage of flammable liquid/gas, don't disassemble, crush, puncture short external contacts, expose to temperature above 60 °C (140 °F), sunshine or like, expose to extremely low air pressure or dispose of in fire or water. Replace only with specified batteries. The symbol indicating 'separate collection' for all batteries and accumulators shall be the crossed-out wheeled bin shown below:

According to EU Batteries Directive 2013/56/EU, any batteries, accumulators and button cells containing more than 0.0005 % mercury, more than 0.002% cadmium or more than 0.004 % lead, shall be marked with the chemical symbol for the metal concerned: Hg, Cd or Pb respectively. Please Refer to the symbol below:

Hg


Cd

Pb

Download and set up the Application

- Go to the "App Store" or "Google Play", and search "SenssunLife" to download the Application.
- Or in another most convenient way by scanning the QR code below

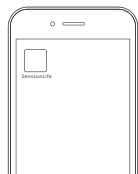MovingLife

Connect the scale with your phone

Open the App and turn on the Bluetooth of your phone, sign up your account and follow the guideline to set up your profile and connection step by step. Being aware of it, the scale has to be on when you start to connect with it. All measurements you taken can be saved in history for detail review anytime anywhere by getting access to your own account.


Weight unit conversion

You will find the unit conversion button at the bottom of scale, press it to convert the weight unit, you will see it changed on LCD.

Start Measurement

For the first measurement, keep the scale connected with the App, then put it on a solid and flat floor, and step on it for the measurement. Once the weight is stable and flashing, fat measurement will be proceeded. You will review your weight and body fat display on the screen for several time, then the scale will be off automatically.

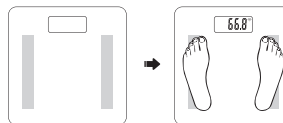- Do not put the scale on a carpet or other soft surface, it will lead to a wrong weight measurement.
- The first measurement has to be taken with Bluetooth connection, the first weight data will be reference for next auto user recognition.
- To make a proper measurement of body fat, please keep your bare feet on the electrodes of scale. The scale platform should be clean and dry.
- If you want to review all measurement data, open the App, and check them in detail interface.

Warning indications

Low Fat% Indication
The Fat% is too low. Please lead a more nutritional diet and take good care of yourself.

High Fat% Indication
The Fat% is too high. Please watch your diet and do more exercise.

Calibration
Error exists, please re-test to get the correct result.

Over-load Indication
The weighing subject on the platform exceeds the maximum capability of scale. Please step off to avoid damage.

Low Battery Indication
The battery power is running low, please replace with a new battery.

WARRANTY

- This product is covered by the SENSSUN Product and Warranty for 12 months from the date of its original purchase.
- If manufacturing defects occur, please contact us immediately to launch a warranty claim via:
Email: **support@camry.com.cn**
We will instruct you on how to return the defective unit back to us for repair or replacement.
- No return will be accepted without seller authorization.

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