
RENPHO

USER MANUAL

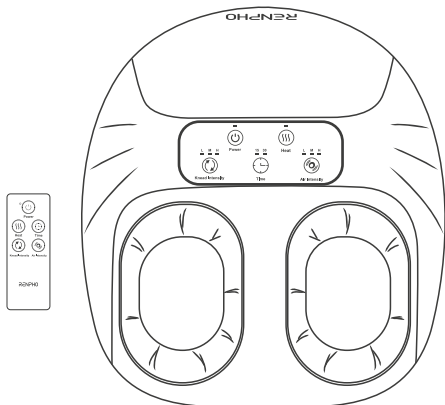SHIATSU FOOT MASSAGER

Model: RP-FM069R

CONTENTS

ENGLISH

01 ~ 11

DEUTSCH

12 ~ 22

FRANÇAIS

23 ~ 33

ITALIANO

34 ~ 44

ESPAÑOL

45 ~ 55

Thank you for your purchase of RENPHO Foot Massager. The massager offers a powerful massage. While enjoying massager, you should use it with care and safety. To ensure safety and avoid any property damage, please read this manual carefully and use as indicated. Non-compliance may cause serious injuries or severe property damage.

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. **WARNING** – To reduce the risk of burns, fire, electric shock, or injury to persons.
 - a) This appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
 - b) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
 - c) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
 - d) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 - e) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
 - f) Do not carry this appliance by supply cord or use cord as a handle.

- g) Keep the cord away from heated surfaces.
- h) Never operate the appliance with the air openings blocked.
Keep the air openings free of lint, hair, and the like.
- i) Never drop or insert any object into any opening.
- j) Do not use outdoors.
- k) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- l) To disconnect, turn all controls to the off position, then remove plug from outlet.
- 4. **DANGER:** Do not use while bathing or in a shower.
- 5. **DANGER:** Do not place or store appliance where it can fall or be pulled into a tub or sink.
- 6. **DANGER:** Do not place in or drop into water or other liquid.
- 7. **WARNING:** Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
- 8. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- 9. **WARNING** - Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- 10. **WARNING** - Do not stand on or in appliance. Use only while seated.
- 11. **DANGER** – To reduce the risk of electric shock. The appliance should never be left unattended when plugged in. Unplug from the outlet when not in use, and before cleaning, putting on or taking off parts.
- 12. **WARNING:** Do not carry this appliance by supply cord or use cord as a handle.
- 13. **WARNING:** Keep the cord away from heated surfaces.
- 14. **DANGER:** Do not reach for an appliance that has fallen into water. Unplug immediately.
- 15. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- 16. The switch or controlled unit must not be allowed to get wet

during washing and that during drying the cord must be positioned to ensure that water does not flow into the switch or control.

17. **SERVICING OF DOUBLE-INSULATED PRODUCTS.** In a double-insulated appliance, two systems of insulation are provided instead of grounding. No grounding means is provided on a double-insulated product, nor should a means for grounding be added to the appliance. Servicing a double-insulated product requires extreme care and knowledge of the system, and should be done only by qualified service personnel. Replacement parts for a double-insulated appliance must be identical to the parts they replace.
18. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.

About the Battery

1. Do not expose the appliance or battery to excessive temperatures.
2. Replace all batteries of a set at the same time.
3. Different types of batteries or new and used batteries are not to be mixed.
4. Exhausted batteries are to be removed from the appliance and safely disposed of.
5. If the appliance is to be stored unused for a long period, the batteries should be removed.
6. Do not use modified or damaged batteries.
7. Be aware of the risk of terminals of the battery-operated appliance or battery being short-circuited by metal objects.
8. This appliance contains non-rechargeable batteries, these batteries are not to be recharged.
9. Always insert batteries correctly with regard to the polarities (+ and-) marked on the battery.
10. Do not short-circuit batteries.
11. Keep batteries out of the reach of children.

12. Do not charge batteries.
13. Do not force discharge batteries.
14. Do not weld or solder directly to batteries.
15. Do not dismantle batteries.
16. Do not deform batteries.
17. Do not dispose of batteries in fire.
18. Do not allow children to replace batteries without adult supervision.
19. Store unused batteries in their original packaging away from metal objects. If already unpacked, do not mix or jumble batteries.

WARNING

DANGER – To reduce the risk of electric shock.

1. Always unplug the device after use and before cleaning.
2. Never use pins or other metal fasteners with the device.
3. Never unzip the device.

WARNING – To reduce the risk of fire, electric shock or injury to persons.

1. Please consult your doctor before use if you are a user of the following types.

- Children under 9 years old
- Pregnant or postpartum women
- Persons with a pacemaker or implanted medical devices
- Persons suffering from diabetes or other sensory impairment
- Persons with open wounds, discolored areas or any part of the body that is swollen, burned, inflamed or where skin eruptions or sores are present
- Persons with increased or reduced sensitivity to heat
- Persons suffering from osteoporosis, weak bones, abnormal or curved spine or other acute pains
- Persons with injured massage areas
- Persons having received cardiac or visceral surgery
- Persons easy to get sprained or injured
- Persons feeling any discomfort (e.g. inflamed skin, redness, itching, etc.)

2. Do not use the device outdoors.
3. Do not use the device if the covers are detached.
4. Do not use the device with wet hands or feet.
5. Do not allow pets or animals to lie or chew on the device and its power cord.
6. Never reach or touch the device that has fallen into the water. Unplug the device immediately.
7. Position the power cord properly to avoid cable-related injuries.

SAVE THESE INSTRUCTIONS

ABOUT THE RENPHO FOOT MASSAGER

1. Specifications

Model: RP-FM069R

Rated Voltage: AC 120V 60HZ (US&CA)

AC 220-240V 50/60Hz (EU&UK&AU)

Rated Power: 40W

Heating Area: Bottom

Time: 15mins/30mins

Material: ABS + Lycra fabric

Heating Temperature: Below 50°C/122°F

2. What's in the box

Foot Massager x 1

AAA Battery x 2

Remote control x 1

3. Product Introduction


Control panel

Heating on bottom

Detachable
foot cover
with zipper

Remote control

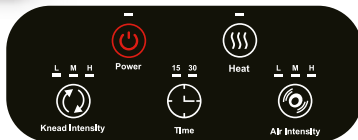4. Function

Power: Press the Power button to turn on/off.

Heat: Press to enable / disable heat function.

Knead Intensity: Press to adjust knead intensity among low, medium, high or off.

Timer: Press to select a 15 / 30 min auto-off timer.

Air Intensity: Press to adjust air pressure intensity: low, medium, high or off.

USING THE RENPHO FOOT MASSAGER

Enjoying a foot massage for 15 or 30 minutes once or twice daily can relieve muscle pain, improve blood circulation. And we recommend you limit massaging to 15-30minutes. Longer massage time may cause mild pain.

1. Connect the plug to the power source. The Power indicator will flash and the massager enters standby mode.
2. Please remove the protective film from the touch panel before use. Press the Power button. The knead and air intensity will be at default low level (L). The default auto shut-off will be 15 minutes. The heat function is turned on by default.
3. Press the Time button to change the auto shut-off to 30 minutes. If you press the Time button again during the 15-minute timer, the 30-minute timer will restart without counting the previous time in.
4. Press the Heat button to turn on/off the heat function.

TIPS: Use your thumb to press buttons to get better touch sensitivity.

5. The foot massager fits up to US Men's size 14 (EU Men's size 48).
6. To remove and wash the foot cover:

PRODUCT CARE AND MAINTENANCE

1. It is prohibited for non-professional personnel to repair the massager or replace parts.
2. Do not use or store it in wet environments or high temperatures.
3. Unplug the power plug from the outlet after using.
4. If high temperature or foul odor is experienced while using the product, stop using immediately and unplug the device.
5. Please keep the product clean. You can use a soft cotton cloth and neutral detergent to clean it.
6. The product possesses a temperature-protection device; It will automatically shut-off when the temperature exceeds the safe limit.

TROUBLESHOOTING

Common problems	Cause	Solutions
The appliance is not working.	1. Disconnected power or loose plug, causing poor contact. 2. The power switch is off or the timer has been reset. 3. The motor automatically shuts off due to the device being on too long, activating overheating protection.	1. Re-plug in the appliance. 2. Turn the power switch on again. 3. Use the appliance after the motor cools down.
The massage intensity is too high or too low.	1. Your foot is positioned at the wrong angle. 2. Suitable intensity level is not selected.	1. Adjust the contact angle to get the moderate intensity. 2. The appliance has 3 intensity levels. Choose the mode you prefer. We suggest starting with the low intensity for a gentle introduction.
No heat or insufficient heat.	The heating function is manufactured according to international safety standards, which means it may take a bit longer, around 7-10 minutes, before you start feeling the heat.	If you press the foot cover to see the light, the heating is turned normally. The rate of heating is related to the environment. If you are in a cold area, the heating time may take longer.
The touch panel is not working properly.	1. The plastic film on the touch panel is not removed. 2. Fingers are wet. 3. Insufficient touch area.	1. Remove the plastic film before use. 2. Dry your fingers before touching the button. 3. Use your thumb to press buttons to get better touch sensitivity.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

TEL: +1(844) 417 0149 (US&CA)

Monday-Friday 9:00AM-4:30PM

1800 MY RENPHO (1800 69 7367) (AU ONLY)

Monday-Friday 9:00AM-6:00PM (AEST)

+44(0)7434 666088 (UK ONLY)

Monday-Friday 9:00AM-5:00PM (GMT)

Email: support@renpho.com (US&CA)

support-au@renpho.com (AU)

support-eu@renpho.com (UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a retail store or other collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2014/35/EU, Directive 2011/65/EU, The Electrical Equipment (Safety) Regulations 2016, The Electromagnetic Compatibility Regulations 2016, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:
<https://renpho.eu/pages/compliance-data>

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