

# waterpik® waterflosser™

## QUICK START GUIDE

WP-660/662/663/670/672/  
673/674/675/676

Register your product and  
receive over \$50 in savings.  
[waterpik.com/registration](http://waterpik.com/registration)

Visit [waterpik.com/HowTo660](http://waterpik.com/HowTo660)  
to watch How To Video


### 1 Unit Prep



**A.** Fill reservoir with warm water and replace it firmly on the base.




**B.** Insert tip until it clicks into place.

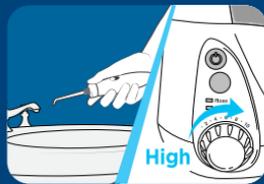

**C.** Turn pressure to high, point tip into sink, turn power and water switches ON until water flows. Turn switches OFF.


### 2 Water Floss

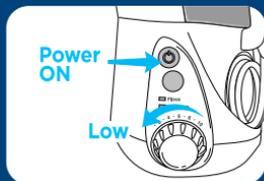

**A.** Start on low pressure and increase to your liking. Turn power ON.



**B.** Turn water ON. Close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink.

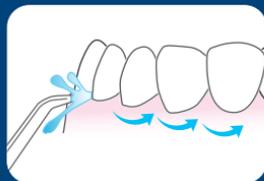

**C.** Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

**TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.**

# waterpik® waterflosser®

## GUIDE DE DÉMARRAGE RAPIDE

WP-660/662/663/670/672/  
673/674/675/676

Enregistrez votre produit  
et bénéficiez de plus de  
50 \$ d'économies.  
[waterpik.com/registration](http://waterpik.com/registration)  
Veuillez visiter le site Web  
[waterpik.com/HowTo660](http://waterpik.com/HowTo660) pour  
visionner la vidéo d'instructions

### 1 Préparation de l'appareil



**A.** Remplissez le réservoir avec de l'eau tiède et remettez-le fermement en place sur son socle.




**B.** Insérez l'embout jusqu'à ce qu'il s'emboîte.



**C.** Réglez la pression à un niveau élevé, orientez l'embout vers l'évier, mettez l'appareil sous tension et réglez l'interrupteur d'alimentation en eau sur ON (marche) jusqu'à ce que l'eau s'écoule. Mettez les interrupteurs en position OFF (arrêt).


### 2 Utilisation de l'hydropulseur

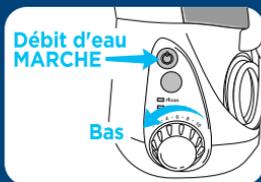

**A.** Démarrez sur le réglage le plus bas et augmentez la pression selon vos préférences. Allumez l'appareil.



**B.** Réglez l'interrupteur d'alimentation en eau sur la position ON (marche). Gardez vos lèvres suffisamment fermées pour empêcher les éclaboussures, tout en laissant l'eau s'écouler de votre bouche dans le lavabo.

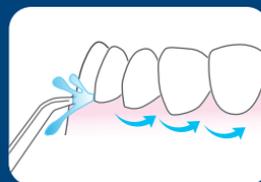

**C.** Dirigez le jet d'eau vers la ligne gingivale à un angle de 90°. Suivez la ligne gingivale et faites de courtes pauses entre les dents.

**POUR ÉVITER LES ÉCLABOUSSURES : GARDEZ L'EMBOUIT DANS LA BOUCHE LORSQUE L'APPAREIL EST EN MARCHÉ.**