
SonicBrush Slim20

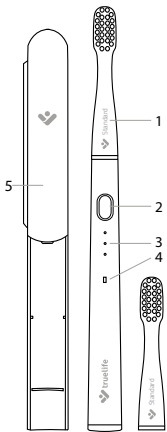EN – INSTRUCTIONS

PACKAGE CONTENTS

- TrueLife SonicBrush Slim20 brush
- Travel case
- 2x Standard heads
- USB-C cable
- Instructions
- Safety instructions

PRODUCT DESCRIPTION

1. Interchangeable head
2. ON / OFF button and change of modes
3. Modes indicator
4. Charging indicator
5. Travel case

TOOTHBRUSHING

1. For maximum efficiency, divide the cleaning into 4 sections and devote 30 seconds to each of them. After 30 seconds, the brush will take a short pause and alert you to move to the next section. After 2 minutes, the brush turns itself off. Cleanse at least twice a day, in the morning and in the evening.
2. Firmly fit the brush head onto the plastic shaft with the bristles facing the front of the handle.
3. Moisten the head, apply toothpaste and rub it on the teeth.
4. Switch on the brush with the ON / OFF button. To select another mode, press the button again within 5 seconds of power on and select one of the three modes. If you would like to end cleaning earlier, press the ON / OFF button.
5. In each section, start by cleaning the outside of the teeth, then the inside of the teeth, and finally clean the surfaces. Hold the brush at an angle so that the bristles gently touch the teeth and the edges of the gums at the same time. When cleaning, do not press the brush and do not make circular movements like with a manual brush, but move slowly from one tooth to another.
6. Rinse the brush head and body under running water and then wipe them with a cloth or towel.

CLEANING MODES

SLOW – For cleaning sensitive teeth and gums

MEDIUM – For general cleaning and plaque removal

FAST – For removing stubborn dental plaque and whitening teeth

REPLACEMENT OF HEADS

We recommend changing the heads every 3 months. After this period of use, wear occurs where the frayed and bent fibers have a reduced ability to remove dental plaque and bacteria accumulate on them.

CHARGING

1. When the brush is almost discharged, the charging indicator will flash red quickly.
2. Open the rubber cover on the bottom of the brush and connect the USB-C cable to the charging USB-C connector, the other end to the adapter and the adapter to the electrical socket. You can also charge the brush from a computer or power bank. Charge the sonic brush using the included USB-C cable. Do not use fast chargers for charging.
3. The indicator flashes red while charging. Once the brush is charged, the charging indicator lights up green. Charging time is 4 hours.

SPECIFICATIONS

Battery life	30 days (2x daily 2 minutes)
--------------	------------------------------

Charging	USB-C cable
Number of moves	70,000 movements/min
Modes	SLOW, MEDIUM, FAST
Timer	2 min (30s interval)
Memory	Yes
Additional heads	Slim-series
Brush dimensions	16 x 16 x 20.3 cm
Brush weight	33 g
Waterproof	IPX7
Input voltage	5 V – 1 A
Battery	320 mAh, Li-ion
Charging time	4 hours

MANUFACTURER'S RESPONSIBILITY

TrueLife products are warranted against manufacturing defects. In the event of such defects under warranty, elem6 s.r.o. shall repair or replace the defective part or product with a comparable part or product at its option. Rebuilding, altering or otherwise tampering with the product shall be deemed to exceed the scope of this warranty. This is the case even if the product is lost or damaged as a result of a defect, negligence or other error in the product. This warranty does not cover damage resulting from use, accident or normal wear and tear. elem6 s.r.o. assumes no liability for any accident, injury, death, loss or other claim related to or arising out of the use of this product. In no event shall elem6 s.r.o. be liable for incidental or consequential damages related to or arising out of the use of this product or any of its components. The return or replacement of parts and/or products may be subject to shipping, handling, replacement and/or restocking charges.

SAFETY WHEN WORKING WITH THE DEVICE

1. Read the safety information leaflet carefully.
2. Read the current version of the manual, which is available for download at: www.truelife.eu.
3. Printing errors and changes to information are subject to change.

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