



SHARPER IMAGE®

AIR COMPRESSION BOOTS

Please read all instructions and warnings prior to use. Keep this manual as it contains important information for future reference.

WHAT YOU GET

- a. Air Compression Boots
- b. Control Unit
- c. Air Hose Extension
- d. Power Adapter with DC Jack


a.



b.



c.

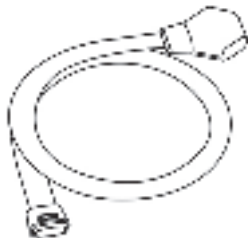

d.



• Accessory type and quantity may vary depending on the model. •

• The images shown in this manual are for illustrative purposes only and may be subject to change. •

THE FUTURE OF RECOVERY

01

Plug the barrel jack on the end of the power cord into the input port on the control unit.

Plug the power adapter into a wall outlet.

Note: Do not use the power adapter to charge or power other devices.

Hose Release
Button



02

Locate the foot or feet you want to use for your massage. The foot is/are labeled L (left leg) and R (right leg).

Find a comfortable place to sit or recline with your legs extended for the duration of your massage and remove any footwear (wearing socks is recommended). Slide on each boot over the corresponding leg and then zip the boot up all the way to the top of the zipper track.

Once the zipper is in place and fully up, the zipper takes down fully fits against the zipper track. The tabs should point towards the base of the boots. Failure to properly lock the zipper may cause the boots to open during inflation, reducing compression.

IMPORTANT: Never attempt to stand or walk while wearing the boots. The boots are intended for massage only. They are not designed for stability or mobility. Using the boots for purposes other than their intended use may damage the boots or result in personal injury.

03


IMPORTANT: Make sure the hose release is free of lint, debris or when connecting each hose, press each connector into each connection port slowly and firmly until it clicks securely into place.

To use both ports at the same time, plug the end of the air hose extension with a single connector into the air hose connection port on the control unit. Next, connect each end of the hose to the connector (one port on the other end of the air hose extension).


To use a single port, connect the boot air hose directly to the control unit connector. Using the air hose extension with using the air hose extension with a single port, seal the unused connection port with the rubber plug located on the bottom of the connector (bottom of the top side features the Sharper Image logo).

The air hose can be disconnected by pressing each hose release button while gently pulling the connector away from the connection port. Always pull the plastic connector gently. Never pull or push using great force.


Air Hose Extension (Top View)

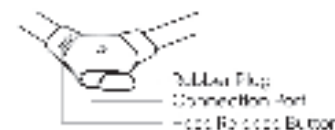

Air Hose Extension (Bottom View)

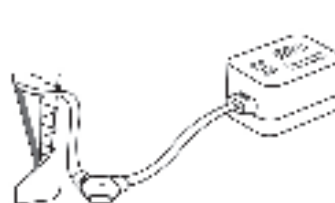

Single-Foot Usage

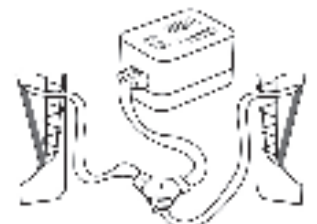

Both-Feet Usage

To power on the massage, press and hold the power button. To power off the massage, press and hold the power button again.

Quick press the timer button after any massage starts, pressing the timer button will increase the massage time by fixed 10 mins, 20 mins, return to initial time.

Press the mode button to cycle through the 3 massage modes.

- Full Relaxation – all zones warmness and then relaxation
- Good all – zones each compress sequentially and then all zones decompress
- Sequence – zones each compress and decompress sequentially
- Zones – select up to 8 zones for targeted pressure massage
- Manual – select up to 8 zones and manually adjust compressions


In zones mode or manual mode, make your compression force choices using the arrow and adjust buttons.

After all mode selections have been made, press and hold the select button to start your massage.

Quick press the manual select button to cycle through the 5 compression levels.

Quick press the power button during a massage to begin decompress on. The massage will automatically power off after 5 minutes and the body will be fully decompressed.

Quick press the decompress button and hold the power button to turn off the device.

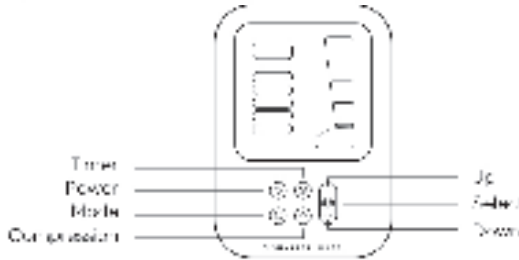

NOTE: The images shown here are for illustrative purposes only and may be subject to change, the actual descriptions of control panel and their locations may vary depending on the model.

CARE INSTRUCTIONS

- Spot clean only. Wipe up and wipe with a damp cloth.
- Avoid use where a plasma is visible.
- Avoid touching before connecting to power source.
- Do not submerge in water.
- Do not use detergent or soap when cleaning.
- Hang the body in the flat for long-term storage.
- Never touch or wear the body.
- Do not clean or directly bend any of the needed hoses.

SAFETY WARNINGS

PLEASE READ BEFORE USE

- Read the instruction book carefully before use.
- As with any massage, bruising may occur. Notify a medical level and check out your use of this massage machine immediately if you feel any unusual physical discomforts.
- If you experience numbness, tingling, or any other sensation, stop your massage immediately. Seek a doctor and medical advice if these symptoms persist.
- Do not fall asleep or fall asleep during the massage.
- Do not use attachments or any other device that may be included.
- Do not leave this massage unattended when it is plugged into a power source.
- Do not use this massage machine for medical purposes or treatments.
- Keep away from sources of water and moisture or fumes.
- Do not touch the power cable or parts with wet hands.
- Do not touch this item by the power cable, hose, or attachments.
- Do not use it for long.
- Thoroughly inspect product before every use. If the product appears worn, frayed, splintered, cracked, or broken in any manner, discontinue use and discard immediately.
- This product is not intended to treat, cure, or prevent a medical condition.
- Do not use this item if you are pregnant or have an injury or illness. Consult your physician if you have a medical condition that may prevent safe use of this item.

WARNING: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FCC STATEMENT

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

SHARP IMAGE NAME AND LOGOS (INCLUDING THE 77 LOGO, TOMORROW'S TOMORROW, AND POWERBOOST) ARE TRADEMARKS OR REGISTERED TRADEMARKS. DISTRIBUTED BY AND © 2022 MERCHSOURCE, LLC. ALL RIGHTS RESERVED.

PRINTED IN MEXICO
SC:12-16-2024 / Elbrd: J.M
OBPN:1018195
SHNMDL: 1018195_B
ADD: 1018196

8 1 4 1 0 0 3 0 6

