PROJECT E
BEAUTY

LightAura

LED FACE MASK

USER GUIDE

Model Number: PE021

MULTILINGUAL INSTRUCTIONS

MULTIPLE LANGUAGES USER MANUAL

For other languages (French, Spanish, Italian, German),
scan the following QR code

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WARNING

- Discontinue use if you experience any discomfort or persistent skin redness lasting more than two hours. It could indicate light sensitivity or another condition. Consult your healthcare provider for advice.
- Never operate the device if the cord or plug is damaged, if it is not functioning properly, or if it has been dropped, damaged, or submerged in water. Please return the device to Project E Beauty for examination and repair.
- Ensure the device is securely positioned to prevent accidental falls or being pulled into a tub or sink. Maintain a safe distance from water sources at all times, even during bathing or showering. Should the device accidentally come into contact with water, promptly unplug it to prevent any potential hazards.
- Keep the device away from flammable aerosol products.
- Use only the provided AC power supplies or chargers specifically designed for your LED mask. Using incompatible power sources may cause damage to the device.
- If the mask starts to overheat, turn it off immediately.

CONTRADICTIONS:

DO NOT USE the device if any of the following apply to you

1. Pregnant - Do not use if you are pregnant or suspect that you might be. The device has not been tested on pregnant individuals, and the effects, if any, are unknown.
2. Migraine - Avoid use if you have experienced a migraine within the last two years, as the lights may increase eye sensitivity.
3. Seizures - Do not use if you have a history of seizures triggered by light.
4. Lesions - Refrain from using this device over any suspicious or cancerous lesions. Doing so may delay proper medical treatment.

Failure to use and maintain the LED mask in accordance with the instructions outlined in this manual will invalidate the warranty.

BENEFITS

Each color of LED light serves a specific function and can be easily selected using the control box.

Red Light (630nm): The Anti Ager

Collagen synthesis reduces fine lines and wrinkles for a more elastic and youthful complexion

Blue Light (460nm): The Acne Savior

Blue light creates free radicals that destroy acne-causing bacteria. This process prevents breakouts and regulates oil production

Green Light (520nm): The Brightener

Fades dark spots, combats hyperpigmentation, and balances skin tone by significantly reducing melanin production

Yellow Light (630+520nm): The Healer

Stimulates the production of red blood cells to calm sensitive skin, reduces redness and remove toxins from the skin

Purple Light (460+630nm): The Multi-Tasker

Combines the benefits of blue and red light with the added benefit of targeting uneven skin tone and inflammation

Cyan Light (460+520nm): The Calmer (Blue and Green)


Purifies, soothes and heals the skin. Reduces appearance of capillaries. Beneficial for the treatment of rosacea and other skin conditions accompanied by redness

White Light (~mix): The Balancer

White light integrates all wavelengths on visible spectrum and therefore penetrates deepest into the skin. Powerful combination of light improves overall skin health for total rejuvenation

WHAT'S IN THE BOX?

- Main unit x 1
- Control Box x 1
- Head Strap x 1
- Power Cord x 1
- User Manual x 1

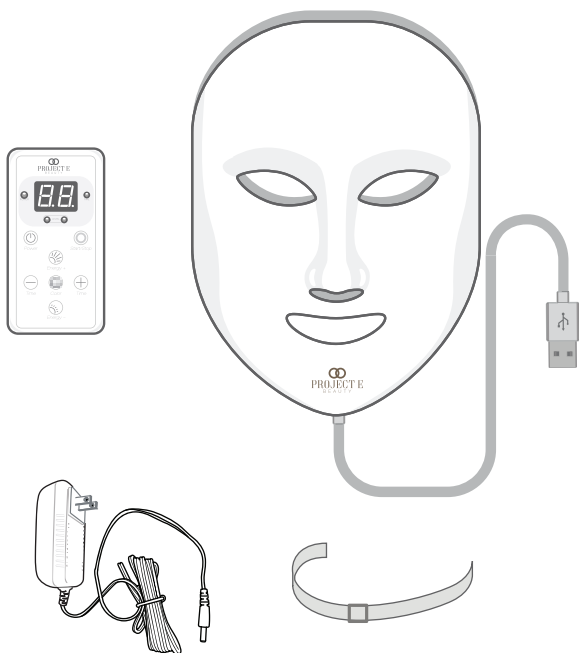LIGHT SENSITIVITY TEST

Recommended before use

To determine if you are sensitive to light, please follow the instructions below. It is advised to test a small, less sensitive area of skin, such as your forearm.

Press the "Power" button to turn on the mask and then press the "Start/Stop" button to activate the mask LEDs. Use the "Color" button to switch to different colors.

Place the device on an area of your forearm for a 3-minute treatment.

Wait for 24 hours. You should NOT experience any adverse reactions during this time.

A negative reaction to light may include pain, discomfort in the treatment areas, or any redness/dryness of the skin lasting more than two hours. This reaction can be attributed to medication, supplements, or a potential skin issue.

Eye Safety

The design of the lights and the mask intentionally excludes the eye area for safety reasons.


For a comfortable experience, we recommend starting with the low energy mode and gradually increasing the energy level as you become accustomed to the treatment.

Please note that it is NOT recommended to directly shine the light source into the eyes or the eyes of others for an extended period of time as this can cause sensitivity.

HOW TO USE

Follow these steps for proper operation during a treatment:

1. Insert the head strap into each side of the LightAura mask to ensure comfortable wear on the face. Adjust the strap for a customized fit, ensuring the mask is not worn too tightly against the face.

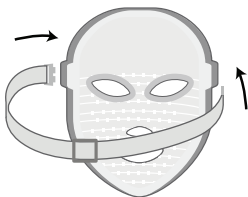2. Connect the USB cord to the jack on the control box. If the USB plug doesn't fit, attempt to insert it in the opposite direction carefully.

3. If needed, disinfect the surface of the mask.
Ensure your face is clean before using the mask.


4. Turn on the unit by pressing the Power button. The green indicator LED should light up.

5. Activate the mask LEDs by pressing the Start/Stop button.

6. Use the color button on the control box to choose your preferred color setting for the mask.

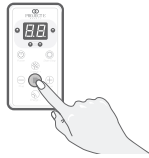7. Adjust the power level using the Energy +/- buttons, with the maximum level being five (5).

8. Set the desired duration for the mask using the duration +/- buttons, with a maximum duration of 60 minutes.

9. Place the mask over your face and secure it with the strap. Adjust as necessary for comfort.

10. Avoid looking directly into the LEDs. While the powerful bulbs are safe for skin use, they can cause sensitivity to the eyes if stared at for an extended period of time.

It is recommended not to use the mask for more than 20 minutes per session, with the suggested operating time being between 10 and 15 minutes. Start with one 10-minute treatment, four times a week, unless otherwise directed by a doctor. Once you have achieved your initial goals, you may reduce usage to twice a week.

The mask will automatically shut off once it reaches the set time limit. To turn the device off before the chosen time setting, use the ON/OFF switch on the control box. If you have sensitive skin, start with shorter periods of time and gradually increase. Always prioritize your comfort and never exceed your own comfort level.

MAINTENANCE

1. When unplugging the device, ensure a safe removal by firmly holding the plug and gently disconnecting it from the AC outlet.
2. After each use, diligently cleanse the mask surface with a soft cloth to maintain hygiene.
3. Promptly dry the surface if it encounters any moisture, chemicals, or water to prevent damage.
4. Maintain a clean and sterile environment for the device at all times.
5. Store the device in a clean and dry location to preserve its functionality.

TECHNICAL SPECIFICATIONS

Model no.	PE021
Red Light Photon Wavelength	630 nm $\pm$ 5
Blue Light Photon Wavelength	460 nm $\pm$ 5
Green Light Photon Wavelength	520 nm $\pm$ 5
Purple Light Photon Wavelength	460+630 nm $\pm$ 5
Cyan Light Photon Wavelength	460+520 nm $\pm$ 5
Yellow Light Photon Wavelength	630+520 nm $\pm$ 5
White Light Photon Wavelength	460+630+520nm $\pm$ 5
Output Power	Max 17.8W
Voltage	100V - 240V, 60Hz/50Hz

1 YEAR WARRANTY

ACTIVATE YOUR 1 YEAR FREE WARRANTY NOW.

The product comes with a full one year warranty against manufacturer defects. It does not include damage caused by misuse, accident, abuse, neglect or other external causes. The warranty is void if an unstable voltage or inappropriate power supply has been used.

All we need is your name, email and order ID. That's it!

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