

MEDTECH[®]
Caring for life

MPV 1 Manipol Massager



MEDTECH MPV 1 Manipol Massager User Manual

[Home](#) » [Medtech](#) » MEDTECH MPV 1 Manipol Massager User Manual 

Contents

- 1 MEDTECH MPV 1 Manipol Massager
- 2 Product Information
- 3 Product Usage Instructions
- 4 Frequently Asked Questions
- 5 INTRODUCTION
- 6 Massage principle
- 7 Features
- 8 Usage
- 9 Product Part Introduction
- 10 How to use
- 11 Massage Heads
- 12 Repairing and maintenance
- 13 Technical parameters
- 14 Warranty conditions
- 15 WARRANTY
- 16 Customer Service
- 17 Documents / Resources
 - 17.1 References

MEDTECH[®] 
Caring for life



Product Information

Specifications

- **Brand:** Manipal
- **Model:** Massager MPV1
- **Size:** 105 mm (L) x 136 mm (H)
- **Colour:** Frosting
- **Warranty:** Refer to warranty conditions

Product Usage Instructions

Powering On and Off

- To turn on the Manipol Massager MPV1, press the ON/OFF button.
- Press the same button again to turn it off.

Adjusting Speed and Modes

- Use the Speed + and Speed – buttons to adjust the massage intensity. Navigate through different modes by pressing the Modes button.

Protective Mesh Cover

- Ensure the protective mesh cover is in place before using the massager to prevent any accidents or damage.

Frequently Asked Questions

- **Q: What is the size of the Manipol Massager MPV1?**
 - **A:** The size of the massager is 105 mm (L) x 136 mm (H).
- **Q: How do I adjust the speed of the massage?**
 - **A:** You can adjust the speed using the Speed + and Speed – buttons on the device.
- **Q: Is there a warranty for this product?**
 - **A:** Please refer to the warranty conditions provided in the user manual.

INTRODUCTION

- Please read the operation book before using it, especially the part of attention. Keep it properly. Thank you.

Massage principle

- Massage is a traditional Indian medicine external treatment, is through a variety of means of action to the surface of the particular site, regulating the body's physiological, and pathological conditions, achieving the effect of therapy and health care as a method. Effects of massage on the principle, mainly reflected in the change of biological information system function, adjustment and correction of abnormal cutting position three.
- It is the main treatment method for certain diseases or auxiliary measures, use very extensively.
- According to the theory of traditional Indian medicine, especially the main and collateral channels theory emphasizes on the human body through the main and collateral channels, acupoints and visceral exists between the organic immanent connection. Visceral disease can be reflected on the surface of the body surface through the main and collateral channels, collaterals, and acupoint massage to stimulate, and dredge main and collateral channels, lubricating joints, promoting qi and blood running through the main and collateral channels, can also be disease "information" to convey to the diseased organs, strengthen the internal organs function, strengthen the body's resistance, thus play a role in health care.
- From the modern medical knowledge of anatomy. Stimulation of human skin, muscle, joint neural, vascular and lymphatic etc, to promote the local blood circulation, improve the new supersedes the old., promote the body's natural resistance to disease, promote the absorption of the inflammatory exudate, relieve muscle spasm and pain.
- After thousands of years of medical practice, it proves massage has a curative effect on soft tissue injury occurring inflammation, adhesion, spasm, contractures, hyperplasia, degeneration, fibrosis, vascular nerve compression and many other pathologies, is an important means of soft tissue injury rehabilitation.

Features

1. Comfortable and convenient with handle, massage various parts of the body.
2. The massage strength is adjustable.
3. By continuously rotating the massage, vibration waves will be into surface parts, and they will decompose and the consumption of fat, to achieve the effect of weight loss.
4. The machine has three changeable massage heads, different massage head generates different massage effects, are convenient to use, and the users can choose according to different requirements.


Usage



Arm



Shoulder

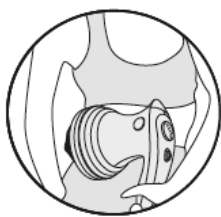

Back



Leg



Foot



Warning

When using the product, it must be covered with a protective cloth cover, lest in clamp injury.

Product Part Introduction



How to use

1. Press the "⏻" key to turn the massager on or off.
2. Press the "🖱️" key to use the auto massage mode. There are a total of three massage modes to choose from.
3. Press the "📊" key to increase or decrease the massage intensity level. There are a total of five intensity levels to choose from.
4. The massager will automatically turn off after 15 minutes.
5. After use, you must turn off the product by pressing the off key.

Massage Heads



**BALL
MESSAGE HEAD**



**SCRAPING
HEAD**




**WAVY
MASSAGER HEAD**



**FROSTING
MASSAGER HEAD**

Notice

1. After using and before cleaning, you must unplug.
2. If the machine operation is informal, please cut off power immediately, as soon as possible to contact the dealer.
3. If the supply cord is damaged, you must use a special cord or repair it by the manufacturer.
4. Pregnant women, skin disease, heart disease and other diseases, inflammation of the skin, do not use it, please
5. The children or the disabled can not use the machine.
6. Do not use it over 15 minutes in the same position.
7. Do not use this machine in damp places, such as a bathroom.
8. It is suggested that using it an hour after the meal.
9. If replacing the massage head, please make sure the machine has stopped.

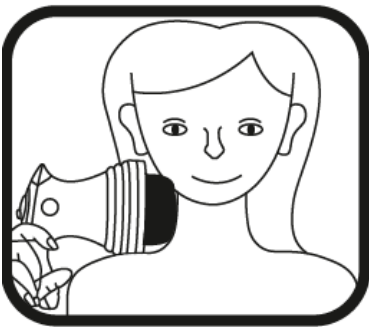

- Use Protective clothing correctly to protect the hair
- If use it wrongly, the hair will be tangled.



Repairing and maintenance

1. After using, place the product in a clean place, with no direct sunlight.
2. Do not pull the cord or pull the product with the cord.
3. Do not make water into the product.
4. Please use a soft cloth to clean the machine.

Technical parameters

1. **Voltage:** AC 220V

2. **Frequency:** 50Hz
3. **Rated power:** 25W

Warranty conditions

- Manipol Massager MPV1 is made under strict quality control & utmost care is taken to ensure long trouble-free service.
- This unit is warranted for 12 months from the date of purchase.
- We undertake to repair or replace at our discretion any part which was originally defective, provided the unit was not tampered with or mishandled.
- The warranty is only against manufacturing defects & shall not be applicable in case of damage due to accident, acts of nature, loss of components, intentional damage, and improper installation or if it was opened or serviced by authorized persons.
- The warranty is void if the set is not operated according to instructions given in the "Instruction manual booklet" If any repair work is carried out by persons other than authorized company personnel, defects are caused by improper or reckless use, which shall be determined by company personnel.
- Warranty does not cover accessories external to the equipment, defects caused due to causes beyond control like lightning, abnormal voltage, acts of God or while in transit to a service centre or purchasers/premises.

WARRANTY

MANIPOL MESSENGER MPV1

WARRANTY CARD

- Subject to Mumbai Jurisdiction only
- **No.**.....
- **Model**.....
- **Serial No.**.....
- **Invoice No. Date:**.....
- **Stamp and Signature of Dealer**.....
- **Name and Address of Dealer**.....

Customer Service

- **Marketed by: MEDTECH LIFE PVT. LTD.**
- **Regd. address: B6, Byculla Service Industries,**
- **Off D. K. Marg, Byculla, Mumbai-400 027, Maharashtra, INDIA. or Consumer complaint, Please contact Consumer Car**
- **Officer at Regd.**
- **Address C/O +91 720 8088 721**
- support@medtechlife.com
- www.medtechlife.com
- **AR01017-00-00**

Documents / Resources



[MEDTECH MPV 1 Manipol Massager](#) [pdf] User Manual
MPV 1, MPV 1 Manipol Massager, Manipol Massager, Massager

References

- [User Manual](#)

[Manuals+](#), [Privacy Policy](#)

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.