



Omnipod 5 Insulet Provided Controller User Guide

[Home](#) » [omnipod](#) » Omnipod 5 Insulet Provided Controller User Guide 

Contents

- [1 Omnipod 5 Insulet Provided Controller](#)
- [2 Omnipod 5 Onboarding](#)
- [3 Step 1 – Creating an Omnipod® ID](#)
- [4 Step 2 – Reading and Validating Data Privacy Consent](#)
- [5 Step 3 – Linking your Omnipod Account with a Glooko® Account](#)
- [6 Prepare for your Training Day](#)
- [7 Training Day Checklist](#)
- [8 Additional Resources](#)
- [9 FAQ](#)
- [10 Documents / Resources](#)
 - [10.1 References](#)
- [11 Related Posts](#)



Omnipod 5 Insulet Provided Controller



Specifications


- Compatible with Dexcom G6, Dexcom G7, and FreeStyle Libre 2 Plus sensors
- Sensors are sold separately and require a separate prescription

Onboarding Step-by-Step Guide



Thank you for choosing the Omnipod® 5 Automated Insulin Delivery System, integrated with the leading sensor brands.*

Kickstart your journey with our Step-by-Step Onboarding Guide for Omnipod 5.

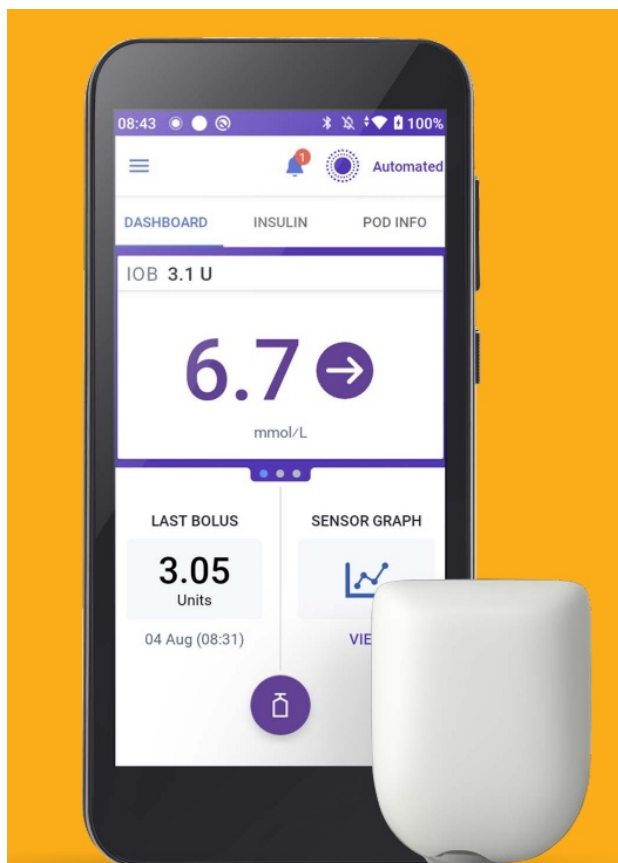

Omnipod 5 Onboarding

Before you get started on Omnipod 5, you must complete your Omnipod 5 Onboarding online before your Omnipod 5 product training.

During the Onboarding, you will create an Omnipod ID and complete the consent screens. You will also be provided with information about how your personal data is processed.

When you activate the Controller for the first time, you must enter your Omnipod ID and password.

Step 1 – Creating an Omnipod® ID

After your order is processed by Insulet, you will receive a “Complete Your Omnipod® 5 Onboarding Now” email. Open the email and select Start Omnipod® 5 Onboarding and log in with your or your dependent’s existing Omnipod ID.


If you did not receive an email:

1. Go to www.omnipod.com/setup or scan this QR code:
2. Select your country.



If you do not have an Omnipod ID

3a. Select Create Omnipod® ID.

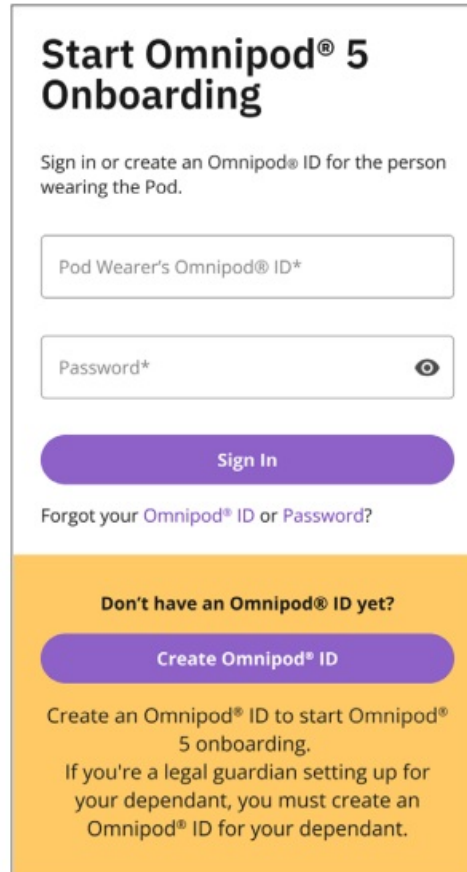

Start Omnipod® 5 Onboarding

Sign in or create an Omnipod® ID for the person wearing the Pod.

Pod Wearer's Omnipod® ID*

Password*

Sign In


[Forgot your Omnipod® ID or Password?](#)

Don't have an Omnipod® ID yet?

Create Omnipod® ID

Create an Omnipod® ID to start Omnipod® 5 onboarding.
If you're a legal guardian setting up for your dependant, you must create an Omnipod® ID for your dependant.


1. Fill in the form with your information, or the details of the dependant if you are acting as a parent or legal guardian. You will receive an email from Insulet to complete setting up your account.
2. Open the "Omnipod® ID set up nearly complete" email. Make sure you check your Junk or Spam folder if you do not see the email.
3. Select Set Up Omnipod® ID in the email. The link is valid for 24 hours.
4. Follow the on-screen instructions to review your information and set up your ID and password.
5. Follow the on-screen instructions to set up two-factor authentication by email (required) or SMS text message (optional).
6. Enter the verification code sent by email or SMS text message to complete account set up.
7. Log in with your new Omnipod ID and password.
8. Follow the on-screen instructions to verify your account if logging in from a different device.

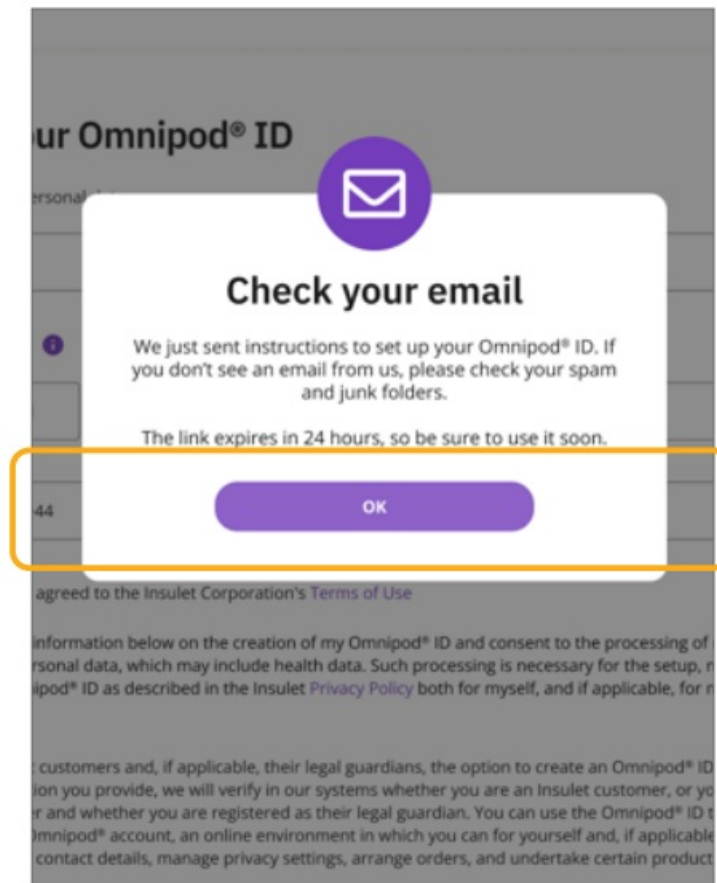

OR

If you already have an Omnipod ID

3b. Log in with your existing Omnipod ID and password.

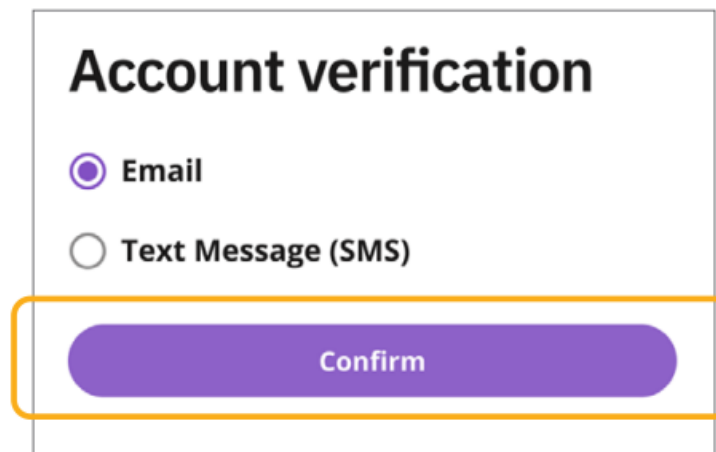

Parents and Legal Guardians

Make sure you create the Omnipod ID on behalf of the customer in your care. Select I am a legal guardian for a dependent who will wear Omnipod® 5 at the top of the Create Omnipod® ID form.



The Omnipod ID:

- should be unique
- should be at least 6 characters long
- should not contain special characters (e.g. !#£%&*~@)
- should not contain blank spaces

Passwords

- should be at least 8 characters long
- should include upper case, lower case, and number.
- should not include your (or the customer's) first name, last name, or Omnipod ID
- should only contain the following special characters (!#\$%+-<>@_)

Step 2 – Reading and Validating Data Privacy Consent

At Insulet, the safety and security of our Users and products is paramount in everything we do. We are dedicated to making the lives of people with diabetes easier and simplifying diabetes management. Insulet respects the privacy of every one of our customers and is committed to the protection of their personal information. We have dedicated teams that are focused on keeping customer information safe from unauthorised access.

After setting up your account, you must review and consent to the following data privacy policies:

1. Omnipod 5 Terms & Conditions – Required
2. Omnipod 5 Consents – Each type of consent must be agreed to individually:
 - Product Use – Required
 - Data Privacy Introduction – Required
 - Product Research, Development and Improvement – OptionalSelect Skip and Continue to opt-out
If you select Agree and Continue, a few optional questions display

Step 3 – Linking your Omnipod Account with a Glooko® Account

Glooko is the Omnipod 5 data management platform that enables you to:

- See your glucose and insulin data
- Share your data with your healthcare provider to support informed system adjustments
 - We recommend that you link your Omnipod ID to your Glooko account. If you do not have a Glooko account you can create one during setup by following these steps
 - Ask your healthcare provider for their clinic's ProConnect code to share your diabetes data


ProConnect Code:



Link a Glooko Account

After consenting to the data policies, the Omnipod 5 website prompts you to link your Glooko account.

1. Select Link on Omnipod 5
2. Select Continue to allow Omnipod 5 to send you to Glooko to log in or create a Glooko account
3. Within Glooko:
 - Select Sign Up for Glooko if you or the customer do not already have a Glooko account
Follow the on-screen instructions to create a Glooko account
 - Select Log In if you or the customer already have a Glooko account

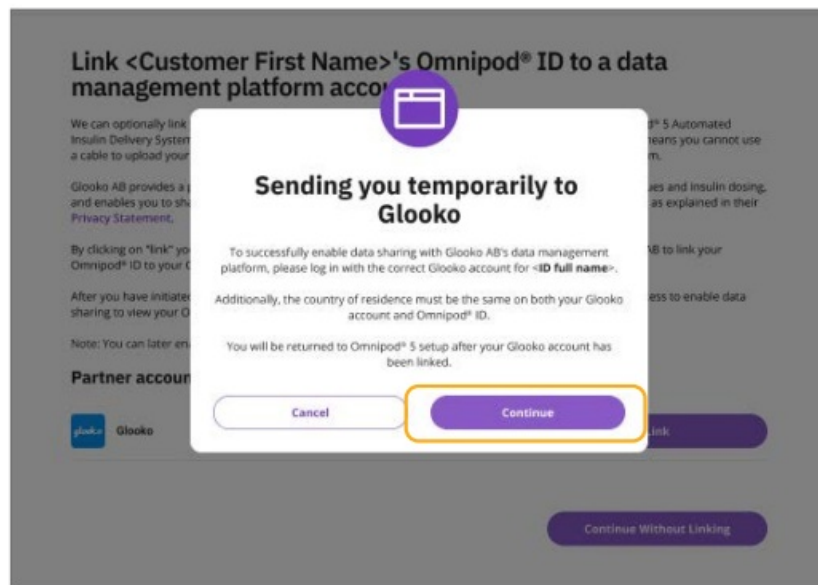


Share Glooko Data with your Healthcare Provider

After you create an account and log in, Glooko prompts you to share your Omnipod 5 data with your medical team.

1. In the Glooko app, enter the ProConnect Code your healthcare provider provided.
2. Select Share Data.
3. Select the You share data with Insulet checkbox.
4. Select Continue. You have completed setting up Glooko, but you must return to Omnipod 5 to finish sharing your data.
5. Select Return to Omnipod 5.
6. Select Agree on the Data Sharing with Glooko consent.
7. Select Continue.

Omnipod 5 sends you a confirmation email that your onboarding is complete. Once you begin using the Omnipod 5 System, Omnipod 5 will share your data with your healthcare provider through Glooko.

Congratulations on completing Omnipod® 5 Onboarding.

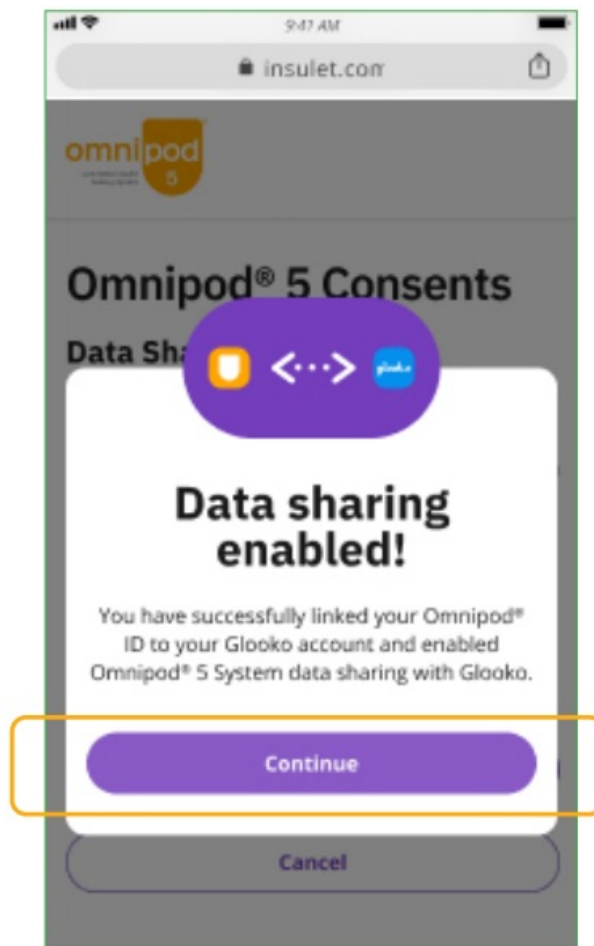

Prepare for your Training Day

In preparation for starting on Omnipod 5, please follow guidance from your healthcare provider regarding any changes to your current therapy (including any insulin therapy adjustments). You must be trained by your healthcare provider and/or Insulet Clinical team before you start on Omnipod 5.

Omnipod 5 Starter Kit

- If you are receiving your training at home, we will send you the Omnipod 5 Starter Kit and box(es) of Omnipod 5 Pods. You will also need a vial of the rapid acting insulin† prescribed by your healthcare provider.
OR
- If you are being trained in the hospital, your Omnipod 5 Starter Kit and box(es) of Omnipod 5 Pods will be there. Remember to take a vial of the rapid acting insulin† if you are using this already.

If you are expecting the delivery of your Omnipod 5 Starter Kit and Pods, and have not received these within 3 days of your scheduled training, please contact Customer Care on 0800 011 6132 or +44 20 3887 1709 if calling from abroad.



Sensors*

Dexcom Sensor

- Please come to training wearing an active Dexcom G6 or Dexcom G7 Sensor using the Dexcom app on a compatible smartphone. Also ensure that your Dexcom receiver is switched off.†

FreeStyle Libre 2 Plus Sensor

- Please ensure your healthcare provider has provided you with a prescription for FreeStyle Libre 2 Plus Sensors.
- If you are currently using a FreeStyle Libre Sensor, continue to wear this sensor when you attend your Omnipod 5 training.
- Please bring a new, unopened FreeStyle Libre 2 Plus Sensor with you to the Omnipod 5 training.

Insulin

Remember to bring a vial of rapid-acting insulin‡ to your training.

Sensors are sold separately and require a separate prescription.

†The Dexcom G6 sensor must be used with the Dexcom G6 mobile app. The Dexcom G6 receiver is not compatible.

The Dexcom G7 sensor must be used with the Dexcom G7 app. The Dexcom G7 receiver is not compatible.

‡ NovoLog®/NovoRapid®, Humalog®, Trurapi®/Truvelog/Insulin aspart Sanofi®, Kirsty®, and Admelog®/Insulin lispro Sanofi® are compatible with the Omnipod 5 System for use up to 72 hours (3 days).

Training Day Checklist

Checklist

- Have you created your Omnipod ID and password? It is important you remember your Omnipod ID and password as you will use this to log into the Omnipod 5 Controller during your training.
- Have you completed your onboarding?
- Have you accepted all mandatory consent where we provide you with information on the processing of your

personal data?

- (Optional) Did you complete linking your or your dependent's Omnipod ID with the Glooko account?
- Did you see the 'Onboarding completed!' screen and did you get the confirmation email?
- Do you have a vial of rapid-acting insulin* for your training?
- Are you wearing an active Dexcom Sensor using the Dexcom app on a compatible smartphone and ensured your Dexcom receiver is switched off?

OR

- Do you have a FreeStyle Libre 2 Plus unopened Sensor ready to be activated at your training?

Omnipod ID

- Omnipod ID:
- Password:

Glooko account

- Email (username):
- Password:

Dexcom/ FreeStyle Libre 2 Plus User ID

- Username/email address:
- Password:
- ProConnect Code:*

Additional Resources

To be fully prepared for your Omnipod 5 training, we encourage you to watch the 'How-To Videos' prior to your product training.

These and other additional online resources can be found at: [Omnipod.com/omnipod5resources](https://omnipod.com/omnipod5resources)



Pod shown without necessary adhesive.

If you have any questions or concerns regarding Omnipod 5 not answered by the online resources, please contact the Omnipod team on:
0800 011 6132* or +44 20 3887 1709 if calling from abroad.



If you have any questions about your treatment, please contact your diabetes team.

©2025 Insulet Corporation. Omnipod, the Omnipod logo, and Simplify Life are trademarks or registered trademarks of Insulet Corporation in the United States of America and other various jurisdictions. All rights reserved. Dexcom, Dexcom G6 and Dexcom G7 are registered trademarks of Dexcom, Inc. and used with permission. The sensor housing, FreeStyle, Libre, and related brand marks are marks of Abbott and used with permission. Glooko is a trademark of Glooko, Inc. and used with permission. All other trademarks are the property of their respective owners. The use of third-party trademarks does not constitute an endorsement or imply a relationship or other affiliation. Insulet International Limited 1 King Street, 5th Floor, Hammersmith, London W6 9HR. INS-OHS-01-2025-00163 V1

FAQ

How do I link my Glooko account with Omnipod 5?

After consenting to data policies, select “Link” on Omnipod 5 and continue to log in or create a Glooko account. Share data with your healthcare provider by entering the ProConnect Code provided and following the on-screen instructions.

Documents / Resources



[Omnipod 5 Insulet Provided Controller](#) [pdf] User Guide
5 Insulet Provided Controller, 5 Insulet Provided Controller, Provided Controller, Controller

References

- [Setup | Omnipod](#)
- [User Manual](#)

[Manuals+](#). [Privacy Policy](#)

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.