

ORACURA OC150 LITE SMART WATER FLOSSER



# ORACURA OC150 LITE Smart Water Flosser Operational Manual

[Home](#) » [ORACURA](#) » ORACURA OC150 LITE Smart Water Flosser Operational Manual 

## Contents

- 1 ORACURA OC150 LITE Smart Water Flosser
- 2 Product Specification
- 3 Configuration
- 4 Steps to Use The Water Flosser®
- 5 3 Operation Modes
- 6 Maintenance
- 7 Troubleshooting
- 8 Register Your Warranty
- 9 Complementary Product Offers
- 10 DISPOSAL
- 11 Safety Guidelines
- 12 ORACURA® Product Range
- 13 FREQUENTLY ASKED QUESTIONS
- 14 References
- 15 Related Posts

# ORACURA®

ORACURA OC150 LITE Smart Water Flosser




## Product Specification



## Accessories Included



## Configuration

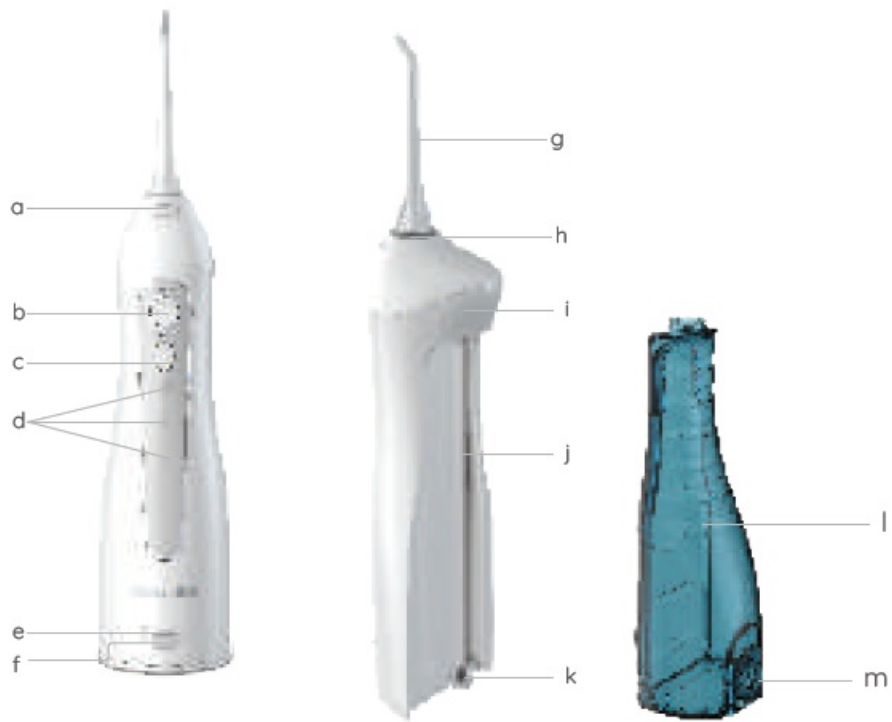

- a] Nozzle Release Button
- b] On/Off Button
- c] Mode Button
- d] Mode Indicators
- e] Charging Point
- f] Charging Point Rubber
- g] Nozzle
- h] Identification Ring
- i] Handle
- j] Water Tank Suction Pipe
- k] Water Tank Pipe Filter
- l] Water Tank
- m] Water Tank Cap

### **Charging The Water Flosser®**



1. Take out the charging point rubber and insert the water flosser® plug securely into the charging point socket. Before inserting the plug, ensure that there are no traces of liquid inside the socket.
2. Connect the charging adaptor to the power outlet.
3. The “CH” indicator blinks when the water flosser® is charging and it will turn off once it is fully charged.
4. After charging, cover the charging point again with charging point rubber. If the “CH” indicator blinks during or after water flossing, it means that the water flosser® battery needs to be recharged. Full recharging will take about 4 hours.

### **Steps to Use The Water Flosser®**

1. Attach the nozzle to the water flosser® handle.
2. Open the water tank cap and fill the tank with water. OR Slide the Water tank down to take it off the water flosser®. Fill the water tank with water and connect it to the water flosser® again.
3. Select the desired mode by pressing the “Mode” button.
4. Place the tip in the mouth. Keep your mouth slightly open, to let the water flow out and to avoid any water splash.
5. Press the “On/O” button to switch on the water flosser®. Start with gently cleaning the back teeth and work gradually towards the front teeth.

- Direct the stream at a 90-degree angle to your gumline.
  - Glide the tip along the gumline and pause briefly between teeth.
  - Clean the inside and outside of both the upper and lower teeth.
  - Clean between the teeth and around any braces, crowns, and bridgework thoroughly.
6. Turn on the device by pressing the “On/O” button again.
  7. Remove the nozzle by pressing the nozzle release button.

**Note:** Avoid using water warmer than 40°C

### Scan for Live Demo



You can schedule a live product demonstration by scanning the above QR code.

## 3 Operation Modes

You can select the desired operation mode by pressing the mode button on your water flosser® handle.

- **Normal Mode:** Higher water pressure ensures advanced cleaning. Ideal for everyday use.
- **Soft Mode:** The low water pressure allows for gentle cleaning of the mouth. Ideal for first-time users, children (age 8+ years), the elderly, and people with gum problems.
- **Pulse Mode:** The flow of water comes out in short bursts, which is ideal for cleaning & massaging the gums.

### Steps After Using The Water Flosser®

To prevent bacteria, make sure there are no water traces inside the water tank. If you're not going to use the flosser for a long time, be sure to wipe and dry the same before preserving it.

1. Open the water tank cap and dispose of the water from the tank.
2. Press the On/Off button to turn on the water flosser® and drain out the remaining water.
3. Press the On/Off button again to turn off the water flosser®, and close the water tank cap.
4. Wipe the flosser with a dry cloth.

## Maintenance

Use only water or neutral washing agents to clean the water flosser®. Never use caustic or abrasive agents (For

eg. vinegar or limescale removers). Do not use water/ cleaning fluids over 40°C. Keep the Water Flosser® in a cool place and avoid sunlight.

### **Main Unit**



1. To remove dirt, wipe the unit with a damp cloth (using liquid soap/ water) and then wipe it again with a dry cloth.
2. Do not submerge the main unit in water.

### **Nozzles**



1. Wash with water and wipe with a soft cloth.
2. Replace the nozzle after every 3-4 months.

### **Water Tank Suction Pipe**



1. Wash with water and wipe with a soft cloth.
2. Do not bend, pull or twist.

### **Water Tank & Tank Cap**



1. Wash the tank with water and wipe with a soft cloth.
2. Keep the tank dry and away from sunlight when not in use for a long period.
3. Wipe o any grime or dirt with a cloth on the tank cap.

### **Troubleshooting**

#### **Your ORACURA® Water Flosser®**



Problems	Reasons	Solutions
It does not work.	You have just purchased the product or it has been left unused for more than 3 months.	Charge the product for at least 4 hours.
It only works for a few minutes even after charging.	The battery has reached the end of its life cycle. Charging time was insufficient	Contact the service centre for support. Charge the product for at least 4 hours.
Water Pressure is very low	There is dirt stuck in the nozzle. The nozzle is deformed. There is dirt stuck in the filter	Replace the nozzle. Clean the filter. Fill up the tank
Water is not coming out of the product.	The tank is empty The product is tilted too much while in use	Fill up the tank. Use the product in the vertical position.
	There must be air caught in the pipe.	Remove the tip, close the open mouth with your thumb and shake the water flosser®.



**Please Note:** The nozzle tips should be replaced after every 3-4 months.

### Register Your Warranty

ORACURA's OC150 Smart Water Flosser® is covered under 18 + 6 months warranty period. **18-months** of limited warranty is applicable for material, manufacturing, or operational defects existing at the time of delivery or for the occurrence of any operational error during regular use. Any mishandling of the product or alteration/repair done by an unauthorized person will void the warranty. Consumer is required to submit the original purchase invoice as proof of purchase date and if requested, the entire product, to support a warranty claim.



For the extended 6 months warranty, register your product by scanning the QR code or visiting the link given below. <https://oracura.in/pages/register>



If you face any problems while using our product or have any complaints regarding the same, we advise you to visit our return policy page. <https://oracura.in/pages/return-policy>

## Complementary Product Offers



Go Beyond Manual Brushing with Our Sonic Electric Toothbrushes. If you have purchased an ORACURA® product, then this is your chance to explore various offers complementing your purchased product. Scan the QR code for further details.

Visit: <https://oracura.in/collections/sonic-electric-toothbrush>



## DISPOSAL


The rechargeable battery inside your water flosser® is easy to remove for recycling. At the end of the water flosser®'s life and prior to disposal, please remove the rechargeable battery from the handle by following the instructions below. During the below operation, be sure to protect your eyes, hands, fingers, and the surface on which you work.

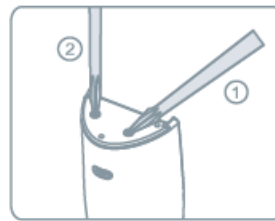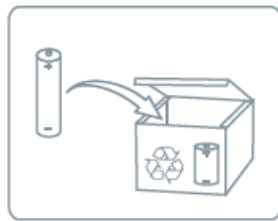


Ref. Fig. 1



Ref. Fig. 2

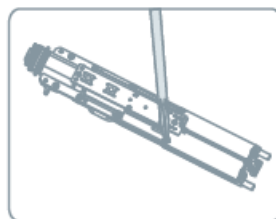

Ref. Fig. 3

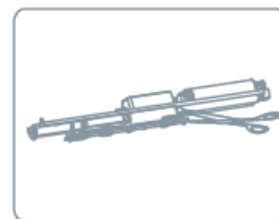

Ref. Fig. 4

1. Remove 2 rubber feet given at the bottom of the handle. (Ref. Fig. 1)
2. Insert a screwdriver into the slot located at the bottom of the handle and turn counter-clockwise to release the bottom cap. (Ref. Fig. 1)
3. Slide the water tank down on the water flosser® and unscrew the two screws beside the water tank suction pipe. (Ref. Fig. 2)

4. Remove the top cover and put the inner component unit out, then unscrew all the inner screws. (Ref. Fig. 3)
5. Remove the battery wire and take the batteries out from the inside. (Ref. Fig. 4)

If you have recycling questions, please contact your local waste management office

## **Safety Guidelines**

Read All Instructions Carefully Before Using The Water Flosser®

### **CAUTION**

- Toddlers or infants should not be allowed to use the product. Keep it out of the reach of children.
- Never try to alter, disassemble or repair the product.
- Do not place or store the product where it can fall or be pulled into a tub or sink.
- Do not use the product under wet conditions. (For eg. bathing)
- Be careful while using the product. The strong force of water spray may hurt sensitive organs of your body

### **WARNINGS**

- This product can be used by children aged 8 years and above. People with reduced physical, sensory or mental capabilities, lack of experience, and knowledge can also use the product if given proper supervision or instruction.
- Use the product only with the power supply unit that is provided.
- Make sure to remove the battery before scrapping the product. The battery should be disposed of safely.
- Please refer to the “disposal” section for any concerns regarding the removal of the battery.
- The product must be disconnected from the supply while removing the battery.
- Do not fill the tank with water warmer than 40°C.
- Never use the adaptor if it has a damaged cord or plug.
- Keep the cord away from heated surfaces.
- Do not use the product without a nozzle.
- Do not submerge the product in water.
- Use the product only for oral cleaning purposes.
- Do not only use mouthwash or any antibacterial rinses in the water tank.
- The concentrated chemical may damage the product by causing the water tank to crack.
- People who experience tooth or gum pain may not be able to use this product and should consult their dentist.
- Do not press the nozzle too firmly against your teeth or gums.
- You may experience gum bleeding during the initial stage of using the flosser. However, with regular usage, the product shouldn't hurt anymore

## **ORACURA® Product Range**

ORACURA® also offers various products such as:

- Water Flossers
- Sonic Electric Toothbrushes

- Combo Packs
- Spare Parts
- Accessories

For the complete product range, Visit: <https://oracura.in/>

### Keep Your Product Up & Running, Always!

Accessories & Spare Parts To Buy with OC150



### Service Station:

ORACURA® Solutions Pvt. Ltd. 6th Floor, Unit 1, Plot No 233/2 Vidhani Kalpavruksha, MIDC Industrial Area, T.T.C. Industrial Area, Mumbai – Pune Hwy, Shiravane, Nerul, Navi Mumbai, Maharashtra, 400706

- Visit Us: [www.oracura.in](http://www.oracura.in)
- Email: [support@oracura.in](mailto:support@oracura.in)



## FREQUENTLY ASKED QUESTIONS

What brand manufactures the ORACURA OC150 LITE Smart Water Flosser?

The ORACURA OC150 LITE Smart Water Flosser is manufactured by ORACURA.

What is the power source for the ORACURA OC150 LITE Smart Water Flosser?

The power source for the OC150 LITE Smart Water Flosser is battery-powered.

What special feature does the ORACURA OC150 LITE Smart Water Flosser offer?

The OC150 LITE Smart Water Flosser features a rechargeable battery for convenient use.

What are the primary product benefits of the ORACURA OC150 LITE Smart Water Flosser?

The main product benefit of the OC150 LITE Smart Water Flosser is its ability to effectively remove plaque for improved oral hygiene.

How many units are included in the ORACURA OC150 LITE Smart Water Flosser package?

The package contains one unit of the OC150 LITE Smart Water Flosser.

What is the price of the ORACURA OC150 LITE Smart Water Flosser?

The price of the OC150 LITE Smart Water Flosser is \$43.00.

What are the product dimensions of the ORACURA OC150 LITE Smart Water Flosser?

The product dimensions of the OC150 LITE Smart Water Flosser are 2.91 x 2.24 x 10.63 inches.

How much does the ORACURA OC150 LITE Smart Water Flosser weigh?

The weight of the OC150 LITE Smart Water Flosser is 14.11 ounces.

What is the item model number for the ORACURA OC150 LITE Smart Water Flosser?

The item model number for the OC150 LITE Smart Water Flosser is OC150.

Does the ORACURA OC150 LITE Smart Water Flosser require batteries? If so, how many?

The OC150 LITE Smart Water Flosser requires 1 Lithium Ion battery for operation.

What should I do if my ORACURA OC150 LITE Smart Water Flosser is not turning on?

Firstly, ensure that the battery is properly inserted and charged. If using rechargeable batteries, try charging them fully. If the issue persists, check for any visible damage to the device or battery compartment.

How can I fix it if the water flow from my ORACURA OC150 LITE Smart Water Flosser is weak or inconsistent?

Check if the water tank is properly filled and seated. Inspect the nozzle for any blockages or debris, and clean it thoroughly. If the problem continues, try adjusting the pressure settings.

What steps should I take if my ORACURA OC150 LITE Smart Water Flosser is leaking water?

Ensure that all connections are securely tightened, including the nozzle and water tank. Check for any cracks or damage to the water tank. Clean the sealing surfaces and ensure they are free from debris.

How can I troubleshoot it if the pressure settings on my ORACURA OC150 LITE Smart Water Flosser are not functioning properly?

Ensure that the device is properly charged or has fresh batteries if applicable. Try resetting the device according to the manufacturer's instructions.

What should I do if the battery life of my ORACURA OC150 LITE Smart Water Flosser seems shorter than expected?

Ensure that the device is fully charged before initial use and recharge it according to the manufacturer's instructions. If the battery still does not hold a charge, it may need replacement.

**DOWNLOAD THE PDF LINK:** [ORACURA OC150 LITE Smart Water Flosser Operational Manual](#)

---

## References

- [User Manual](#)

## **Manuals+. Privacy Policy**

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.