



Sharper Image Facial Rejuvenator User Manual

[Home](#) » [Sharper Image](#) » Sharper Image Facial Rejuvenator User Manual 

Contents


1 SHARPER
IMAGE®

2 Related Posts

SHARPER IMAGE®

PROFESSIONAL FACIAL REJUVENATOR Item No. 206823

Thank you for purchasing the Sharper Image Professional Facial Rejuvenator. Please take a moment to read this guide and store it for future reference.

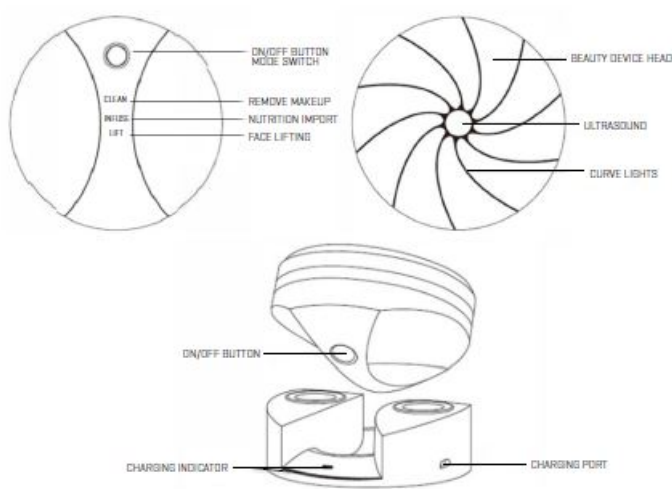

EFFECTS


Using facial texture positioning, the device brightens, lifts and tightens the skin while removing wrinkles.

INSTRUCTIONS

Place the beauty device into the charging base and make sure the on/off button is in the same direction as the charging indicator. During the charging process, the charging indicator will flash red. After it's charged, the charging indicator will turn green.

Note: If the on/off button and the charging indicator are not in the same direction, the charging indicator won't turn red. This means it's unable to be charged.

DEEP CLEANSING (Removes makeup)

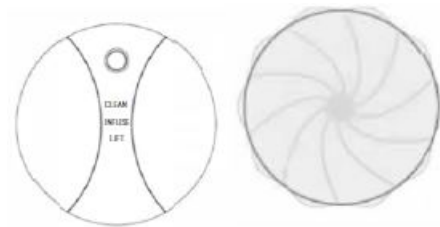

STEP 1

Fix the makeup removing cotton on the device head and apply makeup removing lotion/oil on the cotton.


STEP 2

Long press the on/off button and choose CLEAN mode. Ultrasound function is now activated (green light turns on).

STEP 3

Hold the device parallel to the cheek (as a starting point), then gently rotate from the outside of the face to the chin, then from the bottom up to the forehead (repeat 3 times).

ION-IMPORT

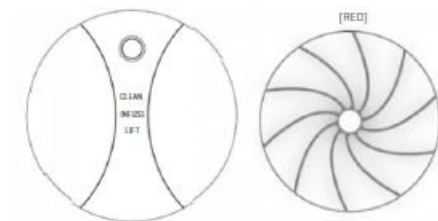

STEP 1

Press the on/off button to select INFUSE mode.


STEP 2

The device activates negative ion (red light turns on).

STEP 3

Apply toner evenly on the face and hold the device parallel to the chin (as a starting point), then gently lift the face from the inside to outside.

LIFTING

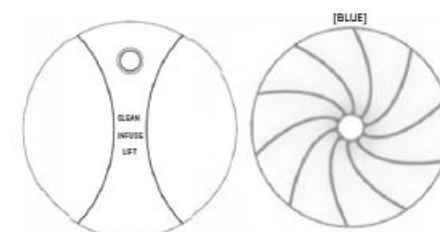

STEP 1

Press the on/off button to select LIFT mode.


STEP 2

The device activates EMS function (blue light turns on).

STEP 3

Apply serums evenly on the face and rotate the device from the inside to the outside, then from the bottom up to the forehead.

ACCESSORIES

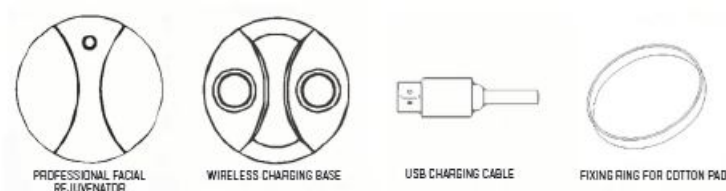

FUNCTIONS

1. Ultrasound is a high-level sound wave technology that penetrates deep beneath the dermal layer. In the skin, it increases the temperature of the surrounding tissue, by increasing oxygenation of the cells, stimulating the regeneration of collagen protein and aiding in the removal of waste products.
2. Green light (520nm wavelength): Brightens complexion, soothes skin, reduces oil secretion and balances water and oil.
3. Red light (630nm wavelength): Bio-active light opens pores for better absorption, improving active cells and brightening skin tone. Provides oxidation while reducing the effects of aging.
4. Blue light (470nm wavelength): Fights acne and repairs skin while leaving no scars.
5. Ion-: Negative ions (with assistance of high frequency vibration) effectively promote the absorption of nutrients, improving skin from the inside out. Delivers a deep moisturizing effect.
6. EMS (Electrical Muscle Stimulation): Through the periodic change of micro-currents, the skin of the face becomes more contoured.

SPECIFICATIONS

Size	2.6 X 1.6 inches
Weight	0.27 lbs.
Material	ABS Plastic and Stainless Steel

WARRANTY / CUSTOMER SERVICE

Sharper Image branded items purchased from SharperImage.com include a 1-year limited replacement warranty. If you have any questions not covered in this guide, please call our Customer Service department at 1 (877) 210-3449. Customer Service agents are available Monday through Friday, 9:00 a.m. to 6:00 p.m. ET.

SHARPER IMAGE®

The SHARPER IMAGE® name and logo are registered trademarks.

©Sharper Image. All Rights Reserved.

[Sharper-Image-Facial-Rejuvenator-Manual-Optimized](#)

[Sharper-Image-Facial-Rejuvenator-Manual-Original](#)

Manuals+.