



Snailax SL489 Vibration Foot Massager User Manual

[Home](#) » [SNAILAX](#) » Snailax SL489 Vibration Foot Massager User Manual 

Contents

- [1 Snailax SL489 Vibration Foot Massager](#)
- [2 Safety Instructions](#)
- [3 Features](#)
- [4 TECHNICAL DATA](#)
- [5 Instructions for use](#)
- [6 Battery Charging](#)
- [7 Troubleshooting](#)
- [8 Care and Maintenance](#)
- [9 Warranty](#)
- [10 FREQUENTLY ASKED QUESTIONS](#)
- [11 VIDEO – PRODUCT OVERVIEW](#)
- [12 References](#)
- [13 Related Posts](#)



Snailax SL489 Vibration Foot Massager



Safety Instructions

READ ALL INSTRUCTIONS BEFORE USE. SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning.
- Do not wet, do not use pins, and never remove the cover.
- This item is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
- This appliance should never be left unattended when plugged in.
- Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures, and direct sunlight.
- Do not use any damaged cables, plugs or loose sockets.
- In case of malfunctions, disconnect immediately from the mains.
- Do not use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Misuse or incorrect use excludes any liability for damage.
- Do not use this device while driving.
- Do not use it while you sleep.
- Do not use with infants, a helpless person or a person insensitive to heat.
- Overlong application can result in skin burns.
- Do not use this product for medical treatment.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not exceed 15 minutes at a time.

If you suffer from any of the medical conditions listed below, consult with a doctor before using this relaxation massage device:

- Suffer from a heart condition
- have a sensory disorder or pacemaker
- suffer from osteoporosis or diabetes
- have a malignant tumor
- have neurological disorder or perceptual disorder on the head
- have fractured or head surgery of your massaging area
- head skin disease and allergic dermatitis
- suffering from acute diseases or are undergoing medical treatment we advise you to consult with your doctor before using this massager.
- The packaging material supplied should not be used as a toy.

Failure to observe the above instructions may constitute misuse of the product and can cause serious injury or burns

Features

Thank you for purchasing the SNAILAX Hand Massager With normal care and proper treatment, it will provide years of reliable service

FEATURES

1. Compression Massage with soothing heat function.
2. 3 levels of massage intensity for options
3. Optional Vibration massage
4. Built-in 10-minute timer to automatic shutdown
5. Built-in rechargeable Lithium battery
6. Easy to operate the Control buttons

CONTENTS

- Hand Massager
- Charging Adapter




TECHNICAL DATA

- **Model:** SL-489
- **Weight:** 2.4 Lbs
- **Adapter Input:** 100-240V~ 50/60 Hz

- **Rated Voltage:** 5.0VDC 2.0A
- **Rated Power:** 8W
- **Default Automatic Runtime:** 10 minutes
- **Battery Spec:** 3.7V Lithium battery, 2200mAH

Instructions for use

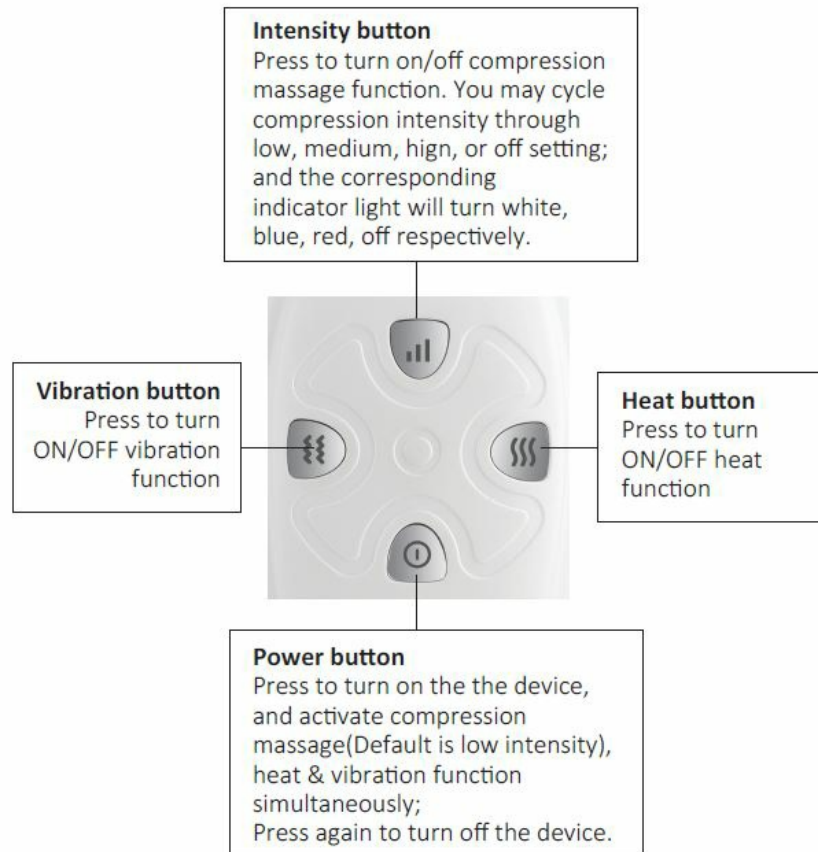

Note

This device will switch off automatically after the 10-minute timer runs out. You should allow your body some rest periods. To avoid straining your muscles too much, we recommend not exceeding a continuous massage of 10 Minutes per time and not exceeding 3 times one day

Battery Charging

1. When the battery is low, the power indicator light will start flashing.
2. Connect the charging adapter cable to the corresponding port on the massager



3. When charging, the power indicator light will start flashing.
4. When full charge, the power indicator will stop flashing.

5. Normally it takes about 3 hours to fully charge.
6. It is safe to use while charging.
7. After fully charged, the device can work continuously for around 1.5 hours

Note: before using this product for the first time, please fully charge it. If the device is not in use for a long time, the product should be charged once a month to prevent the battery life from shortening.

Battery Disposal Warning:

The unit contains a Lithium Ion battery.

Please contact your local recycling center to find out how to properly dispose of the product.

Troubleshooting

Trouble	Reason/Solution
Can not start the massager	Battery Low, Please charge it before using
The massager stopped abnormally	Battery runs out of power, charge the battery; Self-timed 10 mins is out, restart the product
Massager slight noise	It is normal working noise of the internal mechanism

Care and Maintenance

Care and cleaning instructions

1. Always unplug the unit before cleaning it.
2. Before cleaning, allow the unit to cool.
3. Dampen a cloth in water or a 3%–5% mild detergent solution.
4. Wipe the dirty areas with a wet cloth.
5. Do not use a sharp brush, sharp tool, gasoline or diesel to clean it.
6. Do not use a wet cloth to clean the control panel, lest cause failure, only use a dry cloth to wipe the control panel.
7. Never submerge any part of the unit in liquid.
8. Wait for the unit to become fully dry before use.

Maintenance

1. Please keep it into a box and keep in a safe, dry and cool place.
2. Do not keep it in any place with high temperature, near fire or direct sunlight, or it may faded or damaged
3. Do not contact benzene, dilient or any other corrosive liquid.

Warranty

If you have any issue about the product, please feel free to contact us by sending email to support@snailax.com
We will strive to provide the best service possible within 24 hours.

- **30 days unconditionally Return** Snailax product can be returned to receive a full refund for any reason within 30 days. Please contact our customers service (support@snailax.com), our staff will contact you within 24 hours.
- **90 days refund/replace** Snailax product can be returned / replaced within 90 days if the product break down in the period of proper use.
- **12 months Warranty** If the product break down within 12 months in the period of proper use, customers can enjoy the relevant warranty to get them replaced.
- **Attention!** No warranty will be given to any force majeure and man-made causes to a defective product, such as improper care, personal tear down and willful damage, etc.

Extend Warranty for Free

- Enter the following URL or scan the QR code below to find Snailax facebook page and like it, enter "Warranty" to messenger to extend your warranty from 1 year to 3 years.



<https://www.facebook.com/snailax> OR 2) Send message "Warranty" and email us support@snailax.com to extend your warranty from 1 year to 3 years.

Have question ?

- **Tel:** 734-709-6982
- Monday-Friday 9:00AM-4:30PM
- **Email:** support@snailax.com

FCC STATEMENT

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant

to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio / TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this unit is not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

2750 Carpenter Road, Suite 1 B, Ann Arbor, 48108 Tel:734-709-6982 support@snailax.com www.snailax.com

FREQUENTLY ASKED QUESTIONS

What is the brand of the hand massager?

The brand is Snailax.

Can you provide the specific model number for the Snailax hand massager?

The model number is SL489.

What is the color of the Snailax SL489 Hand Massager?

The color is White.

What are the product dimensions of the Snailax SL489 Hand Massager?

The product dimensions are 8.27 x 5.59 x 10.43 inches.

What is the intended use for the Snailax SL489 Hand Massager?

It is designed for use on Hands.

What is the power source for the Snailax SL489 Hand Massager?

The power source is Battery Powered.

How much does the Snailax SL489 Hand Massager weigh?

The weight is 2.91 Pounds.

What functions does the Snailax SL489 Hand Massager offer for hand relaxation?

It combines compression, optional vibration, and heat soothing functions to provide a comprehensive massage experience targeting each acupuncture point of the fingers, palms, and wrists.

How does the vibration feature in the Snailax SL489 Hand Massager benefit users?

The optional vibration helps soothe sore muscles and soft tissues, effectively alleviating finger numbness, arthritis, joint pain, or carpal tunnel pain.

What is the purpose of the soothing heat therapy in the Snailax SL489 Hand Massager?

The gentle heat therapy aims to eliminate fatigue and coldness while promoting blood circulation.

How many massage modes does the Snailax SL489 Hand Massager offer, and what are they?

It offers 3 massage modes to meet different massage needs.

Is the Snailax SL489 Hand Massager wireless, and how long does the battery last?

Yes, it is wireless and rechargeable. The built-in lithium battery can last 90 minutes after a full charge.

Where can the Snailax SL489 Hand Massager be conveniently used?

It can be used at home, in the office, or during traveling.

Is there a timer feature in the Snailax SL489 Hand Massager, and what is the duration?

Yes, there is a built-in 10-minute timer to ensure the optimal massage duration.

Is the Snailax SL489 Hand Massager suitable for both left and right hands?

Yes, it is designed for both left and right hands and accommodates a full range of hand sizes.

VIDEO — PRODUCT OVERVIEW



DOWNLOAD THE PDF LINK: [Snailax SL489 Vibration Foot Massager User Manual Vibration-Foot-Massager-User-Manual.mp4](#)

REFERENCE: [Snailax SL489 Vibration Foot Massager User Manual-Device.Report](#)

References

- [User Manual](#)