

Contents [[hide](#)]

- [1 Songhui CS01 Smart Body Fat Scale](#)
- [2 Specifications](#)
- [3 Display instructions](#)
- [4 Directions for use](#)
- [5 APP download and measure](#)
- [6 FAQs](#)
- [7 Documents / Resources](#)
 - [7.1 References](#)

Songhui

Songhui CS01 Smart Body Fat Scale



Specifications

- **Dimensions:** 270mm x 125mm x 30mm
- **Weight:** 105g
- **Size:** 280*280mm
- **Weight:** 1.3 kg
- **Battery:**2 AAAA
- **Division value:** 0.1 kg
- **Measurement Range:** 0.2-180kg
- **Display:** LED
- **Material:** ABS, Safety tempered glass

Information

The Smart Body Fat Scale is a cutting-edge device designed to accurately measure various body metrics to help you track your health and fitness progress.

Display instructions

1. Lo(low battery)
2. O-Ld(Overweight prompt)
3. err(test error)
4. Unable to boot or display(replace the battery)

Safety precautions


- Be sure to read the safety precautions before use.
- This product should not be used by those equipped with a pacemaker, artificial lung, or wearing an electrocardiograph.
- Do not use this product on slippery surfaces

Directions for use

- The correct way to measure
- Please put the product on a hard, flat floor; do not put it on the carpet or any other soft ground, or it will cause inaccurate measurement data. •
- Clear before each measurement.

Reset method: After lightly stepping on the scale, the screen number will automatically

return to 0.0lb or 0.00lb

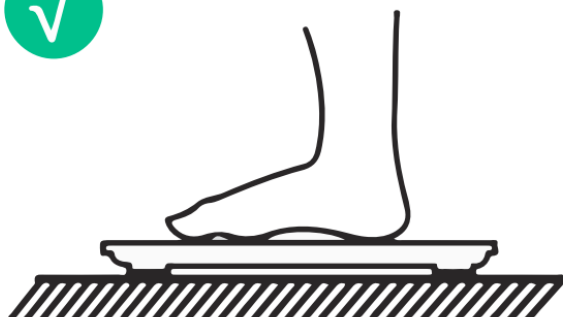

Place the scale on level ground



Do not place on carpets or soft floors

The correct way to stand



APP download and measure


Step 1: APP Download

- Scan the QR code to download the “UREVO” app.
- iOS can be downloaded through the App Store, and Android can be downloaded through Google Play.



Step 2: How to connect the scale for the first time

- Click “Add Device”, enter the search device page, and click “weight scale” to search.h;
- Activate the scale, search for the scale device, and connect it.
- You can start taking weight data.

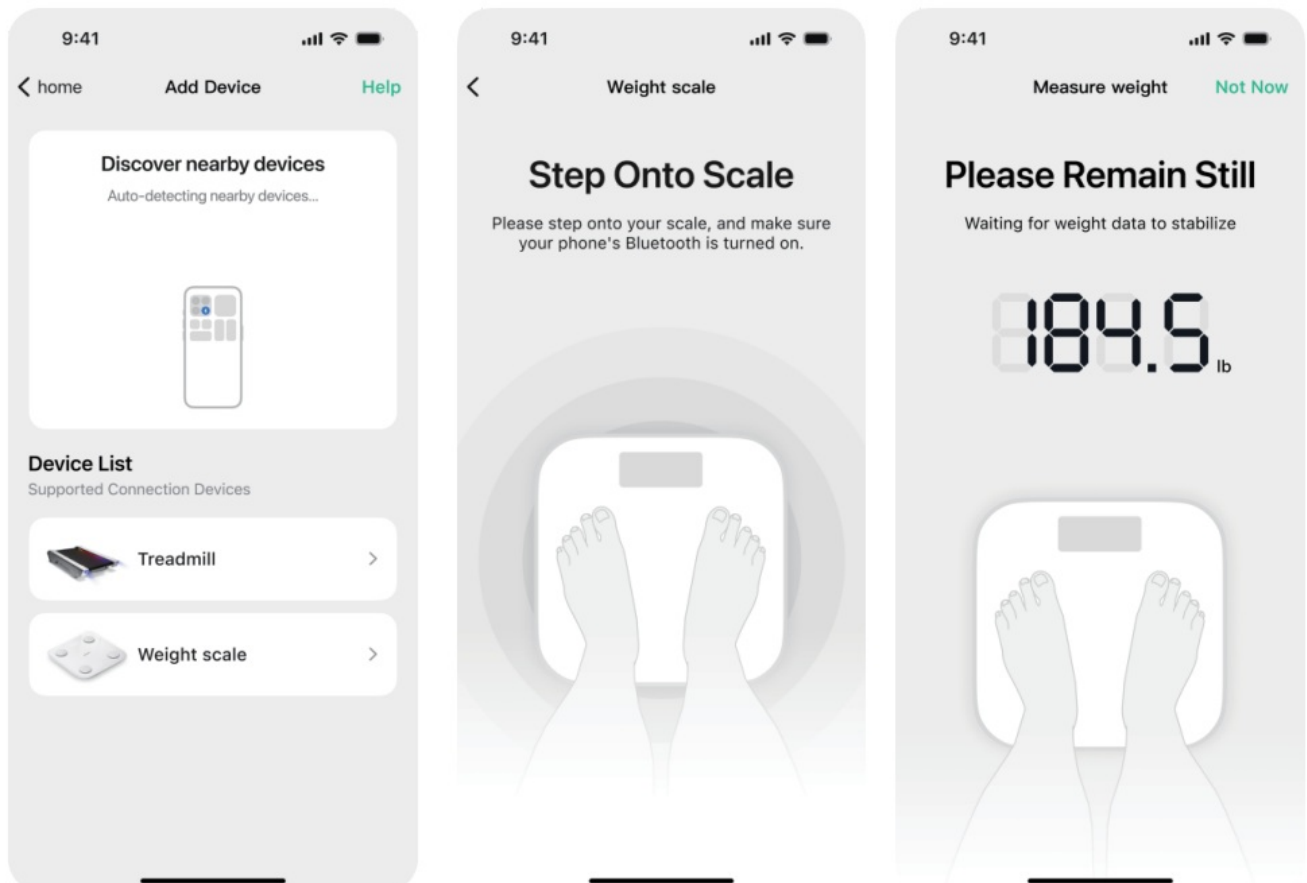


Step 3: View data. Go to the Scale home page to view measurements and more.



Units Switch: Press the unit key at the back of the scale to choose the weighing unit you would like to use.

Press here to switch units kg/lb

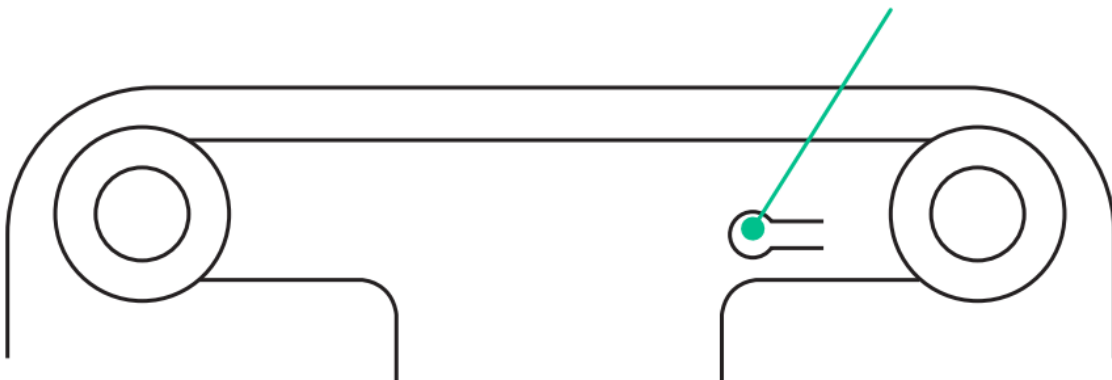

FCC Caution

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This device may not cause harmful interference, and
2. This device must accept any interference received, including interference that may cause undesired operation. Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, our stand to Part 15 of the FOC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. His equipment generates, uses, and can radiate radio frequency energy and, if not installed and used according to the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help. The device has been evaluated to meet the general Ki exposure. The device can be used in portable exposure conditions without restriction.

FAQs

Q: How often should I use the Smart Body Fat Scale?

A: It is recommended to use the scale regularly, such as once a week, at the same time of day for consistent measurements.

Q: Can multiple users use the scale?

A: Yes, many smart scales have the capability to store data for multiple users. Check the user manual for instructions on setting up profiles.

Documents / Resources



[Songhui CS01 Smart Body Fat Scale \[pdf\]](#) Instruction Manual
CS01 Smart Body Fat Scale, CS01, Smart Body Fat Scale, Body Fat Scale, Fat Scale

References

- [User Manual](#)

Body Fat Scale, CS01, CS01 Smart Body Fat Scale, Fat Scale, Smart Body Fat Scale,

Songhui Songhui

Leave a comment

Your email address will not be published. Required fields are marked *

Comment *

Name

Email

Website

☐ Save my name, email, and website in this browser for the next time I comment.

Post Comment

Search:

e.g. whirlpool wrf535swhz

Search

[Manuals+](#) | [Upload](#) | [Deep Search](#) | [Privacy Policy](#) | [@manuals.plus](#) | [YouTube](#)

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.