




## TaiDoc TD-1242BT Thermometer User Guide

[Home](#) » [TaiDoc](#) » TaiDoc TD-1242BT Thermometer User Guide 

*TaiDoc TD-1242BT Thermometer User Guide*

# ™

## Thermometer Quick Start Guide

### TaiDoc Clever TD-1242BT


The enclosed thermometer is provided for your use during your participation in the care management program to record your body temperature. The following information explains the process for recording your body temperature and viewing trend data.

Please see the enclosed manufacturer's guide for complete instructions relating to the TaiDoc Clever TD-1242BT thermometer, including important information regarding when you should consult your physician prior to using this device.

#### Record Your Vitals

1. On the tablet, select the thermometer icon.
2. Aim the device at the center of the forehead 1 to 2.5 inches away from the skin.
3. Press and release the scan/power button to take a measurement. A double beep sound indicates a reading was taken.


4. A number will display on the device. Look at the tablet. The device will wirelessly connect with the tablet and record your data automatically.
5. Select YES to confirm.

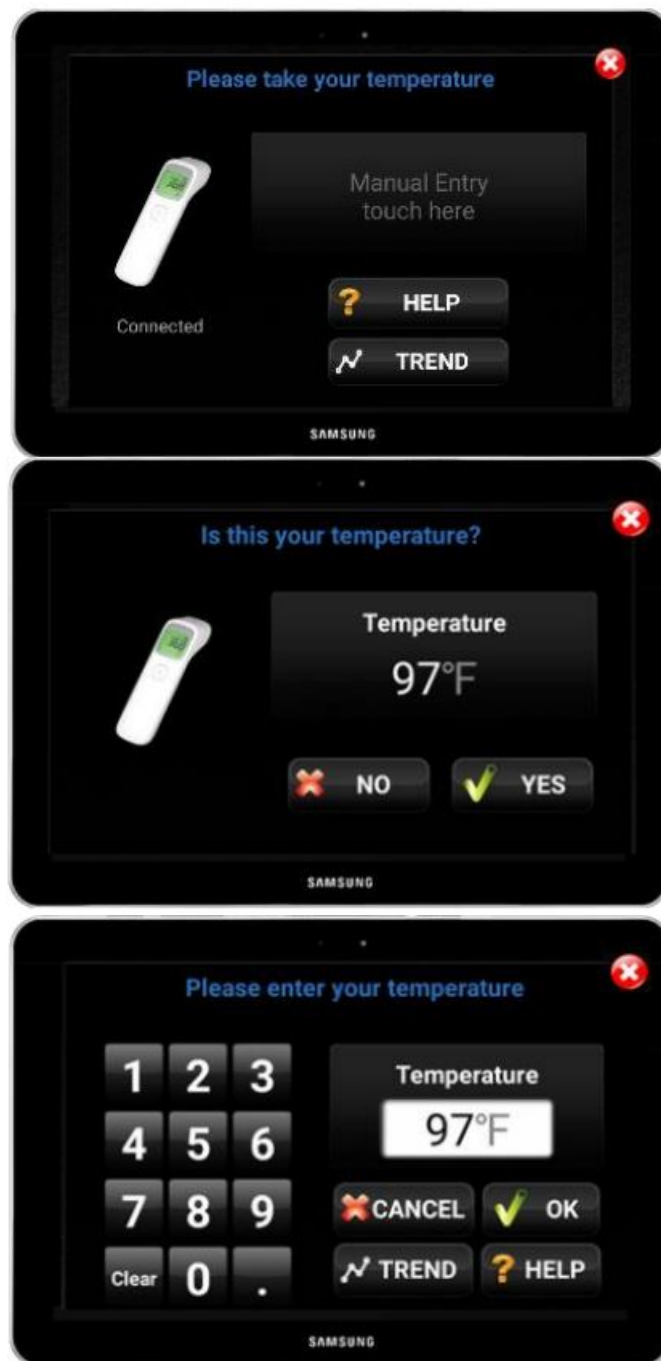

If the device does not transmit to the tablet, move the device closer to the tablet. The thermometer and the tablet should be within eyesight of each other. To record your data manually, select Manual Entry. In the field provided, record your temperature.

Note: The tablet shows the temperature unit of measurement ( or ) that is applied in the Care Team Portal, which may differ from the device unit of measurement.

#### **View Trend Data**

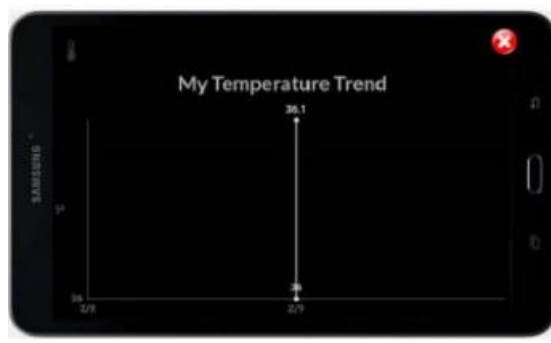

To view the changes in your biometric data over time, select the TREND button on the entry screen. Biometric data is available for the devices used for your program.

### Reminders

If you have an emergency, call 911. The program is not for emergencies.

Do not turn off or unplug the tablet while you are participating in the program.

Check your tablet to ensure it is plugged into the wall and the cord is securely plugged into the tablet.

If you forget your password, select HELP and then RESET. Your Care Team will be notified and assist you as needed.

DO NOT throw away the health kit box your tablet and devices came in. The box will be used to package the tablet and devices to be returned at the end of your remote care program.

### Tips for Accurate Readings

Do not take a reading while moving or talking.

Wait 30 minutes after exercising, bathing, or being in extreme temperatures before taking a reading.

The thermometer turns off automatically after 30 seconds.

**DO NOT USE THIS PROGRAM FOR EMERGENCY OR URGENT CARE. IN AN EMERGENCY, CALL 911 OR GO TO THE NEAREST EMERGENCY ROOM.** The program and program nurses and other representatives cannot diagnose your problems or suggest treatment. This program is not a substitute for your doctor's care. Your health information is kept confidential in accordance with the law. Any information and any device provided through this program is for your information only. Participation in the program and use of the device is voluntary, and your ability to use the device may differ depending on the nature of your medical condition. The program and device do not provide real-time monitoring or act as an emergency alert system, and are not intended for use when time-critical care is required.

### Contents

[1 Documents /](#)

[Resources](#)

[2 Related Posts](#)

### Documents / Resources



## [TaiDoc TD-1242BT Thermometer](#) [pdf] User Guide

### TD-1242BT, Thermometer