

Medline MDS86100RF

Medline MDS86100RF Toilet Safety Rail Instruction Manual

Model: MDS86100RF

1. INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your Medline MDS86100RF Toilet Safety Rail. This product is designed to offer stability and support for individuals requiring assistance when sitting down or standing up from a toilet. Please read all instructions carefully before assembly and use.



Image 1.1: Medline Toilet Safety Rail installed on a standard toilet.

Safety Information

- **Maximum Weight Capacity:** Do not exceed the maximum weight capacity of 250 lbs (113 kg).
- Ensure all components are securely fastened before each use.
- Do not use if any parts are damaged or missing. Contact Medline customer support for replacement parts.
- This product is intended for use as a toilet assist device only. Do not use for purposes other than its intended use.
- Consult a healthcare professional if you have questions regarding the suitability of this product for your specific needs.

2. COMPONENTS

Verify that all the following components are present before beginning assembly:

- Toilet Safety Rail Frame (Aluminum)
- Padded Armrests (pre-attached)
- Adjustable Mounting Bracket
- Mounting Knobs (2)
- Height Adjustment Pins (pre-installed)
- Slip-Resistant Rubber Tips (pre-installed on legs)



Height adjustable from 27.5" to 31.5"

Image 2.1: Detail of the comfortable foam armrests.

3. SETUP AND INSTALLATION

The Medline Toilet Safety Rail is designed for easy installation onto most standard toilets. No special tools are required.

1. **Prepare the Toilet:** Lift the toilet seat and lid. Locate the two bolts that secure the toilet seat to the toilet bowl.

2. **Position the Bracket:** Place the adjustable mounting bracket of the safety rail over the toilet seat bolts. Ensure the bracket is centered and flush against the toilet bowl.
3. **Secure the Rail:** Re-attach the toilet seat and lid, aligning the holes with the bracket and toilet bowl. Use the provided mounting knobs to tighten the bracket securely to the toilet seat bolts. Hand-tighten only; do not overtighten.
4. **Adjust Width:** The width between the arms can be adjusted. Loosen the knobs on the horizontal bar of the frame, slide the arms to the desired width (between 15-18 inches), and re-tighten the knobs.
5. **Verify Stability:** Gently test the rail by applying downward pressure to ensure it is stable and does not wobble. If unstable, re-check all connections and tighten the mounting knobs.



Image 3.1: The adjustable mounting bracket secures the rail to the toilet seat bolts.



Image 3.2: The bracket is designed to fit under most standard toilet seats.

4. OPERATING INSTRUCTIONS AND ADJUSTMENTS

Using the Safety Rail

The safety rail provides a secure grip for assistance. When sitting down, slowly lower yourself while holding onto the armrests. When standing up, push down on the armrests to assist in rising.

Height Adjustment

The height of the armrests can be adjusted to suit individual user needs, ranging from 27.5 inches to 31.5 inches from the floor.

1. Locate the spring-loaded adjustment pins on each leg of the safety rail.
2. Press the pin inward and slide the leg up or down to the desired height. Ensure the pin fully engages in one of the adjustment holes.
3. Repeat for the other leg, ensuring both sides are adjusted to the same height for stability.



**Easily mount to any toilet
with provided adjustable bracket**

Image 4.1: The spring-loaded pin allows for easy height adjustment.



Image 4.2: Ensure the pin is fully engaged in the desired height setting.

5. MAINTENANCE

Cleaning

The Medline Toilet Safety Rail features a durable aluminum frame with a scratch-resistant powder coating, making it easy to clean.

- Wipe down the frame and armrests regularly with a damp cloth and mild soap.
- Avoid using abrasive cleaners or solvents, as these may damage the finish.
- Ensure the rail is dry before next use.

Inspection

Periodically inspect the safety rail for any signs of wear, damage, or loose components.

- Check the mounting bracket and knobs for tightness.
- Inspect the rubber tips on the legs for wear; replace if damaged to maintain slip resistance.

- Ensure height adjustment pins engage properly.

6. TROUBLESHOOTING

Issue: The safety rail feels wobbly or unstable.

Solution:

- Ensure the mounting knobs securing the bracket to the toilet seat bolts are fully tightened.
- Verify that both legs are adjusted to the same height and the adjustment pins are fully engaged.
- Check that the width adjustment knobs on the horizontal bar are securely tightened.

Issue: The rail does not fit my toilet.

Solution:

- The rail is designed to fit most standard toilets. Ensure the adjustable mounting bracket is correctly positioned over the toilet seat bolts.
- Confirm that the width between the arms is adjusted appropriately for your toilet and user comfort.

7. SPECIFICATIONS

| Feature | Detail |
|--------------------------------|-------------------------------------|
| Model Number | MDS86100RF |
| Brand | Medline |
| Material | Aluminum |
| Color | Silver/Gray |
| Minimum Height | 27.5 Inches |
| Maximum Height | 31.5 Inches |
| Width Between Arms | 15-18 Inches (Adjustable) |
| Handle Width | 19-22 Inches (Adjustable) |
| Maximum Weight Capacity | 250 Pounds |
| Product Dimensions (L x W x H) | 27.5-31.5"L x 19-22"W x 27.5-31.5"H |
| Item Weight | 4.39 pounds |
| UPC | 080196755403 |

250
pounds

Front to back
27.5–31.5"

27.5–31.5"

19–22"

**Sturdy aluminum frame
is scratch-resistant**

Image 7.1: Overview of product dimensions and weight capacity.

8. WARRANTY INFORMATION




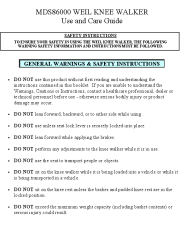

Medline products are manufactured to high-quality standards. For specific warranty details, please refer to the warranty card included with your product or visit the official Medline website. Keep your purchase receipt as proof of purchase.


9. CUSTOMER SUPPORT

If you have any questions, require assistance, or need to order replacement parts, please contact Medline customer support:


- **Website:** Visit the official Medline website for support resources.
- **Phone:** Refer to the contact information provided on the product packaging or Medline's website.

Related Documents - MDS86100RF

| | |
|---|---|
|  | <p>Medline Foldable Toilet Safety Frame User Instructions and Warranty</p> <p>User instructions, warranty information, and safety guidelines for the Medline Foldable Toilet Safety Frame. Includes assembly, cleaning, and important safety warnings.</p> |
|  | <p>Medline Basic Steel Rollator with 6\"</p> <p>Comprehensive user instructions, assembly guide, operating procedures, maintenance tips, and warranty information for the Medline Basic Steel Rollator with 6" Wheels (Model 100RA).</p> |
|  | <p>Medline Freedom Rollator User Manual and Safety Guide</p> <p>Comprehensive user manual for the Medline Freedom Rollator, including safety instructions, assembly guide, operating procedures, maintenance tips, and warranty information. Features include adjustable height, 250 lbs weight capacity, and lightweight aluminum construction.</p> |
|  | <p>MDS86000 Weil Knee Walker Use and Care Guide</p> <p>This guide provides essential safety instructions, fitting procedures, operating methods, maintenance tips, and warranty information for the Medline MDS86000 Weil Knee Walker, a mobility aid designed for individuals recovering from leg injuries or surgery.</p> |
|  | <p>Medline MDS86000G4 Four Wheeled Kneewalker Assembly and User Guide</p> <p>Comprehensive assembly and user guide for the Medline MDS86000G4 Four Wheeled Kneewalker. Includes parts identification, detailed step-by-step assembly instructions, safety guidelines, adjustment procedures for knee pad height, hand brakes, and handle height, as well as instructions on how to fold the device.</p> |



**STEP-BY-STEP
ASSEMBLY GUIDE.**




Trigger walker with 5" (12.7 cm) wheels, Black
MDS86410TRIG


[Medline MDS86410TRIG Trigger Walker Assembly Guide](#)

Step-by-step assembly instructions for the Medline MDS86410TRIG folding trigger walker with 5-inch wheels and glide brakes. Includes setup, adjustments, safety guidelines, and maintenance tips.

Documents - Medline – MDS86100RF



Durable medical equipment product guide
Innovative products and support



[\[pdf\] User Manual Guide Catalog](#)

Durable medical equipment product guide Medline 14 dic 2022 — G90 214 8 Follow your rollator s manual on how to adjust the handle as Ideal design for patients with dexterity limitations LIT474 CAT DME 19W3805227 medline media catalog Docs MKT |||

Durable medical equipment product guide Innovative products and support Elevate care with medical ... 113.4 kg G5-200X1 Aluminum Foldable Toilet Safety Rails with Padded Armrests 300 lb 136.1 kg **MDS86100RF** Aluminum Foldable Toilet Safety Rails 250 lb 113.4 kg G3-0300H Aluminum Toilet Safet...

lang:en **score:34** filesize: 6.56 M page_count: 64 document date: 2022-12-05


FOLDABLE TOILET SAFETY FRAME
CAGNE D'APPUI DE SÉCURITÉ
PLIABLE POUR TOILETTE
BARANDILLA DE SEGURIDAD
PLIABLE PARA INODORO




[Medline Foldable Toilet Safety Frame User Instructions and Warranty](#)

User instructions, warranty information, and safety guidelines for the Medline Foldable Toilet Safety Frame. Includes assembly, cleaning, and important safety warnings.

lang:fr **score:31** filesize: 396.01 K page_count: 6 document date: 2021-06-29



Durable medical equipment product guide
Empowering health, enhancing lives.



[\[pdf\] Guide Catalog](#)

Durable medical equipment product guide mobility1st ca Elevate care with that goes the distance High quality engineering long lasting performance and a comprehensive unique offering makes Medline your CAT DME 2024 09 |||

Durable medical equipment product guide Empowering health, enhancing lives. Elevate care with medi ... Toilet Safety Rails Easy-to-clean aluminum frame mounts securely onto bowl Item No. Description **MDS86100RF** Aluminum Foldable Toilet Safety Rails G30300 Aluminum Toilet Safety Rails with Padded A...

lang:en **score:26** filesize: 7.44 M page_count: 72 document date: 2024-03-28



[\[pdf\]](#) Guide Catalog

Durable medical equipment product guide Feb 8 2024 — Features Compact and lightweight it moves with the user effortlessly stabilizes easily locking casters Quiet motor provides 24 7 continuous CAT

DME020824 medline media catalog Docs MKT WP |||

Durable medical equipment product guide Empowering health, enhancing lives.

Elevate care with medi ... Toilet Safety Rails Easy-to-clean aluminum frame mounts securely onto bowl Item No. Description **MDS86100RF** Aluminum Foldable Toilet Safety Rails G30300 Aluminum Toilet Safety Rails with Padded A...

lang:en score:25 filesize: 61.04 M page_count: 72 document date: 2024-02-01



[\[pdf\]](#) Guide Catalog

Durable medical equipment product guide Jan 13 2025 — using the console buttons or remote control

Transfer resident to a walker wheelchair Attach components seat motor All In 4 Easy Steps CAT DME

Mobile medline media catalog Docs MKT WP |||

Durable medical equipment product guide Empowering health, enhancing lives.

Elevate care with medi ... Toilet Safety Rails Easy-to-clean aluminum frame mounts securely onto bowl Item No. Description **MDS86100RF** Aluminum Foldable Toilet Safety Rails G30300 Aluminum Toilet Safety Rails with Padded A...

lang:en score:25 filesize: 12.57 M page_count: 80 document date: 2025-01-09

[\[pdf\]](#) Guide Catalog

Durable medical equipment product guide Mobility 1st*** Compatible with Medline Bariatric Homecare

Bed MDR107004 only Overbed Tables Standard overbed tables various features and material

optionsCAT DME mobility1st ca 2024 09 CAT DME srsItid AfmBOOpM5eAcheBzGoCEi4L84UlWbaYFhB

bTm8c2POH1mlkAs0SkPp4 optionsmobility1st bTm8c2POH1mlkAs0SkPp4mobility1st AfmBOooJErLpm

Qz53J3qtCBXL7kLkJ9tUDA7HwSu538vW m2y8wELg0 m2y8wELg0mobility1st

AfmBOOpYepzXMNtgCs0mcNSK15cyvmzMSOQ9ebWCykHcHi2FZu2KNAnK

AfmBOOpYepzXMNtgCs0mcNSK15cyvmzMSOQ9ebWCykHcHi2FZu2KNAnKmobility1st AfmBOor0fDe

E2LleAaS688LUjKdN1t 6xgwrD8pplkLrIGDQ VGuvE8 AfmBOoppZxn8EzyX

4SVjgFk76wXLil3j6lpNnPJlyyZwLI9NvvllKoymobility1st VGuvE8mobility1st AfmBOoqg0t

oDwBsHTB9RJphbORsZUFB4ivfBCXCwSAOw43gMicK4ZJK optionsv1mobility1st

oDwBsHTB9RJphbORsZUFB4ivfBCXCwSAOw43gMicK4ZJKmobility1st

AfmBOoqbCwAjjbtqLA8JkBTfFXVmFfKoCFIBicqe5bR1zP8Lx NTnPg NTnPgmobility1st

AfmBOoqT6v9kQi0BTXUBStHsZH5dLvIpFtF0wBVg0qnIVSjMawJqZYM AfmBOoqRcvHmFVS3hS

TIDInAEjNWEDAwb3aGHRR0 GQCe30yp2SykUImobility1st

AfmBOoqT6v9kQi0BTXUBStHsZH5dLvIpFtF0wBVg0qnIVSjMawJqZYMmobility1st

AfmBOoquaFF67AoRdXlv 1hzikYMI9x rNFf fQvHekcehG GZurr62i optionsv1v1mobility1st

GZurr62imobility1st AfmBOorB tagJJBWwUzLCaEcWVIEMaC7BIAiHTWS5XjMa vaGdE54Lz

vaGdE54Lzmobility1st AfmBOoom2bqyCmo03PV2QbiCo3tns2Y y8bG3fu8jv2eW7EFGenRkM

AfmBOoqSIU5VohCiw xB0 xKae2EDqeT32komkchLYWzUwbwzNAI8Jf8mobility1st


optionsv1v1v1mobility1st


AfmBOoqjnIQBSplH6CVuepuPzrKUsfKTopsHep6kM9u12PuxCoiYTKlimobility1st
optionsv1v1v1v1mobility1st
AfmBOoq0EEHxvStG75ViJKziBEAQB64dbddrBQC90CSOTDq17F8DWEB2mobility1st
AfmBOoqetqwpe92p4Is h0RgT3G514jKgv6BS3PEi 7DAecnrJ7FPiBmobility1st AfmBOooNEecKvW9
nSqMv6ESb MDiejV4Q7ZFOu4yHTU2u4YvpVwqaoSmobility1st
AfmBOoq7LWfHqBbitLxMJEAAdSVA9WJ0BNBDNcA qbPbpf5VMzlsgGkTrmobility1st
AfmBOooKX3tW2crBCJXPEIHfbXbHdRlyVh2y50ky2WWi02yD1IC04IR2mobility1st
y8bG3fu8jv2eW7EFGenRkMmobility1st AfmBOoqXs53w riZeem wn7No9f9SkINDzQCUCRotJ68KPhd
wQdGdbO AfmBOorHkBmveff8Q9xJXUb6zfDurVGgVjYcMWpWdsGj4pvgyeiHvpWnmobility1st
optionsv1v1v1v1v1mobility1st wQdGdbOmobility1st AfmBOookzuE0EjbuP3r83KDz23OAx
vRz92XtLLvoiUDXcf cPDxpWAC cPDxpWACmobility1st AfmBOopxr3IPT
nH2tFVKXhYzoeq5O6QpHnk011FLnHXQqRTCaPRr5hL AfmBOoqV2Uz
VmalWOszDjllBnGfwm8zOk1LOhCWAHUyaENAh3Gigpmobility1st
nH2tFVKXhYzoeq5O6QpHnk011FLnHXQqRTCaPRr5hLmobility1st
AfmBOoojAF7RwFP0YEVsjd6zPXlrd0dsw8tnug80G44kaQvPgnt pZvv
AfmBOoqyu5VCmZK2jmi5WY5BlsYjhahrTlflwTjiOmUqHl8r6WI 04cUumobility1st
optionsv1v1v1v1v1v1mobility1st pZvmmobility1st AfmBOopyQmwYmpQfgHUQNUqTPNPZJvbNK2Pp
AqndOjCK8uQSowMzyj1 AfmBOor G37m69gUsgqWkt66nGh H1PSXndQJOaL
7ODvGbOYmpWYf8mobility1st AfmBOoru
yBNEbk9l7jdelQmk935yMLeXMPBICpkolqB62toklbCwbhymobility1st AfmBOoqz
AaWm4PCR7b3NOzbBGHNICqls1GrCBEBhTDCmQjqOchS0 ZKmobility1st 1stDurable optionsG2
101KRX1 Bath Chair Back 1 ea 57 G2 201KRX1 without Shower 400 lb 181 4 kg 13 75 20 5 12 E0240 Y
1 v1v1v1v1v1v1mobility1st AfmBOoo
XGuB3wa4pOteDGtK0zuEA7otR1d00bH7dpOfa7cxKvZzuNHgmobility1st
AqndOjCK8uQSowMzyj1mobility1st
AfmBOoqwcL2mZUwEgfvnOdxDZju28luyT2yS275TrKYIoSbORQZ5mtYN
AfmBOopWrcbXjDKpM3bsitYXhlwQ05m3fG ytxw3Za tN8fsPRrkDdE2mobility1st
AfmBOoqwcL2mZUwEgfvnOdxDZju28luyT2yS275TrKYIoSbORQZ5mtYNmobility1st AfmBOorD Z3w2ll9f
uMNckontU8ey l4gDgZkmneSWQiv44 Og1mCxD
AfmBOoou98BIP2EjTKlqxR0npwPOW9YJZsIWbNKHiscd7M p7wzin08ymobility1st
1 v1v1v1v1v1v1v1mobility1st Og1mCxDbmobility1st AfmBOooULr5
0oJESfSZvDKvsSup4JBYUFQXaipL9M9XgSe90iRKjym8
0oJESfSZvDKvsSup4JBYUFQXaipL9M9XgSe90iRKjym8mobility1st AfmBOoqSch BPPSYCvLvR2Ut1xjP
z2KP5rfX7rOrWELuvCXi0fRuJgV AfmBOooKpTGF0CoUrf7Z nYbeg8b5Ak0R2KII1YixNY1
2inLFBsIEQrmobility1st AfmBOorBBdWhWkRI
9JdhIO5qGL8hSSSTk4KpV85w0sKI0sxcaUylYg8mobility1st AfmBOopPxB3Z633TUF3g2E
Wzcc1fJJHxZD1zWVMY3g9zgHMXtXm20HFmobility1st AfmBOopT8AkPlnGC2WmG7JAZ CP
iOYojVYQrWvqq5P3nr30acJqmjumobility1st z2KP5rfX7rOrWELuvCXi0fRuJgVmobility1st
AfmBOopzNObmRA88JyD ZqPrNxpzLSwyF9SQV5fRGjD3cNAjL6JZz4HE
1 v1v1v1v1v1v1v1v1mobility1st AfmBOorfj0ytQUHKBMRG3jO
dtkXb4EVk5zOp4ZKL40HP5ot4Wjia2V6mobility1st
ZqPrNxpzLSwyF9SQV5fRGjD3cNAjL6JZz4HEmobility1st AfmBOooVumTOvliE9TncN KL
qDS8nzkPvvRQ2l3GgMHl5v3NuNd7lGG AfmBOoonpR0E8f mMvzZR0bRZsHZpE3JF
dto2MkVz6HpQwmZ3dHxRzmobility1st AfmBOopzzW7X5Nmm 5QW1wsAF2te98ewK9S6DmON1WINT
A0S60ncjvQmobility1st 1 v1v1v1v1v1v1v1v1mobility1st
qDS8nzkPvvRQ2l3GgMHl5v3NuNd7lGGmobility1st AfmBOorluVR5j9zHB
oVkd0TlODu9FZ8ilMddsp44csimBCHNw88b0Hm
AfmBOopUd0wgVz357jEXzD6c0xYo0HZAxrldofx1ONyMGCpjkmTookqpmobility1st AfmBOoqhJd
IJHzgwpCo4KK88sefQudaW6AtJftCKGklk10CmGJx r2mobility1st AfmBOookf
jk1qtvtAyyq5C6VrMfKKhaWN WgNEvZYdM 2eCHitFoe51mobility1st ||| ||| Durable medical equipment
product guide Empowering health, enhancing lives. Elevate care with medical equipment that goes the
distance High-quality engineering, long-lasting performance, and a comprehensive and unique product

[illegible]

[illegible]

Durable medical equipment product guide Empowering health, enhancing lives.
Elevate care with medi ... Toilet Safety Rails Easy-to-clean aluminum frame mounts
securely onto bowl Item No. Description **MDS86100RF** Aluminum Foldable Toilet
Safety Rails G30300 Aluminum Toilet Safety Rails with Padded A...
lang:en score:24 filesize: 7.41 M page count: 72 document date: 2024-03-28

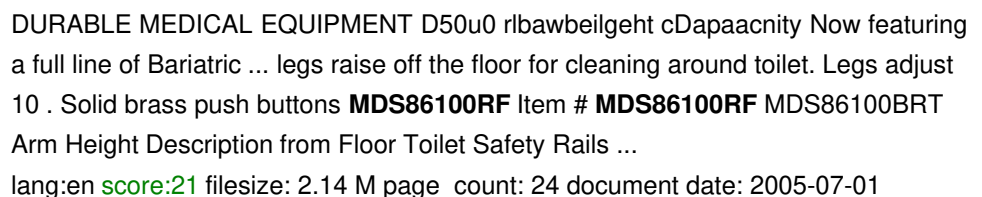

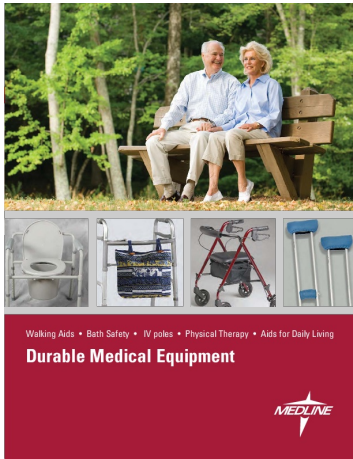


[\[pdf\]](#) Accessories Catalog

38051 Medline DME ps page 1 28 @ Normalize 2 Catalog 2005PatientAidsBroch Jim patientaids medline 80 literature dme

Walking Aids Bath Safety IV poles Physical Therapy Aids for Daily Living Durable Medical Equipme ... ly to toilet. Rail legs raise off the floor for cleaning around toilet. Legs adjust 10 . Item # **MDS86100RF** MDS86100BRT Description Arm Height From Floor Depth Toilet Safety Rails 28 to 32 17...

lang:en **score:21** filesize: 2 M page_count: 28 document date: 2008-06-06

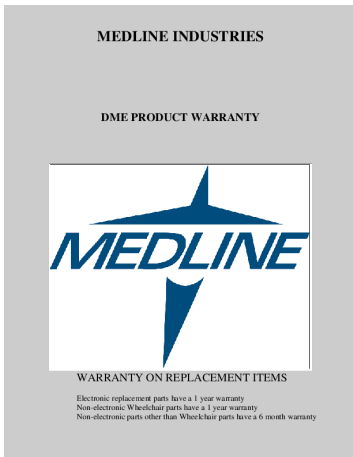

[\[pdf\]](#) Declaration of Conformity Warranty

Microsoft Word MEDLINE Warranties csr version 07 14 2010 doc shopton INDUSTRIESMEDLINE 100020119 Warrantyassets unilogcorp 187 ITEM DOC Warranty |||

MEDLINE INDUSTRIES DME PRODUCT WARRANTY WARRANTY ON REPLACEMENT ITEMS Electronic replacement parts h ... er. BATH PRODUCTS THIS WARRANTY COVERS THE MEDLINE AND GUARDIAN PRODUCT LINES MDS80315 MDS80316 **MDS86100RF** MDS80314 6 MONTH WARRANTY: Handles Wearable parts Grips Screws LIMITED LIFETIME W...

lang:en **score:21** filesize: 1.15 M page_count: 21 document date: 2010-09-16




[\[pdf\]](#) Guide Catalog

Durable medical equipment product guide 5 days ago — Rollators provide users with both a mobility aid and seat to rest on when tired Most rollators allow handle adjustment for proper fit comfort To CAT DME srsItd AfrmBOooNI7uSRz BlmNUug9KlZ9uR s1hUkFyeSnkbnvkBTzMX3N0ISC mobility1st ca 2024 09 |||

Durable medical equipment product guide Empowering health, enhancing lives. Elevate care with medi ... Toilet Safety Rails Easy-to-clean aluminum frame mounts securely onto bowl Item No. Description **MDS86100RF** Aluminum Foldable Toilet Safety Rails G30300 Aluminum Toilet Safety Rails with Padded A...

lang:en **score:21** filesize: 7.41 M page_count: 72 document date: 2024-03-28

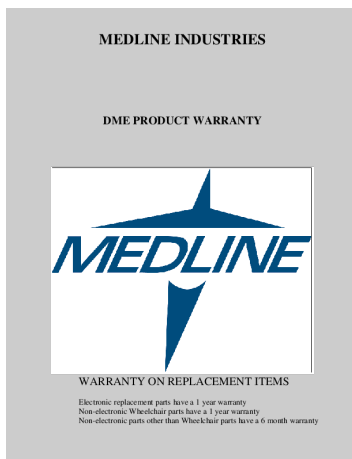

[\[pdf\]](#) Declaration of Conformity Warranty

Microsoft Word MEDLINE Warranties csr version 07 14 2010 doc shopton INDUSTRIES prod

1831725212 i sears s d mp tc 10048557

MEDLINE INDUSTRIES DME PRODUCT WARRANTY WARRANTY ON REPLACEMENT ITEMS Electronic replacement parts h ... er. BATH PRODUCTS THIS WARRANTY COVERS THE MEDLINE AND GUARDIAN PRODUCT LINES MDS80315 MDS80316 **MDS86100RF** MDS80314 6 MONTH WARRANTY: Handles Wearable parts Grips Screws LIMITED LIFETIME W...

lang:en score:21 filesize: 1.15 M page_count: 21 document date: 2010-09-16




[\[pdf\]](#)

39489 Medline Fall Prevention MKT209341LIT410R ps page 7 @ Preflight 2 Layout 1 Jim medline

product literature MedSurg

FALL PREVENTION ALARMS WEARABLES BATH SAFETY LIFTS B E D / B E D S I D E A I D S F O R D A I ... 10 25cm . Item No. Description Arm Height From Floor Depth Distance Between Arms Outside **MDS86100RF** Toilet Safety Rails 28 to 32 171/2 171/2 to 201/2 71cm x 81cm 45cm 45cm x 52cm ...

lang:en score:21 filesize: 5.23 M page_count: 16 document date: 2009-09-14

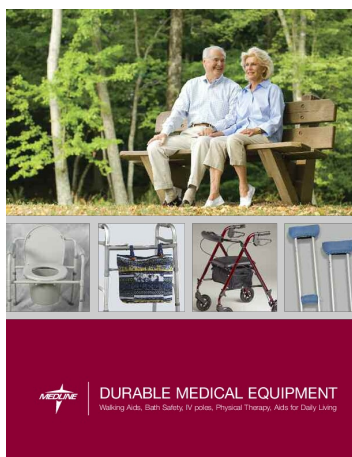

[\[pdf\]](#) Accessories

2005PatientAidsBroch 2 PatientAids medline product literature DME

DURABLE MEDICAL EQUIPMENT Walking Aids, Bath Safety, IV poles, Physical Therapy, Aids for Daily Livi ... floor for cleaning around toilet. Legs adjust 10 25cm . 250lb 113kg weight capacity. Item # **MDS86100RF** MDS86100BRT Description Arm Height From Floor Toilet Safety Rails 28 to 32 71cm to ...

lang:en score:20 filesize: 1.09 M page_count: 28 document date: 2010-01-08



[pdf]

38797 Medline Fall prevention MKT208354LIT410R qxd Layout 1 Jim LIT410R MKT208354 cf1 bettymills

product more info

FALL PREVENTION ALARMS WEARABLES BATH SAFETY LIFTS B E D / B E D S I
D E A I D S F O R D A I ... 10 25cm . Item No. Description Arm Height From Floor
Depth Distance Between Arms Outside **MDS86100RF** Toilet Safety Rails 28 to 32
171/2 171/2 to 201/2 71cm x 81cm 45cm 45cm x 52cm ...

lang:en score:19 filesize: 1.92 M page_count: 16 document date: 2009-01-26



[pdf] Guide Catalog

LIT474 CAT DME 19W6094872 medline media catalog Docs MKT MedEmporium |||

Durable medical equipment product guide Empowering health, enhancing lives.

Elevate care with medi ... Toilet Safety Rails Easy-to-clean aluminum frame mounts
securely onto bowl Item No. Description **MDS86100RF** Aluminum Foldable Toilet
Safety Rails G30300 Aluminum Toilet Safety Rails with Padded A...

lang:en score:16 filesize: 48.72 M page_count: 80 document date: 2025-01-02

Bath Safety

Elevated Toilet Seats

- Makes using the toilet easier for those with trouble bending
- Available in a variety of locking mechanisms and mounting methods



MDS80326
Extra Wide with Legs
600 lb 272 kg Capacity
Steel legs go to the floor and soft foam armrests provide maximum support in getting on and off the toilet.



G30300
Locking with Padded Arms
Padded armrests provide getting on and off toilet locking device secures and holds toilet in locked position.



G30300
Locking with Molded Arms
Arms molded to locking on and off the toilet. Locking device secures the toilet to the toilet.



MDS80326
Universal Locking
Simple to use universal locking device. Steel legs go to the floor and soft foam armrests provide maximum support in getting on and off the toilet.



MDS80326
Hinged Roll-On
Roll-on design. Hinged seat and soft foam armrests provide maximum support in getting on and off the toilet. Hinged seat and soft foam armrests provide maximum support in getting on and off the toilet.



G30300
Basic Roll-On Seat
Simple design. Hinged seat and soft foam armrests provide maximum support in getting on and off the toilet. Hinged seat and soft foam armrests provide maximum support in getting on and off the toilet.

| Item Number | Description | Weight Capacity (lb/kg) | Dimensions (in/cm) | Qty |
|-------------|--------------------------|-------------------------|--------------------|------|
| MDS80326 | Extra Wide with Legs | 600/272 | 45/50 | Each |
| G30300 | Locking with Padded Arms | 300/136 | 45/50 | Each |
| G30300 | Locking with Molded Arms | 300/136 | 45/50 | Each |
| MDS80326 | Universal Locking | 300/136 | 45/50 | Each |
| MDS80326 | Hinged Roll-On | 300/136 | 45/50 | Each |
| G30300 | Basic Roll-On Seat | 300/136 | 45/50 | Each |

Click here to view detailed product descriptions and prices. Call 800-397-5899 to order today!

[pdf] Catalog

Medline Bath Safety Catalog Product Toilet Frame Rails G30300H MDS86100RF Active Forever product
catalog medline safety rails activeforever |||

Bath Safety Elevated Toilet Seats M akes using the toilet easier for those with trouble
bending A vailable in a variety of locking mechanisms and mounting methods
MDS80326 Extra Wide with Legs 600 lb 272 kg Capacity Steel legs go to the floor
and soft foam armrests provide maximum support in ...

lang:en score:15 filesize: 87.21 K page_count: 1 document date: 2018-07-03



MEDICARE
Part A Interim
Part B Carrier
DMEPOS Regional Center

June 6, 2018

David Jacobs, President
Medline Industries, Inc.
One Medline Place
Mundelein, IL 60060-4455

Re: Toilet Safety Rails (Model MDS86100RF)

Dear Mr. Jacobs:

This letter is in response to your recent inquiry for coding verification of the above listed product(s) manufactured by your company. The Statistical Analysis Practice Medical Equipment Regional Center (SAPMERC) has reviewed the documentation and information submitted for HCPCS Coding. The SAPMERC conducts reviews of products to determine the correct HCPCS code(s) of DMEPOS product(s) for Medicare billing.

It is our determination that the Medline HCPCS code(s) to bill the four Durable Medical Equipment Regional Centers (DMERCs) is/are:

E0243 Toilet rail, each.

This HCPCS coding decision applies to submitted product(s) as presented to and reviewed by the SAPMERC. Any modifications to the product(s) could change the HCPCS code and would need to be reviewed for coding verification. The assignment of a HCPCS code to the product(s) listed in no way is construed as an approval or endorsement of the product(s) by SAPMERC or Medicare, nor does it imply or guarantee claim reimbursement or coverage. For questions regarding claim coverage or reimbursement please contact your regional DMERC.

Should you disagree with this coding decision, a re-review of the product(s) can be requested. The SAPMERC will provide a re-review if the request is made within 45 days of the date of this letter and additional documentation is provided supporting the request. If a request for a re-review is made after 45 days, the request is treated as a new Coding Verification Review and a complete application must be submitted along with the additional documentation supporting the request.

Palmetto CBA
Statutory Authority: Florida Health Care Insurance Code, § 400.01
One Office Blvd., Suite 100, Gainesville, FL 32601-1944
A CMS Contracted Intermediary and Carrier

[Click here to view detailed product descriptions and prices.](#) Call 800-397-5899 to order today!

[\[pdf\]](#)

HCPCS Coding for Medline Toilet Safety Rail Product Toilet Frame Rails G30300H MDS86100RF Active

Forever product hcpcs letter activeforever |||

Click here to view detailed product descriptions and prices. Call 800-397-5899 to order today ...

lang:en score:14 filesize: 139.67 K page_count: 1 document date: 2018-07-03