

[Manuals.plus](#) /

> [Sonic Alert](#) /

> Sonic Alert SBP100 Portable Vibrating Alarm Clock User Manual

Sonic Alert RA49238

Sonic Alert SBP100 Portable Vibrating Alarm Clock User Manual

Model: RA49238 | Brand: Sonic Alert

INTRODUCTION

The Sonic Alert SBP100 Portable Vibrating Alarm Clock is designed to provide a reliable wake-up solution for individuals who are heavy sleepers, deaf, or hard of hearing, especially while traveling. This compact device features a powerful built-in vibrator, a loud audible alarm, and a clear digital display. Its portable design allows for use in various settings, ensuring you wake up on time for important events.



Image: The Sonic Alert SBP100 Portable Vibrating Alarm Clock in white, showing its compact design and digital display.

PACKAGE CONTENTS

- Sonic Shaker™ Travel Alarm Clock (White)
- 3 AAA batteries (included)
- Travel Case (may vary by package)

SETUP

1. Battery Installation

The SBP100 alarm clock is powered by 2 AA batteries for the vibrator and 1 AAA battery for the clock display. All required batteries are included.

1. Locate the battery compartment on the underside of the alarm clock.
2. Open the battery compartment cover.
3. Insert 2 AA batteries and 1 AAA battery, ensuring correct polarity (+/-).
4. Close the battery compartment cover securely.

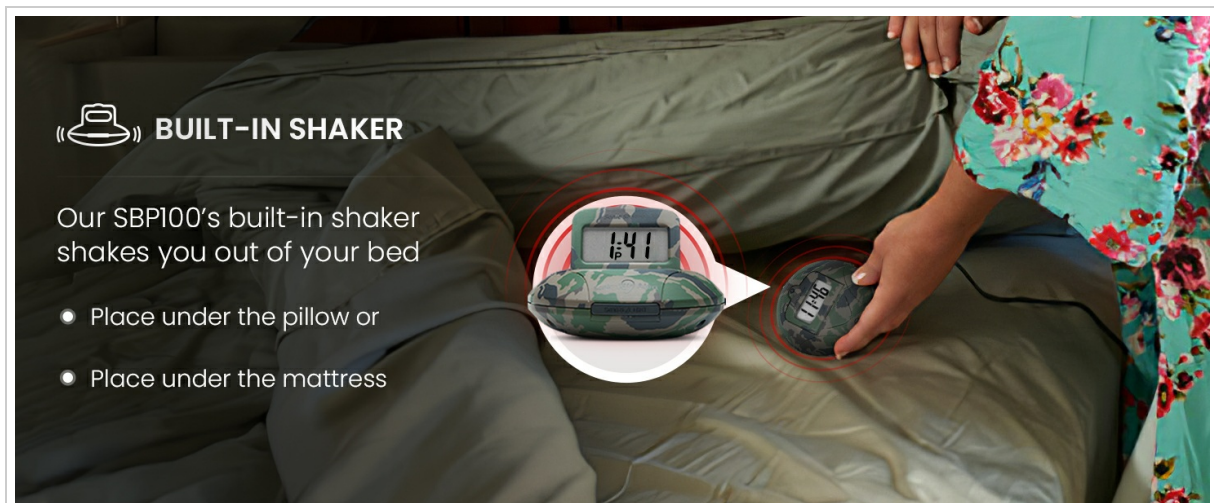


Image: The battery compartment of the Sonic Alert SBP100, illustrating where to insert the AA and AAA batteries.

2. Initial Placement

The SBP100 is designed for versatile placement to ensure effective waking.

- **Under Pillow:** Place the alarm clock under your pillow for direct vibration. Use the included pillow strap to secure it.
- **Under Mattress:** For stronger vibration, place the alarm clock between your mattress and box spring.
- **Bedside Table:** The alarm can also be placed on a bedside table, utilizing its loud buzzer.
- **Portable Use:** The compact design allows it to be clipped to clothing, a backpack, or carried in a pocket.

Sonic Alert

Specialized Hearing Solutions

THE PERFECT TRAVEL ALARM CLOCK



Built-In Shaker



12/24 Hour Format



Pop-Up Display



Battery Operated



Snooze/Light Button



Portable



Image: The Sonic Alert SBP100 positioned under a pillow, secured by its strap, for silent vibrating alarm.



Clip to Pillow



Clip to Pocket



Clip to Backpack



Clip to Clipboard

Image: Various ways to carry and secure the SBP100, including clipping it to a pillow, placing it in a pocket, attaching it to a backpack, or using it on a clipboard.

OPERATING INSTRUCTIONS

1. Setting Time and Alarm

The time and alarm controls are located on the top of the display unit when it is popped up.

Easy To Use Features

Enjoy a hassle-free experience with our portable alarm clock



Sonic Alert
Specialized Hearing Solutions

Sonic Shaker
11:46
Sonic Alert

-  Easy to use and operate
-  Pop-up Display with Light
-  Easy to Set Function Keys

Image: A detailed view of the SBP100's control buttons for setting time and alarm, located above the digital display.

1. **To Set Time:** Press the **TIME SET** button. Use the **HOUR** and **MINUTE** buttons to adjust the current time.
2. **To Set Alarm:** Press the **ALARM SET/OFF** button. Use the **HOUR** and **MINUTE** buttons to set your desired alarm time.
3. **To Turn Alarm On/Off:** Use the **ON/OFF** switch located on the right side of the clock. Slide it to **ON** to activate the alarm or **OFF** to disable it.

2. Alarm Modes

The SBP100 offers three alarm modes: vibration only, buzzer only, or both. The mode switch is located on the left side of the alarm clock.

Built-in Super Shaker

The bed shaker is a bed vibrating unit that vibrates when the alarm goes off to ensure that everyone wakes up on time.



Image: Illustration of the three alarm modes available on the SBP100: vibration, audible alarm, or a combination of both.

- **VIB:** Activates the vibrator only. Ideal for silent waking or for individuals who are deaf or hard of hearing.
- **BUZZ:** Activates the loud 90dB pulsating buzzer only.
- **VIB BUZZ:** Activates both the vibrator and the loud buzzer simultaneously for maximum alert.

3. Snooze Function and Display Light

The SNOOZE/LIGHT button is located on the front of the alarm clock.

- **Snooze:** When the alarm sounds, press the **SNOOZE/LIGHT** button to activate a 4-minute snooze period. This can be reactivated for up to one hour from the original alarm setting.
- **Display Light:** Press the **SNOOZE/LIGHT** button to illuminate the display for a few seconds, allowing you to view the time in low-light conditions. This manually lit display helps conserve battery life.

4. Test Function

A test button is located on the top of the alarm clock, allowing you to check the vibrator, buzzer, and battery condition.

1. Ensure the alarm is set to **ON**.
2. Press the **TEST** button. The alarm will activate according to the selected mode (VIB, BUZZ, or VIB

BUZZ).

3. If the alarm does not activate or is weak, check the battery condition.

MAINTENANCE

1. Battery Replacement

When the display becomes dim or the alarm functions weaken, it is time to replace the batteries. Refer to the "Battery Installation" section under Setup for instructions.

2. Cleaning

Wipe the alarm clock with a soft, dry cloth. Do not use liquid cleaners or abrasive materials.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Alarm does not sound or vibrate.	<ul style="list-style-type: none">• Alarm is set to OFF.• Incorrect alarm mode selected.• Dead or low batteries.	<ul style="list-style-type: none">• Ensure the Alarm ON/OFF switch is set to ON.• Verify the mode switch (VIB, BUZZ, VIB BUZZ) is set correctly.• Replace all batteries with fresh ones. Use the TEST button to check.
Display is dim or blank.	Low or dead AAA battery.	Replace the AAA battery.
Alarm is not loud enough.	Alarm mode set to VIB only.	Change the mode switch to BUZZ or VIB BUZZ.

SPECIFICATIONS

- **Model Number:** RA49238
- **Product Dimensions:** 4.71 x 9.42 x 7.07 inches
- **Item Weight:** 7.1 ounces
- **Power Source:** 2 AA batteries (for vibrator), 1 AAA battery (for clock) - all included
- **Display:** 0.75" Pop-up digital display with convenience light
- **Alarm Volume:** Up to 90dB (buzzer)
- **Snooze Duration:** 4 minutes (continuous for up to 1 hour)

WARRANTY AND SUPPORT

Warranty Information

Sonic Alert provides a 1-year warranty for this product. The warranty covers manufacturing defects and ensures the product functions as intended under normal use. Please note the following conditions:

- The warranty is valid for products purchased from authorized sellers.
- The warranty may not cover products received through state or governmental programs.

- International product purchases may not be covered.
- Within the warranty period, your product will typically be replaced once.
- Proof of purchase and authorization for return are required for warranty claims.

Customer Support

For technical assistance, warranty claims, or any questions regarding your Sonic Alert SBP100, please contact Sonic Alert customer service. Contact details can typically be found on the official Sonic Alert website or within the product packaging.

For more information, you may visit the [Sonic Alert Store on Amazon](#).

PRODUCT VIDEO OVERVIEW



Video: An instructional video demonstrating the features and operation of the Sonic Alert SBP100 Portable Vibrating Alarm Clock.