

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [TRANQUILITY](#) /

› [Tranquility Slimline Original Adult Disposable Brief User Manual - X-Small \(Model CSPK-2166\)](#)

TRANQUILITY CSPK-2166

Tranquility Slimline Original Adult Disposable Brief User Manual

Model: CSPK-2166 | Size: Adult X-Small

PRODUCT OVERVIEW

The Tranquility Slimline Original Adult Disposable Briefs are designed to provide superior absorbency and comfort with a slimmer profile, making them suitable for everyday use. These briefs feature a cloth-like back sheet and micro-hook closure tabs for a secure and comfortable fit. They are latex-free and engineered with advanced technologies to promote skin health and prevent leaks.



The packaging for Tranquility Slimline Original Adult Disposable Briefs, X-Small size, indicating a quantity of 10 and heavy protection.

Key Features:

- **Premium Slim Fit:** A classic full-fit brief with a discreet profile for versatile use.
- **Peach Mat Core:** Designed to combat odor, promote dry skin, and aid in pH neutralization, ideal for sensitive skin.
- **Kufguard Technology:** Inner leg cuffs and leg elastics channel fluid into absorbent zones, preventing leaks.
- **Dual Wetness Indicator:** Changes from yellow to blue when moisture is present, signaling when a change is needed.
- **Micro-Hook Closure Tabs:** Securely fasten the brief and allow for refastening.
- **Latex-Free:** Suitable for individuals with latex sensitivities.

APPLICATION INSTRUCTIONS

Proper application ensures maximum comfort and leak protection. Follow these steps for a secure fit:

1. **Prepare the Brief:** Lay the brief flat and gently fold it lengthwise to activate the leg gathers.
2. **Position the Brief:** With the user lying on their side, slide the brief between the legs from front to back, ensuring the Kufguards® inner leg cuffs are upright and positioned around the groin area.
3. **Secure the Tabs:** Bring the front panel of the brief up. For the first fitting, grip both the blue target tab and the patterned tab together and place them on the front of the brief. Ensure the brief fits snugly with the tape tabs placed on the hips.
4. **Adjust for Fit:** Check for any gaps around the legs and adjust the tabs as needed. The refastenable tabs allow for multiple adjustments to achieve a comfortable and secure fit.

Refastenable Tape Tab System




Step 1

When fitting (or applying) the brief for the first time, grip both the Blue Target Tab and the Patterned Tab together and place on the front of brief for a proper fit.



Step 2

To remove or open the brief, grip the Patterned Tab alone, pulling away from Blue Target Tab.

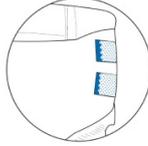

Step 3

To replace or close the brief, return the Patterned Tab to the Blue Target Tab.

Good Fit
Should fit like underwear



Brief should fit snugly with the placement of the tape tabs **on the hips** and no gaps around the legs.

A three-step diagram illustrating how to apply and adjust the refastenable tape tabs for a proper fit, showing how to grip and place the tabs.



Briefs Size Guide

Size	Waist/Hip	Weight	Women's Pant Size	Men's Pant Size
Youth	16" - 22"	28 - 42 lbs	N/A	N/A
X-Small	18" - 26"	42 - 90 lbs	0 - 3	18 - 26
Small	24" - 32"	85 - 100 lbs	4 - 6	24 - 32
Medium	32" - 44"	100 - 150 lbs	8 - 12	32 - 44
Large	44" - 58"	150 - 200 lbs	14 - 18	44 - 58
X-Large	56" - 64"	200 - 250 lbs	18 - 24	56 - 64
2XL	60" - 80"	250+ lbs	26 - 28	60+
3XL	64" - 96"	250+ lbs	28+	64+
4-5XL	70" - 108"	250+ lbs	28+	70+

Illustrations demonstrating a good fit for the brief, with tape tabs on the hips and no gaps around the legs, contrasted with examples of briefs that are too big or too small.

USAGE GUIDELINES

The Tranquility SlimLine Original Briefs are designed for effective containment and skin protection. Observe the following for optimal performance:

- **Wetness Indicator:** The brief features a Dual Wetness Indicator that turns from yellow to blue when moisture is present. This visual cue helps determine when a change is necessary.
- **Change Frequency:** For optimal effectiveness and skin health, it is recommended to change the brief every 8 hours, or sooner if the wetness indicator signals saturation.
- **Absorbent Capacity:** These briefs have a capacity of 10.9 oz / 322 mL.



An inside view of the brief highlighting key features: Soft, Cloth-like Outer Layer, Micro-hook Closure Tabs, Dual Wetness Indicator, Gentle Leg Gathers, Kufguards Inner Leg Cuffs, and Latex-Free.



An image showing the brief and its stated absorbent capacity of 10.9 oz / 322 mL, equivalent to 1 to 1.5 cups of normal urine void.

DISPOSAL AND CARE

Proper disposal and hygiene are important for user comfort and environmental considerations.

- **Disposal:** The briefs are designed for easy, mess-free removal. When finished with wear, rip apart the side seams (micro-hook tear-away seams) and pull the protective brief away from the body. Roll the soiled brief inward and secure with the tabs before discarding in a waste receptacle. Do not flush.
- **Storage:** Store unused briefs in a cool, dry place away from direct sunlight and moisture.
- **Skin Care:** Regularly clean and dry the skin to maintain hygiene and prevent irritation.

TROUBLESHOOTING

If you encounter issues with your Tranquility Slimline Original Briefs, consider the following:

- **Leaks:**
 - **Poor Fit:** Ensure the brief is snug around the waist and legs. Refer to the 'Application Instructions' and 'Size Guide' to confirm the correct size and proper tab placement. Gaps around the legs can lead to leaks.
 - **Kufguards Position:** Verify that the Kufguards® inner leg cuffs are upright and properly positioned to contain fluid.

- **Capacity Exceeded:** If leaks occur frequently before the wetness indicator changes, the brief's capacity may be insufficient for the user's needs. Consider a brief with higher absorbency or more frequent changes.

- **Discomfort or Bunching:**

- **Incorrect Size:** A brief that is too large or too small can cause bunching or pinching. Consult the 'Size Guide' to ensure the appropriate size is being used.
- **Improper Application:** Re-apply the brief, ensuring it is smoothed out and the tabs are secured evenly without over-tightening.

- **Tabs Tearing:**

- **Over-tightening:** Avoid pulling the tabs too tightly during application. The brief should be snug but not strained.
- **Rough Handling:** Apply and adjust tabs gently.

PRODUCT SPECIFICATIONS

Tranquility Briefs Size Guide

Size	Waist/Hip	Weight	Women's Pant Size	Men's Pant Size
Youth	16" - 22"	28 - 42 lbs	N/A	N/A
X-Small	18" - 26"	42 - 90 lbs	0 - 3	18 - 26
Small	24" - 32"	85 - 100 lbs	4 - 6	24 - 32
Medium	32" - 44"	100 - 150 lbs	8 - 12	32 - 44
Large	44" - 58"	150 - 200 lbs	14 - 18	44 - 58
X-Large	56" - 64"	200 - 250 lbs	18 - 24	56 - 64
2XL	60" - 80"	250+ lbs	26 - 28	60+
3XL	64" - 96"	250+ lbs	28+	64+
4-5XL	70" - 108"	250+ lbs	28+	70+

Absorbent Capacity Under Pressure
10.9 oz. / 322 mL

Normal urine void is 1 to 1½ cups

A table detailing the size guide for Tranquility briefs, including waist/hip measurements, weight ranges, and women's/men's pant sizes for various brief sizes from Youth to 4-5XL.

General Product Details:

- **Model Number:** CSPK-2166
- **Product Dimensions:** 21 x 15.75 x 12 inches
- **Item Weight:** 1.21 Pounds (for the X-Small pack)
- **Department:** Unisex Adult
- **Manufacturer:** Principle Business Enterprises
- **Incontinence Protector Type:** Adult Diaper
- **Age Range (Description):** Adult

WARRANTY INFORMATION

Specific warranty details are not provided with this product. For any concerns regarding manufacturing defects or product performance, please retain your proof of purchase and contact the manufacturer directly. Tranquility products are designed and manufactured with quality in mind.

CUSTOMER SUPPORT

For further assistance, questions, or feedback regarding your Tranquility Slimline Original Adult Disposable Briefs, please contact Principle Business Enterprises, the manufacturer of Tranquility products. You can typically find contact information on the product packaging or by visiting the official Tranquility website.

Always refer to the latest information provided by the manufacturer for the most accurate and up-to-date support resources.