

Mölnlycke 0925438

Mölnlycke Mefix Self-Adhesive Fabric Tape User Manual

Model: 0925438

1. INTRODUCTION

Mefix self-adhesive fabric tape is designed for secure fixation of dressings, gauze, tubes, and cannulas. Its porous, nonwoven fabric allows air and water vapor to pass through, promoting skin breathability. The stretchable material conforms to body contours, ensuring comfort and unrestricted movement. Mefix is non-absorbent, preventing fluid strike-through, and provides gentle yet secure fixation.

Key Features:

- Designed to secure gauze and wound dressings to the user's body.
- Porous and non-absorbing fabric minimizes fluid strike-through.
- Skin-friendly adhesive layer allows skin to breathe.
- Latex-free.

2. SAFETY INFORMATION

Important Considerations:

- Avoid using on broken skin or open wounds.
- Change dressings regularly to maintain a clean wound environment.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

3. SETUP

Unpacking and Preparation:

Mefix tape comes in a roll. The release paper is pre-measured in 4-inch (10cm) intervals and features a perforated split backing for easy removal and application. The tape can be easily cut to the desired shape and size.



Image: Mefix Self-Adhesive Fabric Tape roll. This image displays the product as a roll of white fabric tape with a branded core, partially unrolled.



Image: Mefix Self-Adhesive Fabric Tape product box. The box features the Mefix logo and multilingual descriptions of the product.

4. APPLICATION INSTRUCTIONS

Step-by-Step Guide:

Before applying, ensure the skin around the wound area is clean and dry.

4.1. Prepare the Dressing

Clean the wound area thoroughly. Apply the primary dressing (e.g., gauze) to the wound as instructed by a healthcare professional.

4.2. Cut the Tape

Cut a piece of Mefix tape to the appropriate size, ensuring it is large enough to cover the primary dressing and adhere securely to the surrounding skin.

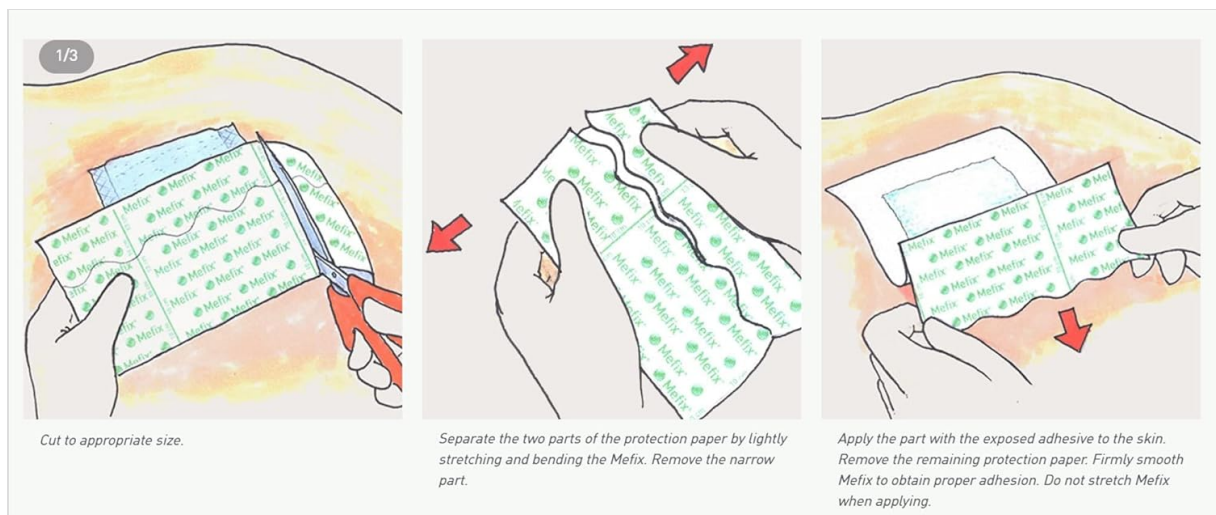


Image: Mefix Application Steps. This image illustrates the process of cutting the tape, peeling the backing, and applying it to a dressing.

4.3. Prepare for Application

Separate the two parts of the protection paper by lightly stretching and bending the Mefix tape. Remove the narrow part of the release paper.

4.4. Apply the Tape

Apply the part of the tape with the exposed adhesive to the skin, securing the dressing. Do not stretch Mefix tape when applying to avoid tension on the skin.

4.5. Secure the Dressing

Remove the remaining protection paper. Firmly smooth Mefix to obtain proper adhesion, ensuring there are no wrinkles or air bubbles.

5. CARE AND STORAGE

Mefix tape is designed for single use. Store the unused roll in a cool, dry place, away from direct sunlight and extreme temperatures. Ensure the packaging is sealed to protect the tape from contamination.

6. TROUBLESHOOTING

Common Issues and Solutions:

- **Tape not adhering properly:** Ensure skin is clean, dry, and free of oils or lotions before application. Avoid stretching the tape during application.
- **Skin irritation:** If redness, itching, or irritation occurs, remove the tape immediately and consult a healthcare professional. Ensure the tape is not applied too tightly.
- **Difficulty removing backing paper:** The release paper is perforated. Gently stretch and bend the tape to help separate the backing.

7. SPECIFICATIONS

Feature	Detail
Brand	Mölnlycke
Model Number	0925438
Material	Nonwoven fabric, paper backing
Size	2 inches (5cm) wide x 11 yards (10m) long
Product Dimensions (Roll)	3 x 2.9 x 2.2 inches
Weight	2.82 ounces
Latex Free	Yes
Porous	Yes
Stretchable	Yes
First Available Date	September 18, 2006

8. WARRANTY AND SUPPORT

For product inquiries, technical support, or to report any issues, please contact Mölnlycke customer service. Specific warranty information for this product is typically provided with the purchase documentation or can be obtained directly from the manufacturer.

Manufacturer Contact:

Please refer to the product packaging or Mölnlycke's official website for the most current contact information.