

Welch Allyn 72200

Welch Allyn #72200 3.5V NiCad Rechargeable Battery User Manual

INTRODUCTION

This manual provides essential instructions for the proper use, installation, and maintenance of the Welch Allyn #72200 3.5V NiCad Rechargeable Battery. Please read this manual thoroughly before using the product to ensure safe and efficient operation.

SAFETY INFORMATION

Observe the following safety precautions to prevent injury or damage to the product:

- Use only with compatible Welch Allyn devices as specified in the "Product Overview" section.
- Do not attempt to disassemble, crush, or puncture the battery.
- Do not expose the battery to extreme temperatures, direct sunlight, or fire.
- Keep the battery away from water or other liquids.
- Dispose of the battery according to local regulations for NiCad batteries. Do not dispose of in household waste.
- If the battery shows signs of damage, leakage, or overheating, discontinue use immediately.

PRODUCT OVERVIEW

The Welch Allyn #72200 is a 3.5V Nickel-Cadmium (NiCad) rechargeable battery designed for use with specific Welch Allyn diagnostic instrument handles. It is color-coded black for identification.



An image showing the Welch Allyn #72200 3.5V NiCad rechargeable battery next to its blue and white retail packaging. The battery itself is white with black text, indicating 'Welch Allyn', 'Ni-Cd', '3.5V', and 'Made in China'. The packaging also displays '#72200' and '3.5V RECHARGEABLE BATTERY'.

Key Features:

- **Battery Type:** Nickel-Cadmium (NiCad)
- **Voltage:** 3.5V
- **Model Number:** 72200
- **Compatibility:** Designed for Welch Allyn handles including models 71000, 71010, 71015, 71020, 71050, 71051, 71054, 71055, and 71670.
- **Color Coding:** Black

SETUP

1. **Unpacking:** Carefully remove the battery from its packaging. Inspect for any visible damage.
2. **Initial Charge:** For optimal performance, fully charge the battery before its first use. Refer to your Welch Allyn device's manual for specific charging instructions.

3. Installation:

- Ensure your Welch Allyn instrument handle is turned off.
- Locate the battery compartment on your instrument handle.
- Insert the #72200 battery, ensuring correct polarity (+ and - terminals align). The battery is designed to fit snugly.
- Secure the battery compartment cover.

OPERATING INSTRUCTIONS

Once installed and fully charged, the Welch Allyn #72200 battery powers your compatible Welch Allyn diagnostic instrument handle. The operation of the instrument itself is detailed in its respective user manual.

Charging:

- Always use the original Welch Allyn charger or a compatible charging base designed for your instrument handle.
- Connect the instrument handle with the installed battery to the charger.
- Monitor the charging indicator on your instrument handle or charger. Charging times may vary.
- Avoid overcharging for extended periods, although modern chargers typically have overcharge protection.

MAINTENANCE

- **Cleaning:** Keep the battery terminals clean and free from dust or corrosion. Use a dry, soft cloth for cleaning.
- **Storage:** If storing the battery for an extended period, ensure it is partially charged (around 50%) and stored in a cool, dry place. Avoid storing fully charged or fully depleted.
- **Battery Life:** NiCad batteries can develop a "memory effect" if repeatedly partially discharged and recharged. To mitigate this, occasionally allow the battery to fully discharge before recharging it completely.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Battery not charging	Incorrect installation; faulty charger; battery at end of life; dirty terminals.	Re-install battery correctly; check charger functionality; clean battery and charger terminals; replace battery if old.
Short operating time after full charge	Battery memory effect; battery at end of life.	Perform a full discharge/recharge cycle; replace battery.
Instrument not powering on	Battery not charged; incorrect installation; instrument fault.	Charge battery; check installation; consult instrument manual for troubleshooting.

If issues persist, contact Welch Allyn customer support or a qualified service technician.

SPECIFICATIONS

Attribute	Detail
Brand	Welch Allyn
Model Number	72200
Battery Type	NiCad (Nickel-Cadmium)
Voltage	3.5V
Item Weight	3.2 ounces
Package Dimensions	3.62 x 1.93 x 1.34 inches
Recommended Use	Diagnostic Instruments
First Available Date	December 18, 2006

WARRANTY AND SUPPORT

Specific warranty information for this Welch Allyn #72200 battery is typically provided with the original purchase documentation or can be obtained directly from Welch Allyn or your authorized dealer. For technical support or service inquiries, please refer to the contact information provided by Welch Allyn or your product supplier.

Note: This product is listed as a "Replacement (Non-OEM)" in some contexts. Ensure compatibility and warranty implications with your specific Welch Allyn device.


© 2023 Welch Allyn. All rights reserved. Information subject to change without notice.

Related Documents - 72200

 Upute za uporabu	Welch Allyn Connex Spot Monitor Upute za uporabu Ovaj dokument pruža detaljne upute za uporabu Welch Allyn Connex Spot Monitor uređaja, uključujući postavljanje, rad, upravljanje podacima o pacijentu, alarme i rješavanje problema.
 Upute za uporabu	Welch Allyn HSCRIBE: Uživatelská příručka pro analýzu Holterova monitorování Tato uživatelská příručka poskytuje podrobné informace o systému Welch Allyn HSCRIBE pro analýzu Holterova monitorování, včetně instalace, provozu, údržby a řešení problémů.

CP 150 12-lead resting electrocardiograph  Directions for use © 2015 Welch Allyn	Welch Allyn CP 150 12-Lead Resting Electrocardiograph: Directions for Use This manual provides comprehensive instructions for the Welch Allyn CP 150 12-lead resting electrocardiograph, covering setup, operation, maintenance, troubleshooting, and compliance information for healthcare professionals.
 SureTemp Plus Wall Holder	SureTemp Plus Wall Holder: Installation and Mounting Guide User instructions for mounting the SureTemp Plus Wall Holder by Welch Allyn. Includes steps for installation, label application, and security features, with recommended dimensions and clearances.
 Welch Allyn Connex Spot Monitor	Welch Allyn Connex Spot Monitor Desk Stand Startup Guide Comprehensive startup guide for assembling and setting up the Welch Allyn Connex Spot Monitor Desk Stand (Model 7000-DST) by Hillrom, detailing connection of sensors and accessories.
 Welch Allyn HScript	Welch Allyn HScript Holter Analysis System User Manual Comprehensive user manual for the Welch Allyn HScript Holter Analysis System, detailing its capabilities for ECG data acquisition, analysis, reporting, and system configuration for clinical professionals.

Documents - Welch Allyn – 72200

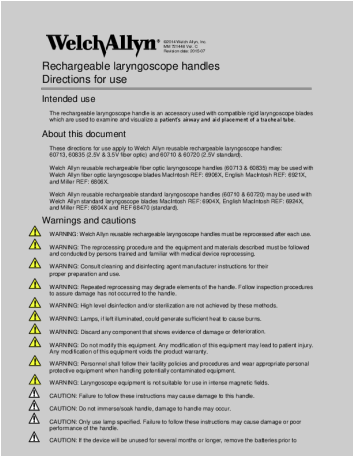

[pdf] User Manual Accessories

untitled Jamie Strong Rechargeable laryngoscope handles Directions for use Hillrom Welch Allyn Laryngoscope Handles User has validated the above instruction as being capable of preparing these re

The user must ensure that WEL60835 20170503040710914 welch allyn rechargeable manual image tigermedical Manuals quickmedical

2014 Welch Allyn, Inc. MM 721448 Ver. C Revision date: 2015-07 Rechargeable laryngoscope handles Dir ... 2.3Wh 1 hr NA 72000, 2.5V 900mAh 2.3Wh 1 hr NA 72000, 2.5V 900mAh 2.3Wh 1 hr HPX060 **72200**, 3.5V 620mAh 2.2Wh hr HPX003 *

Total run time is estimated based on a fully charged batt...

lang:en score:49 filesize: 266.52 K page_count: 8 document date: 2015-07-16



[StarLine Diagnostic Equipment & Supply Catalog](#)

Explore a comprehensive catalog of diagnostic medical equipment and supplies from leading manufacturers including StarLine, American Diagnostic Corporation (ADC), Welch Allyn, SECA, Detecto, and 3M Littmann. Discover a wide range of products such as scales, thermometers, stethoscopes, blood pressure monitors, ECG machines, spirometers, and ophthalmoscopes designed for healthcare professionals.

lang:en score:47 filesize: 2.11 M page_count: 48 document date: 2002-10-30




[\[pdf\] User Manual Instructions](#)

Rechargeable Battery Handles and Chargers General Use Instructions Welch Allyn Inc MaterialNo 720272Rev B ©2013WelchAllyn Printed in U S A WelchAllyn 4341 State Street Road Skaneateles NY 13153 0220 Tel 800 535 6663 welchallyn welch allyn rechargeable battery handles chargers user manual quickmedical

ENGLISH Material No. 728485 DIR No. 80017989 Rev. D 2018 Welch Allyn, Inc. Rechargeable Battery Ha ... s during the two years from the date of manufacture. This warranty does not apply to the use of the **72200** or 72300 battery in other products. Conforms with: Handles IEC/UL/CSA/EN 60601-1 Safety and...

lang:en score:46 filesize: 6.13 M page_count: 4 document date: 2018-09-05

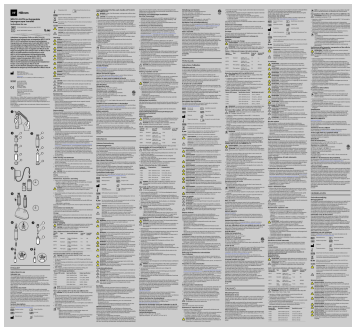

[\[pdf\] Instructions](#)

Rechargeable laryngoscope handles Instructions for use Welch Allyn Inc 80029308LIT hillrom co uk content dam aem us en sap documents LIT 80029 |||

WELCH ALLYN rechargeable laryngoscope handles Instructions for use 775105, 80029308 Ver. B, Revision ... 1 hr NA 2.3Wh 60713 Fiber-optic 72000, 2.5V 900mAh 1 hr 2.3Wh HPX060 60835 Fiber-optic **72200**, 3.5V 620mAh 2.2Wh 3/4 hr HPX003 * Total run time is estimated based on a fully charged ba...


lang:en score:41 filesize: 1.22 M page_count: 4 document date: 2022-08-18



[\[pdf\]](#) Instructions Guide

Welch Allyn Laryngoscopy Repair Help Learn How to Fix It Yourself Welch Rechargeable Fiber Optic Laryngoscope Handle Maintenance Instructions documents cdn ifixit 4TJXVKZvTfWrckrM thumbnailUrl guide images igi 2VRFv5GOvoj3HHSP standard type title |||


Instructions Menu Form No. 608354 Rev. A Form No. 608354 Rev. A Form No. 608354 Rev. A Welch All ... erence Chart Style Model No. Battery Lamp Lamp Cartridge 3.5V Rechargeable 60835 Welch Allyn **72200** 00300 608351-501 The CE mark on this product indicates it has been tested to and conforms ...
lang:en score:39 filesize: 158.24 K page_count: 2 document date: 2020-05-18

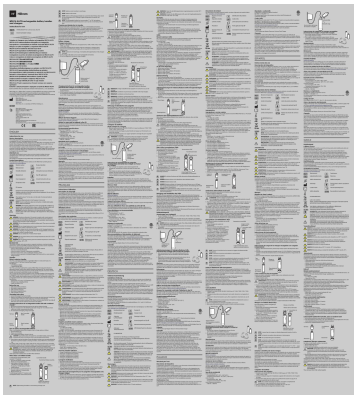

[\[pdf\]](#) Instructions

Instructions for use Welch Allyn® rechargeable battery handles and chargers Allyn Inc 80030374LIT hillrom co uk content dam aem us en sap documents LIT 80030 |||


WELCH ALLYN rechargeable battery handles and chargers Instructions for use 776846, 80030374 Ver. A, ... s during the two years from the date of manufacture. This warranty does not apply to the use of the **72200** or 72300 battery in other products. Welch Allyn reserves the right to make changes without not...
lang:en score:39 filesize: 3.07 M page_count: 4 document date: 2023-06-13

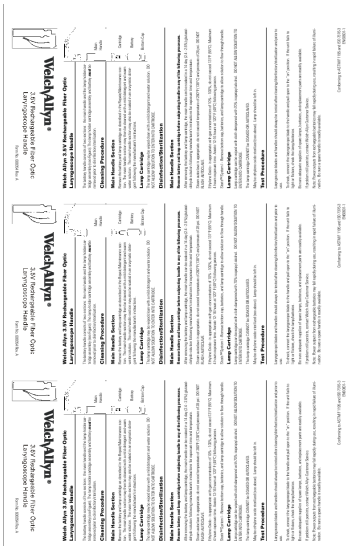

[\[pdf\]](#) Instructions Guide


Welch Allyn Laryngoscopy Repair Help Learn How to Fix It Yourself Welch 3 5V Rechargeable Fiber Optic Laryngoscope Handle 60835 Manual documents cdn ifixit dV5g2ap5Mw5IfXRR thumbnailUrl guide images igi 46HhncjZVlbWUISg standard type title |||

Instructions Menu Form No. 608354 Rev. A Form No. 608354 Rev. A Form No. 608354 Rev. A Welch All ... erence Chart Style Model No. Battery Lamp Lamp Cartridge 3.5V Rechargeable 60835 Welch Allyn **72200** 00300 608351-501 The CE mark on this product indicates it has been tested to and conforms ...
lang:en score:39 filesize: 158.24 K page_count: 2 document date: 2020-05-18

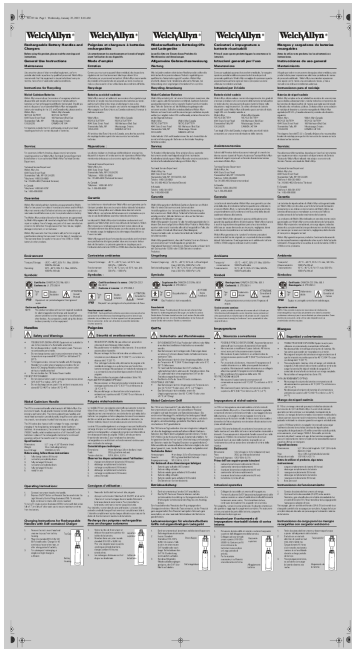


[\[pdf\]](#) Instructions

H PM791327 4 k rev L prn Jewella Tiger Medical Rechargeable Battery Handles and Chargers Before using this product please read the warnings instructions General Use Instructions Maintenance WEL71000 A 20141215055330011 rechargeinstructions 20080116 nickelcadmium handle image tigermedical Manuals welchallyn documents EENT Power Sources Wall Units Nickel Cadmium Handle 791327.fm Page 1 Wednesday, January 22, 2003 8:46 AM Rechargeable Battery Handles and Chargers Befo ... s during the two years from the date of manufacture. This warranty does not apply to the use of the **72200** or 72300 battery in other products. Environment Les clients rsidant en Amrique du Nord doiven... lang:it **score:38** filesize: 343.38 K page_count: 2 document date: 2003-01-22

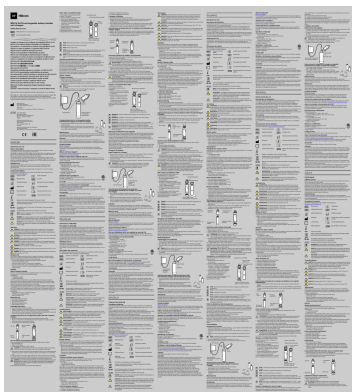

[\[pdf\]](#) Instructions

Instructions for use WELCH ALLYN rechargeable battery handles and chargers Welch Allyn Inc Nov 30 2022 — Poignées et chargeurs avec batteries rechargeables wiederaufladbare Akkugriffe und Ladegeräte Impugnatura e caricabatterie per 80029169LIT hillrom co uk content dam aem us en sap documents LIT 80029 ||| WELCH ALLYN rechargeable battery handles and chargers Instructions for use 774888, 80029169, Ver. B, ... s during the two years from the date of manufacture. This warranty does not apply to the use of the **72200** or 72300 battery in other products. Welch Allyn reserves the right to make changes without not... lang:en **score:38** filesize: 2.7 M page_count: 4 document date: 2022-08-24

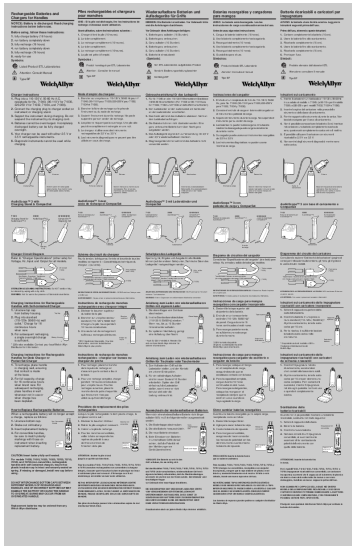


[\[pdf\]](#) User Manual Instructions

Rechargeable Batteries and Chargers for Handles iFixit Piles rechargeables et chargeurs pour manches
AVIS Si la pile est déchargée lire les instructions de recharge ci dessous avant utilisation Avant Welch
Allyn Oscope Charges Instructions documents cdn ifixit BssXOUTEpXC1KfD1 |||

Diagnostic Instruments Rechargeable Batteries and Chargers for Handles NOTICE:
Battery is discharged ... ble Handle 710 Wall Charger UL 70720 2.5v Handle with CE
Wall Charger 71000 3.5v Handle with **72200** 3.5v Rechargeable Self-contained
Charger Battery Black printing 71010 3.5v Rechargeable H...

lang:es **score:38** filesize: 508.08 K page_count: 2 document date: 0000-00-00

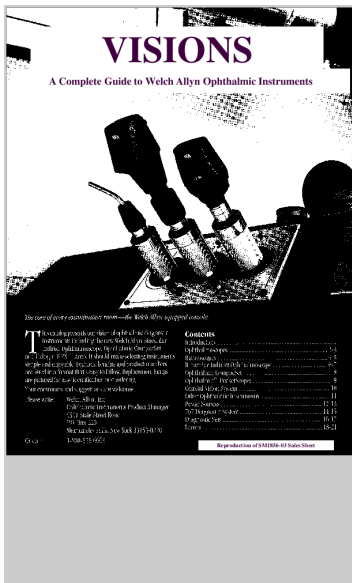

[\[pdf\]](#) Guide Catalog

VISIONS A complete guide SM1856 03 Sales sheet ophthalmoscopes retinoscopes binocular indirect
power source lamp battery compac set pocketscope 767 diagnostic welchallyn 80 medical support
manuals

VISIONS A Complete Guide to Welch Allyn Ophthalmic Instruments Reproduction of
SM1856-03 Sales Sheet ... specify country and voltage No. 71410-optional wall
bracket for use with above desk chargers No. **72200**-NiCad rechargeable battery
color-coded black U.S. Patents Pending 3.5v Convertible Handle Re...

lang:en **score:38** filesize: 2.85 M page_count: 23 document date: 1999-03-24




[\[pdf\]](#) User Manual Instructions

Instructions for use Welch Allyn® 719 Series Lithium Ion Handle Allyn Inc praxisdienst fr WARNUNG
773422 DIR 80027573 Ver A Revision date 2021 01 Serie lithiumion håndtag Lithiumionengriff der Mango
de alimentación ión litio serie 124302 handle with li ion table charger quick user manual out media |||

Welch Allyn 719 Series Lithium Ion Handle Instructions for use 773422, DIR
80027573, Ver. A, Revisio ... acjami przez dwa lata od daty produkcji. Niniejsza
gwarancja nie obejmuje korzystania z akumulatorw **72200** i 72300 w innych
produktach. Firma Welch Allyn zastrzega sobie prawo do wprowadzania zmian w p...


lang:en **score:37** filesize: 1.53 M page_count: 2 document date: 2021-02-01

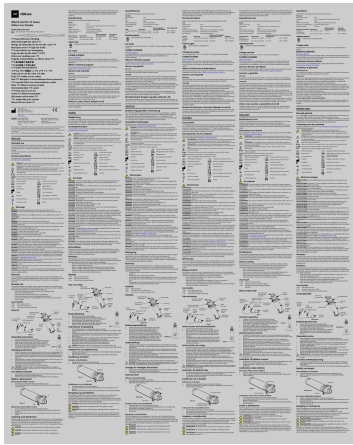

[\[pdf\]](#) Instructions

Instructions for use Welch Allyn 719 Series Lithium Ion Handle Inc 80029317LIT hillrom content dam aem us en sap documents LIT 80029 |||

WELCH ALLYN 719 Series Lithium Ion Handle Instructions for use 775104, 80029317, Ver. B, Revision da ... acjami przez dwa lata od daty produkcji. Niniejsza gwarancja nie obejmuje korzystania z akumulatorów **72200** i 72300 w innych produktach. Firma Welch Allyn zastrzega sobie prawo do wprowadzania zmian w p... lang:en score:36 filesize: 3.98 M page_count: 4 document date: 2022-08-25

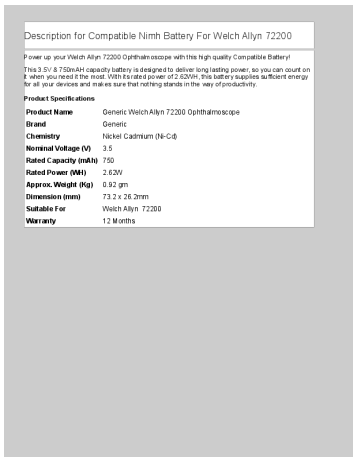


[\[pdf\]](#) Specifications Warranty

Download Compatible Nimh Battery For Welch Allyn 72200 Ophthalmoscope 3 5 750mAH at best price hospitalstore index dispatch attachments getfile attachment id 2766 |||

Description for Compatible Nimh Battery For Welch Allyn **72200** Power up your Welch Allyn **72200** Ophthalmoscope with this high quality Compatible Battery This 3.5V 750mAH capacity battery is designed to deliver long lasting power, so you can count on it when you need it the most. With its rated pow...

lang:de score:36 filesize: 48.99 K page_count: 1 document date: 2023-06-20

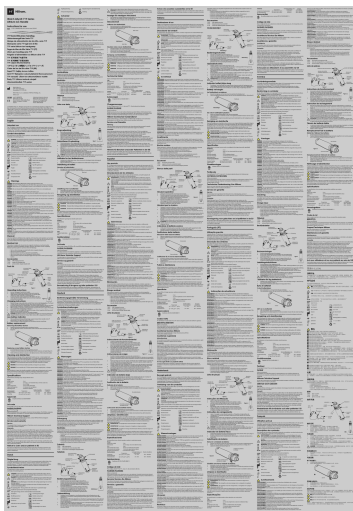

[\[pdf\]](#) User Manual Instructions

Instructions for use Welch Allyn® 719 Series Lithium Ion Handle Allyn Inc 771361 DIR 80026590 Ver B Revision date 2020 08 Serie lithiumion håndtag 80026590LIT emeai welchallyn content dam documents sap LIT 80026 |||

Welch Allyn 719 Series Lithium Ion Handle Instructions for use 771361, DIR 80026590, Ver. B, Revisio ... acjami przez dwa lata od daty produkcji. Niniejsza gwarancja nie obejmuje korzystania z akumulatorów **72200** i 72300 w innych produktach. Firma Welch Allyn zastrzega sobie prawo do wprowadzania zmian w p... lang:en score:33 filesize: 4.29 M page_count: 2 document date: 2020-08-07



[Welch Allyn 719 Series Lithium Ion Handle - Instructions for Use](#)

Comprehensive instructions for the Welch Allyn 719 Series Lithium Ion Handle, covering intended use, operation, charging, battery replacement, cleaning, specifications, safety warnings, and warranty information.

lang:en score:32 filesize: 4.22 M page_count: 2 document date: 2021-02-01



[\[pdf\] User Manual Instructions](#)

Instructions for use Welch Allyn® 719 Series Lithium Ion Handle Allyn Inc Hillrom Serie lithiumion håndtag Lithiumionengriff der Mango de alimentación ión litio serie Manopola agli ioni di 80027573LIT hillrom co uk content dam aem us en sap documents LIT 80027


Welch Allyn 719 Series Lithium Ion Handle Instructions for use 773422, DIR 80027573, Ver. A, Revisio ... acjami przez dwa lata od daty produkcji. Niniejsza gwarancja nie obejmuje korzystania z akumulatorw **72200** i 72300 w innych produktach. Firma Welch Allyn zastrzega sobie prawo do wprowadzania zmian w p...
lang:en score:32 filesize: 1.4 M page_count: 2 document date: 2021-02-01



[\[pdf\] Guide Catalog](#)

Veterinary Product Procedure Guide Vet Catalog Barbara Collins vetguide lib store yahoo net allheart Veterinary Product Procedure Guide For more than three-quarters of a century, Welch Allyn veteri ... 3--Handle sub assembly top section No. 71050-C--Extra bottom section only, 100-130v, 50/60 Hz No. **72200**--Nickel-cadmium rechargeable battery color coded black for Nos. 71000, 71010, 71015, 71020, ...

lang:en score:32 filesize: 1.99 M page_count: 50 document date: 2000-03-21

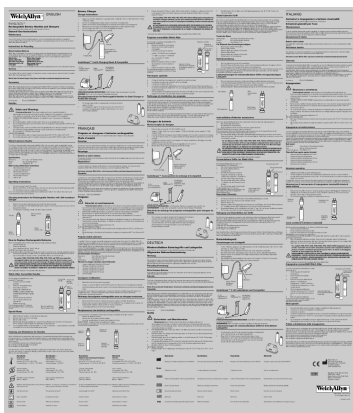


[\[pdf\] Instructions](#)

Rechargeable Battery Handles and Chargers General Use Instructions Wlech Allyn Inc Welch Array 710 Alternation Nickel Cadmium Ability Handle cdn shopify s files 1 1046 1086 |||

ENGLISH Material No. 720272 Rev. B 2013 Welch Allyn, Inc. Rechargeable Battery Handles and Charger ... s during the two years from the date of manufacture. This warranty does not apply to the use of the **72200** or 72300 battery in other products. Conforms with: Handles IEC/UL/CSA/EN 60601-1 Safety and...


lang:en score:31 filesize: 6.26 M page_count: 4 document date: 2013-10-03

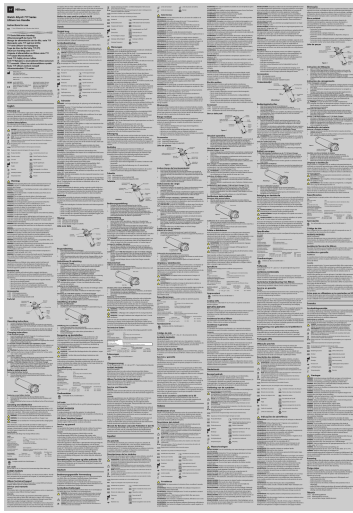

[pdf] User Manual Instructions

Instructions for use Welch Allyn® 719 Series Lithium Ion Handle Allyn Inc Hillrom 28 févr 2023 — Serie lithiumion håndtag Lithiumionengriff der Mango de alimentación ión litio serie Manopola agli ioni di
80025954LIT hillrom co uk content dam aem us en sap documents LIT 80025 |||

Welch Allyn 719 Series Lithium Ion Handle Instructions for use 771362, DIR 80025954, Ver. B, Revisio ... acjami przez dwa lata od daty produkcji. Niniejsza gwarancja nie obejmuje korzystania z akumulatorów **72200** i 72300 w innych produktach. Firma Welch Allyn zastrzega sobie prawo do wprowadzania zmian w p...
lang:en score:31 filesize: 1.53 M page_count: 2 document date: 2020-08-07

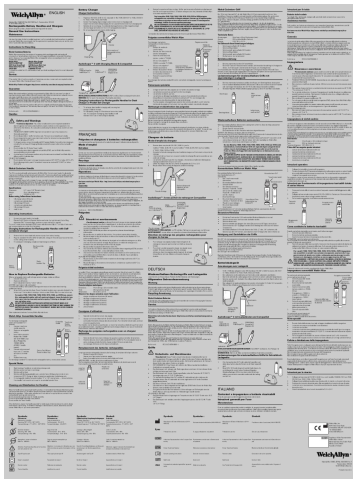

[pdf] Instructions

Directions for use Rechargeable Battery Handles and Chargers ML Welch Allyn Inc 80017989LIT assets hillrom is content

ENGLISH Material No. 728485 DIR No. 80017989 Rev. E, Revision date: 2019-09 2019 Welch Allyn, Inc. ... s during the two years from the date of manufacture. This warranty does not apply to the use of the **72200** or 72300 battery in other products. Conforms with: Handles IEC/UL/CSA/EN 60601-1 Safety and...
lang:en score:30 filesize: 6.1 M page_count: 4 document date: 2019-09-13



AMPOLETAS DE MEDICINA
WELCH ALLYN

Código	Descripción	
WA0380	AMPOLLETA DIAGNOSTICO WELCH ALLYN 03.800 3.5V ORIGINAL	
MWA-8300	AMPOLLETA RETINOSCOPIO WELCH ALLYN 08300 3.5V ORIGINAL	
WA-9500	AMPOLLETA DIAGNOSTICO WELCH ALLYN 08500 3.5V ORIGINAL HEADLIGHT	
03421	AMPOLLETA WELCH ALLYN 08.500 3.5V BH601 HEADLIGHT ZINGG	
WA-1900	AMPOLLETA COLPISCOPIO WELCH ALLYN 01900 ORIGINAL	
72200-ZINGG	BATERIA WELCH ALLYN 72200 3.5V	
WA-9500U	AMPOLLETA PASELON WELCH ALLYN 09500U 16.50W ORIGINAL	

ILUMINA LIMITADA
Carmen Covarrubias 680-A Ñuñoa
Teléfono: 222371651 / 222379124

www.zingg.cl
ilumina@zingg.cl
WhatsApp: +56954187903



[pdf]

Andres Callis WELCH ALLYN 08500 3.6V ALTERNATIVA HEADLIGHT DIAGNOSTICO Ilumina LimitadaBH601cdn shopify s files 1 0687 1262 8424 BH601 v 1742503022 |||

AMPOLLETAS DE MEDICINA WELCH ALLYN C digo WA03800 MWA-8300 WA-8500 BH601 WA-9800 **72200**-ZINGG WA-9500U Descripci n AMPOLLETA DIAGNOSTICO WELCH ALLYN 03.800 3.5V ORIGINAL AMPOLLETA RETINOSCOPIO WELCH ALLYN 08300 3.5V ORIGINAL AMPOLLETA DIAGNOSTICO WELCH ALLYN 08500 3.6V ORIGINAL HEADLIGHT AMPOLLET...
lang:tl score:29 filesize: 1.52 M page_count: 1 document date: 2023-03-02



[\[pdf\]](#) Instructions

Download Fiber Optic Laryngoscope Systems Welch Allyn Hillrom 80025926LIT hillrom co uk content dam aem us en sap documents LIT 80025 |||

Welch Allyn rechargeable laryngoscope handles Instructions for use 771384 DIR 80025926 Ver. B Revisi ... 1 hr NA 2.3Wh 60713 Fiber-optic 72000, 2.5V 900mAh 1 hr 2.3Wh HPX060 60835 Fiber-optic **72200**, 3.5V 620mAh 2.2Wh 3/4 hr HPX003 * Total run time is estimated based on a fully charged ba...

lang:en **score:28** filesize: 992.05 K page_count: 0 document date: 0000-00-00




[\[pdf\]](#) Instructions

DownloadDownloadFiber Optic Laryngoscope Systems Welch Allyn HillromFiber

Hillrom80025926LIT80025926LIT80025926LIThillrom co uk content dam hillrom aem us en sap documents LIT 80025 80025926LIThillrom 80025926LIT |||

Welch Allyn rechargeable laryngoscope handles Instructions for use 771384 DIR 80025926 Ver. B Revi ... 1 hr NA 2.3Wh 60713 Fiber-optic 72000, 2.5V 900mAh 1 hr 2.3Wh HPX060 60835 Fiber-optic **72200**, 3.5V 620mAh 2.2Wh 3/4 hr HPX003 * Total run time is estimated based on a fully charged ba...

lang:en **score:28** filesize: 992.05 K page_count: 0 document date: 2023-02-16

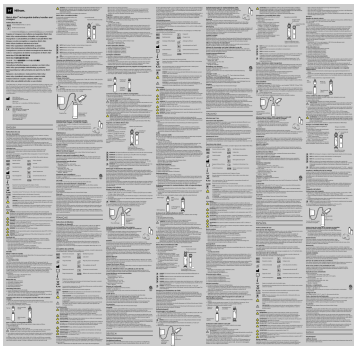


[\[pdf\]](#) Instructions

Download 3 5 V Nickel Cadmium Rechargeable Handles Welch Allyn Hillrom 80027574LIT hillrom co uk content dam aem us en sap documents LIT 80027 |||

Welch Allyn rechargeable battery handles and chargers Instructions for use 773423 DIR 80027574 Rev. ... s during the two years from the date of manufacture. This warranty does not apply to the use of the **72200** or 72300 battery in other products. Welch Allyn reserves the right to make changes without not...

lang:en **score:27** filesize: 2.41 M page_count: 0 document date: 0000-00-00

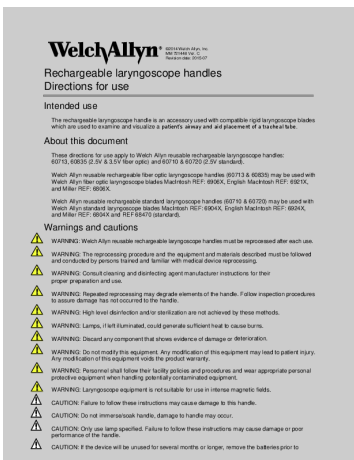


[\[pdf\]](#) Accessories

untitled LIT DFU RECHARGABLE LARYNGOSCOPE HANDLE Jamie Strong Télécharger

Laryngoscopes à fibres optiques Welch Allyn Hillrom 80017985LIT hillrom ca content dam aem us en sap documents 80017 |||

2014 Welch Allyn, Inc. MM 721448 Ver. C Revision date: 2015-07 Rechargeable laryngoscope handles Dir ... 2.3Wh 1 hr NA 72000, 2.5V 900mAh 2.3Wh 1 hr NA 72000, 2.5V 900mAh 2.3Wh 1 hr HPX060 **72200**, 3.5V 620mAh 2.2Wh hr HPX003 * Total run time is estimated based on a fully charged batt...

lang:en **score:27** filesize: 273.98 K page_count: 8 document date: 2015-07-16

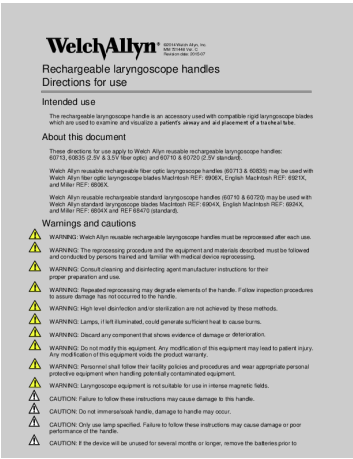


[pdf] Accessories

untitled LIT DFU RECHARGABLE LARYNGOSCOPE HANDLE Jamie Strong Download Fiber Optic Laryngoscope Systems Welch Allyn Hillrom 80017985LIT www.stg.hillrom.com/content/dam/aem/us/en/sap/documents/80017.iii

2014 Welch Allyn, Inc. MM 721448 Ver. C Revision date: 2015-07 Rechargeable laryngoscope handles Dir ... 2.3Wh 1 hr NA 72000, 2.5V 900mAh 2.3Wh 1 hr NA 72000, 2.5V 900mAh 2.3Wh 1 hr HPX060 **72200**, 3.5V 620mAh 2.2Wh hr HPX003 * Total run time is estimated based on a fully charged batt...

lang:en score:27 filesize: 274.21 K page_count: 8 document date: 2015-07-16

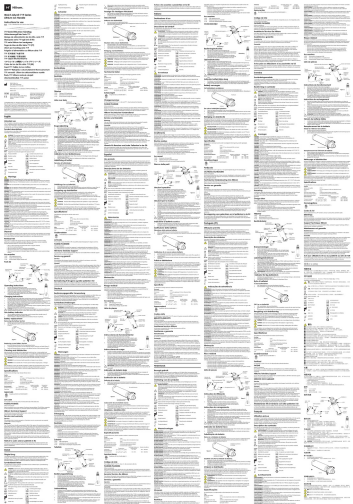

[pdf] Instructions


Instructions for use Welch Allyn® 719 Series Lithium Ion Handle Allyn Inc use719 Welch 3 5V 71900 USB Usegrovemmedical Customer grmein Customerpages specpages

71900 Use.iii Welch Allyn 719 Series Lithium Ion Handle Instructions for use 771361, DIR 80026590, Ver. B, Revis ... jami przez dwa lata od daty produkcji. Niniejsza gwarancja nie obejmuje korzystania z akumulator w **72200** i 72300 w innych produktach. Firma Welch Allyn zastrzega sobie prawo do wprowadzania zmian w p... lang:en score:27 filesize: 4.29 M page_count: 2 document date: 2020-08-07



[pdf] Reference Guide Guide Accessories Catalog

REPLACEMENT PARTS iFixit USE ONLY SPECIFIED CHARGER RECHARGEABLE BATTERY 4 8 VOLTS 1 8AH CATALOG NO 72240 ASM IN USA EXCLUSIVELY FOR WELCH ALLYN NICKEL CADMIUM MUST BE Welch Allyn Vital Signs Replacement Parts documents cdn ifixit EmcixN2I2CdNDIFp.iii REPLACEMENT PARTS Lamps 158 Batteries 164 Blood Pressure Cuff Reference Guide 165 Blood Pr ... lacement Batteries 71960 72000 72100 Current Instruments Lithium Ion Handle 71900 70700, 71500 **72200** 72240 72300 72420 72600 72610 USE ONLY SPECIFIED CHARGER RECHARGEABLE BATTERY 4.8 VOLTS -1.8... lang:en score:27 filesize: 1.2 M page_count: 16 document date: 0000-00-00

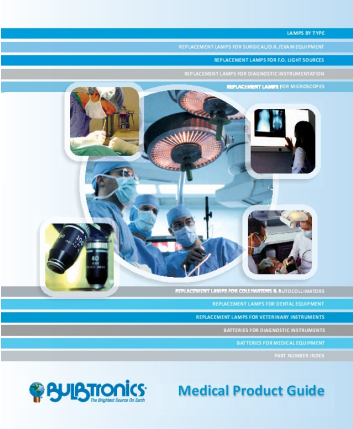

Bulbtronics Medical Product Guide: Replacement Lamps & Batteries


Explore the Bulbtronics Medical Product Guide, a comprehensive catalog of replacement lamps and batteries for medical, dental, and scientific equipment. Find solutions for surgical, diagnostic, microscopy, and more, featuring leading brands like OSRAM, Philips, and GE. lang:en score:27 filesize: 3.99 M page_count: 240 document date: 2010-10-06



[pdf] Accessories Catalog

ama catalog part I cdr basha1 Maximised performance and durability in a compact set that travels easily with you from one patient to the next high diagnostic tools Welch Allyn 1534175452 Diagnostic sets amaproducts au plugin catalogue |||

DIAGNOSTIC SETS Welch Allyn Pocket LED Otosopes, Ophthalmoscopes and Sets
Maximised performance an ... s Rectal Adapter 13820 12810 11110 12820 11710
11720 11730 11735 11630 11620 11740 Macroview 3.5V **72200** 72300 88644 88755
88745 88658 71600 88639 71610 88667 88671 88727 88740 Other Accessories - ...
lang:en score:25 filesize: 1.9 M page_count: 9 document date: 2018-05-14

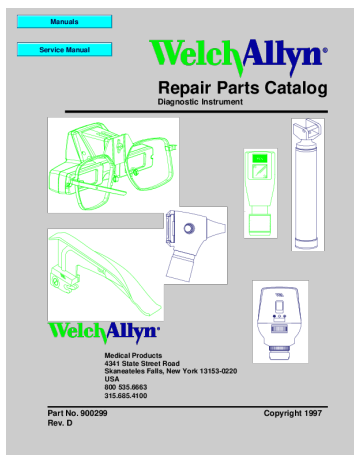

[pdf] User Manual Guide Catalog

Welch Allyn Optometry Repair Help Learn How to Fix It Yourself Welch Diagnostic Instrument Parts
Catalogdocuments cdn ifixit heVFoPhgYAXbRsME thumbnailUrl guide images igi Av1QLhbKjQgJk53L
standard type title CatalogRepair CatalogWelch heVFoPhgYAXbRsMEdocuments ||| ||| Manuals
TechView Menu Service Manual Repair Parts Catalog Diagnostic Instrument Medical Products 4341 State
Street Road Skaneateles Falls, New York 13153-0220 USA 800 535.6663 315.685.4100 Part No. 900299
Rev. D Copyright 1997 Revision Hist
Manuals TechView Menu Service Manual Repair Parts Catalog Diagnostic Instrument
Medical Products 4 ... 4 8 710125-502 9 717126-502 10 710127-503 11 71050 12
717409 13 717410-501 14 717411-1 15 **72200** Description Release Button Bottom
Cap Spring Assy Retainer Handle Tube Assy Top Cap Sleeve To...
lang:en score:25 filesize: 3.07 M page_count: 138 document date: 2020-05-18



[pdf]

Advancing Frontline Care Welch Allyn AJM Produktkatalog 2023 ajmedical se media wysiwyg ajm |||
Hillrom r del av Baxter AdvancFRinOgNTFLrIoNnEtCliAnReE Care
ProduktkaPtraolodguktkatalog Advancin ... ingsbara batteriets design. Vi
rekommenderar att du endast anvnder kta Welch Allyn-batterier 72300 **72200** 71960
719-BATT till nickelkadmiumhandtag med laddare till nickelkadmiumhandtag med
universall...
lang:no score:25 filesize: 3.91 M page_count: 29 document date: 2023-03-17




[pdf]

W57750 qxd emergency medicine Welch Allyn Medical Device Depot s Blog Page 5 Click here for a list of WAEmergencyMedicine meddevicedepot |||

EMERGENCY MEDICINE 6 Welch Allyn. Advancing Frontline Care. You re on the front lines of patient ... r, 1 case 25 rolls ECG Chart Paper for use with PIC 30, PIC 40, PIC 50, 40 mm box of 3 #71960 #72200 7 Replacement Lithium Ion Battery 6 Replacement NiCad Rechargeable Battery black for Des...


lang:en score:24 filesize: 2.25 M page_count: 36 document date: 2010-09-02

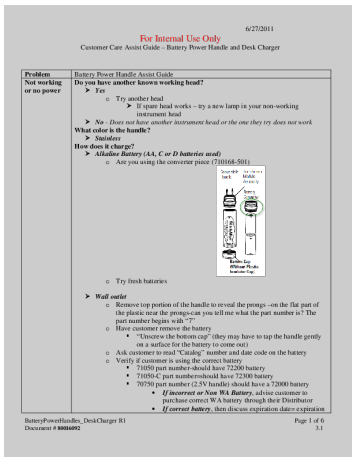

[pdf] Troubleshooting Guide Guide Catalog

Error Code 6 Welch Allyn Family of Companies Battery Power Handles Desk Charger Troubleshooting part number begins with 7 2 April 2001 Switched to Paliney type rheostat be compatible PanOptic a rheostats welchallyn content dam documents upload docs Internal Only Documents Service Article 6/27/2011 For Internal Use Only Customer Care Assist Guide Battery Power Handle and Desk Charger P ... ode on the battery o Verify if customer is using the correct battery 71050 part number-should have 72200 battery 71050-C part number=should have 72300 battery 70750 part number 2.5V handle should...

lang:en score:22 filesize: 112.94 K page_count: 6 document date: 2011-06-27

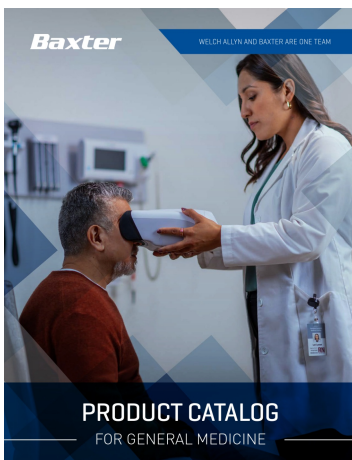

[pdf] User Manual Catalog

PRODUCT CATALOG McKesson Medical Surgical The Welch Allyn Scale Tronix 4802D Pediatric Infant helps enhance daily and routine weighing of newborns infants CONFIGURATIONS XB series Baxter Product General Medicine Catalog 1 mms mckesson content 2025 06 |||

WELCH ALLYN AND BAXTER ARE ONE TEAM PRODUCT CATALOG FOR GENERAL MEDICINE SOLUTIONS FOR PRIMARY CAR ... ger 71670: Welch Allyn 3.5 V rechargeable power handle for universal desk charger, nickel-cadmium #72200 battery 719-WAL: Wall mounting bracket for universal desk charger WELCH ALLYN DIAGNOSTIC SET ...


lang:en score:22 filesize: 4.04 M page_count: 29 document date: 2024-07-15



[pdf] Westone Audiology & Hearing Protection Catalog

Explore Westone Laboratories' comprehensive catalog featuring impression materials, clinical supplies, advanced hearing protection, and consumer audio products for hearing healthcare professionals and individuals.

lang:en score:21 filesize: 9.35 M page_count: 189 document date: 2019-11-04

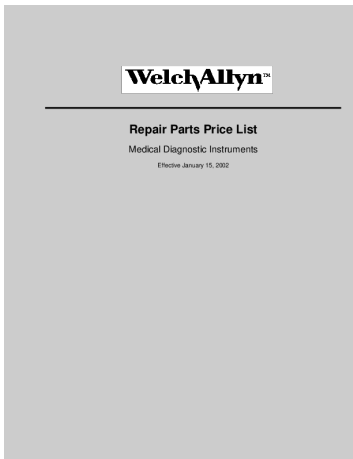


[\[pdf\]](#) Troubleshooting Guide Catalog Price List

Welch Allyn 2002 Repair Parts Price List iFixit Jan 15 · NOTES This 0016 72000 2 5V RECHARGEABLE BATTERY RED 33 35 70720 5V Medical Diagnostic Instruments documents cdn ifixit S2taryEhHQbyblcX

|||

Repair Price Menu TechView Menu Repair Parts Price List Medical Diagnostic Instruments Effective Jan ... 1050 3.5V RECHARGEABLE BOTTOM SECT 104.65 0012 717409 0013 717410-501 0014 717411-1 0015 **72200** BOTTOM CAP INSULATOR BOTTOM CAP ASSY BOTTOM CAP-... lang:en **score:21** filesize: 505.69 K page_count: 52 document date: 0000-00-00

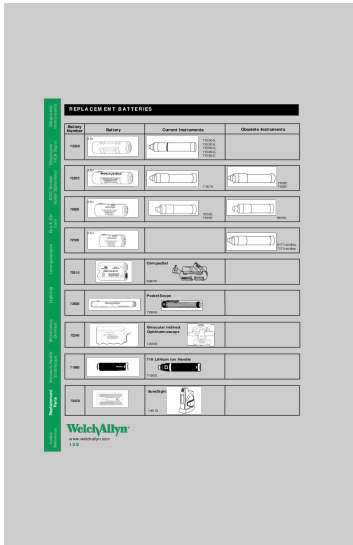


[\[pdf\]](#) Catalog

Master Catalog Medical Instruments Unknown Scheda tecnica articolo Acquista Welch Allyn 72200

BATTERY 3 5V Black Rechargeable Original a prezzi convenienti su Mercateo wa 158 v 394 static mercateo a5 4579243ec47b4292ae7760e4194b12e0 |||

Diagnostic Instruments ECG/Stress/ Diagnostic -- Holter/Spirometry Vital Signs Eye Ear Care Laryngoscopes REPLACEMENT BATTERIES Battery Number 3.5v 72300 Battery **72200** 3.5v 72000 2.5v 72100 2.5v Current Instruments 71000-A, 71020-A, 71000-C, 71020-C, 71055-C 71670 70700 71500 72610 7...

lang:nl **score:21** filesize: 67.55 K page_count: 1 document date: 2006-10-11

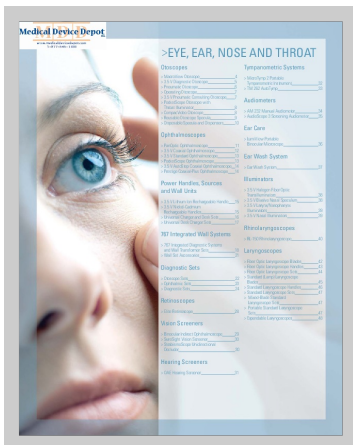

[\[pdf\]](#) User Manual

307109 aa TOC indd Hihosilver EYE EAR NOSE AND THROAT CAN med Healthcare AM 232 Manual Audiometer 34 AudioScope 3 Screening 35 Ear Care LumiView Portable Binocular Microscope 36 Wash System 37 Illuminators 5 V Halogen Fiber Optic Transilluminators 38 Bivalve Nasal Speculum Larynx Nasopharynx 39 WAEENT meddevice depot

EYE, EAR, NOSE AND THROAT Oscopes Tympanometric Systems MacroView Oscope 4 MicroTymp 2 ... y orange for 71000-A/71000-C 710168-501 Battery Converter Ring section only for 71000-A/71000-C **72200** Replacement NiCad Rechargeable Battery black for 71000 and 71670 71140 Universal Desk Charg... lang:en **score:21** filesize: 2.12 M page_count: 46 document date: 2010-09-02



[\[pdf\]](#) Catalog

SM3053 EU r3 MC14993 Product Catalog A4 indd PRODUCT CATALOGUE Welch Allyn Prestige Coaxial Plus Ophthalmoscope with Halogen HPX lamp head only 04900 LED Universal Desk Charger 2 Lithium Ion handles EN emeai welchallyn content dam documents upload docs eu
PRODUCT CATALOGUE EVERYTHING YOU NEED AT THE POINT OF CARE
 MEDICAL L THERMOM BLOOD PR PHYSICAL ... neded that you use only genuine Welch Allyn batteries. 72300 for NiCad Handle w/ Power Transformer **72200** for NiCad Handle w/ Universal Charger 71960 for Lithium Ion Handle Visit www.welchallyn.com ...


lang:en **score:21** filesize: 7.7 M page_count: 44 document date: 2018-12-12



[\[pdf\]](#) User Manual Documentation Catalog

Product Reference Catalogue Braun ThermoScan® PRO 4000 Thermometry 04000 600 Ear Thermometer without rechargeable batteries 53020 0000 Connex VM Software is the backbone that automates vitals documentation process to eliminate manual steps and errors 1163290 henryschein assets Medical
 Product Reference Catalogue Eye, Ear, Nose Throat Get everything you need at the front lines of ... 020-A Rechargeable Handle 3.5 V Nickel-Cadmium Rechargeable Handles: Accessories 72300 710171-501 **72200** 710168-501 Replacement NiCad Rechargeable Battery orange for 71020-A/71020-C Battery Charge...

lang:en **score:21** filesize: 12.38 M page_count: 232 document date: 2012-04-04

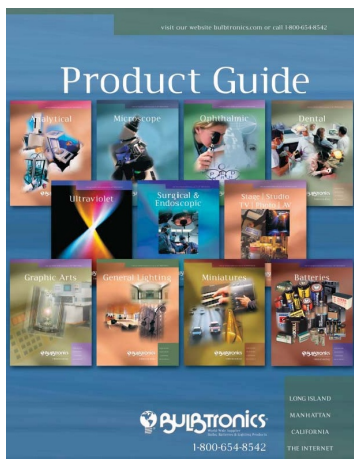

[\[pdf\]](#) Specifications Guide Catalog

000 CVR INS MASTER Catalog Set 1 An Introduction to the Product Guide Bulbtronics C Coiled filament is wound into a helical coil F96T8 SP41 HO 332361 F4T5 CW 10004 20416 237834 F40 41U RS EW ALTO For first time in 27 year history we have gathered information about our entire product line of bulbs batteries and related lighting products one place guide bulbtronics media Documents |||
 An Introduction to the Product Guide For the first time in Bulbtronics 27 year history, we have ga ... n 76600 Manufacturer Part # 11429 12034 X-04.99.623 1919-P5020 CJB360 2.4V 360H11 3-6v 72000 72100 **72200** 72600 Disposable 6 Pack 76600 07600 Bulb Ordering Code CJB114 CJB120 X0199623 1919P5020 CJB36...

lang:en **score:20** filesize: 2.78 M page_count: 216 document date: 2010-12-23

Carley

1503 West 228th Street, Terrace, BC H0T 0J0 Canada Office: (604) 525-8474 info@carleybattery.com

Rechargeable Batteries for Diagnostic Power Handles

Batteries also available without Carley Logo and Reorder number

7000 [®] 3.0W NiCd Rechargeable Battery [mm]	7000 [®] 3.0W NiCd Rechargeable Battery [mm]	7000 [®] 3.0W NiCd Rechargeable Battery [mm]	7000 [®] 3.0W NiCd Rechargeable Battery [mm]
CARLEY # C4-0400-00	CARLEY # C4-0400-00	CARLEY # C4-0400-00	CARLEY # C4-0400-00
60700 7000	7000A 7100D	7000A 7100D-A	7100D 7200T
60700 7000	7000B 7100D	7000B 7100D-B	7200D 7200T
60700 7000	7000C 7100D	7000C 7100D-C	
60700 7000	7000E 7200T	7000E 7200T	
70700 7000	7200T 7200T		

Quality Manufactured and Color Coded Low Cost - Phase Savings on all Your Customers
Packaged 1 per Box For Optimal Performance Charge for 18 Hours Overnight

Let Carley Lamps be your one-stop source for all of your medical replacement lamp and battery needs.
Log onto www.CarleyBattery.com and create an online account today or call (604) 525-8474
and let one of our trained experts get you the right lamp or battery for your equipment.

© 2010 Carley Battery Ltd. All rights reserved. Carley, Carley logo, and Carley Lamp are trademarks of Carley Battery Ltd. in Canada and other countries. Other names may be trademarks of their respective owners. The information herein is provided as a guide only. Please refer to the actual product specifications for detailed information. Carley Battery Ltd. does not warrant the accuracy or completeness of the information herein. Carley Battery Ltd. is not responsible for errors or omissions in this document. Carley Battery Ltd. reserves the right to change the information herein at any time without notice.

[\[pdf\]](#) Datasheet