

Geemarc Amplipower60

Geemarc Amplipower60 Ultra Amplified Corded Telephone User Manual

Model: Amplipower60

1. INTRODUCTION

The Geemarc Amplipower60 is an ultra-amplified corded telephone designed to provide clear and loud communication, especially beneficial for individuals with hearing impairments. This manual provides detailed instructions for setting up, operating, and maintaining your telephone to ensure optimal performance.



Figure 1: Geemarc AmpliPOWER60 Ultra Amplified Corded Telephone

This image shows the main view of the Geemarc AmpliPOWER60 telephone, highlighting its large keypad and ergonomic design.

2. PRODUCT FEATURES

The AmpliPOWER60 telephone is equipped with several features to enhance usability and accessibility:

- **Amplified Landline Phone:** Specifically designed to help hard of hearing individuals make and hear phone calls easier, and prevent missing incoming calls.
- **Loud Ringer & Strobe Notification:** Features a loud ringer and an extra bright flashing incoming call strobe for visual alerts.
- **Extra-Large Keypad:** Equipped with easy-to-read numbers for convenient dialing.
- **Quick Dialing:** Includes 3 color-coded buttons for quick access to emergency services or family members, and 9 quick dial keys for frequently used contacts.
- **Extra Loud & Hearing Aid Compatible:** Provides up to 65dB of extra clear TrueSound, up to 15dB for the speakerphone, and up to 95dB for the ringer. It is also hearing aid compatible.

AMPLIFIED LANDLINE PHONE



Helps hard of hearing individuals make & hear phone calls easier



Prevents missing incoming calls at home



Easy to use



Figure 2: Key Benefits of the Amplified Landline Phone

This image illustrates the primary advantages of the Amplipower60, including ease of use, assistance for the hard of hearing, and preventing missed calls.



Figure 3: Sound Amplification and Compatibility

This image details the sound amplification capabilities of the phone, including up to 65dB TrueSound, 15dB speakerphone, and 95dB ringer, along with its hearing aid compatibility.

EXTRA LOUD & HEARING AID COMPATIBLE

- Up to 65dB of extra clear TrueSound
- Up to 15dB of sound for the speakerphone
- Up to 95dB of sound for the ringer



Figure 4: Dialing Convenience

This image highlights the large numbers and quick dial buttons (M1-M9) that make dialing easy and convenient for users.

3. SETUP GUIDE

Follow these steps to set up your Geemarc Amplipower60 telephone:

1. **Unpack Contents:** Carefully remove all components from the packaging. Ensure you have the telephone base, handset, coiled handset cord, line cord, and power adapter.
2. **Install Batteries:** The phone requires 4 AAA batteries (not included) for certain functions like memory retention during power outages. Locate the battery compartment on the underside of the phone, insert the batteries according to the polarity markings, and close the compartment.
3. **Connect Handset:** Plug one end of the coiled handset cord into the handset and the other end into the jack labeled 'HANDSET' on the side or back of the telephone base.
4. **Connect Line Cord:** Plug one end of the line cord into the jack labeled 'LINE' on the telephone base and the other end into your telephone wall jack.
5. **Connect Power Adapter:** Plug the small end of the power adapter into the 'DC IN' jack on the telephone base and the other end into a standard electrical outlet.
6. **Initial Test:** Lift the handset. You should hear a dial tone. If not, recheck all connections.

4. OPERATING INSTRUCTIONS

4.1 Making a Call

1. Lift the handset or press the **SPKR** button to activate the speakerphone.
2. Listen for a dial tone.
3. Dial the desired telephone number using the large keypad.
4. To end the call, replace the handset or press the **SPKR** button again if on speakerphone.

4.2 Receiving a Call

When the phone rings and the strobe light flashes:

- Lift the handset to answer the call.
- Alternatively, press the **SPKR** button to answer on speakerphone.

4.3 Using Amplification

The Amplipower60 is designed to provide significant sound amplification. The amplification is typically active by default or can be adjusted:

- **Handset Volume:** Use the **VOLUME** slider on the side of the phone to adjust the incoming call volume through the handset. This can provide up to 65dB of clear TrueSound.
- **Speakerphone Volume:** When using the speakerphone, the volume can be adjusted to provide up to 15dB of sound.
- **Ringer Volume:** The ringer can be set to provide up to 95dB of sound. Adjust the ringer volume switch on the side or back of the phone.
- **Tone Control:** Use the **TONE** button to adjust the frequency response for clearer sound, especially for certain types of hearing loss.

4.4 Speed Dial / Memory Buttons

The phone features 3 color-coded emergency buttons (M, CS) and 9 quick dial memory buttons (M1-M9).

1. **Programming Memory Buttons:** Refer to the detailed programming instructions in the full manual for setting up emergency and quick dial numbers. Typically, this involves pressing a **PROG** button, then the memory button, dialing the number, and pressing **PROG** again.
2. **Using Memory Buttons:** Lift the handset or press **SPKR**, then press the desired memory button (M, CS, M1-M9) to dial the stored number.



Figure 5: User Experience with Amplified Sound

This image depicts a woman comfortably using the Amplipower60, illustrating the clarity and ease of communication provided by its amplified features.

5. MAINTENANCE

To ensure the longevity and optimal performance of your Geemarc Amplipower60 telephone, follow these maintenance guidelines:

- **Cleaning:** Use a soft, damp cloth to clean the telephone's surface. Do not use abrasive cleaners, waxes, or solvents.
- **Avoid Liquids:** Keep the telephone away from water and other liquids. If liquid spills on the phone, immediately unplug it and allow it to dry completely before reconnecting.
- **Temperature and Humidity:** Store and operate the phone in a dry environment, away from extreme temperatures and high humidity.
- **Battery Replacement:** Replace the AAA batteries periodically (e.g., once a year) to ensure memory retention functions properly during power outages.
- **Cord Care:** Avoid tangling or stretching the handset and line cords excessively, as this can damage the internal wiring.

6. TROUBLESHOOTING

If you experience issues with your Amplipower60 telephone, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
No dial tone	Line cord or power adapter not connected properly.	Ensure all cords are securely plugged into the phone and wall jacks. Check power outlet.
Cannot hear caller / Low volume	Volume setting is too low; amplification not active.	Adjust the VOLUME slider. Ensure the tone control is set appropriately.
Caller cannot hear me / Muffled sound	Handset or speakerphone microphone issue; speaking too far from microphone.	Ensure handset cord is fully inserted. Speak directly into the microphone. Test with another phone if possible.
Static or Hissing sound	Line interference; high amplification setting.	Check phone line for damage. Try adjusting the volume or tone control. If persistent, contact your telephone service provider.
Ringer not working or too quiet	Ringer volume set too low or off.	Adjust the ringer volume switch on the phone.

If the problem persists after trying these solutions, please contact Geemarc customer support or your retailer for further assistance.

7. SPECIFICATIONS

Feature	Detail
Manufacturer	Geemarc
Brand	Geemarc
Model Number	Amplipower60
Item Weight	2.5 pounds
Product Dimensions	8.75 x 7.5 x 1.5 inches
Batteries Required	4 AAA batteries (for memory retention)
Color	Black
Material Type	Other materials
Telephone Type	Corded
Power Source	Corded Electric




Feature	Detail
Dialer Type	Single Keypad
Answering System Type	Digital
Amplification (Handset)	Up to 65dB TrueSound
Amplification (Speakerphone)	Up to 15dB
Ringer Volume	Up to 95dB

8. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the documentation included with your purchase or contact the manufacturer, Geemarc, directly. You may also reach out to the seller, Sonic Alert MI, for assistance.

Keep your purchase receipt as proof of purchase for any warranty claims.

Related Documents - Amplipower60

	<p>Geemarc CL595 Amplified Corded Phone Quick Start Guide</p> <p>A quick start guide for the Geemarc CL595 amplified corded telephone, detailing its features, setup, operation, SOS function, answering machine, and technical specifications. This hearing aid compatible phone is designed for ease of use with large buttons and clear display.</p>
	<p>Geemarc CL595 User Guide</p> <p>Comprehensive user guide for the Geemarc CL595 telephone, covering setup, features, and troubleshooting.</p>
	<p>Geemarc AMPLIDECT COMBI-PHOTO 295 Quick User Guide</p> <p>Get started with your Geemarc AMPLIDECT COMBI-PHOTO 295 amplified cordless telephone. This quick guide covers setup, features like big buttons and hearing aid compatibility, and accessing the answering machine. For full details, visit www.geemarc.com.</p>




[Geemarc CL595 Amplified Corded Phone User Guide](#)

User guide for the Geemarc CL595 amplified corded phone, detailing setup, features, and troubleshooting. Learn how to use its amplified sound, large buttons, and various functions for enhanced communication.



[Geemarc Photophone 110: Forstærket Foto-ID-Telefon Brugervejledning](#)

Denne brugervejledning giver detaljerede instruktioner til Geemarc Photophone 110, en forstærket foto-ID-telefon designet til ældre og brugere med høreapparater. Lær om opsætning, funktioner, fejlfinding og garanti.

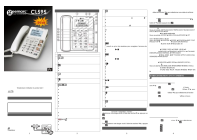

[Geemarc CL595 Corded Phone Quick Start Guide](#)

A quick start guide for the Geemarc CL595 corded phone, covering installation, base overview, menu settings, SOS functions, answering machine operation, and technical specifications.