

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Encore](#) /
- › [Encore Deluxe Battery and Manual Vacuum Erection Device User Manual](#)

Encore 44020-001

Encore Deluxe Battery and Manual Vacuum Erection Device User Manual

Model: 44020-001 | Brand: Encore, Inc.

1. INTRODUCTION


The Encore Deluxe Battery and Manual Vacuum Erection Device (VED) is designed to assist individuals in achieving and maintaining an erection. This system offers a non-invasive approach, utilizing vacuum therapy to create penile rigidity. It includes both battery-powered and manual pump options for user flexibility, along with various accessories to ensure proper fit and comfort. This manual provides detailed instructions for the safe and effective use, setup, maintenance, and troubleshooting of your device.

The American Urology Association recommends vacuum therapy as a treatment option for erectile dysfunction.

2. COMPONENTS OVERVIEW

Your Encore Deluxe VED kit includes the following components:

- Vacuum Cylinder
- Battery-Powered Pump
- Manual Pump
- Seven (7) Constriction Rings of varying sizes
- Ring Loading Applicator
- Lubricant
- Discreet Carrying Case with foam cutouts
- Instruction Booklet and Instructional DVD

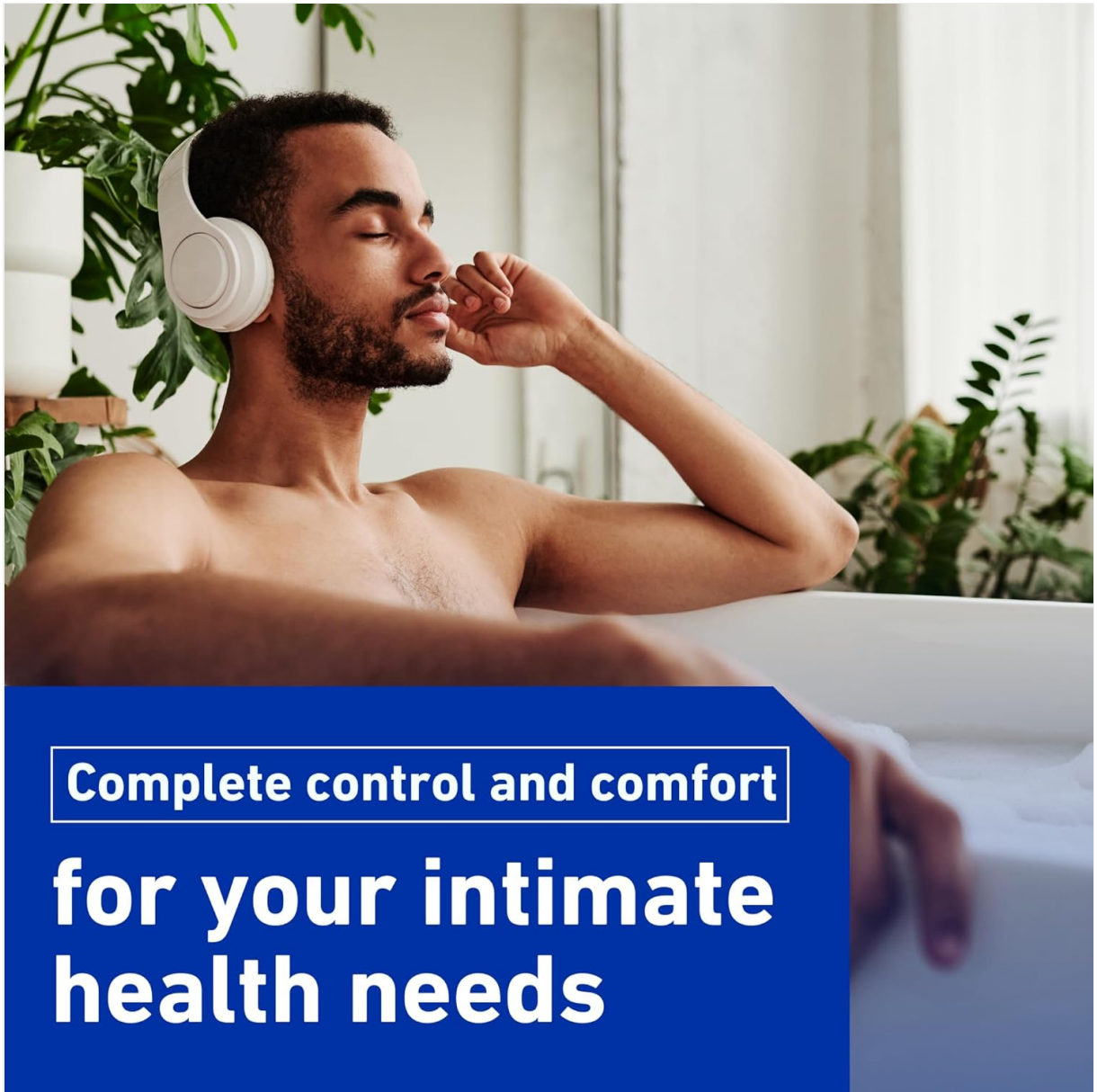

Complete control and comfort

**for your intimate
health needs**

Figure 2.1: All components of the Encore Deluxe Vacuum Erection Device kit.

Seven incrementally different sized rings for personalized comfort and support



Figure 2.2: A selection of the seven different sized constriction rings included for personalized comfort and support.

Features a patented ring ejector for effortless ring removal after use



Figure 2.3: The ring loading applicator and ejector for effortless ring placement and removal.

Comes with a tube of lubricant and zippered carrying case for convenient storage and travel



Figure 2.4: Included lubricant and the discreet carrying case with foam cutouts for organized storage and travel.

3. SETUP

1. **Unpack Components:** Carefully remove all items from the carrying case and inspect them for any damage.
2. **Battery Installation (for Battery Pump):** If using the battery-powered pump, ensure fresh batteries are installed according to the markings inside the battery compartment.
3. **Assemble Pump to Cylinder:** Attach your chosen pump (manual or battery) securely to the open end of the vacuum cylinder. Ensure a tight seal.
4. **Select Constriction Ring:** Choose a constriction ring that provides a comfortable yet firm fit. Use the ring loading applicator to place the selected ring onto the base of the cylinder.
5. **Apply Lubricant:** Apply a generous amount of the provided lubricant to the opening of the cylinder and to the base of the penis to ensure a comfortable seal and ease of use.

4. OPERATING INSTRUCTIONS

Before first use, it is recommended to watch the instructional DVD for a visual guide.

Your browser does not support the video tag.

Video 4.1: Official instructional video demonstrating the use and features of the Encore Deluxe Battery and Manual Vacuum Erection Device.

4.1 General Usage

1. Place the lubricated cylinder over the penis, ensuring a tight seal against the body.
2. Begin pumping to create a vacuum.
3. Once an erection is achieved, carefully slide the constriction ring from the base of the cylinder onto the base of the penis.
4. Remove the cylinder. The constriction ring should remain at the base of the penis to maintain the erection.
5. The constriction ring should not be worn for more than 30 minutes.

4.2 Using the Manual Pump

Attach the manual pump to the cylinder. Squeeze the pump bulb repeatedly to create a vacuum. A release valve is typically present to relieve pressure if needed.

4.3 Using the Battery-Powered Pump

Attach the battery-powered pump to the cylinder. Activate the pump according to its specific controls (e.g., power button). The pump will automatically create a vacuum. Ensure batteries are fresh for optimal performance.

4.4 Removing the Constriction Ring

After use, or if discomfort occurs, use the ring ejector sleeve to safely and effortlessly remove the constriction ring from the base of the penis. Do not pull or force the ring off without using the ejector or appropriate technique.

5. MAINTENANCE AND CARE

1. **Cleaning:** After each use, disassemble the cylinder, pumps, and rings. Wash all components with mild soap and warm water. Rinse thoroughly and allow to air dry completely before storage. Do not use harsh chemicals or abrasive cleaners.
2. **Lubricant:** Replenish lubricant as needed. Only use water-based lubricants compatible with medical devices.
3. **Storage:** Store all components in the provided discreet carrying case in a clean, dry place, away from direct sunlight and extreme temperatures.

6. TROUBLESHOOTING

- **No Vacuum or Weak Vacuum:**

- Check all connections between the pump and cylinder for a secure fit.
- Ensure sufficient lubricant is applied to create a tight seal against the body.
- For battery pump: Check battery levels and replace if necessary.
- Inspect the cylinder and pump for any cracks or damage that could compromise the seal.

- **Difficulty Achieving Erection:**

- Ensure proper technique as described in the operating instructions and instructional video.
- Verify that the chosen constriction ring size is appropriate.

- Consult a healthcare professional if issues persist.
- **Discomfort During Use:**
 - Release the vacuum immediately if pain or severe discomfort occurs.
 - Ensure adequate lubrication.
 - Re-evaluate constriction ring size; a different size may be more comfortable.

7. SPECIFICATIONS

Model Number	44020-001
Manufacturer	Encore, Inc.
Package Dimensions	12.83 x 9.92 x 2.91 inches
Weight	1 Pound
First Available Date	May 4, 2004

8. IMPORTANT SAFETY INFORMATION

Consult a Healthcare Professional: Always consult with a healthcare professional before using any medical device, especially if you have pre-existing medical conditions or are taking medications.

Usage Duration: Do not wear the constriction ring for more than 30 minutes at a time. Prolonged use can lead to injury.

Proper Fit: Ensure the constriction ring is the correct size. An improperly sized ring can cause discomfort or injury.

Lubrication: Always use adequate lubrication to prevent skin irritation and ensure a proper seal.

Discomfort: If you experience pain, numbness, coldness, or any unusual sensation, immediately release the vacuum and remove the constriction ring. If symptoms persist, seek medical attention.

Hygiene: Clean all components thoroughly after each use to prevent infection.

Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This device is not a dietary supplement.

9. WARRANTY AND SUPPORT

For warranty information, technical support, or to purchase replacement parts, please refer to the contact information provided in the original product packaging or visit the manufacturer's official website. Keep your purchase receipt for warranty claims.