

[Manuals.plus](#) /

› [Drive Medical](#) /

› Drive Medical Aluminum Rollator Walker Fold Up and Removable Back Support, Padded Seat, 6" Wheels, Red

Drive Medical R726RD

Drive Medical Aluminum Rollator Walker User Manual

Model: R726RD

PRODUCT OVERVIEW

The Drive Medical Aluminum Rollator Walker is designed to provide stable and comfortable mobility support. This 4-wheeled adult walker features a durable aluminum frame, a padded seat and backrest, and 6-inch non-marring wheels suitable for both indoor and outdoor use. It includes a convenient under-seat storage pouch and adjustable ergonomic handles for user comfort.



Image: The Drive Medical Aluminum Rollator Walker, showcasing its red frame, black padded seat, and four wheels. This image provides a general view of the product.

SETUP AND ASSEMBLY

The rollator comes mostly pre-assembled. Follow these steps to complete the setup and prepare your rollator for use:

1. **Unpacking:** Carefully remove all contents from the carton.
2. **Unfolding:** Ensure the rollator is fully opened. The seat should be in the down position and the folding bar under the seat should be locked.
3. **Backrest Installation:** Open the frame and slide the padded backrest into the frame brackets. Ensure the snap buttons protrude through the holes to secure it. The backrest is hinged and can be folded up or down for storage.
4. **Handle Height Adjustment:** Insert the hand grip posts into the designated holes. Adjust the handle height to your preference by removing the tightening knobs, setting the handles to the desired height, and then re-inserting and tightening the knobs securely. Ensure both handles are set to the same height for balance and stability. The handle height is adjustable between 32 and 37 inches.

assembly instructions

Do not attempt to adjust or operate Rollator without carefully reading all instructions in this manual.

1. First remove all contents from carton.
2. Open frame and slide padded backrest into frame brackets making sure snap bottoms protrude through holes.
3. Remove tightening knobs from packaging and insert handle tube into frame at desired height lining up holes appropriately. Remove tightening knobs from packaging.



WEIGHT CAPACITY 300 LBS (136 KG)

safety precautions

1. Please make sure the Rollator is fully opened, seat is down and the folding bar under seat is locked in the down position.
2. Be sure that both knobs are fully tightened and are at equal heights. Test Rollator for stability before applying full weight to hand grips.
3. Test the loop locks before operating the Rollator. When wheels are locked (see illustrations for locking and releasing locks below), lock bar on rear wheels should prevent wheels from rolling. If the wheels roll while brakes are locked, please take unit back to dealer for adjustment. When locks are released, the wheels should roll freely.
4. Always lock loop locks before sitting and do not use seat when unit is on an incline or uneven ground.
5. Replace tires when significantly worn or damaged.

Important Safety Notice

Do not have anyone push you while you are seated on the Rollator. This is a walking aid only and is not to be used as a transportation device.

loop lock operation

Always engage and disengage wheel locks on both sides simultaneously.

1. To Slow Down or Come to a Temporary Stop – Pull up slowly on the loop locks with four fingers until unit slows down or comes to a temporary stop.
2. To Temporarily Lock Wheels in Position – Press down on both loop locks with palms until you hear a “click” that will tell you wheels are locked.
3. To Release Wheel Locks – Pull up on both loop locks with palms of hand to original position.



Image: An illustrative diagram detailing the assembly process, including the placement of hand grips, backrest, and the folding hinge for portability. Key components like the orthopedic hand grips, loop locks, backrest, wide comfortable seat, and hanging bag for storage are labeled.



Image: A visual representation of the rollator's adjustable handle height, showing how the handles can be raised or lowered to achieve a comfortable and ergonomic reach for the user.

OPERATING INSTRUCTIONS

Proper operation of your rollator ensures safety and comfort. Familiarize yourself with the braking system and seating function.

Using the Brakes (Loop Lock Operation)

1. **To Slow Down or Come to a Temporary Stop:** Pull up slowly on the loop locks with four fingers until the unit slows down or comes to a temporary stop.
2. **To Temporarily Lock Wheels in Position (Parking Brake):** Press down on both loop locks with your palms until you hear a "click". This indicates that the wheels are locked. Always engage the parking brake when sitting on the rollator.
3. **To Release Wheels:** Pull up on both loop locks with your palms to return them to their original position.

Important Safety Note: Always ensure both brakes are fully tightened and apply equal pressure. Test the rollator for stability before applying full weight to the hand grips. When wheels are locked, they should prevent the rollator from rolling. If the wheels roll while brakes are locked, professional adjustment may be required.

Using the Seat

The rollator features a padded seat for resting. Before sitting, always ensure the rollator is on a flat, even surface and the parking brakes are fully engaged to prevent movement. Do not push the rollator while you are seated on it, as it is a walking aid and not intended as a transportation device.

Folding and Storage

The rollator features a side-to-side folding design for easy storage and transport. To fold, lift the seat and then pull up on the strap located underneath. The rollator will collapse, maintaining a standing position for convenient storage.



MAINTENANCE

Regular maintenance ensures the longevity and safe operation of your Drive Medical Rollator Walker.

- **Periodic Inspection:** Drive Medical recommends periodic visual inspection of the product to ensure all parts and hardware are secure, in good working order, and free from wear, tear, fraying, or looseness. Check for any obstructions or impediments to normal, safe operation.
- **Brake System Adjustment:** The hand brake system is factory preset. However, it may require adjustment over time due to normal wear. To adjust, hold the brake firmly in place. To tighten, turn the nut clockwise. To loosen, turn the nut counter-clockwise.
- **Wheel Replacement:** Replace tires when they are significantly worn or damaged. Non-marring wheels are designed to protect floors.
- **Cleaning:** Clean the frame and seat with a damp cloth and mild detergent. Avoid abrasive cleaners.

instructions to adjust the hand brake system


Important Safety Notice

The hand brake system fitted on this rollator complies with the manufacturer's quality standards. The hand brake system has been preset (adjusted) to the best performance settings. To better suit you personally, please follow the instructions set out below. You may need to adjust the hand brake systems from time to time as they will wear out with normal use.

Hold firmly in place do not turn.

To tighten the hand brake adjustment. "turn the nut clockwise"

To loosen the hand brake adjustment. "turn the nut counter clockwise"

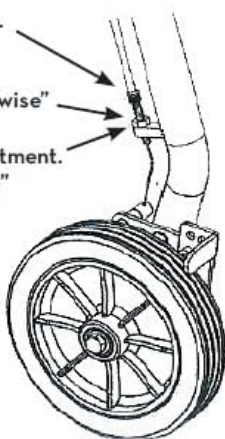

maintenance

Drive recommends periodic visual inspection of this product by the user, to ensure that all parts and hardware are secure, that components are in good working order and not worn, torn, frayed, or loose, and that there are no obstructions or impediments to normal, safe operation. If any part or mechanism appears to be unsafe or compromised, do not use the product. If you experience problems, please contact the provider that sold you the product or visit www.drivemedical.com/help for more info.

lifetime limited warranty

Your Drive brand product is warranted to be free of defects in materials and workmanship for the lifetime of the original consumer purchaser.

This device was built to exacting standards and carefully inspected prior to shipment. This Lifetime Limited Warranty is an expression of our confidence in the materials and workmanship of our products and our assurance to the consumer of years of dependable service.

This warranty does not cover device failure due to owner misuse or negligence, or normal wear and tear. The warranty does not extend to non-durable components, such as rubber accessories, casters, and grips, which are subject to normal wear and need periodic replacement.

If you have a question about your Drive device or this warranty, please contact an authorized Drive dealer.

⚠ WARNING: This product can expose you to chemicals, including bisphenol A (BPA), which is known to the State of California to cause birth defects or other reproduction harm. For more information go to www.P65Warnings.ca.gov.



© 2024 Medical Depot, Inc. All rights reserved.

Drive is a trademark of Medical Depot, Inc.

Port Washington N.Y. 11050 USA Made in China

Image: A detailed diagram showing the hand brake mechanism and instructions for tightening or loosening the brake adjustment nut to ensure proper brake function.

TROUBLESHOOTING

This section addresses common issues you might encounter with your rollator.

- **Brakes Not Engaging/Holding:**
 - Ensure the brake cables are not snagged or damaged.
 - Adjust the hand brake system as described in the Maintenance section.
 - If adjustment does not resolve the issue, or if cables are damaged, contact a qualified technician or the provider for repair or replacement.
- **Rollator Feels Unstable/Wobbly:**
 - Verify that all tightening knobs and screws are securely fastened, especially after assembly or adjustment.
 - Ensure the handle heights are set equally on both sides.
 - Check wheels for any obstructions (e.g., hair, debris) that might impede smooth rolling.
- **Difficulty Folding/Unfolding:**
 - Ensure the seat is lifted before attempting to pull the folding strap.
 - Check for any obstructions or debris in the folding mechanism.

SPECIFICATIONS

Feature	Detail
Model Number	R726RD
Material	Aluminum
Color	Red
Item Weight	13 pounds (5.9 kg)
Product Dimensions (L x W x H)	31 x 24.2 x 7.6 inches (folded)
Adjustable Handle Height	32 - 37 inches
Weight Capacity	300 pounds (136 kg)
Wheel Size	6 inches (non-marring)
Seat Features	Padded seat and removable hinged backrest
Storage	Under-seat pouch



Image: A diagram illustrating the key dimensions of the rollator, including height range, seat dimensions, and overall footprint, along with its 300 lbs weight capacity.

WARRANTY INFORMATION

Your Drive Medical Rollator Walker is backed by a **Limited Lifetime Warranty** from Drive Medical. This warranty covers defects in materials and workmanship for the lifetime of the original consumer purchaser.

This device is built to exacting standards and carefully inspected prior to shipment. The Limited Lifetime Warranty is an expression of confidence in the materials and workmanship of our products and our assurance to the consumer of years of dependable service.

Exclusions: This warranty does not cover device failure due to owner misuse or negligence, or normal wear and tear. The warranty does not extend to non-durable components, such as rubber accessories, casters, and grips, which are subject to normal wear and need periodic replacement.

SUPPORT AND CONTACT

For further assistance, questions, or to address any issues not covered in this manual, please utilize the following resources:

- **Online Support:** Visit the official Drive Medical website for additional information and support resources:
www.drivemedical.com/help
- **Product Documentation:** A PDF version of the user manual may be available for download: [User Manual \(PDF\)](#)
- **Authorized Dealer:** If you have questions regarding your Drive Medical device or its warranty, please contact an authorized Drive Medical dealer.

Proposition 65 Warning: This product can expose you to chemicals, including bisphenol A (BPA), which is known to the State of California to cause birth defects or other reproductive harm. For more information go to www.P65Warnings.ca.gov.

© 2024 Drive Medical. All rights reserved.

Drive is a trademark of Medical Depot, Inc. Port Washington N.Y. 11050 USA. Made in China.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question

Open the live page

<https://manuals.plus/asin/B002VWK3F2>

More documents and user manuals: <https://manuals.plus/>