



Manuals.plus /

- › REDMON Since 1883 /
- › REDMON Precision Digital Pet Scale, Small - User Manual

## REDMON Since 1883 7455

# REDMON Precision Digital Pet Scale

Model: 7455  
Brand: REDMON Since 1883

## INTRODUCTION

---

The REDMON Precision Digital Pet Scale is designed to accurately measure the weight of small pets. This scale is an essential tool for monitoring your pet's health, managing their weight, and ensuring correct dosages for weight-dependent products. It features precision digital technology for clear, easy-to-read weight displays and includes functions for unit conversion, tare, and memory.

## SAFETY INFORMATION

---

Please read all safety instructions before using the scale.

- Ensure the scale is placed on a stable, flat surface to prevent tipping and ensure accurate readings.
- Keep the scale away from water or excessive moisture to prevent damage to electronic components.
- Do not overload the scale. The maximum weight capacity is 55 pounds (25 kg).
- This product contains chemicals known to the State of California to cause cancer and birth defects or other reproductive harm. Handle with care and wash hands after use.
- Batteries are not included. Use two AAA batteries. Ensure correct polarity when inserting batteries. Remove batteries if the scale will not be used for an extended period.

## PACKAGE CONTENTS

---

Verify that all items are present in the package:

- REDMON Precision Digital Pet Scale (Model 7455)
- User Manual (this document)

*Note: Two AAA batteries are required for operation and are not included.*

## SETUP

---

### 1. Battery Installation:

Locate the battery compartment on the underside of the scale. Open the compartment cover. Insert two AAA batteries, ensuring the correct positive (+) and negative (-) orientation. Close the battery compartment cover securely.

### 2. Placement:

Place the scale on a hard, flat, and stable surface. Uneven surfaces can affect accuracy. Avoid placing it on carpets or soft flooring.



Image: Top view of the REDMON Precision Digital Pet Scale, showing the flat weighing surface and control panel.

## OPERATING INSTRUCTIONS

Familiarize yourself with the control buttons: **UNIT**, **HOLD**, and **On/Off TARE**.

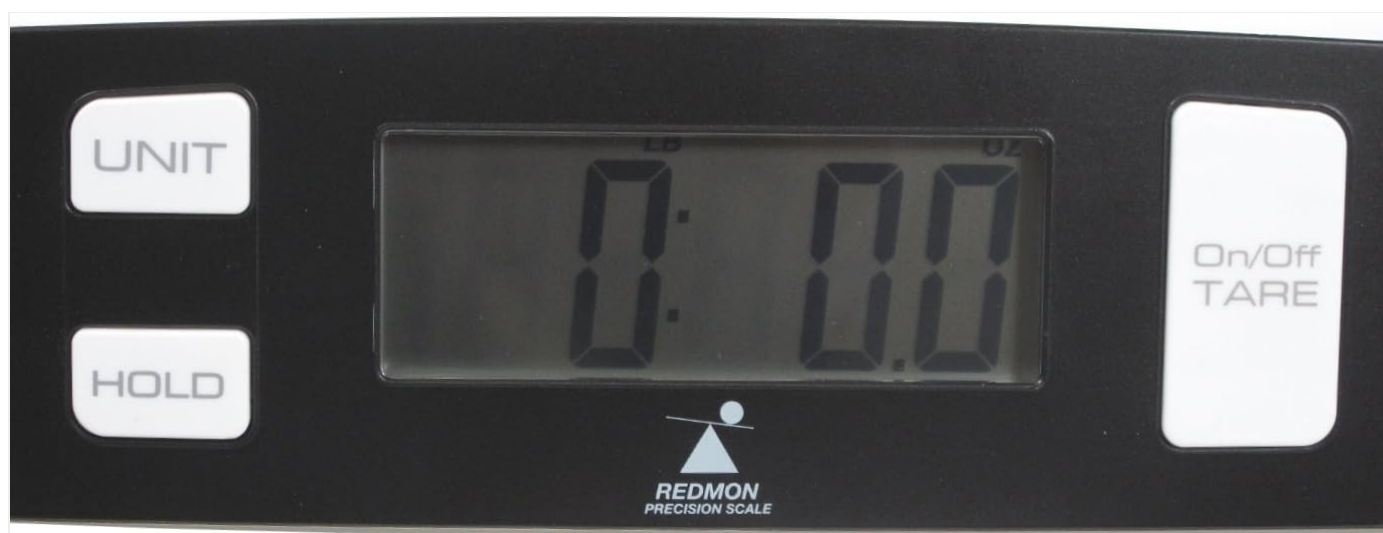


Image: Close-up view of the scale's digital display and control buttons (UNIT, HOLD, On/Off TARE).

### 1. Turning On/Off:

Press the **On/Off TARE** button to turn the scale on. The display will show "0.00" when ready. To turn off the scale, press and hold the **On/Off TARE** button for a few seconds until the display turns off. The scale also features an auto-off function to conserve battery power if left inactive.

### 2. Changing Units (Pounds/Kilograms):

With the scale on, press the **UNIT** button to toggle between pounds (lb) and kilograms (kg). The selected unit will be

indicated on the display.

### 3. Tare Function:

The tare function allows you to subtract the weight of a container (e.g., a pet carrier or blanket) so that only the pet's weight is measured.

- Place the empty container on the scale.
- Wait for the weight to stabilize.
- Press the **On/Off TARE** button. The display will reset to "0.00".
- Now, place your pet inside the container. The display will show only your pet's weight.

### 4. Hold Function:

If your pet is restless, the **HOLD** function can help capture a stable weight reading.

- Place your pet on the scale.
- Press the **HOLD** button. The scale will attempt to lock onto a stable weight reading.
- Once a stable reading is obtained, it will be displayed and held, even if the pet moves or is removed.
- Press **HOLD** again to release the reading and return to live weighing.

### 5. Memory Feature (Previous Weight):

The scale can recall the last measured weight, which is useful for tracking your pet's progress over time.

- After a weight has been measured and displayed, the scale automatically stores it.
- To view the previous weight, turn the scale on and then press the **HOLD** button twice quickly. The display will briefly show the last recorded weight.

### 6. Weighing Your Pet:

Ensure the scale is on a flat surface and displaying "0.00". Gently place your pet on the center of the scale. Wait for the weight to stabilize on the display. If your pet is moving, use the **HOLD** function.



Image: A small dog comfortably lying on the REDMON Precision Digital Pet Scale during a weighing session.

- **Cleaning:** Wipe the scale surface with a soft, damp cloth. Do not use abrasive cleaners or immerse the scale in water.
- **Storage:** Store the scale in a dry place at room temperature. Its slim design allows for easy storage.
- **Battery Replacement:** Replace batteries when the low battery indicator appears on the display.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Scale does not turn on.	Batteries are dead or incorrectly installed.	Check battery polarity. Replace with new AAA batteries.
Inaccurate readings.	Scale is on an uneven surface; weight exceeds capacity; pet is moving excessively.	Place scale on a hard, flat surface. Ensure pet's weight is within 55 lb capacity. Use the HOLD function for moving pets.
Display shows "Err" or "OL".	Overload error.	Remove weight immediately. The scale's capacity is 55 lbs (25 kg).
Display shows "Lo".	Low battery indicator.	Replace both AAA batteries.

For further assistance, please refer to the [Troubleshooting Guide \(PDF\)](#) or contact customer support.

## SPECIFICATIONS

Feature	Detail
Model Number	7455
Brand	REDMON Since 1883
Maximum Weight Capacity	55 lbs (25 kg)
Minimum Weight Capacity	3.5 oz (0.10 kg)
Measurement Units	Pounds (lb) / Kilograms (kg)
Display Type	Digital LCD
Power Source	2 x AAA Batteries (not included)
Dimensions (L x W x H)	25 x 16 x 1.75 inches
Item Weight	5.5 pounds
Material	Plastic
Special Features	Delta Weight Calculation, Previous Weight Memory, Low Battery Indicator, Tare Function, Hold Function

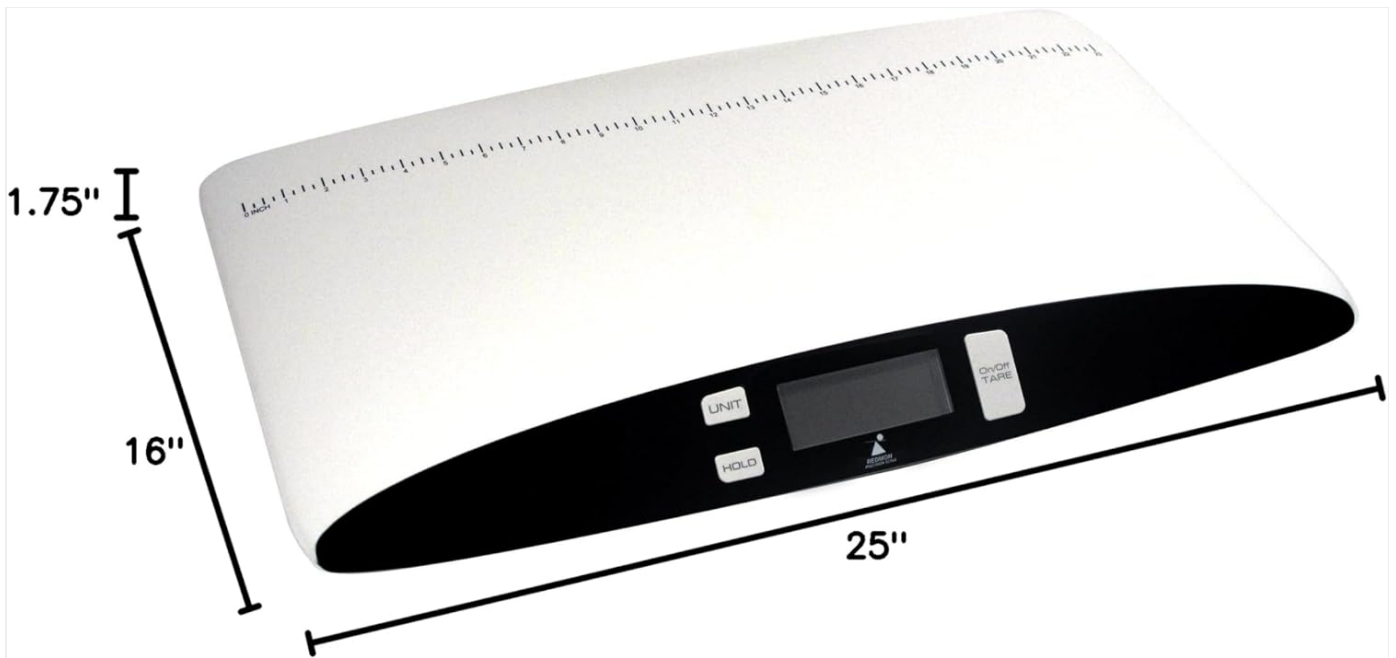


Image: Diagram showing the dimensions of the REDMON Precision Digital Pet Scale.

## WARRANTY AND SUPPORT

---

**Warranty:** The REDMON Precision Digital Pet Scale comes with a 1 Year Limited Warranty. Please retain your proof of purchase for warranty claims.

**Support:** For additional information, detailed instructions, or technical support, please refer to the official user guide or troubleshooting guide:

- [User Guide \(PDF\)](#)
- [Troubleshooting Guide \(PDF\)](#)

You may also visit the [REDMON Since 1883 brand page](#) for more product information.

