

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

> [Keuco](#) /

> [Keuco SK Griff WC Plan Care 34903 Right-Hand Grab Bar with Wireless Release - User Manual](#)

Keuco 34903011838

Keuco SK Griff WC Plan Care 34903 Right-Hand Grab Bar with Wireless Release

Model: 34903011838 | Brand: Keuco

[Introduction](#) [Safety Instructions](#) [Package Contents](#) [Setup & Installation](#) [Operating Instructions](#)
[Maintenance](#) [Troubleshooting](#) [Specifications](#) [Warranty & Support](#)

1. INTRODUCTION

This manual provides detailed instructions for the installation, operation, and maintenance of your Keuco SK Griff WC Plan Care 34903 Right-Hand Grab Bar with Wireless Release. Please read this manual thoroughly before installation and use to ensure proper function and safety. Keep this manual for future reference.

2. SAFETY INSTRUCTIONS

- **Professional Installation Recommended:** Installation should be performed by a qualified professional to ensure secure mounting and compliance with local building codes.
- **Weight Capacity:** Do not exceed the maximum weight capacity specified for this grab bar. Overloading can lead to failure and injury.
- **Secure Mounting:** Ensure the grab bar is securely fastened to a structural support (e.g., wall studs, blocking) capable of bearing the intended load. Do not mount to drywall alone.
- **Regular Inspection:** Periodically check all mounting screws and connections for tightness. Retighten if necessary.
- **Cleaning:** Use only mild, non-abrasive cleaners. Harsh chemicals can damage the finish.
- **Wireless Release:** Ensure the wireless release mechanism is tested regularly for proper function.

3. PACKAGE CONTENTS

Please verify that all components are present before beginning installation:

- Keuco SK Griff WC Plan Care 34903 Right-Hand Grab Bar

- Mounting Hardware (screws, anchors - specific type may vary based on wall material)
- Wireless Release Module (integrated or separate, depending on model variant)
- Installation Template (if included)

If any parts are missing or damaged, please contact Keuco customer support.

4. SETUP & INSTALLATION

Proper installation is crucial for the safety and functionality of the grab bar. Refer to the technical drawings for precise measurements and mounting points.

4.1. Tools Required

- Drill
- Drill bits appropriate for wall material
- Level
- Measuring tape
- Pencil
- Screwdriver
- Stud finder (for drywall installation)

4.2. Installation Steps

1. **Determine Mounting Location:** Identify the desired height and position for the grab bar. Ensure it is accessible and provides adequate support. Use a stud finder to locate wall studs if mounting to drywall.
2. **Mark Drilling Points:** Hold the grab bar or use the provided template against the wall. Use a level to ensure it is perfectly horizontal. Mark the drilling points with a pencil.
3. **Drill Pilot Holes:** Drill pilot holes at the marked locations. The size of the drill bit will depend on your wall material and the type of anchors/screws being used.
4. **Insert Anchors (if necessary):** If mounting into hollow walls or masonry, insert appropriate wall anchors into the pilot holes.
5. **Mount the Grab Bar:** Align the grab bar with the drilled holes and secure it using the provided screws. Tighten all screws firmly, but do not overtighten.
6. **Test Security:** Gently apply pressure to the grab bar to ensure it is securely mounted and does not wobble.
7. **Wireless Release Setup:** If your model includes a separate wireless release module, follow its specific instructions for pairing and placement. For integrated units, ensure the release button functions correctly.

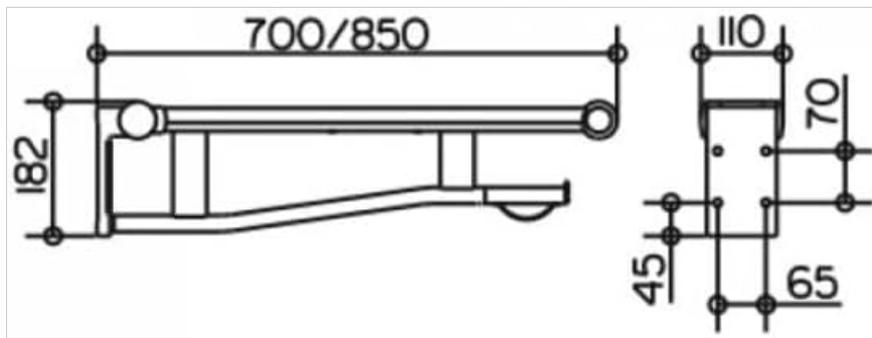


Figure 1: Technical drawing illustrating the overall dimensions and profile of the Keuco SK Griff WC Plan Care grab bar. Key measurements include a variable length of 700mm or 850mm, a total height of 182mm, and a depth of 110mm from the wall, with a 70mm grip height.

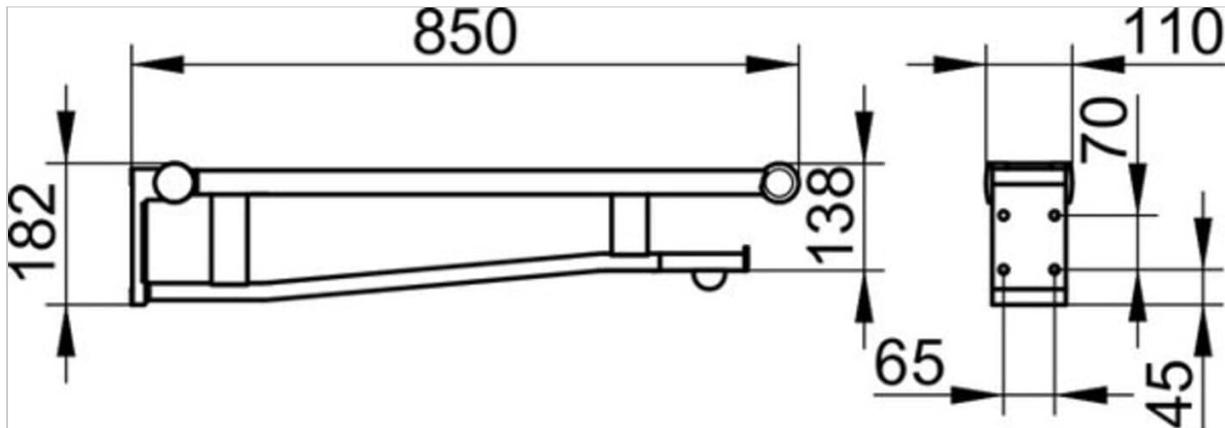


Figure 2: Detailed technical drawing of the Keuco SK Griff WC Plan Care grab bar, specifically showing the 850mm length variant. It highlights a total height of 182mm, a grip height of 138mm from the bottom edge, and a wall projection of 110mm with a 70mm grip depth.

5. OPERATING INSTRUCTIONS

The Keuco SK Griff WC Plan Care grab bar is designed to provide stable support. The integrated wireless release mechanism offers additional functionality.

5.1. Using the Grab Bar for Support

- Grasp the grab bar firmly with your hand.
- Apply pressure as needed to assist with sitting, standing, or maintaining balance.
- Ensure your grip is secure before transferring weight.

5.2. Using the Wireless Release

The wireless release function is typically used for specific care scenarios, allowing for remote activation of a call system or other integrated functions. Refer to the specific instructions provided with your wireless system for pairing and operation.

- Locate the wireless release button on the grab bar.
- Press the button firmly to activate the connected function.
- Ensure the connected receiver is within range and properly configured.

6. MAINTENANCE

Regular cleaning and inspection will prolong the life and maintain the appearance of your grab bar.

6.1. Cleaning

- Clean the grab bar regularly with a soft cloth and a mild, non-abrasive cleaning solution.
- Avoid using harsh chemicals, abrasive cleaners, scouring pads, or solvents, as these can damage the chrome finish and plastic components.
- Rinse thoroughly with clean water and dry with a soft cloth to prevent water spots.

6.2. Inspection

- Periodically check all mounting screws for tightness. Retighten if any looseness is detected.
- Inspect the grab bar for any signs of damage, corrosion, or wear. If damage is found, discontinue use and contact a qualified professional for assessment or replacement.
- Test the wireless release function periodically to ensure it is operating correctly.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Grab bar feels loose or wobbly	Mounting screws are loose; improper installation; wall structure compromised.	Tighten all mounting screws. If problem persists, consult a professional to inspect the installation and wall structure. Do not use if unstable.
Wireless release not functioning	Low battery in wireless module; module out of range; pairing issue; receiver malfunction.	Check and replace batteries in the wireless module. Ensure the module is within range of its receiver. Re-pair the module with the receiver according to its specific instructions. Consult the receiver's manual for troubleshooting.
Surface finish appears dull or stained	Residue from cleaning products; hard water spots; lack of regular cleaning.	Clean with a mild, non-abrasive cleaner and a soft cloth. Rinse thoroughly and dry immediately. For stubborn spots, a specialized chrome cleaner may be used sparingly.

8. SPECIFICATIONS

Product Dimensions: 10.31 x 36.69 x 4.41 inches (26.19 x 93.20 x 11.20 cm)

Item Weight: 12.91 pounds (5.86 kg)

Manufacturer: KEUCO

ASIN: B005E7O1RM

Item Model Number: 34903011838

Batteries Required: No (for the grab bar itself; wireless release module may require batteries)

Material: Chrome/Light Grey finish (specific materials may vary)

Type: Right-Hand Grab Bar

Feature: Integrated Wireless Release

9. WARRANTY & SUPPORT

9.1. Warranty Information

Keuco products are manufactured to high-quality standards. For specific warranty terms and conditions, please refer to the warranty card included with your product or visit the official Keuco website. Keep your proof of purchase for warranty claims.

9.2. Customer Support

If you have any questions, require technical assistance, or need to report a problem, please contact Keuco customer support through their official website or the contact information provided with your product documentation.

Website: www.keuco.com (Please check for regional variations)

Service Hotline: Refer to your regional Keuco website for local contact numbers.