

[Manuals.plus](#) /

> [Drive Medical](#) /

> Drive Medical 4-Wheel Rollator Walker with Seat (Model 10257RD-1) Instruction Manual

## Drive Medical 10257RD-1

# Drive Medical 4-Wheel Rollator Walker with Seat (Model 10257RD-1) Instruction Manual

## INTRODUCTION

The Drive Medical 4-Wheel Rollator Walker with Seat is engineered to offer reliable support, comfortable mobility, and user-friendly features for individuals seeking enhanced stability. Its robust steel frame ensures a steady walking experience both indoors and outdoors. Key features include a padded seat for resting, a removable hinged backrest for support, smooth-rolling 7.5-inch caster wheels for various terrains, and convenient hand brakes for speed control and secure locking. An under-seat pouch provides practical storage for personal items, promoting daily independence.

This rollator is designed for durability with a 350-lb weight capacity and a long-lasting steel frame, making it a dependable companion for daily support, rehabilitation, and maintaining an active lifestyle.



## 4 Wheel Rollator

Get the confidence you need to stay on the go indoors and outdoors with the 4 Wheel Rollator

1. **Large 7.5-inch caster wheels** handle bumps and diverse terrain with ease
2. **Under-seat storage pouch** transports belongings securely
3. **12-inch padded seat** lets you rest on the go

## SETUP AND ASSEMBLY

---

The Drive Medical 4-Wheel Rollator Walker (Model 10257RD-1) is designed for straightforward assembly. Follow these steps to prepare your rollator for use:

1. **Unfold the Rollator:** Carefully pull the frame apart until it locks securely into its open position. Ensure all locking mechanisms are engaged.
2. **Adjust Handle Height:** Locate the easy-turn levers on the handle posts. Loosen the levers, adjust the handles to a comfortable height (typically between 31-35 inches), and then securely tighten the levers. The handles should be at wrist height when standing upright with arms relaxed.
3. **Position the Backrest:** The rollator includes a removable hinged backrest. Position it up for back support when seated, or fold it down for compact storage.
4. **Check Brakes:** Before first use, ensure the hand brakes are functioning correctly. Squeeze the brake levers to slow down or stop. Push the levers down to lock the brakes for stationary seating.



Image: Adjusting the handle height for personalized comfort and proper posture.

## Assembly Video Guide

Video: This video provides a detailed review of the Drive Medical Rolling Walker, including insights into its assembly process and features. (Source: Stuff with Leon)

## OPERATING INSTRUCTIONS

---

Using your Drive Medical Rollator Walker is intuitive and designed for ease of use. Always ensure the rollator is fully unfolded and stable before use.

1. **Walking:** Place your hands on the ergonomic handles. Maintain an upright posture and walk naturally, allowing the rollator to roll smoothly in front of you.
2. **Braking:** To slow down or stop, gently squeeze both hand brake levers simultaneously. For temporary stops, a light squeeze is sufficient.

3. **Locking Brakes for Seating:** When you wish to sit, push both hand brake levers down until they click into the locked position. This prevents the rollator from moving while you sit or stand. Always ensure brakes are locked before attempting to sit or stand.
4. **Sitting:** With the brakes locked, turn around so your back is to the seat. Slowly lower yourself onto the padded seat. Use the backrest for support if desired.
5. **Standing Up:** Ensure the brakes are still locked. Push off from the seat and handles to stand up. Once standing and stable, release the brake levers by pulling them up to unlock the wheels.

## The Stability You Need to Stay Independent

The sturdy and easy-to-use 4 Wheel Rollator from Drive DeVilbiss Healthcare gives you support and comfort you can count on, all at an affordable price

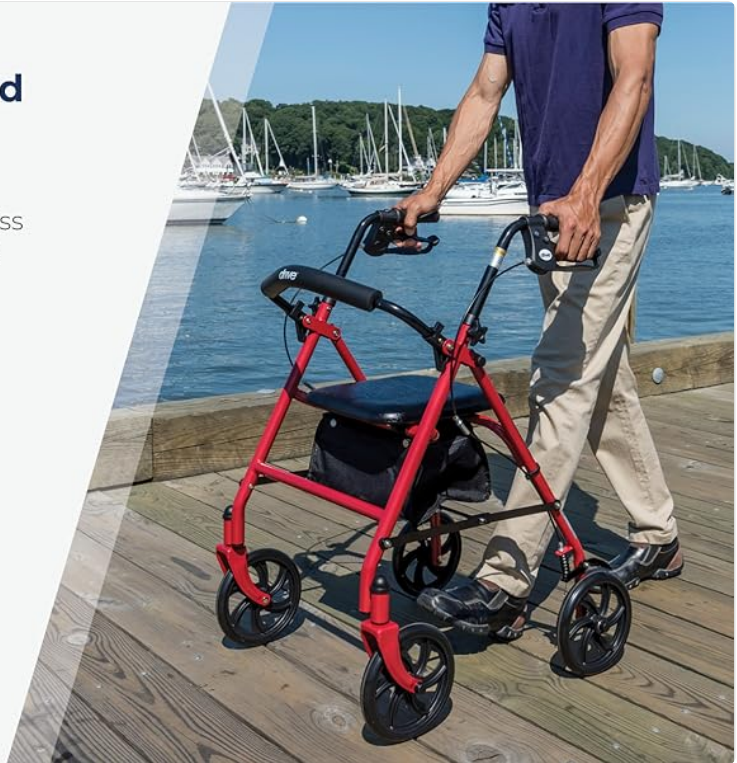


Image: The rollator provides stability for walking on various surfaces.

### Operation Video Guide

Video: A physical therapist reviews the Drive Medical Rollator, demonstrating its features and proper usage for mobility and safety. (Source: Dr. Chuck)

Video: This video showcases the sturdy design and features of the 4-Wheel Walker with Seat, highlighting its ease of use and storage. (Source: CJ's Finds)

### FEATURES

- **Durable Steel Frame:** Provides stable walking support and improved balance, built for long-term use.
- **Padded Seat with Backrest:** A comfortable padded seat allows for resting during longer outings, complemented by a removable hinged backrest for additional support.
- **Smooth 7.5-inch Caster Wheels:** Designed for easy navigation across various indoor and outdoor surfaces.
- **Secure Hand Brakes:** Convenient hand brakes offer easy speed control and can be locked for stationary seating.
- **Height Adjustable Handles:** Ergonomic handles can be adjusted to a user's preferred height (31-35 inches) for optimal comfort and posture.
- **Under-Seat Storage Pouch:** A discreet and convenient pouch located under the seat provides storage for personal belongings.

- **Compact Folding Design:** The rollator folds compactly for easy storage and transport, ideal for active lifestyles.

## Designed for Safety & Convenience



### Sturdy Frame

Steel construction supports up to 300 pounds



### Adjustable Handles

Easy-turn levers let you quickly adjust the handles between 31 and 35 inches for a comfortable reach



### Easy-to-Reach Brakes

Deluxe hand brake can be pressed when needed or pushed down to lock rollator into place



### Compact Fold

Folds easily for compact storage at home or in the trunk of your car

Image: Visual representation of the rollator's main features for safety and convenience.

## Adjustable, ergonomic handles & hand breaks



Image: The padded seat and backrest offer a comfortable resting spot.



**Rest comfortably with  
fold-down backrest**



Image: Ergonomic hand brakes provide secure stopping and parking capabilities.

**Easy to fold, store  
& transport**



**Weights only 20 lbs.**

Image: The rollator's compact folding design allows for easy storage and transport.

## MAINTENANCE

To ensure the longevity and safe operation of your Drive Medical Rollator Walker, regular maintenance is recommended:

- **Cleaning:** Wipe down the frame and seat regularly with a damp cloth and mild soap. Avoid abrasive cleaners.
- **Wheel Inspection:** Periodically check the wheels for wear and tear. Ensure they rotate freely and are free of debris.
- **Brake Check:** Test the hand brakes before each use to confirm they engage properly and hold the rollator securely when locked. Adjust if necessary (refer to assembly instructions for adjustment points).
- **Fastener Tightness:** Regularly inspect all screws, nuts, and bolts to ensure they are tight. Tighten any loose fasteners to maintain structural integrity.

## TROUBLESHOOTING

If you encounter any issues with your rollator, refer to the following common troubleshooting tips:

- **Rollator Does Not Roll Smoothly:** Check wheels for debris (hair, lint, etc.) and remove if present. Ensure wheels are not damaged.
- **Brakes Not Engaging/Locking:** Inspect the brake cables for any kinks or damage. Adjust the brake tension using the adjustment mechanism near the wheels or handles. If issues persist, consult a qualified technician.
- **Unstable Frame:** Ensure all frame components are fully extended and locked into place. Check all fasteners for tightness and tighten if loose.

For persistent issues or concerns not addressed here, please contact Drive Medical customer support.

## SPECIFICATIONS

<b>Model Number</b>	10257RD-1
<b>Product Dimensions</b>	25.5"D x 31"W x 23.5"H
<b>Item Weight</b>	17.9 Pounds
<b>Material</b>	Alloy Steel
<b>Color</b>	Red
<b>Manufacturer</b>	Drive Medical

## WARRANTY AND SUPPORT

For specific warranty information, please refer to the documentation included with your product at the time of purchase. Drive Medical products typically come with a manufacturer's warranty covering defects in materials and workmanship.

For customer support, replacement parts, or further assistance, please visit the official Drive Medical website or contact their customer service department directly. Contact information can usually be found on the product packaging or on the manufacturer's website.