## PHONAK 098-0016

# PHONAK C&C Line Tool Set Instruction Manual

Model: 098-0016

## INTRODUCTION

This instruction manual provides detailed guidance for the proper use and maintenance of the PHONAK C&C Line Tool Set. This set is designed to assist in the routine cleaning of hearing aid components, specifically vents, tubing, and for wax removal, ensuring optimal performance and hygiene of your hearing device.

## PRODUCT OVERVIEW

The PHONAK C&C Line Tool Set includes the following components:


- **Vent/Tubing Cleaning Lines:** Two flexible lines for clearing blockages in hearing aid vents and tubing.
- **Wax Removal Tool:** A specialized tool for carefully removing wax buildup.
- **Cleaning Brush:** A small brush for general cleaning of surfaces and openings.

## SETUP

Before beginning any cleaning procedure, ensure your hands are clean and dry. Work in a well-lit area on a clean, flat surface to prevent loss of small components. Identify the specific parts of your hearing aid that require cleaning.

## OPERATING INSTRUCTIONS

Follow these steps for effective cleaning of your hearing aid components using the PHONAK C&C Line Tool Set.

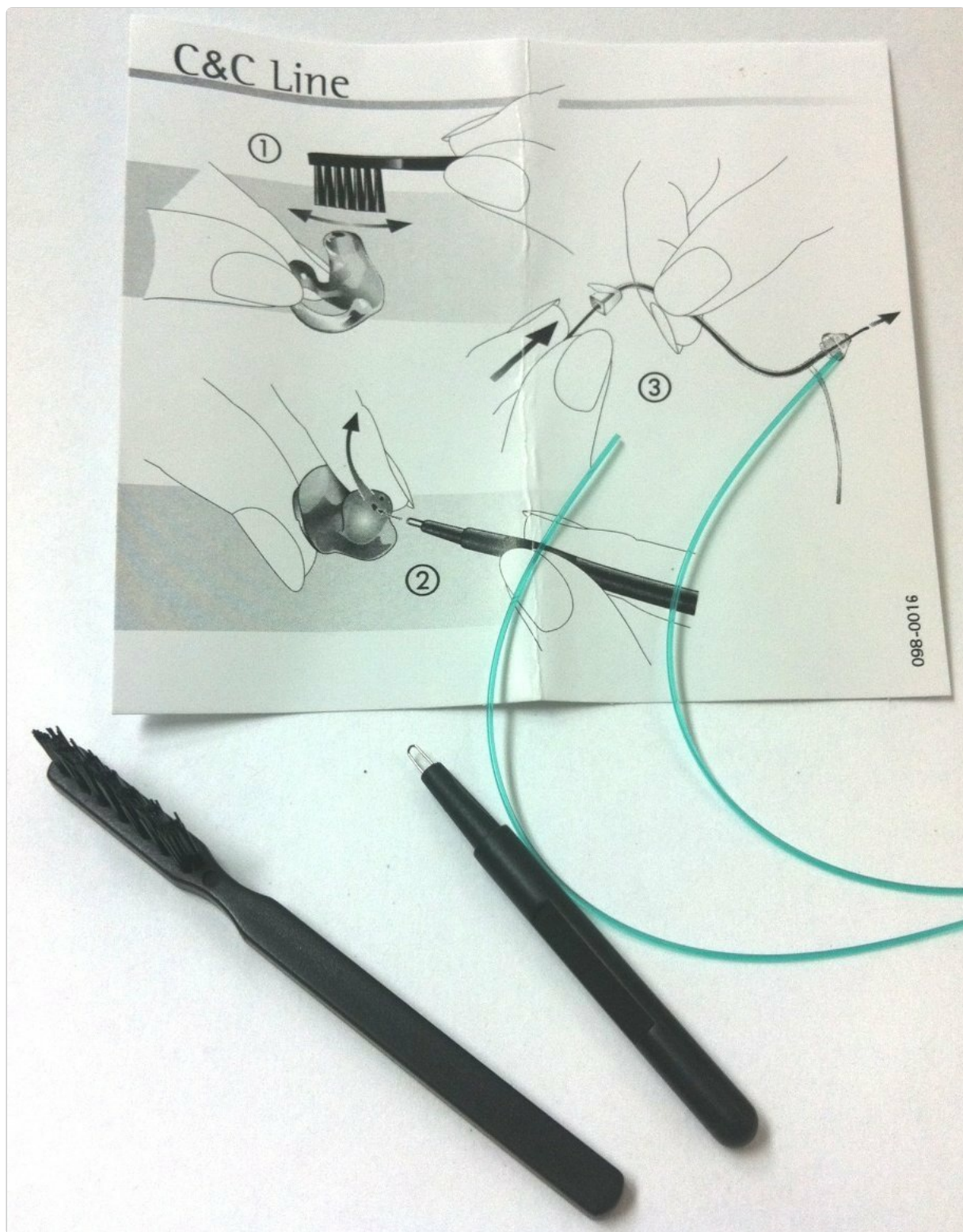*Image: The PHONAK C&C Line Tool Set, including a cleaning brush, a wax removal tool, and two flexible cleaning lines, laid out with an instruction diagram illustrating the cleaning process.*

**1. Step 1: Cleaning the Earmold/Shell Surface (Using Cleaning Brush)**

Gently brush the surface of the earmold or hearing aid shell to remove any visible debris or earwax. Pay particular attention to the sound outlet and microphone openings. Refer to diagram point 1.

**2. Step 2: Clearing Vents (Using Wax Removal Tool)**

If your hearing aid has a vent, carefully insert the pointed end of the wax removal tool into the vent opening. Gently push through to clear any blockages. Avoid excessive force. Refer to diagram point 2.

**3. Step 3: Cleaning Tubing (Using Vent/Tubing Cleaning Lines)**

For hearing aids with tubing, thread one of the flexible cleaning lines through the tubing, starting from one end and pushing it all the way through until it emerges from the other end. This will dislodge any moisture or

wax buildup. Repeat if necessary. Refer to diagram point 3.

MAINTENANCE

After each use, clean the tools to ensure their longevity and hygiene. Wipe the cleaning brush and wax removal tool with a dry, clean cloth. The cleaning lines can be rinsed with warm water and allowed to air dry completely before storage. Store all tools in a clean, dry place.

TROUBLESHOOTING

- **Difficulty clearing tubing:** Ensure the cleaning line is fully inserted and gently rotated. If blockage persists, consult your hearing care professional.
- **Tools appear damaged:** Discontinue use of damaged tools to prevent harm to your hearing aid. Replace the tool set as needed.
- **Hearing aid still sounds muffled after cleaning:** The issue may not be related to external blockages. Consult your hearing care professional for further diagnosis.

SPECIFICATIONS

- **Model Number:** 098-0016
- **UPC:** 610373682521
- **Manufacturer:** PHONAK
- **Package Dimensions:** Approximately 3 x 1.5 x 0.02 inches
- **Item Weight:** Approximately 0.64 ounces

WARRANTY AND SUPPORT

Due to hygiene reasons, this product is not eligible for return or exchange once opened. For any questions regarding the use of this tool set or your PHONAK hearing aid, please contact your authorized PHONAK dealer or hearing care professional.

Related Documents - 098-0016

--	--

Bedienungsanleitung	Klapphilfe	Kalkulator funktionell
Instrucciones de uso	Manueli di istruzioni	Hesbik dila xiz xing
Mode d'emploi	Brükeren leiding	取携演算機
Istruzioni per l'uso	Instrukcja obsługi	사용 설명서
Gebrauchsanweisung	Brükeren leiding	使用说明书
Manual de instrucciones	Procytywanie instrukcji	نشرة إلكترونية
Gebruiksaanwijzing	Spis treści	מדריך

The Multiset set includes necessary tools for wax removal and cleaning of earmolds and hearing systems. It contains a cleaning brush with an integrated magnet system for easy battery removal.

**PRECAUTIONS FOR USE:** External use. Keep out of reach of children.

If you experience problems, contact a hearing aid specialist. All serious incidents related to the device should be the subject of a notification to the manufacturer and the competent authorities in the Member State in which the user is located.

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
A comprehensive user guide for the Phonak Multitool, detailing its use for cleaning hearing aids and earmolds, wax removal, and easy battery replacement with its integrated magnet system. Includes essential safety precautions.

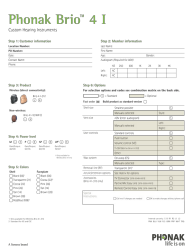	<p><a href="#">Phonak RIC Hearing Aid Maintenance Guide for Marvel, Paradise, and Lumity</a></p> <p>A comprehensive guide from Phonak on how to clean, maintain, and troubleshoot your Marvel, Paradise, and Lumity RIC hearing aids, including filter replacement and moisture management. Learn about ear domes, filter types, drying procedures, and service options.</p>
	<p><a href="#">Easy Line Remote 3.0 User Guide: Control Your Hearing Aids with Your Smartphone</a></p> <p>Learn how to use the Easy Line Remote 3.0 app by Sonova/Phonak to control your hearing aids, adjust settings, access remote support, and manage your hearing diary. This guide covers installation, pairing, features, and troubleshooting.</p>
	<p><a href="#">Phonak RemoteControl User Guide - Enhance Your Hearing Experience</a></p> <p>Comprehensive user guide for the Phonak RemoteControl, detailing setup, operation, troubleshooting, safety, and warranty information for managing hearing aid programs and volume.</p>
	<p><a href="#">Phonak Cassia Hearing Aids: Product Information and Features</a></p> <p>Explore the Phonak Cassia range of hearing instruments, detailing key features like SoundRecover, UltraZoom, and WhistleBlock, along with technical specifications and fitting options for mild to profound hearing loss.</p>
	<p><a href="#">Phonak Brio 4 I Custom Hearing Instruments Order Form and Configuration Guide</a></p> <p>Detailed order form and configuration guide for Phonak Brio 4 I custom hearing instruments, including customer and member information, product options, power levels, colors, and a comprehensive combination matrix for Brio 4 I-10 NW O and Brio 4 I-312 models.</p>