

## Accu-Chek 102 Lancets

# Accu-Chek FastClix Lancets User Manual

For use with Accu-Chek FastClix Lancing Device

[Introduction](#)

[Setup](#)

[Operation](#)

[Maintenance](#)

[Troubleshooting](#)

[Specifications](#)

[Warranty & Support](#)

## 1. INTRODUCTION

The Accu-Chek FastClix Lancets are designed for use with the Accu-Chek FastClix lancing device to obtain a small blood sample for blood glucose testing. These lancets are pre-set in a drum, eliminating the need to handle individual needles, which enhances safety and convenience. Each drum contains 6 sterile lancets, and this package includes 102 lancets, provided in 17 drums.



Image 1.1: Accu-Chek FastClix Lancets product packaging. The box displays 'Accu-Chek FastClix 102 Lancets' and an image of a lancet drum.



**SAME PRODUCT, NEW LOOK**  
**Packaging received may vary**

Image 1.2: Comparison of previous and current Accu-Chek FastClix Lancets packaging. The image indicates 'Same Product, New Look'.

## 2. SETUP

### 2.1 Loading the Lancet Drum

To prepare your Accu-Chek FastClix lanceting device for use, follow these steps to load a new lancet drum:

1. Remove the cap from your Accu-Chek FastClix lanceting device.
2. Insert a new drum of 6 pre-set FastClix lancets into the empty device. Ensure it clicks securely into place.
3. Replace the cap onto the lanceting device.



(lancing device  
not included)

**Compatible with all models of  
Accu-Chek® FastClix lancing devices**

Image 2.1: A hand inserting a lancet drum into the Accu-Chek FastClix lancing device. The lancing device is not included with the lancets.

- 1** Insert drum of pre-set FastClix lancets into empty device



- 2** Place cap over drum and click into place



- 3** Adjust lancing depth setting by twisting cap



## How to load FastClix lancet drum into device

Image 2.2: Step-by-step visual guide showing how to insert the lancet drum and adjust the lancing depth by twisting the cap.

### 2.2 Adjusting Lancing Depth

The Accu-Chek FastClix lancing device offers 11 customizable depth settings to ensure comfortable and effective blood sampling. To adjust the depth:

- Twist the cap of the lancing device to select your desired penetration depth. Start with a shallower setting and adjust as needed for an adequate blood sample.



## Anatomy of FastClix device

Image 2.3: Diagram illustrating the anatomy of the FastClix device, highlighting the customizable depth setting (1-11) and the lancet drum.

### 3. OPERATION

#### 3.1 Performing a Blood Glucose Test

Once the lancet drum is loaded and the depth is set, you can perform a blood glucose test:

1. **Prepare Your Hands:** Wash your hands thoroughly with soap and warm water, then dry them completely.
2. **Prepare the Lancing Device:** Ensure the lancing device is primed. The Accu-Chek FastClix device typically primes with a single click.
3. **Obtain Blood Sample:** Place the lancing device firmly against the side of your fingertip. Press the release button to activate the lancet and obtain a blood drop.
4. **Apply Blood to Test Strip:** Gently squeeze your finger to form a round drop of blood. Apply the blood sample to the designated area of your test strip (used with a compatible blood glucose meter).



# ACCU-CHEK®

EXPERIENCE WHAT'S POSSIBLE



## FastClix Lancing Device

Checking your blood sugar is  
SIMPLY EASY – 1 click is all it takes



Fast &  
Easy



Precision  
Guided



Multi-Depth  
Settings



AST  
Adaptable

Image 3.1: A visual guide detailing five simple steps for blood glucose testing, including hand washing, inserting a test strip, collecting a blood sample, applying the sample, and syncing results.



**Designed for comfort,  
easy to use**

## 4. MAINTENANCE

### 4.1 Cleaning the Lancing Device

Regular cleaning of your Accu-Chek FastClix lancing device is important to maintain hygiene and proper function. Refer to your lancing device's specific instruction manual for detailed cleaning procedures.

### 4.2 Disposing of Used Lancet Drums

After all 6 lancets in a drum have been used, the drum should be safely disposed of. Do not attempt to reuse lancets or lancet drums. Place used drums in a puncture-resistant sharps container immediately after use. Follow local regulations for proper disposal of sharps waste.

## 5. TROUBLESHOOTING

### 5.1 Insufficient Blood Sample

- **Adjust Depth Setting:** If the blood sample is too small, increase the lancing depth setting on your FastClix device.
- **Warm Hands:** Ensure your hands are warm before testing, as this can improve blood flow.
- **Proper Technique:** Press the lancing device firmly against your fingertip.

### 5.2 Lancet Drum Not Advancing

- Ensure the drum is correctly inserted and clicked into place.
- Verify that the lancing device is properly primed.
- If the issue persists, the lancing device may require service or replacement.

## 6. SPECIFICATIONS

Product Name	Accu-Chek FastClix Diabetes Lancets
Model Number	102 Lancets
Quantity	102 Lancets (17 drums of 6 lancets each)
Compatibility	Accu-Chek FastClix lancing device only
Dimensions	2.4 x 2.3 x 1.7 inches (packaging)
Weight	1.13 ounces (packaging)
UPC	365702288102







## 7. WARRANTY & SUPPORT

For questions regarding the Accu-Chek FastClix Lancets or the Accu-Chek FastClix lancing device, please contact Accu-Chek customer support directly. Specific warranty information for lancets is typically not provided as they are disposable medical supplies. Always refer to the documentation included with your lancing device for its warranty details.




For further assistance, visit the official Accu-Chek website or consult your healthcare professional.

## Related Documents - 102 Lancets

	<p><a href="#">Accu-Chek Instant Blood Glucose Meter User Guide and Support</a></p> <p>Comprehensive guide to using the Accu-Chek Instant blood glucose meter, including setup, testing, connecting with the mySugr app, and customer support information.</p>
	<p><a href="#">Accu-Chek Instant Quick Reference Guide: Blood Glucose Monitoring</a></p> <p>A quick reference guide for the Accu-Chek Instant blood glucose meter. Learn how to perform blood glucose tests, set up the device, change batteries, and use the Accu-Chek FastClix lancing device for diabetes management.</p>
	<p><a href="#">Accu-Chek Guide User's Manual: Blood Glucose Monitoring</a></p> <p>Explore the Accu-Chek Guide Blood Glucose Monitoring System with its user-friendly manual. Learn about meter setup, testing procedures, data management, and maintenance for effective diabetes care.</p>
	<p><a href="#">Accu-Chek Guide Blood Glucose Monitoring System: User Guide and Test Strip Information</a></p> <p>Learn how to use the Accu-Chek Guide, Guide Link, and Guide Me blood glucose monitoring systems and test strips for effective diabetes management. Includes instructions, safety information, and performance data.</p>
	<p><a href="#">Accu-Chek Instant Blood Glucose Meter &amp; mySugr App User Guide</a></p> <p>Comprehensive guide to the Accu-Chek Instant blood glucose meter, covering setup, blood sugar testing, connecting to the mySugr app for data logging, and using the mySugr Bolus Calculator. Includes support information for UK and Ireland users.</p>
	<p><a href="#">Accu-Chek Instant Blood Glucose Meter: User Guide and mySugr App Integration</a></p> <p>Comprehensive guide to setting up and using the Accu-Chek Instant blood glucose meter, including connecting with the mySugr app for diabetes data logging and management. Learn about testing, app features, and the bolus calculator.</p>

[pdf] User Manual Guide

accu chek elearning and training content pptx au sites g files iut441 f Accu Chek e Learning Content See Details Source What is diabetes Diabetes Australia Available at diabetesaustralia about what Accessed June 2021 1 2 NEW Simply Easy Pack product 3 Guide Me Getting to know the meter Deal Mobile Manual Pdf \$FFX KHN 0RELOH |||

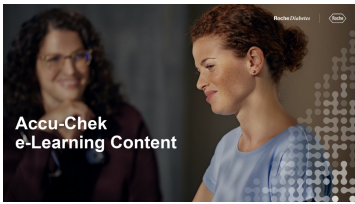

Accu-Chek e-Learning Content What will I learn from this module 1 What is Diabetes 2 NEW Accu-Chek Simply Easy Diabetes Pack product training 3 NEW Accu-Chek Guide Me product training: Getting to know the meter Getting to know the lancing device How to do perform a test Transferring Data ...


lang:en score:41 filesize: 3.34 M page\_count: 36 document date: 0000-00-00

[pdf] Instructions

pi 04807391001 02 EN ES PT indd V1 R1 black 2006 04 qEN Instructions for use The Accu Chek Multiclix lancing device is an easy to that has been developed allow virtually pain free collection of capillary blood from a multiclix accu chek ec file fid 12776 |||

ACCU-CHEK Multiclix LANCING DEVICE qEN Instructions for use The Accu-Chek Multiclix lancing device ... C-DmsC sh1e 22k l43aM0nlCaucnletcitsel it xs: IR:aREnFEcFe0t04d47ru77m7886s69: 835000011 17 drums **102 lancets** : REF 0 4778642001 34 drums 2300C4 lancets : REF 0 4778677001 No t2aCll pack sizes wil...

lang:en score:35 filesize: 2.77 M page\_count: 2 document date: 2006-05-11

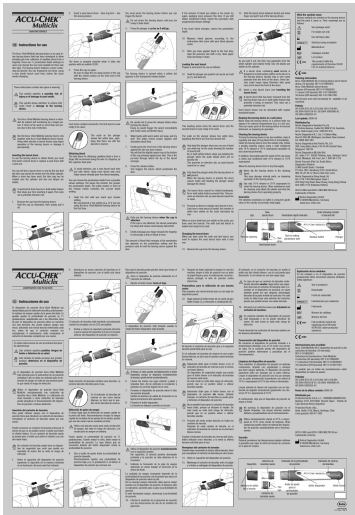


[pdf] Decleration of Conformity Catalog

Sterk Tatjana DCMG~Mannheim MDR IVDR DoC Accu Chek Fastclix Lancets Rev 02 Zertifikate de accu chek sites g files papvje231 2024 08 |||

DocuSign Envelope ID: F48F532C-DEDC-4925-AFEF-EAB0B97110FD EU-Konformittserklrung EU declaration of ... n 09770526 6 lancets 05981310 24 lancets 05208459 24 lancets 05353858 24 lancets 05208475 **102 lancets** 1/4 Roche Diabetes Care GmbH Sandhofer Strae 116; D-68305 Mannheim; Telefon 49-621-759...

lang:de score:32 filesize: 512.6 K page\_count: 4 document date: 2024-02-22

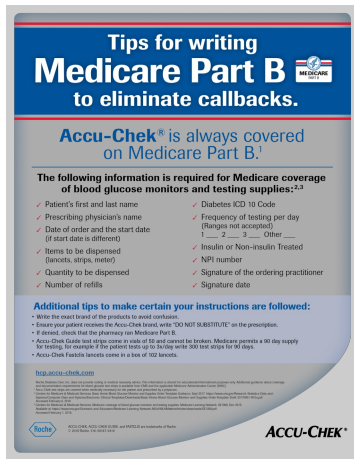

## [pdf] Instructions Guide

Medicare Part B PART Accu Chek Tips for writing 1 test strips are covered when medically necessary the patient and prescribed by a physician 2 Centers Medicaid Services Supplies Order Template Draft 20170901 R10a Accessed February 6 2018 3 Writing to Eliminate Callbacks hcp accu chek sites g files iut1121 f

Tips for writing Medicare Part B MEDICARE PART B to eliminate callbacks. Accu-Chek is always covered ... tests up to 3x/day write 300 test strips for 90 days. Accu-Chek FastClix lancets come in a box of **102 lancets**. [hcp.accu-chek.com](http://hcp.accu-chek.com) Roche Diabetes Care, Inc. does not provide coding or medical necessit...

lang:en score:29 filesize: 192.57 K page\_count: 1 document date: 2019-03-08



## [pdf]

The world s first 1 Click lancing device with a drum 5496 FastClix Fact Sheet Sept20 rochediabetes sites g files iut1376 f microsities en gb resources product information |||

Accu-Chek FastClix The world s first 1-Click lancing device with a drum 2.45 for **102 Lancets\*** Key points: No clinical waste subject to following local guidance Prime and release in 1-Click action - easy to teach 11 individual depth settings Clixmotion technology provides gentle lancing and h...

lang:en score:19 filesize: 452.22 K page\_count: 1 document date: 2020-09-03