

**Chicco 00004904000000**

# Chicco Phisio Clean Nasal Aspirator Instruction Manual

MODEL: 00004904000000

## 1. Product Overview

The Chicco Phisio Clean Nasal Aspirator is designed to gently and effectively remove excess nasal secretions from infants and young children. This manual aspirator helps clear nasal passages, promoting easier breathing and comfort for the child. It features a soft, ergonomic design for ease of use and is made from trusted, latex-free materials.



Image 1.1: The Chicco Phisio Clean Nasal Aspirator, designed for gentle nasal clearing.

## Eliminates excess nasal secretion



Image 1.2: The aspirator's function in eliminating excess nasal secretion.

## 2. Safety Information

Please read and understand all safety instructions before using the product. Keep this manual for future reference.

- Always keep the aspirator and its components out of reach of children when not in use.
- Do not insert the tip too deeply into the child's nostril.
- Never use the aspirator if any part is damaged or missing.
- Do not sterilize in a microwave. Follow cleaning instructions carefully.
- This product is intended for nasal aspiration only. Do not use for other purposes.
- If nasal congestion persists or worsens, consult a healthcare professional.

## 3. Components

The Chicco Phisio Clean Nasal Aspirator typically includes the following parts:

- Aspirator body with mouthpiece
- Soft silicone nasal tip

- Hygienic filter (internal)
- Collection chamber

## Trusted material

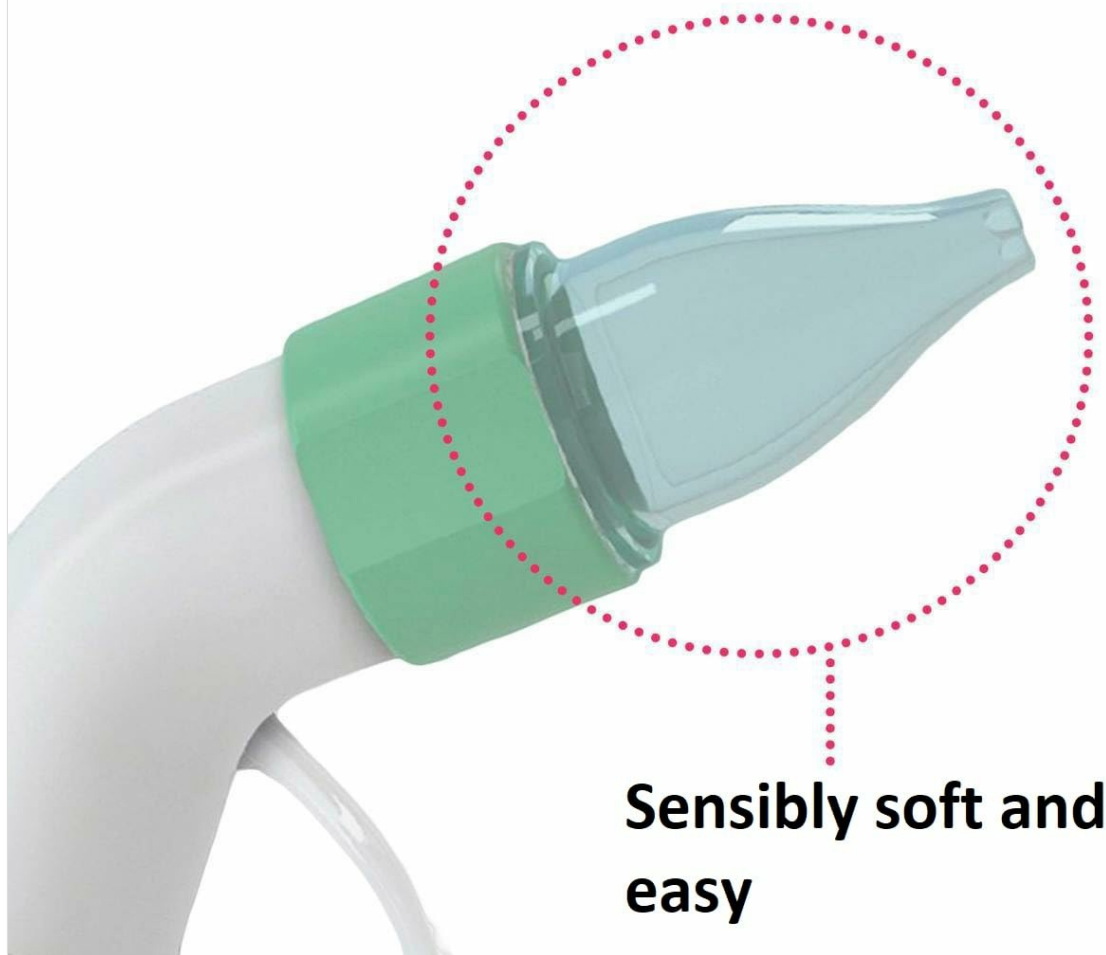


Image 3.1: Close-up of the soft silicone tip, highlighting its gentle material.

## 4. Setup

Before first use and after each subsequent use, ensure all parts are clean and dry.

1. Ensure the hygienic filter is correctly placed inside the collection chamber.
2. Attach the soft silicone nasal tip securely to the collection chamber.
3. Connect the collection chamber assembly to the aspirator body.

## 5. Operating Instructions

Follow these steps for effective and safe nasal aspiration:

1. Position the child comfortably, ideally lying down with their head slightly elevated.
2. Gently insert the soft silicone tip into one of the child's nostrils.
3. Place the mouthpiece of the aspirator in your mouth.
4. Slowly and gently suck through the mouthpiece to create suction. Nasal secretions will be drawn into

the collection chamber.

5. Repeat the process for the other nostril as needed.
6. After use, disassemble and clean all parts immediately.



Image 5.1: Illustration of a parent using the nasal aspirator on an infant.



Image 5.2: A close-up view of the aspirator being used on a baby's nostril.

## 6. Cleaning and Maintenance

Proper cleaning ensures hygiene and extends the life of your aspirator.

- **Disassembly:** Separate all components: mouthpiece, aspirator body, collection chamber, hygienic filter, and silicone tip.
- **Washing:** Wash all parts thoroughly with warm water and mild soap. Ensure all secretions are removed.
- **Rinsing:** Rinse all components under running water to remove any soap residue.
- **Drying:** Allow all parts to air dry completely before reassembling or storing. Ensure no moisture remains to prevent bacterial growth.
- **Filter:** The hygienic filter should be checked regularly and replaced if it appears soiled or damaged.
- **Product Care Instruction:** Hand Wash Only.

## 7. Troubleshooting

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Low or no suction	Improper assembly; blocked tip or filter; insufficient seal.	Ensure all parts are securely connected. Check the silicone tip and hygienic filter for blockages and clean if necessary. Ensure a good seal around the child's nostril.
Secretions not collecting	Filter is saturated or incorrectly placed.	Check the hygienic filter. Replace if saturated or damaged. Ensure it is correctly positioned.
Discomfort for child	Tip inserted too deeply or suction is too strong.	Ensure the tip is gently placed at the entrance of the nostril. Apply gentle, controlled suction. Ensure the child is calm and comfortable.

## 8. Specifications

- **Model Name:** Phisio Clean
- **Item Model Number:** 00004904000000
- **Brand:** Chicco
- **Material:** Silicone (Latex Free)
- **Age Range:** Baby (0m+)
- **Product Dimensions:** 16.54 x 18.9 x 4.72 inches
- **Item Weight:** 2.05 ounces
- **Water Resistance Level:** Waterproof
- **Manufacturer:** Chicco
- **Origin:** Made in Italy



Image 8.1: Indication that the product is Made in Italy.

## 9. Warranty and Support

For warranty information, please refer to the product packaging or the official Chicco website. Chicco products are manufactured to high-quality standards, and specific warranty details may vary by region. For further assistance, product registration, or to explore other Chicco products, please visit the official Chicco website or contact customer support.



# Chicco Digital Catalogue

## Full Range

### Juvenile



**Scan me**

Car seat, walker, highchair,  
bed, boppy etc.

### Nursery



**Scan me**

Cosmetic, feeding  
bottle, soother etc.

### Toys



**Scan me**

Live it.  
*Love it.*



Image 9.1: QR codes linking to Chicco's digital catalogues for various product categories.

[Chicco Juvenile Catalogue](#)

[Chicco Nursery Catalogue](#)

[Chicco Toys Catalogue](#)