



[Manuals.plus](#) /

› [Braun](#) /

› Braun No Touch and Forehead Thermometer - User Manual

Braun BNT300US

Braun No Touch and Forehead Thermometer Instruction Manual

Model: BNT300US | **Brand:** Braun

PRODUCT OVERVIEW

This Braun No Touch and Forehead Thermometer provides clinically accurate and reliable temperature readings. It offers both non-contact and traditional forehead measurement options, making it versatile for adults, babies, toddlers, and children. The device features an on-screen positioning system and guidance light to ensure proper aiming, along with a color-coded display for easy fever interpretation.



Image: Front view of the Braun No Touch and Forehead Thermometer, displaying a temperature reading.

Key Features

- **Non-Contact Measurement:** Take clinically accurate readings from up to 2 inches away, ideal for not disturbing a sleeping child.
- **Easy to Use:** Features an on-screen positioning system and guidance light for accurate aiming.
- **Color-Coded Display:** Temperature readings are displayed with a color-coded system (green for normal, yellow for elevated, red for high) for quick understanding.
- **Versatile Use:** Can be used as both a non-contact and traditional forehead thermometer for the entire family.
- **Pediatrician Recommended:** Braun thermometers are widely recommended by pediatricians for their reliability and accuracy.

BRAUN



FSA/HSA
Eligible

Image: A hand holding the Braun No Touch and Forehead Thermometer, highlighting its eligibility for FSA/HSA.

SETUP AND FIRST USE

Battery Installation

The thermometer requires 2 AA batteries, which are included with the product. To install or replace batteries, slide open the battery compartment cover located on the back of the device. Insert the batteries according to the polarity indicators (+/-) and close the cover securely.



Image: Back view of the thermometer with the battery compartment open, showing two AA batteries inside.

Powering On

Press the power button (marked with a circle and vertical line) located at the bottom of the thermometer to turn on the device. The display will light up, indicating it is ready for use.

OPERATING YOUR THERMOMETER

Taking a Temperature (No-Touch Mode)

Ensure the protective cap is removed from the sensor. Position the thermometer 0-2 cm (0-0.8 inches) from the center of the forehead, between the eyebrows. Press the temperature button (blue button with a thermometer icon). The on-screen positioning system and guidance light will assist in proper aiming. A beep will indicate the measurement is complete, and the temperature will be displayed.



No touch and traditional forehead thermometer in one: same professional accuracy in both modes

Image: A child's temperature being taken using the no-touch mode of the Braun thermometer.

Taking a Temperature (Forehead Touch Mode)

For traditional forehead measurement, gently place the thermometer directly on the forehead. Press the temperature button. The reading will appear on the display after a few seconds.



Image: A mother taking her child's temperature by gently touching the thermometer to the child's forehead.

Color-Coded Fever Guidance

The thermometer's display changes color to provide quick and easy interpretation of the temperature reading:

Color	Meaning
Green	Normal temperature
Yellow	Elevated temperature
Red	High temperature (fever)



Color coded temperature guidance to track progression of fever

Image: The thermometer display showing examples of green, yellow, and red color codes for temperature guidance.

Mute Function

The thermometer features a mute function, allowing you to take readings without sound. This is particularly useful for not disturbing a sleeping child. Locate the volume switch on the side of the device to toggle sound on or off.



Will not wake a sleeping baby

Image: A mother taking a sleeping baby's temperature with the Braun thermometer, demonstrating its silent operation.

Temperature Unit Selection

The device can display temperatures in Fahrenheit (°F) or Celsius (°C). Refer to the full user guide for specific instructions on how to switch between these units.

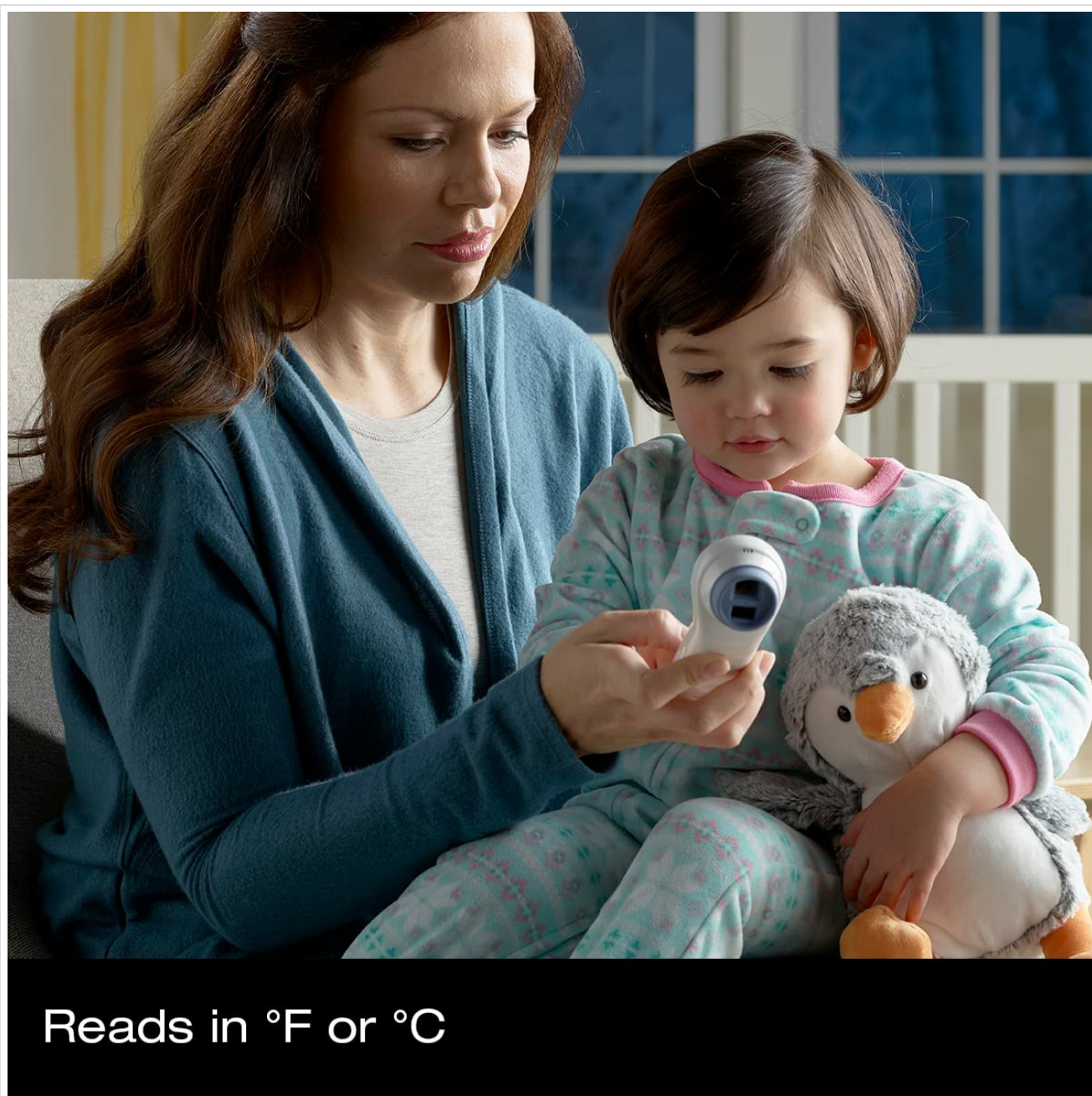


Image: The thermometer display showing a temperature reading in Fahrenheit.

Official Product Videos

There are no official seller-created videos available for this product in the provided data.

CARE AND MAINTENANCE

Cleaning

To clean the thermometer, wipe the probe and body with an alcohol swab or a cotton swab moistened with alcohol (70% Isopropyl). Do not immerse the thermometer in water or any other cleaning solutions. Ensure the thermometer is completely dry before storing.

Storage

Store the thermometer in a dry place, away from direct sunlight, dust, and extreme temperatures. Keep it out of reach of children.

TROUBLESHOOTING

Common Issues

- **"Err" message:** This usually indicates an improper measurement. Ensure the thermometer is positioned correctly and held steady during measurement. Make sure the sensor is clean.
- **Inaccurate readings:** Ensure the sensor is clean and free from debris. Avoid taking measurements immediately after bathing, exercising, or in environments with significant temperature changes. Allow the thermometer to acclimate to room temperature for at least 30 minutes if it has been stored in a very hot or cold environment.
- **No power:** Check battery installation and replace batteries if necessary. Ensure the battery compartment is securely closed.

For detailed troubleshooting steps and error codes, please refer to the complete User Guide.

PRODUCT SPECIFICATIONS

Specification	Detail
Model Number	BNT300US
Product Dimensions	5.32 x 2.76 x 8.07 inches; 7.68 ounces
Batteries	2 AA batteries (included)
Manufacturer Recommended Age	12 months and up
Included Components	Braun Digital No-Touch Forehead Thermometer
Special Feature	Non-Contact
Display Type	Digital
Material	Plastic

Product dimensions (inches)

1.8" length
1.8" width
5.9" height



Image: The Braun thermometer with its key dimensions (length, width, height) indicated.

WARRANTY AND SUPPORT

For comprehensive warranty information, product registration, and customer support, please visit the official Braun Healthcare website or refer to the full User Guide provided with your product.

[Download Full User Guide \(PDF\)](#)