

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Mega Motion](#) /

› [Easy Comfort Superior 3 Position Heavy Duty Big Lift Chair 500 lb Capacity Chaise Lounge Recliner - Ocean BI User Manual](#)

Mega Motion MM-AS-7001-Ocean

Mega Motion Superior Lift Chair User Manual

Model: MM-AS-7001-Ocean | Brand: Mega Motion

IMPORTANT SAFETY INFORMATION

MAXIMUM WEIGHT CAPACITY: 500 lbs. DO NOT EXCEED.

WARNING: This chair has a power recline and lift mechanism. Always keep hands, fingers, and other body parts away from moving parts when operating the chair. Children should be supervised at all times when using this chair.

PRODUCT OVERVIEW

The Mega Motion Superior Lift Chair is a versatile three-position recliner designed for comfort and assistance. It offers closed, partial recline, and full recline positions, along with a power lift function to aid in standing. Featuring a split back design and a chaise seat, it provides comprehensive support. Convenient side pockets are integrated for storage of personal items or the remote control. The chair is equipped with a heavy-duty steel lift mechanism, supporting up to 500 pounds, and includes a built-in battery backup system.

Key Features

- Three-position recline: Closed, Partial Recline, Full Recline
- Power lift mechanism for easy standing assistance
- Split back design for enhanced comfort
- Chaise seat for head-to-toe support
- Convenient side storage pockets
- Heavy-duty steel frame with 500 lb. weight capacity
- Integrated USB charging port
- Built-in battery backup (requires two 9-volt batteries)



This image shows the Mega Motion Superior Lift Chair in its Ocean Blue color, positioned slightly lifted from the floor, with the remote control visible on the side.

COMPONENTS

Your Mega Motion Superior Lift Chair package includes the following main components:

- Lift Chair Unit
- Remote Control

- Power Cord

SETUP INSTRUCTIONS

While specific assembly steps may vary, general setup involves connecting the chair's back to the base, attaching the power cord, and ensuring all connections are secure. Refer to the detailed assembly guide included in your packaging for precise instructions.

1. Carefully remove all components from the packaging.
2. Position the chair base in your desired location, ensuring adequate space for recline and lift functions.
3. Connect the backrest to the chair base according to the alignment guides. Ensure it clicks securely into place.
4. Connect the remote control cable to the corresponding port on the chair.
5. Plug the power cord into the chair's power input, then into a grounded electrical outlet.
6. Install two 9-volt batteries into the battery backup compartment (usually located on the power transformer) to enable emergency operation during a power outage.
7. Test all functions of the chair using the remote control.



This diagram illustrates the key dimensions of the lift chair, including height, depth, and width, to help with placement and fit.

OPERATING INSTRUCTIONS

Remote Control

The lift chair is operated using a wired remote control. The remote typically features two main buttons or a rocker switch:

- **UP Button:** Press and hold to raise the chair to the lift position, assisting you in standing.
- **DOWN Button:** Press and hold to recline the chair through partial and full recline positions, and to return it to the upright, closed position.

Release the button at any point to stop the chair in your desired position.

Chair Positions

- **Closed/Upright Position:** The standard seating position.
- **Partial Recline:** A comfortable, slightly reclined position for relaxation.
- **Full Recline:** A fully extended position, ideal for napping or extended relaxation.
- **Lift Position:** The chair tilts forward, raising the seat to assist the user in standing up.



The lift chair is shown in its standard, upright seating position, ready for use.



This view demonstrates the lift chair in a fully reclined position, with the leg rest extended for full body support.



The chair is shown in its maximum lift position, tilted forward to assist a user in standing up from the chair.

MAINTENANCE

Cleaning

- For general cleaning, use a soft, damp cloth with mild soap and water.
- Avoid harsh chemicals or abrasive cleaners, as they may damage the fabric.
- Spot clean spills immediately to prevent staining.
- Vacuum regularly to remove dust and debris from the fabric.

General Care

- Ensure the power cord is not pinched or tangled.
- Periodically check all connections to ensure they are secure.
- Keep the area around the chair clear to prevent obstruction of moving parts.
- Replace 9-volt batteries in the backup system annually or as needed to ensure functionality during power outages.

TROUBLESHOOTING

If you encounter issues with your lift chair, refer to the table below for common problems and their solutions.

Problem	Possible Cause	Solution
Chair does not move.	Power cord unplugged or loose connection.	Check that the power cord is securely plugged into the wall outlet and the chair.
Remote control is unresponsive.	Remote cable disconnected or damaged.	Ensure the remote control cable is firmly connected to the chair. Inspect for visible damage.
Chair only works intermittently.	Loose electrical connection or motor issue.	Check all power connections. If problem persists, contact customer support.
Battery backup not functioning during power outage.	Dead or missing 9-volt batteries.	Replace the two 9-volt batteries in the battery backup compartment.
Unusual noises during operation.	Mechanism requires lubrication or has a loose component.	Do not attempt to lubricate or repair yourself. Contact customer support for assistance.

If the problem persists after attempting these solutions, please contact Mega Motion customer support.

PRODUCT SPECIFICATIONS

Feature	Detail
Brand	Mega Motion
Model Number	MM-AS-7001-Ocean
Color	Blue
Product Dimensions (D x W x H)	36"D x 39.5"W x 43.5"H
Back Style	Split Back
Special Feature	USB Port
Unit Count	1 Count
Recommended Uses	Lounge
Maximum Weight Recommendation	500 Pounds
Pattern	Solid
Age Range (Description)	Adult
Included Components	Ottoman (as per source data)
Arm Style	Armless
Surface Recommendation	Hard floor
Indoor/Outdoor Usage	Indoor
Seat Height	20 Inches
Fill Material	Polyester
Leg Style	Straight Leg
Reclining Position Count	3
UPC	627554159260
Form Factor	Recliner
Item Weight	143 Pounds
Date First Available	January 1, 2010

WARRANTY AND SUPPORT

Warranty Information

The Mega Motion Superior Lift Chair is backed by a **1-year in-home service warranty**. This warranty covers manufacturing defects and issues arising from normal use within the specified period. Please retain your proof of purchase for warranty claims.


Battery Backup


The chair features a built-in battery backup system designed to allow limited operation during a power outage. This system requires two 9-volt batteries (not included) to function. Ensure these batteries are installed and replaced periodically to maintain emergency functionality.

Customer Support

For technical assistance, warranty claims, or to order replacement parts, please contact Mega Motion customer support. Refer to the contact information provided in your product packaging or visit the official Mega Motion website.

Related Documents - MM-AS-7001-Ocean

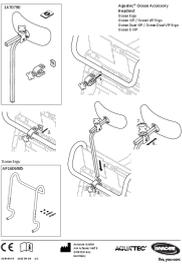	<p>Mega Motion iNFI ATTAIN MM-3603 Infinite Position Power Lift Chair</p> <p>Discover the Mega Motion iNFI ATTAIN MM-3603, a dual-motor infinite position layflat power recline and lift chair. Designed for comfort and mobility, it features USB charging, ergonomic support, and is an FDA Class II Medical Device.</p>
	<p>Mega Motion Trio Grande MM-5001 Power Recline/Lift Chair</p> <p>Discover the Mega Motion Trio Grande MM-5001, a single motor 3-position lift and recline chair featuring breathable superior soft fabric, USB charging port, and whisper-quiet operation for enhanced comfort and mobility.</p>
	<p>Aquatec Ocean Accessory Headrest Mounting Instructions</p> <p>Mounting guide for the Aquatec Ocean Accessory Headrest, compatible with various Ocean shower chair models including Ergo, VIP, Dual VIP, and E-VIP. Provides step-by-step assembly instructions and product details.</p>

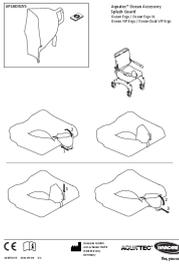


[Aquatec Ocean Splash Guard Mounting Instructions - Invacare](#)

Official mounting instructions for the Invacare Aquatec Ocean Splash Guard accessory, designed for use with the Aquatec Ocean Ergo, Ocean Ergo XL, Ocean VIP Ergo, and Ocean Dual VIP Ergo shower chairs.



[Invacare Aquatec Ocean Ergo Accessories and Adjustments Guide](#)

Explore the Invacare Aquatec Ocean Ergo Family of shower chairs and commodes. This guide details available accessories, customization options, and compatibility for enhanced patient positioning, safety, and toileting.

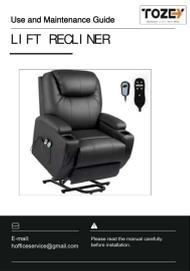

[TOZLY Lift Recliner Use and Maintenance Guide](#)

A comprehensive guide for the use and maintenance of the TOZLY Lift Recliner, covering introduction, assembly, installation, regular maintenance, and warranty information.