

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [NRS Healthcare](#) /
- › [NRS Healthcare Nuvo Corner Shower Stool \(Model N05382\) User Manual](#)

NRS Healthcare N05382

NRS Healthcare Nuvo Corner Shower Stool User Manual

Model: N05382

INTRODUCTION

The NRS Healthcare Nuvo Corner Shower Stool (Model N05382) is designed to provide support and comfort during showering, especially in smaller shower areas. Its corner-fitting design maximizes available space, while the padded, height-adjustable seat enhances user safety and independence. This manual provides essential information for the safe assembly, operation, and maintenance of your shower stool.

SAFETY INFORMATION

Please read all instructions carefully before using the shower stool. Failure to follow these instructions may result in injury or damage to the product.

- **Maximum User Weight:** Do not exceed the maximum weight limit of 127 kg (20 stone).
- **Stability:** Ensure all legs are firmly on the ground and the stool is stable before use. The slightly splayed frame and non-slip feet are designed for stability.
- **Regular Checks:** Periodically inspect the stool for any signs of wear, damage, or loose components. Do not use if damaged.
- **Cleaning:** Clean the stool regularly as described in the Maintenance section to prevent build-up and ensure hygiene.
- **Adjustments:** Only adjust the height when the stool is not in use and ensure all adjustment pins are securely engaged before sitting.
- **Placement:** Position the stool securely in the corner of the shower area.

PACKAGE CONTENTS

The package for the NRS Healthcare Nuvo Corner Shower Stool (N05382) typically includes:

- 1 x Nuvo Corner Shower Stool frame (pre-assembled)
- 1 x Padded seat

SETUP AND ASSEMBLY

The NRS Healthcare Nuvo Corner Shower Stool comes largely pre-assembled, requiring minimal setup.

1. **Unpack the Stool:** Carefully remove all components from the packaging.
2. **Attach the Seat:** The padded seat is designed to clip onto the stool frame. Align the clips on the underside of the seat with the corresponding points on the top of the stool frame and press down firmly until it clicks into place.
3. **Verify Stability:** Before initial use, gently rock the stool to ensure it is stable and all parts are securely connected.

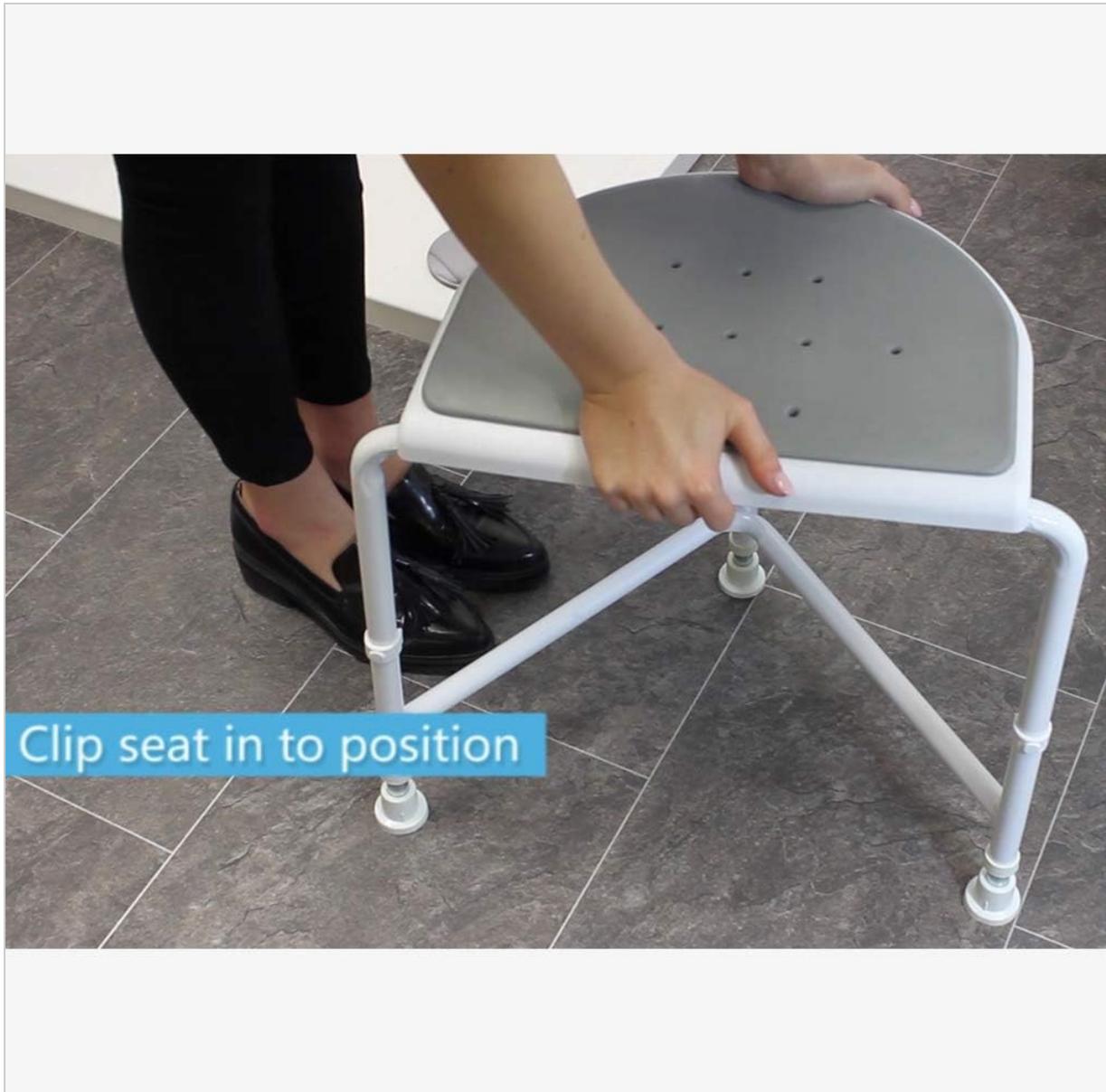


Image: Attaching the padded seat to the stool frame. The seat clips securely into position.

OPERATING INSTRUCTIONS

Using the Shower Stool

1. **Placement:** Place the Nuvo Corner Shower Stool in the corner of your shower area. Ensure the non-slip feet are firmly on the shower floor.
2. **Sitting:** Carefully approach the stool and sit down, ensuring your weight is evenly distributed. The padded seat provides comfort, and the drainage holes prevent water accumulation.
3. **Standing:** When finished, use the stool for support as you carefully stand up.



Image: The Nuvo Corner Shower Stool positioned in a shower corner.



Image: The shower stool in use, providing a comfortable and safe showering experience.

Adjusting the Height

The stool's height can be adjusted to suit individual user requirements. The seat height is adjustable between 45 cm and 60 cm (17.75" – 23.5").

1. **Empty the Stool:** Ensure no one is sitting on the stool before attempting to adjust its height.
2. **Locate Adjustment Pins:** Each leg has a spring-loaded adjustment pin. These are concealed to reduce water access.
3. **Adjust Each Leg:** Press the adjustment pin on one leg and slide the inner leg tube up or down to the desired height. Ensure the pin fully re-engages into one of the adjustment holes. Repeat this process for all three legs, ensuring they are all set to the same height to maintain stability.
4. **Verify Stability:** After adjustment, ensure all pins are securely locked in place and the stool is level and stable before use.



Image: Adjusting the height of the stool legs. Ensure the adjustment pin is fully engaged.

CLEANING AND MAINTENANCE

Regular cleaning and maintenance will prolong the life of your shower stool and ensure hygienic use.

- **Cleaning the Seat:** The padded seat can be easily clipped off for thorough cleaning. Use a mild detergent and warm water. Rinse thoroughly and allow to air dry.
- **Cleaning the Frame:** Wipe the white powder-coated steel frame with a damp cloth and mild detergent. Avoid abrasive cleaners or scouring pads, which can damage the finish.
- **Drying:** After cleaning, ensure all parts are thoroughly dried to prevent water spots or potential corrosion.
- **Inspection:** Regularly check the non-slip feet for wear or damage. Replace if necessary to maintain stability. Inspect the frame for any signs of rust or structural compromise.

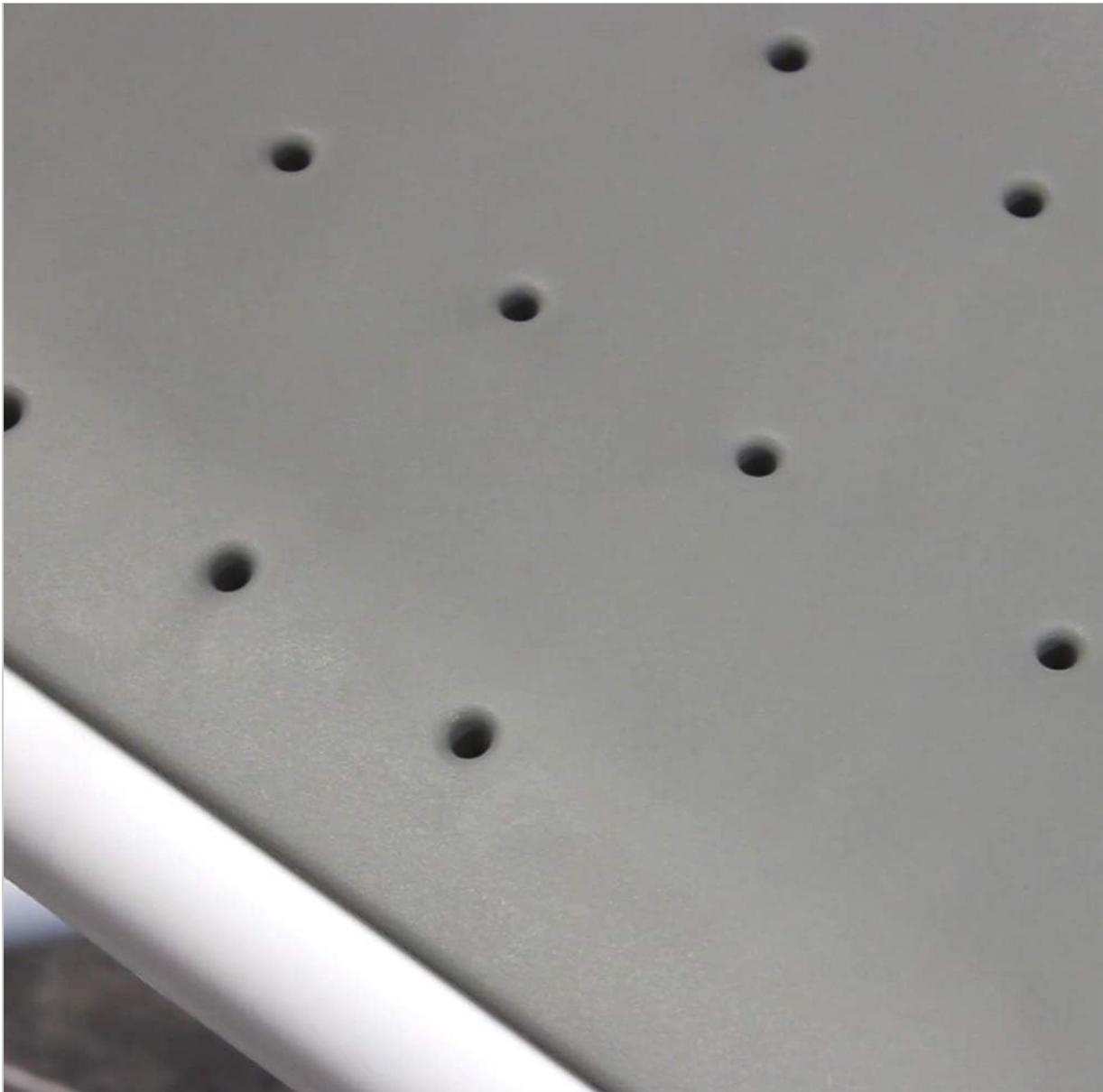


Image: Close-up of the padded seat, showing drainage holes for hygiene and comfort.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Stool feels unstable or wobbly.	Legs not adjusted to the same height. Adjustment pins not fully engaged. Non-slip feet worn or damaged. Uneven shower floor.	Ensure all legs are set to the identical height. Check that all adjustment pins are securely locked into their holes. Inspect and replace worn non-slip feet. Reposition the stool to a more even surface if possible.
Seat does not stay attached.	Seat not clipped on correctly. Clips on seat or frame are damaged.	Press down firmly on the seat until it clicks into place. Inspect clips for damage. If damaged, contact customer support for replacement parts or advice.

Problem	Possible Cause	Solution
Difficulty adjusting leg height.	Pin stuck or mechanism dirty.	Ensure the pin is fully depressed. Clean around the pin and leg tube to remove any debris.

SPECIFICATIONS

Model Number	N05382
Brand	NRS Healthcare
Material	Steel frame, Polyurethane padded seat
Color	White frame, Grey seat
Seat Height Adjustment	450 - 600 mm (17.75 - 23.5 inches)
Seat Width	405 mm (16 inches)
Seat Depth	304 mm (12 inches)
Maximum Footprint (W x D)	560 mm x 390 mm (22 x 15.25 inches)
Product Dimensions (D x W x H)	22.05"D x 21.73"W x 20.24"H (approximate, based on max height)
Item Weight	6.57 Pounds (approximate)
Maximum User Weight Limit	127 kg (20 stone)

WARRANTY AND SUPPORT

For warranty information, replacement parts, or further assistance, please contact NRS Healthcare customer support. Refer to your purchase documentation or the official NRS Healthcare website for contact details. When contacting support, please have your model number (N05382) and date of purchase available.