

## Vive LVA1011

# Vive Raised Toilet Seat Instruction Manual

Model: LVA1011

## 1. PRODUCT OVERVIEW

The Vive Raised Toilet Seat is designed to enhance bathroom safety and accessibility by adding 5 inches of height to your toilet. This elevated riser, complete with comfortable padded handles, makes it significantly easier and safer for individuals with limited mobility, such as the elderly, those recovering from surgery (hip, knee, or back), or handicapped individuals, to sit down and stand up from the toilet without excessive bending or strain.



*Image: The Vive Raised Toilet Seat installed on a standard toilet, clearly indicating the added 5 inches of height.*

## 2. KEY FEATURES

- **Elevated Height:** Adds five inches of height to any standard or elongated toilet bowl, reducing the effort required to sit or stand.
- **Contoured Design:** The seat is ergonomically contoured for safe and comfortable elevated seating.
- **Padded Handles:** Features removable handles with soft foam padding, providing a secure, non-slip grip for added stability.
- **Tool-Free Installation:** Easily installs on most standard or elongated toilet bowls without the need for special tools.
- **Secure Attachment:** Utilizes an adjustable knob and rear wings for a stable and secure fit to the toilet.
- **Durable & Portable:** Constructed from lightweight, durable composite material that supports up to 250 pounds. It is also easy to transport.
- **Easy to Clean:** Can be wiped down with regular household disinfectants or sanitizing wipes, and does not interfere with normal toilet cleaning.

# Easily Sit and Stand

Effortlessly maneuver with 5 extra inches



*Image: A user demonstrating the ease of sitting and standing from the toilet with the added height of the Vive Raised Toilet Seat.*

# Perfect For:

## Elderly, Post-Surgery, Disabled



*Image: A close-up view highlighting the comfortable, padded handles that provide additional support.*

### 3. SETUP

#### 3.1. Unpacking and Component Check

Carefully unpack all components from the packaging. Ensure you have the main raised toilet seat unit and the two removable padded armrests. Inspect all parts for any signs of damage before proceeding with assembly.

# Contoured Seat

## Premium comfort and support



*Image: The main unit of the Vive Raised Toilet Seat with its contoured design and attachment points for handles.*

### 3.2. Attaching the Padded Handles

The padded handles are designed for easy, tool-free attachment. Locate the metal sockets on either side of the raised toilet seat. Depress the brass locking pins on the armrests and slide the armrests into the sockets until the locking pins pop back into place, securing the handles firmly.

*Video: This video demonstrates how to attach the padded handles to the Vive Raised Toilet Seat. It shows the process of depressing the locking pins and sliding the handles into place until they click securely.*

*Video Description: Official Vive Health Direct video demonstrating the simple, tool-free attachment of the padded armrests to the raised toilet seat.*

### 3.3. Installing on the Toilet Bowl

Before installation, remove your existing toilet seat and lid. Position the Vive Raised Toilet Seat over the toilet bowl, ensuring the rear lip of the raised seat fits securely underneath the rear lip of the toilet bowl. Once aligned, securely tighten the adjustable knob at the front of the raised seat. This knob clamps the seat firmly onto the toilet bowl, preventing movement. Verify that the seat is stable and securely fastened before use.

# Easy to Install

Installs in minutes with no tools necessary



*Image: Illustrates the easy, tool-free installation process of the raised toilet seat onto a standard toilet.*

*Video: This video provides a comprehensive overview of the Vive Raised Toilet Seat, including its features and how it attaches to a toilet.*

*Video Description: Official Vive Health Direct video showcasing the features and installation of the Vive Raised Toilet Seat.*

## 4. OPERATING INSTRUCTIONS

Once securely installed, the Vive Raised Toilet Seat is ready for use. The added height and padded handles provide stable support for sitting and standing. Always use the padded handles for support when lowering yourself onto or rising from the seat to maximize stability and reduce strain on joints.



# Sturdy & Stable

## Supports up to 250 pounds



Image: A user demonstrating the stability and support provided by the Vive Raised Toilet Seat during use.

## 5. MAINTENANCE

### 5.1. Cleaning

The Vive Raised Toilet Seat is designed for easy cleaning. It can be wiped down with common household disinfectants or sanitizing wipes. For a more thorough cleaning, the seat can be easily removed from the toilet. The padded handles are also removable, allowing for comprehensive cleaning of all surfaces. Ensure all parts are dry before re-installation.

### 5.2. Storage

When not in use or for travel, the lightweight design and removable handles make the seat easy to store. Store in a clean, dry place away from direct sunlight or extreme temperatures.

## 6. TROUBLESHOOTING

- **Seat feels unstable:** Ensure the adjustable knob at the front is tightened securely. Check that the rear lip of the raised seat is properly seated under the toilet bowl's rear lip.
- **Handles are loose:** Verify that the brass locking pins on the armrests are fully engaged in the metal sockets. Press down firmly on the handles until the pins click into place.

- **Difficulty installing:** Confirm that your toilet is either a standard or elongated size. The seat is designed to fit most common toilet types. Ensure your existing toilet seat and lid have been completely removed.
- **Product deformation/cracking:** While durable, continuous heavy use or exceeding the weight capacity (250 lbs) may lead to wear. Ensure the product is used within its specified weight limit.

## 7. SPECIFICATIONS

Brand	Vive
Model Number	LVA1011
Color	White
Material	Plastic
Style	Universal Size w/ Handles
Shape	Oval
Product Dimensions	18"L x 16.5"W
Item Weight	7.05 Pounds
Weight Capacity	Up to 250 lbs
Operation Mode	Manual
UPC	818323026245



# Comfortable Padded Handles


Soft covered handles provide a secure grip




*Image: Detailed diagram illustrating the dimensions of the Vive Raised Toilet Seat for compatibility reference.*

# Fits Most Toilets

## Standard Round

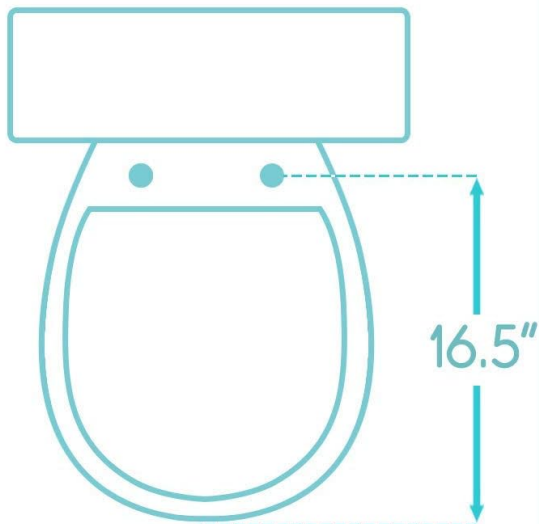

## Elongated

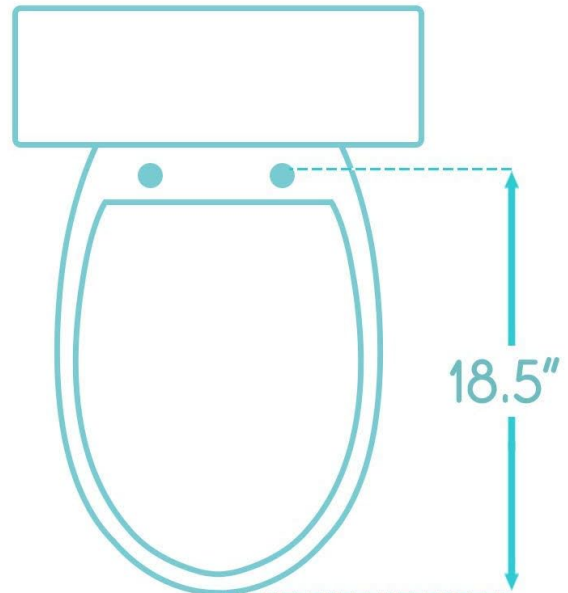

*Image: Diagram showing the compatibility of the Vive Raised Toilet Seat with both standard round and elongated toilet types.*

## 8. WARRANTY & SUPPORT

The Vive Raised Toilet Seat comes with a **Lifetime Guarantee**, ensuring confidence in your purchase. For any questions, concerns, or support needs, please refer to the official Vive Health website or contact their customer service directly.

You can also find additional information and a downloadable user manual (PDF) at the following link [Vive Raised Toilet Seat User Manual](#)