

Invacare T94HCP

Invacare T94HCP Wheelchair Elevating Legrests User Manual

Model: T94HCP

1. INTRODUCTION

This manual provides essential information for the safe and effective use of your Invacare T94HCP Wheelchair Elevating Legrests. These legrests are designed to enhance comfort and support for wheelchair users by allowing leg elevation and providing padded calf support. They feature a swing-away function for easy transfers and tool-free installation for convenience.

2. IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before installation and use. Failure to follow these instructions may result in injury or damage to the product.

- **Compatibility:** These legrests are specifically designed for Invacare wheelchairs. While they may fit other brands, compatibility is not guaranteed. Always verify fitment before use.
- **Secure Attachment:** Ensure legrests are securely attached and locked into place before use. Regularly check all connections for tightness.
- **Weight Capacity:** Do not exceed the weight capacity of your wheelchair or the legrests.
- **Transfers:** Always engage wheelchair brakes and ensure legrests are swung away or removed before attempting transfers to prevent obstruction.
- **Proposition 65 Warning:** This product contains chemicals known to the State of California to cause cancer and birth defects or other reproductive harm.

3. WHAT'S IN THE BOX

Your package should contain the following items:

- Invacare T94HCP Elevating Legrest (1 pair)

4. SETUP AND INSTALLATION

The Invacare T94HCP legrests are designed for tool-free installation, making attachment and removal

straightforward.

1. **Prepare Wheelchair:** Ensure your wheelchair is on a flat, stable surface and the brakes are engaged.
2. **Align Brackets:** Position the legrest's hinge brackets with the corresponding attachment points on your wheelchair's frame. These legrests feature a 1.5-inch spacing between swing-away hinge brackets, specifically for Invacare wheelchairs.
3. **Attach Legrest:** Insert the legrest's pins into the attachment points. The front rigging design allows for easy attachment.
4. **Secure Lock:** Push the legrest firmly until it clicks into place, indicating it is securely locked. Verify that the legrest is stable and does not wobble.
5. **Repeat for Second Legrest:** Follow the same steps to attach the second legrest.

Comfortable and Convenient




With a 1.5-inch spacing between swing away hinge brackets, **these elevating legrests are designed specifically for Invacare wheelchairs.**



Image: Close-up of the legrest hinge bracket, highlighting the 1.5-inch spacing for compatibility with Invacare wheelchairs.

Transfer With Ease

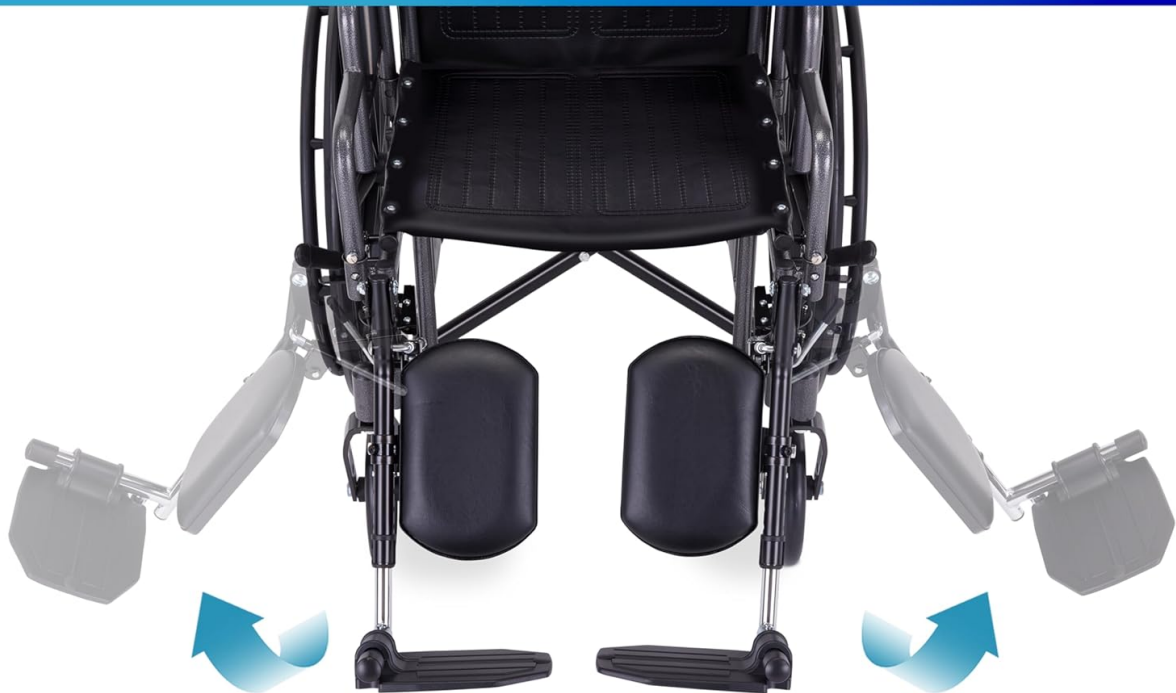

Legrests swing away
clearing the foot space
for easy transfers



Image: Two images showing hands attaching and removing the legrest from a wheelchair, emphasizing the tool-free installation and user-friendly bracket for secure locking.

5. OPERATING INSTRUCTIONS

The Invacare T94HCP legrests offer adjustable features for optimal user comfort and functionality.

5.1 Elevating the Legrests

To elevate or lower the legrests:

1. Locate the adjustment lever on the side of the legrest.
2. Press or pull the lever (depending on design) to release the locking mechanism.
3. Gently raise or lower the legrest to the desired height. The legrests offer 4 adjustment levels.
4. Release the lever to lock the legrest securely in the chosen position. Ensure it is fully engaged before resting your legs.

Durable and Reliable

Easily elevate and lower the legrest to the level of your choice



Elevating the legs can aid in blood circulation

Comfortable calf pads add support



Image: An illustration of the elevating function, showing a person's leg raised on the legrest. Elevating legs can aid in blood circulation, and comfortable calf pads add support.

Adjustable Heights



Image: A diagram showing the legrest at its minimum (16 inches) and extended (19.75 inches) heights, highlighting the 4 adjustment levels and 7.75-inch wide footplate.

5.2 Using the Swing-Away Function

The swing-away feature allows for easy transfers and maneuverability.

1. To swing the legrest away, locate the release lever or button near the attachment point.
2. Activate the release mechanism and pivot the legrest outwards, away from the front of the wheelchair.
3. To return the legrest to its forward position, swing it back until it locks securely into place.

Beneficial and Supportive



Perfectly sized
and non-slip



Sturdy composite
material



Flips up and out of
the way easily



Image: An image demonstrating the swing-away feature of the legrests, showing them pivoted outwards to clear space for easier user transfers to and from the wheelchair.

5.3 Padded Calf Support and Footplate

- The padded calf supports provide comfort and reduce pressure during extended use.
- The composite footplates are designed to accommodate various foot sizes and offer a stable, non-slip surface.
- The footplates can be flipped up and down for enhanced maneuverability or to adjust positioning.

6. MAINTENANCE AND CARE

Regular maintenance ensures the longevity and safe operation of your legrests.

- **Cleaning:** Wipe down the legrests and padded calf supports with a damp cloth and mild soap. Avoid abrasive cleaners or solvents.
- **Inspection:** Periodically inspect all components, including the attachment brackets, adjustment mechanisms, and padding, for any signs of wear, damage, or loose parts.
- **Lubrication:** If any moving parts become stiff, a small amount of silicone-based lubricant can be applied to the pivot points.
- **Storage:** Store the legrests in a clean, dry place when not in use.

7. TROUBLESHOOTING

If you encounter issues with your Invacare T94HCP legrests, refer to the following common solutions:

- **Legrests do not attach securely:** Ensure the legrests are fully aligned with the wheelchair's attachment points and pushed firmly until the locking mechanism engages. Verify your wheelchair is an Invacare model or compatible.
- **Legrests are difficult to elevate/lower or swing away:** Check for any obstructions. Apply a small amount of silicone lubricant to the moving parts if they feel stiff. Ensure the release lever is fully activated.
- **Footplate does not stay in position:** Ensure the footplate is fully flipped down and locked into its operational position.

If problems persist, contact Invacare customer support for assistance.

8. SPECIFICATIONS







Feature	Detail
Model Number	T94HCP
Product Dimensions	23.5 x 7.5 x 23.5 inches
Item Weight	8.5 Pounds
Manufacturer	Invacare
Calf Pads	Padded
Footplate Material	Composite
Installation	Tool-Free

9. WARRANTY AND SUPPORT

For specific warranty information regarding your Invacare T94HCP Wheelchair Elevating Legrests, please refer to the documentation included with your purchase or visit the official Invacare website.

If you require technical assistance, have questions about product functionality, or need to report a defect, please contact Invacare Customer Support:

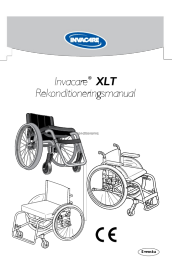- **Invacare Website:** www.invacare.com
- **Customer Service:** Refer to the contact section on the Invacare website for phone numbers and email support.

	<p>Invacare Kuschall Champion Folding Wheelchair - Configuration and Specifications</p> <p>Detailed product information, configuration options, and specifications for the Invacare Kuschall Champion folding wheelchair. Explore frame types, seating adjustments, backrest options, wheel choices, and accessories to customize your mobility solution.</p>
	<p>Invacare Kuschall Champion 2.0 Manual Active Wheelchair - Configuration and Price List</p> <p>Comprehensive configuration guide and price list for the Invacare Kuschall Champion 2.0 manual active wheelchair. Explore a wide range of customizable options for frame, seating, backrests, wheels, accessories, and more.</p>
	<p>Invacare XLT Wheelchair Reconditioning Manual Service & Maintenance Guide</p> <p>Comprehensive reconditioning and maintenance manual for the Invacare XLT manual wheelchair. Covers washing, disinfection, reconditioning procedures, parts, and checklists.</p>
	<p>Invacare Tracer IV Wheelchair Order and Configuration Guide</p> <p>Detailed information on configuring and ordering the Invacare Tracer IV Wheelchair, including frame, seat, arm, back, color, caster, wheel, handrim, and accessory options. Features stock models, pre-configured options, and detailed specifications.</p>
	<p>Invacare Storm 4 X-plore Modulite Powered Wheelchair - Configuration and Price List</p> <p>Detailed configuration options, features, and pricing for the Invacare Storm 4 X-plore Modulite powered wheelchair. Explore seating, drive systems, accessories, and more.</p>
	<p>Invacare Action NG Wheelchair - Features, Specifications, and Options</p> <p>Detailed information on the Invacare Action NG wheelchair, including features, product characteristics, available options, technical specifications, and pricing. This guide covers various configurations, seating options, armrests, footrests, and accessories to customize the wheelchair.</p>