



Manuals.plus /

- › GLO Science /
- › GLO Science Teeth Whitening Gel Vials with Lip Care - Instruction Manual

GLO Science 259120

GLO Science Teeth Whitening Gel Vials with Lip Care - Instruction Manual

Model: 259120

1. PRODUCT OVERVIEW

The GLO Science Teeth Whitening Gel Vials provide a convenient method for teeth whitening. This product includes 6% hydrogen peroxide whitening gel in individual vials and a lip care product. It is designed for use at home or while traveling.



Image 1.1: GLO Science Teeth Whitening Gel Vials and Lip Care.

Key Features:

- **6% Hydrogen Peroxide Formula:** Designed for effective whitening.
- **Sensitivity-Free:** Formulated to minimize tooth and gum sensitivity.
- **Safe for Dental Work:** Can be used on crowns, bridges, and veneers.
- **Portable Vials:** Convenient for on-the-go application.
- **Mint Flavor:** Provides a refreshing taste.
- **Includes Lip Care:** A hydrating lip care product is provided to moisturize and protect lips.

Powerful Whitening, Gentle Formula

6%

Hydrogen Peroxide.
Safer for teeth & gums



No Sensitivity



Vegan &
Cruelty Free



Safe For Crowns,
Bridges & Veneers



With Real
Peppermint Oil



Free From Gluten,
SLS, Latex, Alcohol,
Carbamide Peroxide



Image 1.2: Key features of the GLO Science whitening gel formula.

2. SAFETY INFORMATION

Please read and understand all safety information before using this product.

- Follow device instructions if using with a GLO whitening kit or Advanced GLO BRILLIANT White Smile Teeth Whitening Device.
- Avoid contact with eyes. In case of contact, rinse thoroughly with water.
- This product contains hydrogen peroxide and may cause temporary sensitivity or irritation to teeth and gums in some individuals.
- Do not swallow the gel.
- Rinse mouth thoroughly after use.
- Keep out of reach of children.
- Consult a dentist before use if you have any pre-existing dental conditions, gum disease, or are pregnant or nursing.

3. INGREDIENTS

Teeth Whitening GLO Vial Ingredients:

Distilled Water, Polyethylene Glycol, Hydrogen Peroxide (6%), Potassium Nitrate, Glycerin, Edetate Disodium, Sodium Saccharin, Sodium Fluoride, Sodium Hydroxide, Flavor.

GLO Lip Care Ingredients:

Ricinus Communis (Castor) Seed Oil, Beeswax, Mentha Piperita (Peppermint) Leaf, Stearalkonium Hectorite, Butyrospermum Parkii (Shea) Butter*, Vanilla Planifolia Fruit Extract, Menthol, Alcohol, Aloe Barbadensis Leaf Juice*, Copernicia Cerifera (Carnauba) Wax, Tocopherol, Ammonium Glycyrrhizate, Ubiquinone. *Organic

Includes Nourishing Lip Care

Full of with nourishing ingredients

- **CoQ10** - An essential nutrient which provides anti-oxidant properties. Neutralizes free radicals, with anti-ageing properties
- **Shea butter** - Softens and hydrates dry damaged lips
- **Vitamin E** a powerful antioxidant. Ours is derived from non-GMO sunflower oil for a non-allergenic wholesome ingredient



Image 3.1: GLO Lip Care and its key nourishing ingredients.

4. SETUP AND OPERATING INSTRUCTIONS

4.1 Preparation:

1. Ensure your teeth are clean. Brushing before application is recommended.
2. Apply the provided GLO Lip Care to your lips to moisturize and protect them during the whitening process.



Image 4.1: Applying GLO Lip Care before whitening.

4.2 Application of Whitening Gel:

3. Twist off the cap of a GLO Whitening Gel Vial.
4. Squeeze a small amount of gel onto the precision brush tip.
5. Apply a thin, even layer of the gel directly onto your teeth. Ensure all visible tooth surfaces are covered.
6. The gel is designed to evaporate as it whitens, so no rinsing is required immediately after application unless using with a device that specifies rinsing.
7. If using with a GLO whitening kit and Advanced GLO BRILLIANT White Smile Teeth Whitening Device, follow the specific instructions provided with that device for activation.

Bright Smiles, Delivered.

Brighten and whiten teeth with our clinically proven, power packed GLO vials. Specially-formulated hydrogen peroxide gel whitens without causing pain or sensitivity.

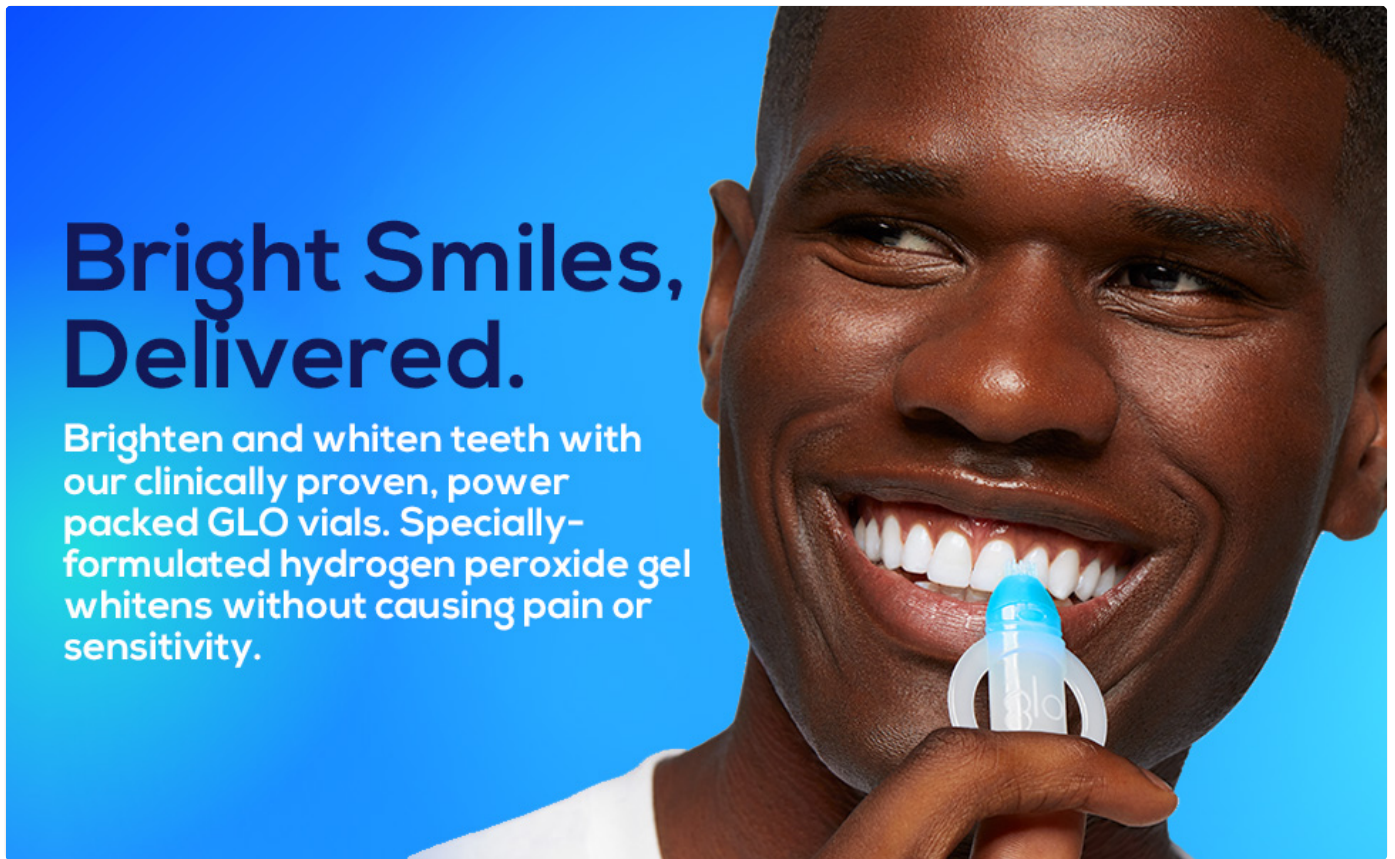


Image 4.2: Applying the GLO Whitening Gel.

Brush on to brighten & whiten



Image 4.3: GLO Vial features and application details.

5. MAINTENANCE AND STORAGE

- Each GLO Vial typically provides 3-4 applications.
- After opening, use the remaining gel in a vial within 7 days for optimal effectiveness.
- Store unused, sealed vials at room temperature. Refrigeration is not required.
- Ensure the cap is securely replaced on opened vials to maintain freshness.

6. TROUBLESHOOTING

If you experience any issues while using the GLO Science Teeth Whitening Gel, consider the following:

- **Increased Tooth Sensitivity:** While the formula is designed to be sensitivity-free, some individuals may experience temporary sensitivity. If sensitivity occurs, reduce the frequency of application or discontinue use and consult a dental professional.

- **Gum Irritation:** Ensure the gel is applied only to the tooth surface and avoid direct contact with gums. If irritation persists, discontinue use.
- **Insufficient Whitening Results:** Ensure consistent application as directed. Results may vary based on individual tooth enamel and the type/severity of stains. For accelerated results, consider using with a GLO whitening device.
- **Gel Consistency:** If the gel appears unusually watery or has changed consistency, discontinue use and contact customer support.

7. SPECIFICATIONS

Feature	Detail
Brand	GLO Science
Model Number	259120
Active Ingredients	Hydrogen Peroxide (6%)
Item Form	Gel
Number of Items	1 (Pack of 10 Gel Vials & 1 Lip Care)
Product Dimensions	8.5 x 3 x 3 inches; 4.66 ounces
UPC	857181002367
ASIN	B00YV0ODF8

8. WARRANTY AND SUPPORT

Specific warranty information for the GLO Science Teeth Whitening Gel Vials is not provided in this manual. For details regarding product warranty or to contact customer support, please refer to the official GLO Science website or contact their customer service directly.

You can visit the GLO Science Store for more information: [GLO Science Store](#)

