

[Manuals.plus](#) /

> [Safety 1st](#) /

> Safety 1st Comfort Check Pacifier Thermometer and Medicine Dispenser User Manual

Safety 1st TH0760300

Safety 1st Comfort Check Pacifier Thermometer and Medicine Dispenser User Manual

1. PRODUCT OVERVIEW

The Safety 1st Comfort Check Pacifier Thermometer and Medicine Dispenser is designed to assist parents in monitoring a baby's temperature and administering liquid medication. This dual-function device combines the familiarity of a pacifier with the utility of a thermometer and a medicine dispenser, aiming to provide comfort while addressing health needs.



Image 1: Components of the Safety 1st Comfort Check Pacifier Thermometer and Medicine Dispenser, showing the pacifier thermometer and the medicine dispenser cup.

2. IMPORTANT SAFETY INFORMATION

Please read all instructions and warnings before using this product. Keep this manual for future reference.

- **WARNING:** This product is not a toy. DO NOT allow children to play with it. When not in use, keep out of reach of children.
- **WARNING:** DO NOT tie the pacifier around a child's neck as it presents a strangulation hazard.
- This product is NOT intended for extended use as a regular pacifier. Use only for temperature checks or medicine administration.
- For indoor home use only.
- Made without BPA.
- Always supervise your child during use.
- Inspect the pacifier nipple before each use. Discard if any signs of damage or weakness are observed.

170109WG

Helping parents every day
Aide les parents tous les jours > www.safety1st.ca

Comfort Check Pacifier Thermometer (Model TH047)
A baby who is feeling sick can be extra fussy. A pacifier thermometer takes your baby's temperature in a soothing way, and allows you to easily keep a watch on fevers.

WARNING: This product is not a toy. **DO NOT** allow children to play with it. When not in use, keep out of reach of children.

FOR INDOOR HOME USE ONLY.
MADE WITHOUT BPA.

⚠ WARNING: DO NOT tie pacifier around child's neck as it presents a strangulation hazard. **NOT** intended for extended use as a pacifier.


Sucette-Thermomètre Confort Check (Modèle TH047)
Un bébé malade peut devenir très difficile. Une sucette-thermomètre prend la température de votre bébé d'une manière apaisante, et vous permet de surveiller facilement la fièvre.

MISE EN GARDE: Ce produit n'est pas un jouet. **NE PAS** laisser les enfants jouer avec. Quand il n'est pas utilisé, gardez hors de la portée des enfants.

POUR USAGE DOMESTIQUE INTÉRIEUR SEULEMENT.
FABRIQUÉ SANS PBA.

⚠ MISE EN GARDE: NE PAS

⚠ MISE EN GARDE. NE PAS ttacher la sucette autour du cou de l'enfant car cela représente un danger d'étranglement. **N'EST PAS** destiné à un usage prolongé comme sucette.

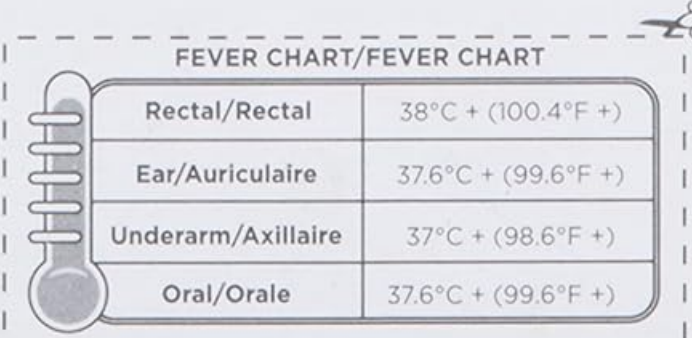

FEVER CHART/FEVER CHART

Rectal/Rectal	38°C + (100.4°F +)
Ear/Auriculaire	37.6°C + (99.6°F +)
Underarm/Axillaire	37°C + (98.6°F +)
Oral/Orale	37.6°C + (99.6°F +)

Use and care instructions enclosed.
Mode d'emploi et d'entretien inclus.



©2017 Dorel Juvenile.
(800) 544-1108 www.safety1st.ca
Made in CHINA. Fabriqué en CHINE.
Styles and colours may vary.
Les styles et couleurs peuvent varier.
Manufactured by/ Fabriqué par:
Dorel Juvenile Group,
A Division of/Une division de Dorel Industries,
2855 Argentia Road, Unit 4, Mississauga, ON L5N 8G6
Distributed by/ Distribué par:
Dorel Juvenile Group, Inc.,
2525 State St., Columbus, IN 47201-7494
Dorel Juvenile Canada,
2855 Argentia Road, Unit 4, Mississauga, ON L5N 8G6
4350-03459 09/16



Image 2: Product packaging displaying important safety warnings and usage instructions.

3. SETUP AND INITIAL USE

3.1 Battery Installation

The pacifier thermometer comes with one unknown type battery included and pre-installed. If replacement is ever needed, refer to the battery compartment on the device. Ensure correct polarity when inserting new batteries.

3.2 Initial Cleaning

Before first use, and after each subsequent use, clean the pacifier thermometer and medicine dispenser thoroughly. Refer to the 'Maintenance and Cleaning' section for detailed instructions.

4. OPERATING INSTRUCTIONS

4.1 Using the Pacifier Thermometer

1. Ensure the pacifier thermometer is clean and dry.
2. Press the power button to turn on the thermometer. A brief display test will show all segments, followed by the last measured temperature or 'Lo' if the ambient temperature is below the measurement range.
3. Place the pacifier in your baby's mouth, ensuring the nipple is fully inserted and the baby is sucking on it.
4. Keep the pacifier in the baby's mouth until the temperature reading stabilizes. This usually takes a few

minutes. The display will show the measured temperature.

5. Read the temperature displayed. *Note: The thermometer displays temperature in Celsius only and cannot be changed to Fahrenheit.*
6. Press the power button again to turn off the thermometer, or it will automatically shut off after a period of inactivity to conserve battery.

4.2 Using the Medicine Dispenser

1. Separate the medicine dispenser cup from the pacifier thermometer.
2. Fill the medicine dispenser cup with the prescribed dose of liquid medication. The cup has clear markings for accurate measurement (e.g., 1/4, 1/2, 3/4, 1 oz).
3. Carefully reattach the pacifier thermometer to the filled medicine dispenser cup, ensuring a secure fit.
4. Place the pacifier in your baby's mouth. As the baby sucks, the medicine will be dispensed through the nipple.
5. Ensure the baby consumes the full dose.
6. After administration, disassemble and clean both parts immediately.



Image 3: The Safety 1st Comfort Check Pacifier Thermometer and Medicine Dispenser, showing the assembled unit ready for use.

5. MAINTENANCE AND CLEANING

Proper cleaning and maintenance are essential for hygiene and product longevity.

- **Cleaning:** After each use, separate the pacifier thermometer from the medicine dispenser cup. Wash both parts thoroughly with warm, soapy water. Rinse well with clean water.
- **Sterilization:** The pacifier nipple can be sterilized by boiling for 5 minutes. Do NOT boil the thermometer unit as it contains electronic components.
- **Drying:** Allow all parts to air dry completely before reassembling or storing.
- **Storage:** Store the product in a clean, dry place, away from direct sunlight and extreme temperatures. Keep out of reach of children.

6. TROUBLESHOOTING

- **No Display/Device Not Turning On:** Check if the battery is correctly installed. If the issue persists, the battery may need replacement.

- **Inaccurate or Inconsistent Readings:** Ensure the pacifier is properly placed in the baby's mouth and held steady for the recommended duration. Ensure the device is clean. If readings remain inconsistent, the thermometer may be faulty.
- **Medicine Not Dispensing:** Check that the medicine dispenser cup is securely attached and that the nipple is not clogged. Ensure the baby is actively sucking.

7. SPECIFICATIONS

- **Model Number:** TH0760300
- **Product Dimensions:** Approximately 4.2 x 2.9 x 0.4 inches
- **Item Weight:** Approximately 0.8 ounces (22.68 g)
- **Battery:** 1 Unknown battery (included)
- **UPC:** 884392597597

7.1 Fever Chart

The following chart provides general guidelines for fever indications based on measurement site:

Measurement Site	Fever Indication
Rectal	38°C + (100.4°F +)
Ear (Auriculaire)	37.6°C + (99.6°F +)
Underarm (Axillaire)	37°C + (98.6°F +)
Oral	37.6°C + (99.6°F +)

Consult a healthcare professional for accurate diagnosis and treatment.

8. WARRANTY AND SUPPORT

For information regarding warranty, returns, or technical support, please refer to the product packaging or contact Safety 1st customer service directly. Contact details can typically be found on the manufacturer's official website or on the product's retail packaging.