

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

Manuals.plus /

› [Juvo](#) /

› Juvo Toilet Aid (Model SATA1) Instruction Manual

Juvo SATA1

Juvo Toilet Aid (Model SATA1) Instruction Manual

Personal Bottom Wiping Tool for Enhanced Independence

1. INTRODUCTION

The Juvo Toilet Aid is a personal hygiene device designed to assist individuals with limited range of motion in performing personal cleaning tasks. This tool provides extended reach and a secure grip for toilet tissue or pre-moistened wipes, promoting self-sufficiency and dignity for users.



Image 1.1: The Juvo Toilet Aid, illustrating its use for both front and back wiping positions.

It is particularly beneficial for seniors, pregnant individuals, those with obesity, or individuals recovering from surgery who may experience difficulty bending, reaching, or twisting.

2. PRODUCT FEATURES

- **18-inch Wand:** Provides extended reach for effective cleaning from either the front or back position.
- **E-Z Squeeze Trigger:** Securely locks toilet tissue or pre-moistened wipes between soft rubberized tips.
- **Dual Release Buttons:** Features two distinct buttons for easy and hygienic release of soiled tissue, one for front use and one for back use.
- **Soft Rubberized Tips:** Designed for comfort and secure tissue retention.
- **Hygienic Cover:** A plastic cover is included to protect the tips after use and during storage.
- **Integrated Caddy:** Comes with a caddy for convenient storage.

MAKES TOILETING EASIER FOR THOSE WITH LIMITED MOBILITY



Image 2.1: Key features of the Juvo Toilet Aid, including its length, grip, and release mechanisms.

3. SETUP AND PREPARATION

The Juvo Toilet Aid requires minimal setup. Upon unboxing, ensure all components are present: the toilet aid wand and the hygienic tip cover.

3.1 Preparing for Use

1. **Gather Tissue:** Fold or ball together 4-5 sheets of standard toilet tissue, or use one flushable wipe.
2. **Insert Tissue:** Place the prepared tissue or wipe between the soft rubberized tips of the toilet aid. Approximately half of the tissue should be sticking out above the tips.
3. **Secure Tissue:** Squeeze the E-Z Squeeze Trigger on the handle to lock the tissue firmly in place. Ensure the tissue is held securely and will not slip during use.

EASY TO USE, EASY TO CLEAN

- 1 Fold or ball together 4-5 sheets of tissue or 1 flushable wipe
- 2 Place between rubber tips. Half of the tissue or wipe should be between the tips/half should be sticking out, above the tips
- 3 E-Z Squeeze Trigger to lock tissue in place
- 4 For **front use** – press Push Button #2 to discharge tissue
- 5 For **back use** – press Push Button #1 to discharge tissue
- 6 Repeat as needed
- 7 Wash tips with warm soapy water
- 8 Store using plastic tip cover

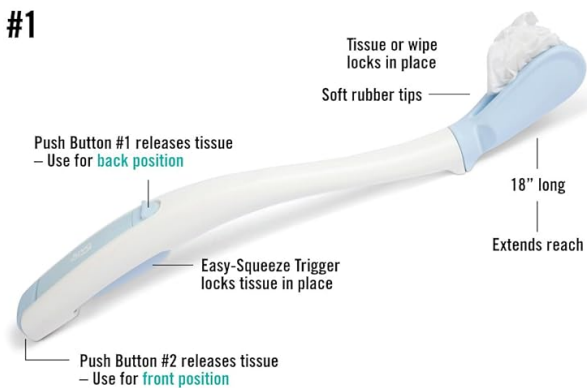


Image 3.1: Visual instructions for loading toilet tissue onto the Juvo Toilet Aid.

4. OPERATING INSTRUCTIONS

The Juvo Toilet Aid can be used from either a front or back approach, depending on individual comfort and mobility needs.

4.1 Using the Toilet Aid

1. **Positioning:** Once tissue is secured, position the toilet aid to reach the desired area.
2. **Cleaning:** Perform cleaning as needed. The extended reach of the wand minimizes strain and discomfort.
3. **Discharging Tissue (Front Use):** For front approach, press **Push Button #2** to release the soiled tissue into the toilet.
4. **Discharging Tissue (Back Use):** For back approach, press **Push Button #1** to release the soiled tissue into the toilet.
5. **Repeat:** Repeat steps 1-4 as necessary with fresh tissue until clean.

A DIGNIFIED SOLUTION... FOR THOSE THAT NEED A LITTLE HELP



Image 4.1: Demonstrations of the Juvo Toilet Aid in various usage scenarios.

5. MAINTENANCE AND CLEANING

Regular cleaning of your Juvo Toilet Aid is essential for hygiene and longevity.

1. **Clean Tips:** After each use, wash the rubberized tips with warm water and soap.
2. **Rinse Thoroughly:** Ensure all soap residue is rinsed off.
3. **Dry:** Allow the tips to air dry or gently pat dry with a clean cloth.
4. **Store:** Place the included plastic tip cover over the cleaned tips when the device is not in use. Store the toilet aid in its caddy or a convenient, clean location, such as a shelf or on top of the toilet tank.

WITHIN REACH

Store product on shelf,
in cabinet or on top of
toilet

Clean tissue holder tips
with warm soapy water

Use plastic cover



Image 5.1: Proper storage of the Juvo Toilet Aid with its protective cover.

6. TROUBLESHOOTING

- **Tissue/Wipe Slips Out:** Ensure the E-Z Squeeze Trigger is fully engaged. If using pre-moistened wipes, try folding them to create more thickness for a secure grip between the rubberized tips.
- **Soiled Tissue Does Not Release:** Verify that you are pressing the correct release button (#1 for back, #2 for front) firmly. If the mechanism feels stiff, ensure no debris is obstructing the release path.
- **Locking Mechanism Wear:** Over extended periods of use, the plastic catch in the locking mechanism may show signs of wear. If the device no longer securely holds tissue, contact Juvo customer support for assistance.

7. SPECIFICATIONS

Model Number	SATA1
Product Dimensions	1.68 x 3.52 x 18.18 inches
Item Weight	6.56 ounces

Manufacturer	Juvo
---------------------	------

8. SAFETY INFORMATION

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is intended as a personal hygiene aid. Always maintain good hygiene practices when using this device.

9. WARRANTY AND SUPPORT

For any questions, concerns, or support regarding your Juvo Toilet Aid, please contact Juvo Products directly. Refer to the product packaging or the official Juvo website for current contact information and warranty details.