

## McKesson 136-74806

# McKesson Convoluted Foam Mattress Topper Overlay Instruction Manual

Model: 136-74806

## INTRODUCTION

The McKesson Convoluted Foam Mattress Topper Overlay is designed to enhance comfort and provide pressure relief. Its innovative convoluted (egg-crate) design helps distribute body weight evenly, reducing pressure points and promoting air circulation. This manual provides instructions for the proper setup, usage, and maintenance of your mattress topper.

## IMPORTANT SAFETY INFORMATION

- **Not Made with Natural Rubber Latex:** This product is manufactured without natural rubber latex.
- **Ventilation:** Upon unpacking, allow the foam topper to air out in a well-ventilated area for a few hours to dissipate any new foam odor.
- **Flammability:** Keep the foam topper away from open flames and high heat sources.
- **Proper Placement:** Ensure the topper is placed securely on the mattress to prevent shifting during use.

## SETUP AND PLACEMENT

1. **Unpack:** Carefully remove the mattress topper from its packaging. The foam may be compressed; allow it to expand fully for several hours, or up to 24-48 hours, before use.
2. **Position:** Place the expanded foam topper directly on top of your existing mattress. Ensure the convoluted (egg-crate) side is facing upwards.
3. **Align:** Align the topper with the edges of your mattress for a secure and comfortable fit.

## Promotes air circulation after surgery



The McKesson Convoluted Foam Mattress Topper Overlay laid out on a chair, demonstrating its full size and the distinctive egg-crate pattern designed for air circulation and pressure distribution.

### USAGE

---

The McKesson Convoluted Foam Mattress Topper Overlay is designed to provide enhanced comfort and support. Its unique honeycomb structure with air perforations helps to:

- **Distribute Weight:** Evenly distributes body weight, relieving pressure from the body.
- **Improve Air Circulation:** Enhances airflow to help avoid pressure ulcers and maintain a comfortable temperature.
- **Reduce Pressure Sores:** Decreases the risk of pressure sores and tissue breakdown by minimizing prolonged pressure on specific areas.



The McKesson Convoluted Foam Mattress Topper Overlay is shown unrolled on a chair, illustrating its design that promotes air circulation, particularly beneficial for post-surgery recovery.

## MAINTENANCE AND CARE

---

Proper care will extend the life and effectiveness of your mattress topper:

- **Cleaning:** The care pads are washable. Follow the instruction: **Hand Wash Only**. Do not machine wash, bleach, iron, or dry clean.
- **Drying:** Air dry thoroughly in a well-ventilated area, away from direct sunlight or heat. Ensure the topper is completely dry before placing it back on the mattress.
- **Inflation:** This foam topper requires no regular inflation, making it easy to maintain.
- **Storage:** If storing, ensure the topper is clean and dry. It can be rolled up for compact storage.





The McKesson Convoluted Foam Mattress Topper Overlay, shown rolled for transport or storage. Its blue convoluted surface is visible.

## TROUBLESHOOTING

- **Odor:** A new foam product may have a slight odor upon unpacking. This is normal and will dissipate within a few hours to a few days in a well-ventilated space.
- **Flattening:** If the topper appears to flatten over time, ensure it is properly supported by your mattress. Rotating the topper periodically may help maintain its shape.
- **Discomfort:** If you experience discomfort, ensure the convoluted side is facing up. Allow sufficient time for your body to adjust to the new support surface.

## SPECIFICATIONS

|                           |                               |
|---------------------------|-------------------------------|
| Brand                     | McKesson                      |
| Model Name                | Mattress Topper Overlay       |
| Item Model Number         | 136-74806                     |
| Color                     | Blue                          |
| Size                      | 72 L X 33 W X 2 H Inch (1 ct) |
| Material                  | Foam                          |
| Special Feature           | Lightweight, Washable         |
| Product Care Instructions | Hand Wash Only                |
| Item Weight               | 2 pounds                      |
| Product Dimensions        | 35 x 5 x 5 inches (packaged)  |



72 inches long

33 inches wide

2 inches thick

The McKesson Convoluted Foam Mattress Topper Overlay is displayed with its key dimensions: 72 inches in length, 33 inches in width, and 2 inches in thickness, providing a clear understanding of its size.

## CONTACT AND SUPPORT





For further assistance or inquiries regarding your McKesson Convoluted Foam Mattress Topper Overlay, please refer to the contact information provided with your purchase or visit the official McKesson website.


## Related Documents - 136-74806



### [McKesson Disposable Vaginal Specula: Instructions for Use and Product Information](#)

Comprehensive guide to McKesson's disposable vaginal specula, covering usage instructions, cautions, indications, sterilization, and product specifications for medical professionals.


|  |  |
|--|--|
|    | <p><a href="#">McKesson Refrigerator/Freezer Thermometer TRIPLE DISPLAY SINGLE PROBE MCK821RFV</a></p> <p>Detailed information and operating instructions for the McKesson Refrigerator/Freezer Thermometer, Triple Display, Single Probe (Model MCK821RFV). Features include high accuracy, min/max memory, alarm settings, and NIST traceability.</p>          |
|   | <p><a href="#">INTEGRA MEDIACLAVE 10/30 Media Sterilizer Quick Start Guide</a></p> <p>Quick start guide for the INTEGRA MEDIACLAVE 10/30 Media Sterilizer, covering setup, operation modes (standard, chocolate agar, autoclave, water bath), dispensing, maintenance, and safety information. Includes manufacturer details and declaration of conformity.</p>  |
|   | <p><a href="#">McKesson Wireless Refrigerator/Freezer Thermometer MCK895WRFT - Triple Display, Single Probe</a></p> <p>User manual and specifications for the McKesson Wireless Refrigerator/Freezer Thermometer, model MCK895WRFT. Details features like triple display, single probe, large digit, wireless connectivity, alarm settings, and error codes.</p> |
|  | <p><a href="#">McKesson Refrigerator/Freezer Thermometer TRIPLE DISPLAY DUAL PROBE MCK821RFV2 User Manual</a></p> <p>User manual and specifications for the McKesson Refrigerator/Freezer Thermometer TRIPLE DISPLAY DUAL PROBE (Model MCK821RFV2). Features include high accuracy, min/max memory, low/high alarm settings, and NIST traceability.</p>          |

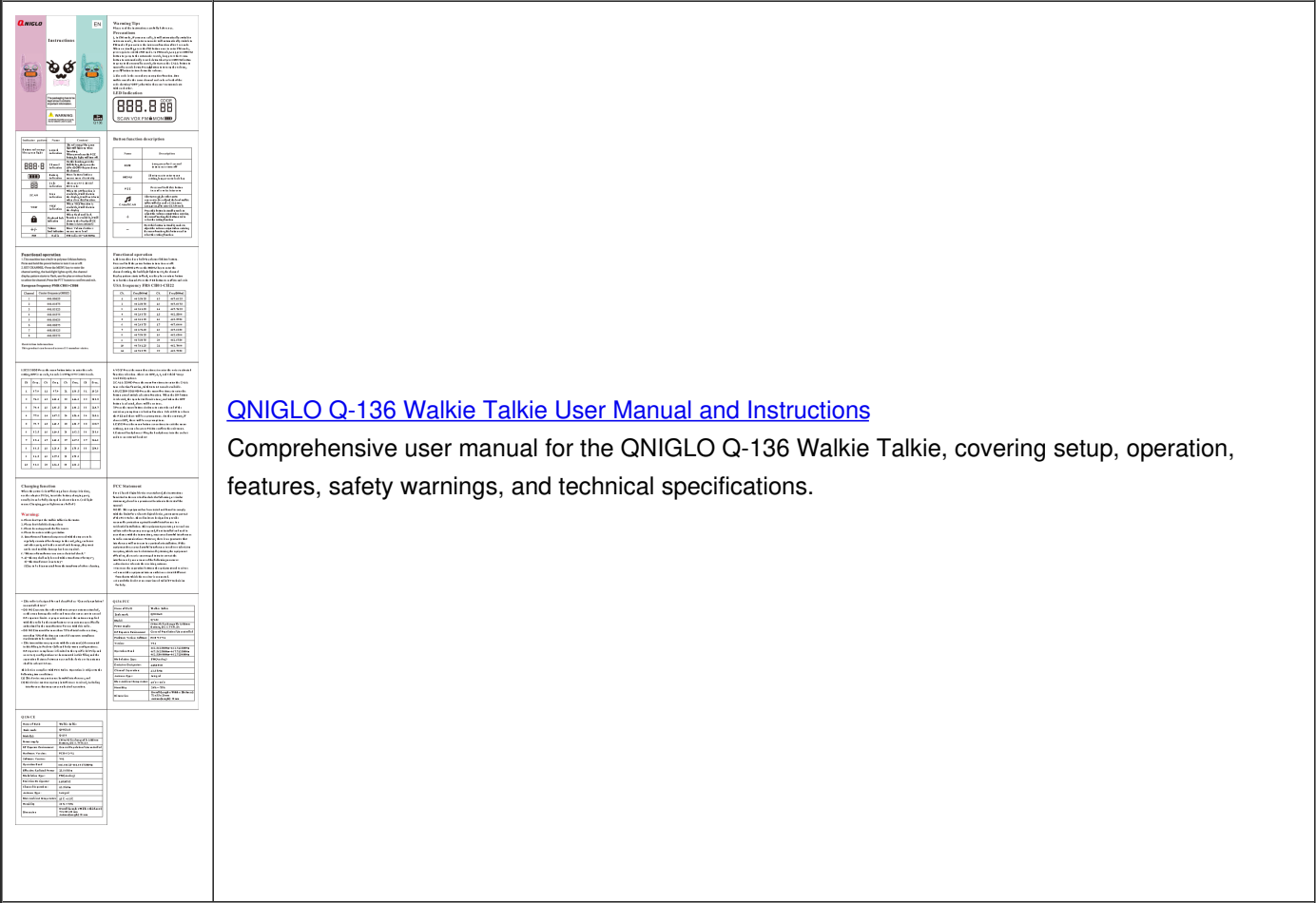

[QNIGLO Q-136 Walkie Talkie User Manual and Instructions](#)

Comprehensive user manual for the QNIGLO Q-136 Walkie Talkie, covering setup, operation, features, safety warnings, and technical specifications.

Documents - McKesson – 136-74806

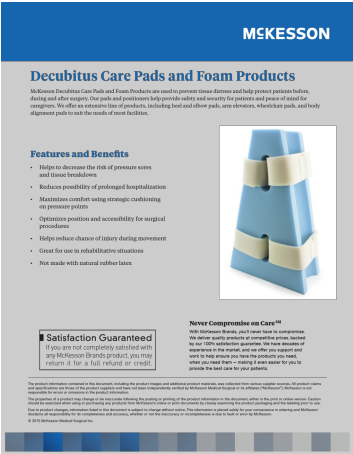

[\[pdf\]](#)

Foam Positioners tear sheet 2015 07 imgcdn mckesson CumulusWeb Click and learn

Decubitus Care Pads and Foam Products McKesson Decubitus Care Pads and Foam Products are used to pre ... l, 8 x 8 x 18 Abduction Pillow Medium, 15 x 22 x 7 Wheelchair Cushion, 18 x 16 x 3 MFR # **136-74806** 136-28656 136-31205 136-28520 136-64185 136-33080 136-21254 136-79631 136-66005 136-75751 ...

lang:en score:25 filesize: 578.27 K page\_count: 2 document date: 2015-08-19